



NOV 02, 2023

OPEN ACCESS



DOI:
dx.doi.org/10.17504/protocols.io.kxygx9opzg8j/v1

Protocol Citation: Sean M Cook, Vera A. Tang, Joanne Lannigan, Jennifer Jones, Joshua A Welsh 2023. Resource 4: rEV Serial Dilution. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.kxygx9opzg8j/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

Created: Jan 30, 2023

Last Modified: Nov 02, 2023

Resource 4: rEV Serial Dilution

Forked from [rEV Acquisition](#)

In 1 collection

Sean M


Cook¹, Vera A. Tang², Joanne Lannigan³, Jennifer Jones¹, Joshua A Welsh⁴

¹Laboratory of Pathology, Translational Nanobiology Section, Centre for Cancer Research, National Institute of Health, National Institutes of Health;

²University of Ottawa; ³Flow Cytometry Support Services;

⁴Translational Nanobiology Section, Laboratory of Pathology, Center for Cancer Research, National Cancer Institute, National Institutes of Health

Translational Nanobiology Section

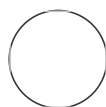

Jennifer Jones

DISCLAIMER






This protocol summarizes key steps for a specific type of method, which is one of a collection of methods and assays used for EV analysis in the NCI Translational Nanobiology Section at the time of submission of this protocol. Appropriate use of this protocol requires careful, cohesive integration with other methods for EV production, isolation, and characterization.

Funders
Acknowledgement:
 NIH
 Grant ID: ZIA BC011502
 NIH
 Grant ID: ZIA BC011503
 NIH
 Grant ID: 4UH3TR002881-03


ABSTRACT



Flow cytometry (FCM) is a common extracellular particles (EPs), including viruses and extracellular vesicles (EVs), characterization method. Frameworks such as MIFlowCyt-EV exist to provide reporting guidelines for metadata, controls, and data reporting. However, tools to optimize FCM for EP analysis in a systematic and quantitative way are lacking. Here, we demonstrate a cohesive set of methods and software tools that optimize FCM settings and facilitate cross-platform comparisons for EP studies. We introduce an automated small particle optimization (SPOT) pipeline to optimize FCM fluorescence and light scatter detector settings for EP analysis and leverage quantitative FCM (qFCM) as a tool to further enable FCM optimization of fluorophore panel selection, laser power, pulse statistics, and window extensions. Finally, we demonstrate the value of qFCM to facilitate standardized cross-platform comparisons, irrespective of instrument configuration, settings, and sensitivity in a cross-platform standardization study utilizing a commercially available EV reference material.

MATERIALS

- 
 Exosome standards fluorescent Millipore Sigma Catalog #SAE0193, referred to as 'rEVs'
- 
 DPBS Thermo Fisher Scientific Catalog #14190144, referred to as 'DPBS'
- 
 Low Protein Binding Collection Tubes (1.5 mL) Thermo Fisher Scientific Catalog #90410, referred to as 'low-binding Eppendorf tubes'
- 
 QbSure Cytek Catalog #B7-10005, referred to as 'QbSure'
- 
 100 nm polystyrene NIST bead Thermo Fisher Scientific Catalog #3100A, referred to as '100 nm NIST bead'

rEV Reconstitution


- 1
 Briefly centrifuge  100 x g, 4°C, 00:05:00 rEVs before opening.

5m
- 2
 Add  100 µL of  4 °C deionized water.

3 Pipet up and down to mix.

Note

Do not vortex

4 Aliquot  20 µL into low-binding Eppendorf tubes.

Equipment

Low Protein Binding Microcentrifuge Tubes

NAME

Microcentrifuge Tubes

TYPE

Thermo Scientific



BRAND

90410

SKU

<https://www.thermofisher.com/order/catalog/product/90410>

LINK




5 Store at  2-4 °C for up to one day (up to 2 weeks at  -80 °C).


rEV Dilutions


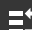
6 Label 10 low-binding Eppendorf tubes:

A	B
Label	Dilution
1	50
2	100
3	200
4	400
5	800

A	B
6	1600
7	3200
8	6400
9	12800
10	25600


7 Add  980 μL of DPBS (preferably  0.1 μm filtered) to the tube 1 from  go to step #6

8 Add  500 μL of DPBS (preferably 0.1 μm filtered) to the remaining tubes.

9 Add  20 μL of the stock reconstituted rEVs to the tube labeled 1 ( go to step #6), pipet up and down to mix.

Note

Do not vortex

10 Pipette  500 μL of the 1:50 dilution in tube 1 into the 1:100 dilution in tube 2, pipet up and down to mix.

Note

Do not vortex

11 Repeat the above step for each of the remaining dilution tubes until the serial dilution has been completed for all dilutions.

- 12 Keep on ice during acquisition.

Flow Cytometer Acquisition

- 13 Use the instrument settings that were established from the previous gain incrementation and scatter calibration resources. On the Aurora, set the window extension to 0. On the CytoFLEX, use the high acquisition mode. Collect for at least 1 minute at a low flow rate with a 30 second recording delay if using plate mode.

Note

To set the window extension to 0 on the Cytex Aurora, navigate to the 'Lasers' tab under the sample acquisition settings and set the window extension to 0.

To change the acquisition mode on the CytoFLEX S, click on the 'Advanced' menu on the top-of-the-screen ribbon. Next click on 'Event Rate Settings' menu, change the acquisition mode to 'High', and press 'Ok'.

- 13.1 Set your instrument to trigger on the most sensitive side scatter channel at the FCM_{PASS} optimal gain output as determined by the FluoSpheres scatter voltration. Set the threshold in the noise such that the event rate is ~1000 events/sec when running DPBS.

Protocol



NAME

Resource 1: Scatter Detector Setting Incrementation for FCM_{PASS}

CREATED BY

Jennifer Jones

PREVIEW

- 13.2 Set the FITC parameter to the FCM_{PASS} optimal gain output as determined by the QbSure fluorescent voltration.

Protocol



NAME

Resource 2: Fluorescence Detector Setting Incrementation for FCMPASS

CREATED BY



Jennifer Jones

PREVIEW

- 14** Collect a DPBS sample in a FACS tube for one minute on low.


Note

This is to ensure the instrument is clean on the side scatter and fluorescent channels before acquiring the rEVs.

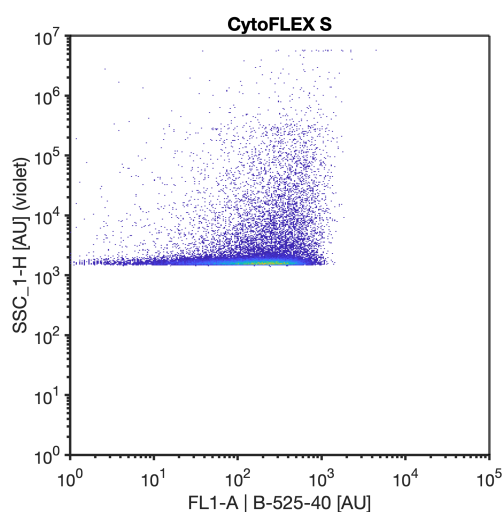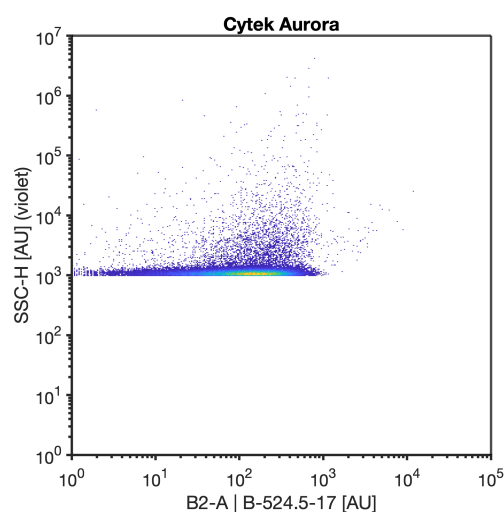

DPBS acquired on the Cytek Aurora and CytoFLEX S.

- 15** Collect all rEVs samples for one minute on low, starting with the highest dilution (least concentrated) and working down to lowest dilution (most concentrated). This will reduce the effect of carry over from high concentrations.

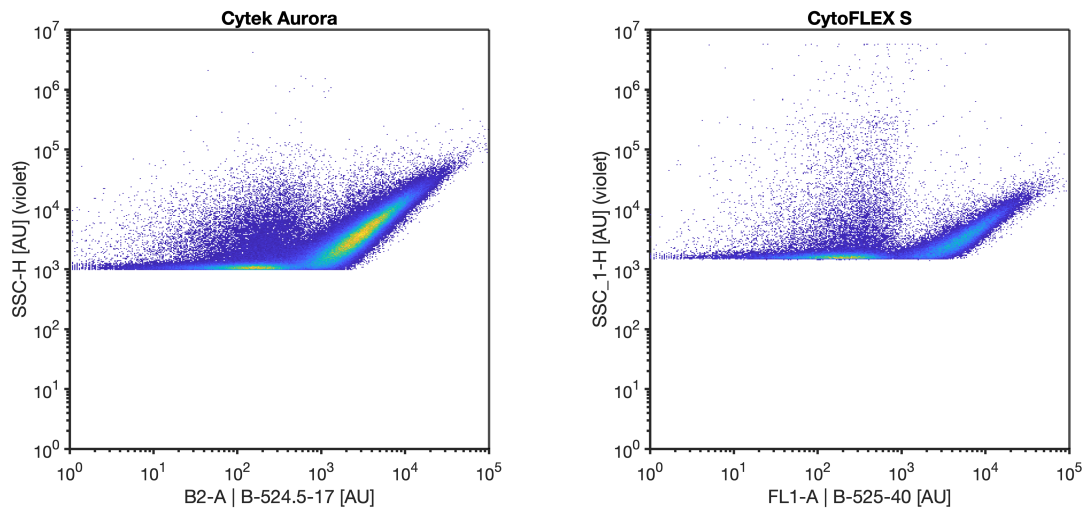


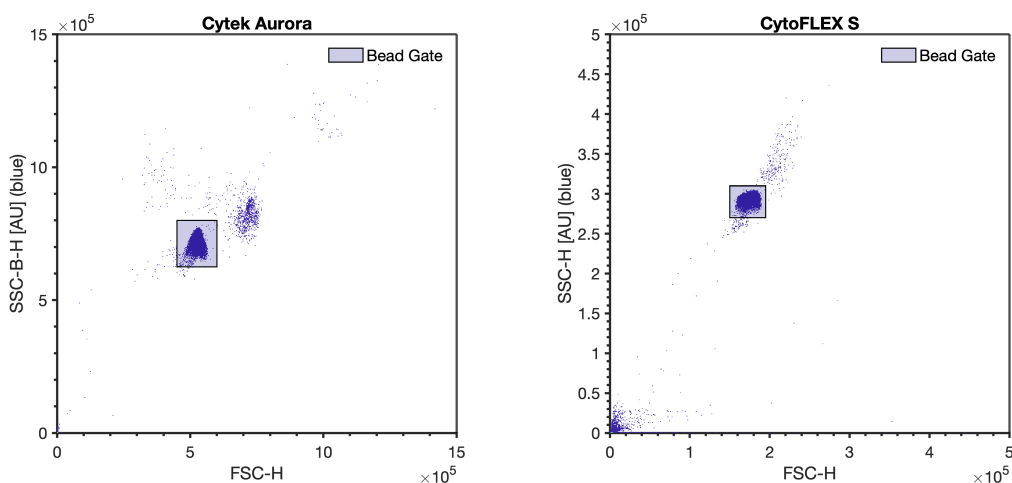
rEVs acquired on the Cytek Aurora and CytoFLEX S.

- 16 Run a DPBS sample in a FACS tube to make sure there are no carryover events; there should be no events in the B1 and B2 channel.

Note

If there is carry over clean instrument with bleach and flush with deionized water as appropriate.

- 17 Label as FACS tube as 'QbSure', and add  500 μL DPBS. Vortex the QbSure beads for 5 sec and add 3 drops QbSure beads into FACS tube. Run QbSure beads on low and collect 10,000 bead events.



QbSure beads gated on blue SSC-H and FSC-H on the Cytex Aurora and CytoFLEX S. The violet SSC settings remain unchanged. In order to limit .fcs file size, a FSC scatter trigger may be used until 10000 QbSure bead events are acquired.

Note

You may need to adjust the B-SSC and FSC detector settings for the bead acquisition. Do not adjust fluorescence detector settings or high sensitivity SSC detector settings. Assuming the QbSure beads were cross calibrated in the previous protocols, this allows for the calibration of the rEV fluorescent data into calibrated units.

- 18** Run the 100 nm polystyrene NIST beads at the same settings as the rEVs (diluted in DPBS to a concentration of 5E6 p/mL) until at least 10000 bead events have been acquired.

Note

Acquiring the 100 nm polystyrene NIST bead allows for the calibration of the rEV side scatter data into calibrated units if the collection half-angle of the instrument has been determined in previous protocols.

- 19** If calibration of data into standard units is desired, the 100 nm polystyrene NIST bead and QbSure beads can be used to calibrate data by following the FCM_{PASS} experiment calibration protocols.

Protocol



NAME

FCMPASS Protocol Collection

CREATED BY

Joshua A Welsh

PREVIEW