

MAY 08, 2023

# Library preparation (dsDNA single indexing, non-UDG, no split)

Forked from [Library preparation \(dsDNA double indexing, non-UDG, 2x split\)](#)

Marcel Keller<sup>1,2</sup>, Christiana L Scheib<sup>1,3</sup>

<sup>1</sup>Institute of Genomics, University of Tartu;

<sup>2</sup>Institute of Forensic Medicine, University of Bern;

<sup>3</sup>University of Cambridge



Marcel Keller

Institute of Forensic Medicine, University of Bern, Institut...

OPEN ACCESS

**DOI:**  
[dx.doi.org/10.17504/protocols.io.n92ldpjexl5b/v1](https://dx.doi.org/10.17504/protocols.io.n92ldpjexl5b/v1)

**Protocol Citation:** Marcel Keller, Christiana L Scheib 2023. Library preparation (dsDNA single indexing, non-UDG, no split). **protocols.io** <https://dx.doi.org/10.17504/protocols.io.n92ldpjexl5b/v1>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working  
We use this protocol and it's working

**Created:** Feb 26, 2023

**Last Modified:** May 08, 2023

**PROTOCOL integer ID:**  
77635

## ABSTRACT

Protocol for the preparation of single indexed double-stranded DNA libraries for Illumina sequencing, optimized for ultra-short ancient DNA molecules, modified from Meyer & Kircher (2010) Cold Spring Harb. Protoc. (doi: [10.1101/pdb.prot5448](https://doi.org/10.1101/pdb.prot5448)). This protocol does not include treatment with UDG (USER) to remove DNA damage in form of deaminated cytosines.

## GUIDELINES

Please read the general guidelines for working in the Ancient DNA protocol collection – University of Tartu, Institute of Genomics.

## MATERIALS

### Reagents:

A	B	C	D	E	F	G
Step	Reagents	Con c.	Unit	Manufactur er	Kit/full description	Product number
End repair	Water			Fisher chemicals	Water for HPLC 1L	10367171
End repair	End Repair Reaction Buffer	10	X	New England Biolabs	NEBNext End Repair Module	E6050L
End repair	End Repair Enzyme Mix	1	X	New England Biolabs	NEBNext End Repair Module	E6050L
Adapter Ligation	Quick Ligation Buffer	5.0	X	New England Biolabs	NEBNext Quick Ligation Module	E6056L
Adapter Ligation	Adapter Mix (2.5 µM)	2.5	µM		*Custom made from Eurofins adapters	n/a

**Keywords:** ancient DNA, aDNA, archeogenetics, archaeogenetics, paleogenetics, palaeogenetics, library preparation

A	B	C	D	E	F	G
Adapter Ligation	Quick T4 Ligase	1.0	U	New England Biolabs	NEBNext Quick Ligation Module	E6056L
Fill in	Water			Fisher chemicals	Water for HPLC 1L	10367171
Fill in	Thermopol Buffer 10X	10	X	New England Biolabs	Bst DNA Polymerase, Large Fragment	M0275
Fill in	dNTP 25mM	25	mM	Thermo Fisher		R1121
Fill in	Bst DNA polymerase	8	U/μl	New England Biolabs	Bst DNA Polymerase, Large Fragment	M0275
PCR	Water			Fisher chemicals	Water for HPLC 1L	10367171
PCR	10x PCR buffer	10	X	Eurogentec	HGS Taq Diamond set	TAQ-I011-1000+
PCR	MgCl <sub>2</sub>	25	mM	Eurogentec	HGS Taq Diamond set	TAQ-I011-1000+
PCR	BSA 20 mg/ml	20	mg/ml	Thermo Fisher		B14
PCR	dNTPs (10 mM)	40	mM	Thermo Fisher		R0191
PCR	HGS Taq Diamond	1	U	Eurogentec	HGS Taq Diamond set	TAQ-I011-1000+
PCR	Universal Primer	10	μM	New England Biolabs	NEBNext Multiplex Oligos for Illumina (Index Primers Set 1-4)	E7335L, E7500L, E7710L, E7730L
PCR	Indexing Primer	10	μM	New England Biolabs	NEBNext® Multiplex Oligos for Illumina® (Index Primers Set 1-4)	E7335L, E7500L, E7710L, E7730L
Purification	PB Buffer			Qiagen	MinElute PCR Purification Kit	19066
Purification	PE Buffer			Qiagen	MinElute PCR Purification Kit	19065
Purification	EB Buffer			Qiagen	MinElute PCR Purification Kit	28006

#### Equipment and consumables:

A	B
Number	Equipment and consumables
1	0.2 ml tube rack
2	1.5 ml tube rack
1	1.5 ml cool block
1	0.2 ml cool block

A	B
	10 µl filter tips
	20 µl filter tips
	100 µl filter tips
	200 µl filter tips
	1000 µl filter tips
[# of samples]×4 (tubes)	0.2 µl PCR strips (8 tubes)
[# of samples]×2+4	1.5 ml tubes
1	5 ml tube
[# of samples]×2	MinElute columns
1	50 ml Falcon (waste)

#### Lab equipment:

Laminar flow hood

Centrifuge (2/1.5 ml)

Heat block

Mini table centrifuge/vortexer

Thermocycler

10 µl pipette

20 µl pipette

100 µl pipette

200 µl pipette

1000 µl pipette

#### Other consumables:

DNA ExitusPlus

Paper towels

#### SAFETY WARNINGS



#### Reagents

*NaOCl (bleach) solution (6%)*

- H290 May be corrosive to metals.
- H314 Causes severe skin burns and eye damage.
- H411 Toxic to aquatic life with long lasting effects.
- EUH206 Warning! Do not use together with other products. May

release dangerous gases (chlorine). Remove from surface after recommended incubation time with water-soaked tissue.



#### *DNA ExitusPlus*

H319 Causes serious eye irritation.



*Guanidinium hydrochloride (GuHCl)* (in PB buffer of Qiagen MinElute kit)

- H302 Harmful if swallowed.
- H332 Harmful if inhaled.
- H315 Causes skin irritation.
- H319 Causes serious eye irritation.



#### *Ethanol*

- H225 Highly flammable liquid and vapor.
- H319 Causes serious eye irritation.



### **Equipment**

#### *UV radiation*

- UV radiation can damage eyes and can be carcinogenic in contact with skin. Do not look directly at unshielded UV radiation. Do not expose unprotected skin to UV radiation.
- UV emitters generate ozone during operation. Use only in ventilated rooms.



## BEFORE START INSTRUCTIONS

### Previous step:

This protocol follows the extract purification protocol.

### Following step:

This protocol ends with the PCR setup. Proceed with the PCR and library purification protocol.

### Equipment and consumables:

A	B
Numbers	Equipment and consumables
1	0.2 ml tube rack
2	1.5 ml tube rack
1	1.5 ml cool block
1	0.2 ml cool block
	10 µl filter tips
	20 µl filter tips
	100 µl filter tips
	200 µl filter tips
	1000 µl filter tips
[# of samples]×4 (tubes)	0.2 µl PCR strips (8 tubes)
[# of samples]×2+4	1.5 ml tubes
1	5 ml tube
[# of samples]×2	MinElute columns
1	50 ml Falcon (waste)

## Preparation

- 1 Turn the hood on full power and open the glass.
- 2 Spray hood and table bench surfaces with DNA Exitus, let sit a minute and wipe down with paper towels.
- 3 Wipe down outside surfaces of reagents/tips with DNA Exitus and place in the hood.
- 4 Label the following tubes:  
5×1.5 ml tubes: ER, AL, FI, EB-1, EB-2  
1×5 ml tube: PCR  
0.2 ml PCR strips: ER, AL, FI, PCR  
Label the 50 ml waste tube, PB tube and PE tube.
- 5 Aliquot EB buffer in tubes EB-1 and EB-2: each [# of samples]×30 µl plus 10%
- 6 Aliquot water for Master Mixes:  
ER: [# of samples]×12.50 µl plus 10%  
FI: [# of samples]×12.20 µl plus 10%  
PCR: [# of samples]×17.00 µl plus 10%
- 7 Aliquot PB buffer to a 50 ml tube: [# of samples]×1000 µl plus 10%
- 8 Prepare PE (wash) buffer by adding ethanol and aliquoting to a 50 ml tube: [# of samples]×1380 µl plus 10%

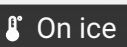


9 Aliquot PE buffer to a 50 ml tube: [# of samples]x1380 µl plus 10%

10 Take DNA extracts out of the freezer to thaw at room temperature. Change gloves.



## Blunt End Repair

11 Use 1.5 ml ER tube to set up the Blunt End Repair Master Mix .



A	B	C	D	E	F
REAGENT	STOCK	FINAL	UNIT	1× VOL (µl)	NOTE
Water [already added]				12.50	
End Repair Reaction Buffer	10	1.000	×	5.00	vortex
End Repair Enzyme Mix	1	0.050	×	2.50	on ice
Master Mix total				20.00	
Template DNA or water				30.00	
REACTION TOTAL				50.00	

Calculate +10% for all Master Mix components.

12 Add 20 µl Master Mix to each tube of the ER strip.



13 Vortex and spin down DNA extracts, add 30 µl of template DNA or water to each tube.



- 14 Mix carefully by resuspending, remove bubbles and spin down.



- 15 Incubate at  20 °C for  00:30:00 .

30m



#### Note

Use free time for preparing Master Mixes for the following steps. Master Mixes should be stored in the fridge.

## MinElute Purification 1

- 16 Take MinElute columns out of the fridge.

- 17 Turn on heat block  37 °C for elution.

- 18 Label columns and tubes with sample ID numbers.

- 19 Add 500 µl PB buffer (binding buffer) to MinElute column. You can use the same tip.




- 20 Add end-repair reaction mix to the PB buffer inside the MinElute columns and mix by resuspending.





21



Spin at  13000 rpm, 00:01:00, discard supernatant, change gloves.

1m


22



Add 690 µl PE buffer (wash buffer), change tip for every sample.

23




Spin at  13000 rpm, 00:01:00, discard supernatant, change gloves.

1m

24



Spin at  13000 rpm, 00:01:00 (dry spin).

1m

25

Put column in new tube, change gloves.

26



Elute in 30 µl EB buffer (elution buffer), change tip for every sample.

27




Incubate at  37 °C for  00:10:00.

10m

28



Spin at  13000 rpm, 00:02:00.

2m

29

Discard the silica column and close the lid.

30

**Note**

Potential stopping point, place end-repaired DNA in freezer until ready for next step.

## Adapter ligation

31

Use 1.5 ml AL tube to set up the Adapter Ligation Master Mix On ice .



A	B	C	D	E	F
REAGENT	STOCK	FINAL	UNIT	1× VOL (μl)	NOTE
Quick Ligation Buffer	5.0	1.00	×	10	vortex
Adapter Mix	2.5	0.25	μM	5	vortex
End Repair Enzyme Mix	1.0	0.10	U	5	on ice
Master Mix total				20.00	
Template DNA or water				30.00	
REACTION TOTAL				50.00	

Calculate +20% for all Master Mix components.

32

Add 20 μl of Master Mix to each tube of the AL strip.



33

Add 30 μl of end-repaired template DNA or water to each tube.



34 Incubate at  20 °C for  00:15:00 .

15m



## MinElute purification 2

35 Take MinElute columns out of the fridge.

36 Turn on heat block  37 °C for elution if not turned on already.

37 Label columns and tubes with sample ID numbers.

38 Add 500 µl PB buffer (binding buffer) to MinElute columns. You can use the same tip.



39 Add adapter ligation reaction mix to the PB buffer inside the MinElute columns and mix by resuspending.



40 Spin at  13000 rpm, 00:01:00 , discard supernatant, change gloves.




41 Add 690 µl PE buffer (wash buffer), change tip for every sample.




42



Spin at  13000 rpm, 00:01:00, discard supernatant, change gloves.

43



Spin at  13000 rpm, 00:01:00 (dry spin).

44

Put column in new tube, change gloves.

45



Elute in 30 µl EB buffer (elution buffer), change tip for every sample.

46




Incubate at  37 °C for  00:10:00.

10m

47



Spin at  13000 rpm, 00:02:00.

48

Discard the silica column, close the lid.

49



#### Note

Potential stopping point, place adapter-ligated DNA in freezer until ready for next step.

## Fill-in reaction

50

Use 1.5 ml FI tube to set up the Adapter Ligation Master Mix  On ice .



A	B	C	D	E	F
REAGENT	STOCK	FINAL	UNIT	1× VOL (μl)	NOTE
Water [already added]				12.2	
Thermopol Buffer	10	1.00	×	5.0	vortex
dNTPs	25	0.40	mM	0.8	vortex
Bst polymerase	8	0.32	U/μl	2.0	on ice
Master Mix total				20	
Template DNA or water				30	
REACTION TOTAL				50	

Calculate +20% for all Master Mix components.

51

Add 20 μl Master Mix to each tube of the FI strip.



52

Add 30 μl of adapter-ligated template DNA or water to each tube.



53

Incubate at  37 °C for  00:30:00 , then at  80 °C for  00:20:00 .

50m



### Note

Use free time for preparing the Master Mix for the following step, Master Mixes should be stored in the fridge

## Library amplification (PCR)

**54** Use 5 ml PCR tube to set up PCR Master Mix  On ice .



A	B	C	D	E	F
REAGENTS	STOCK	FINAL	UNIT	1× VOL (μl)	NOTE
Water [already added]				17	
10x PCR buffer	10	1.00	X	10	vortex
MgCl <sub>2</sub>	25	2.50	mM	10	vortex
BSA	20	1.00	mg/ml	5	vortex
dNTPs 10 mM	40	0.80	mM	2	vortex
HGS Taq Diamond	1	0.02	U	2	on ice
Universal primer 1.0	10	0.2	μM	2	vortex
Master Mix total				48	
Indexing primer	10	0.2	μM	2	vortex
Template DNA or water				50	
REACTION TOTAL				100	

Calculate +10% for all Master Mix components.

**55** Aliquot 48 μl of Master Mix to the tubes of the PCR strip.



**56** Vortex indexes slightly and spin them down with a table centrifuge.

**57** Add 2 μl of the indexing primer (10 μM) to the respective tube with Master Mix.



- 58** Add 50 µl of adapter fill-in reaction mix or water to the respective tube with Master Mix. Mix by resuspending with the pipet.



- 59** Check that the lids are tightly sealed and take the strips to the modern lab for PCR.

