



VERSION 4
NOV 10, 2023

Stereotaxic Injection by Nanoject Protocol V.4


Allen Institute for Brain Science¹

¹Allen Institute

Allen Institute for Brain Science

Tech. support

[Click here to message tech. support](#)

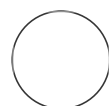

Allen Institute

ABSTRACT

This protocol describes the delivery of a neuronal tracer using the Nanoject II. The surgery uses a stereotaxic system to target specific brain coordinates in the mouse.

Note: Research reported in this publication was supported by the National Institute Of Mental Health of the National Institutes of Health under Award Number U19MH114830. The content is solely the responsibility of the authors and does not necessarily represent the official views of the National Institutes of Health.

ATTACHMENTS

[Standard of Care Protocol.xlsx](#)

[NSBWI-0022_all_procedures_setup_take_down.docx](#)

[AF0111_Stereotaxic_Injection_by_Nanoject.docx](#)

OPEN  ACCESS



DOI:
dx.doi.org/10.17504/protocols.io.bp2l6nr7kgqe/v4

Protocol Citation: Allen Institute for Brain Science 2023. Stereotaxic Injection by Nanoject Protocol.
protocols.io
<https://dx.doi.org/10.17504/protocols.io.bp2l6nr7kgqe/v4> Version created by [Allen Institute](#)

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Nov 10, 2023

Last Modified: Nov 10, 2023

PROTOCOL integer ID:
90789

Keywords: AF0111, surgery,
tracers, injection, nanoject