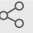




Sep 14, 2022

# 🌐 Gene expression analysis by quantitative Real-Time PCR (qPCR)

 In 1 collection[miquel.vila](#)<sup>1</sup><sup>1</sup>Vall d'Hebron Research Institute1 *Works for me* Share[dx.doi.org/10.17504/protocols.io.ewov1oyz2lr2/v1](https://dx.doi.org/10.17504/protocols.io.ewov1oyz2lr2/v1)

Nuriapenuelas

## ABSTRACT

Quantitative Real-Time PCR (qPCR)

DOI

[dx.doi.org/10.17504/protocols.io.ewov1oyz2lr2/v1](https://dx.doi.org/10.17504/protocols.io.ewov1oyz2lr2/v1)

## PROTOCOL CITATION

miquel.vila 2022. Gene expression analysis by quantitative Real-Time PCR (qPCR). **protocols.io**  
<https://protocols.io/view/gene-expression-analysis-by-quantitative-real-time-cgiztuf6>



## COLLECTIONS

**In vivo reduction of age-dependent neuromelanin accumulation mitigates features of Parkinson's disease**

## LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

## CREATED

Sep 13, 2022

## LAST MODIFIED

Sep 14, 2022

## PROTOCOL INTEGER ID

69945

qPCR plate reagents

## 1 Option 1: TaqMan qPCR

- 1.1 Prepare Taqman Master mix (ul per well):
  - 5ul: Taqman Gene Expression Master Mix (Applied Biosystems, # 4369016)
  - 0.5ul Taqman Gene Expression Assays
- 1.2 Add 5.5ul of Taqman Master mix to each well with an automatic pipette (Multipette E3, Eppendorf)
- 1.3 Add 4.5ul of cDNA per well (10ng cDNA) in technical triplicates to each well of LightCycler 480 Multiwell 384-plate (#4729749001, Roche Diagnostics)

## 2 Option 2: SYBR qPCR

- 2.1 Prepare SYBR Master mix (ul per well):
  - 5ul PowerUp SYBR Green Master Mix (#A25776, Applied Biosystem-ThermoFisher)
  - 1ul Primers (Forward and Reverse, diluted in Nuclease-free water at 5μM)
- 2.2 Add 6ul of SYBR Master mix to each well with an automatic pipette (Multipette E3, Eppendorf)
- 2.3 Add 4ul of cDNA per well (10ng cDNA) in technical triplicates to each well of LightCycler 480 Multiwell 384-plate (#4729749001, Roche Diagnostics)

Thermocycling

Perform PCR using the following cycling conditions in a LightCycler® 480 System (Roche):

- 3
  - 50°C for 2min
  - 95°C for 2min
  - 95°C for 15s + 60°C for 1min (x40cycles)

Data an

- 4 Using the LightCycler Software obtain the threshold cycles (CT) signals from each samples.
- 5 Perform  $\Delta\Delta CT$ -method to analyze the expression data using the endogenous control genes and the reference experimental group.