



JUL 26, 2023

OPEN ACCESS



Protocol Citation: Annika Fendler 2023. DNA extraction from whole blood and buffy coat using the QIAamp DNA mini kit. **protocols.io**

<https://protocols.io/view/dna-extraction-from-whole-blood-and-buffy-coat-usi-cwc8xazw>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working. Tested by Bettina Ergün

Created: Jun 27, 2023

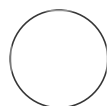
Last Modified: Jul 26, 2023

PROTOCOL integer ID:
84096

DNA extraction from whole blood and buffy coat using the QIAamp DNA mini kit

Annika Fendler¹

¹Charité - University Hospital



Annika Fendler

ABSTRACT

This protocol is for extraction of total RNA from whole blood, or buffy coat or PBMCs using the QIAamp DNA mini kit.

GUIDELINES

This kit is for 400 µL Whole blood or 200 µL Buffy coat

MATERIALS

QIAamp DNA Mini Kit Qiagen Catalog #51304

EtOH

PBS

1.5 and 2 ml LoBind tubes

BEFORE START INSTRUCTIONS

Prepare:

AW1: Add 25ml EtOH to new bottle

AW2: Add 30 ml EtOH to new bottle

RNAse I stocks: 100 mg/ml

Keywords: DNA, whole blood, buffy coat

Preparation of whole blood or buffy coat

- 1 Select whole blood or buffy coat from biobank list and retrieve from storage location. Ensure that sample is correct by scanning the tube barcode and matching it with the biobank list.





1.1

















	Sample ID	Barcode	Sample Type	Volume used	Comment	DNA concentration

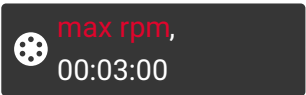
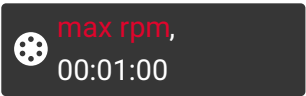


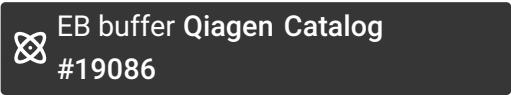
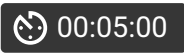

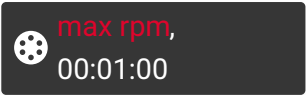

List of samples for DNA extraction

DNA Isolation

22m 30s

- 2 Transfer  200 µL buffy coat or  400 µL whole blood in 1.5 ml collection tubes
- 3 Add  20 µL  40 µL Proteinase K and mix by pipetting up or down

- 4 Add  4 μL  8 μL  RNase A Qiagen Catalog #19101 vortex for  00:00:15 15s
- 5 Add  200 μL  400 μL buffer AL, vortex for  00:00:15 15s
- 6 Incubate  00:10:00 @  56 °C 10m
- 7 Add  200 μL  400 μL ethanol (96–100%) to the sample, and mix again by pulse-vortexing for  00:00:15 . 15s
- 8 transfer to DNA Mini Spin Column
- 9  max x g, 00:01:00 1m
Discard collection tube
- 10 Add  500 μL AW1
- 11  max rpm, 00:01:00 1m
Discard flow-through
- 12 Add  500 μL AW2

- 13  Discard collection tube 3m
- 14 Place column in new 2 ml collection tube
- 15  Discard collection tube 1m
- 16 Place column in new 1.5 ml collection tube
- 17   
- 18 Incubate  @  5m
- 19  1m
- 20 QC: Qubit  and
Tape Station 