



JUN 13, 2023

OPEN  ACCESS

DOI:
dx.doi.org/10.17504/protocols.io.ewov1ob7plr2/v1

Protocol Citation: michela.deleidi, Bianca Marchetti, Federico Bertoli, Carmela Giachino 2023. Luminex-Based Multiplex Cytokine Analysis. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.ewov1ob7plr2/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Mar 31, 2023

Last Modified: Jun 13, 2023

PROTOCOL integer ID:
79821

Keywords: Luminex-Based
Multiplex Cytokine Analysis

Luminex-Based Multiplex Cytokine Analysis

Federico


michela.deleidi^{1,2}, Bianca Marchetti^{3,4}, Bertoli^{1,2},
Carmela
Giachino³

¹Mitochondria and Inflammation in Neurodegenerative Diseases, DZNE, Tübingen-Germany;

²Hertie Institute for Clinical Brain Research, University of Tübingen;

³Neuropharmacology Laboratory, Oasi Research Institute-IRCCS, Troina, Italy;

⁴Biomedical and Biotechnological Sciences, Pharmacology Section, University of Catania-Italy

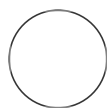

carmela.giachino

ABSTRACT

Multiplexed kit "Mouse High Sensitivity T Cell Magnetic Bead Panel" allows a Simultaneous analyze low levels of cytokine and chemokine biomarker with the High Sensitivity Bead-Based Multiplex Assays using the Luminex technology, in mouse serum, plasma and cell culture samples.

ATTACHMENTS

[k4dub48if.docx](#)



MATERIALS

Materials

- Mouse High Sensitivity T Cell Magnetic Bead Panel" (MHSTCMAG-70K, MILLIPLEX, EMD Millipore Corp., Billerica, MA, USA)
- FLEXMAP 3D® instrument (Luminex Corp., Austin, TX, USA)

Luminex-Based Multiplex Cytokine Analysis

1 Measure Cytokines from murine plasma samples, including IFN- γ , IL-1 β , IL-2, IL-4, IL-6, IL-10, IL-17A, and TNF- α simultaneously using the multiplexed kit "Mouse High Sensitivity T Cell Magnetic Bead Panel" (MHSTCMAG-70K, MILLIPLEX, EMD Millipore Corp., Billerica, MA, USA).

2 Perform the assay according to the manufacturer's instructions with an  Overnight incubation step (17 hours) at  4 °C .



3 Analyze samples as single values and use two controls with corresponding acceptance ranges, provided with the kit, for quality control purposes.



4 Perform the readout on a FLEXMAP 3D® instrument (Luminex Corp., Austin, TX, USA).

5 Acquire data using Luminex xPONENT® software (version 4.3) and determine mean fluorescence intensity (MFI).

6 Use the Bio-Plex Manager™ software (version 6.2) (Bio-Rad, Hercules, CA, USA) for back-calculation of unknown cytokine concentrations.

