




MAY 13, 2023

🌐 NCEM Drop Agar (TM-014)

sandra.crameri¹

¹CSIRO ACDP

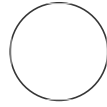

sandra.crameri

ABSTRACT

drop agar

OPEN  ACCESS

DOI:
dx.doi.org/10.17504/protocols.io.rm7vzb4n5vx1/v1

Protocol Citation: sandra.crameri 2023. NCEM Drop Agar (TM-014). **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.rm7vzb4n5vx1/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: May 13, 2023

Last Modified: May 13, 2023

PROTOCOL integer ID:
81826

HEADER

1 SAN:

SPEC No:

OPERATOR:

Agar preparation

- 2 Make [M] 2 Mass / % volume ⊗ Agarose Contributed by users Catalog #A5304 in TC Water. Microwave to dissolve, pour into 🧴 3 mL of the solution into 🧴 5 mL petri dishes and place in fridge to solidify.
- 3 Cut 1 x1 cm squares out of agar with scalpel, lift out onto filter paper

Adsorption & Staining

1m

- 4 Place 🧴 5 µL sample on surface of agar square. Immediately float carbon coated formvar filmed grid face down on the droplet and leave until the droplet diffuses into the agar.
- 5 Stain ⊗ nano-W Contributed by users Catalog ##2018-5ML ⌚ 00:01:00 1m
- 6 Drain & dry using filter paper

