

6



Oct 18, 2021

UW Virology Swift SNAPv2 protocol

Lasata Shrestha¹, Hong Xie¹, Shah A. Mohamed Bakhash¹, Robert J. Livingston¹, Meei-Li Huang¹, Alexander L. Greninger¹, Pavitra Roychoudhury¹

¹Department of Laboratory Medicine and Pathology, University of Washington, Seattle, WA

2	&°
protocol.	
Lasata	Shrestha

Viral whole genome sequencing (WGS) has been instrumental in outbreak investigations, deployment of public health interventions, development as well as evaluation of vaccines and therapeutics. While multiple methods are commercially available for WGS, multiplex amplicon method has proven to be faster, more efficient, scalable, and more cost-effective compared to other methods. Here, we describe the automation of a multiplex amplicon panel for WGS of Severe Acute Respiratory Syndrome Coronavirus-2 (SARS-CoV-2), the causative agent of the current COVID-19 pandemic. The SWIFT Biosciences' primer set amplifies 345 amplicons designed against the SARS-CoV-2 Wuhan-Hu-1 complete genome (NC_045512.2), in a single tube to cover the ~30 kb SARS-CoV-2 genome in less than three hours.

Lasata Shrestha, Hong Xie, Shah A. Mohamed Bakhash, Robert J. Livingston, Meei-Li Huang, Alexander L. Greninger, Pavitra Roychoudhury 2021. UW Virology Swift SNAPv2 protocol. **protocols.io**

https://protocols.io/view/uw-virology-swift-snapv2-protocol-byw4pxgw

whole genome sequencing, SARS-CoV-2, SWIFT Biosciences, PerkinElmer, automation, genome

protocol ,

Oct 10, 2021

Oct 18, 2021



- SuperScript IV Kit <u>Kit content</u>
- RNaseOUT™ Recombinant Ribonuclease Inhibitor Kit details
- SNAP SARS-CoV-2 Kit Kit content Page: 4
- SPRIselect beads https://www.beckman.com/reagents/genomic/cleanup-and-size-selection/size-selection
- 200 Proof Ethanol
- Molecular-grade water
- Qubit®, Nanodrop, or other similar input RNA quantification assay
- qPCR-, electrophoretic-, or fluorometric-based library quantification assay for Illumina® libraries
- Microcentrifuge
- Vortex
- Programmable thermocycler
- Aerosol-resistant tips and pipettes ranging from 1 to 1000 μL
- Pipette tips (e.g., 8-channel or 12-channel), 8-tube strips, an un-skirted 96 well plate, or plate puncher for pre-piercing the foil seal if using single-use UD indexing plates.
- Qubit dsDNA HS Assay Kit
 https://www.thermofisher.com/order/catalog/product/Q32851#/Q32851
- TapeStation DNA ScreenTape & Reagents https://www.agilent.com/en/product/automated-electrophoresis/tapestation-systems/tapestation-dna-screentape-reagents/dna-screentape-analysis-228260
- 80 uL Barrier Sterile 96 Rack Tips https://www.perkinelmer.com/product/80ul-art-sterile-96-rack-10-racks-111624
- 150 uL Barrier Sterile 96 Rack Tips https://www.perkinelmer.com/product/150ul-96-art-tip-box-10-racks-111426
- Polypropylene 384-well microplate, U-Bottom, case of 50
 https://www.perkinelmer.com/product/pp-microplate-384-35ul-u-bottom-50-6008890
- Polypropylene 96-well Microplate, Deep Well U-bottom, 2 mL
 https://www.perkinelmer.com/product/pp-microplate-96-2ml-v-bottom-25-6008880
- HARDSHELL PCR PLATE-96, BLUE/ 50 https://www.perkinelmer.com/product/hardshell-pcr-plate-96-blue-50-6008870
- StorPlate-96V, PP, 96 well, V-bottom, (V), 450μL, 200/box
 https://www.perkinelmer.com/product/storplate-96-v-450-l-200-6008299
- PP RESERVOIR, DW, V, 12 COL, 21mL /25 https://www.perkinelmer.com/product/pp-reservoir-dw-v-12-col-21ml-25-6008700

The critical steps are tagged. Steps 5-108 are optional (for troubleshooting).

Single-strand cDNA synthesis Reagent Set-up

1



Prepare reagents for required number of samples as per the manufacturer's instructions (including 10% overage) linked below.

https://assets.thermofisher.com/TFS-

Assets/LSG/manuals/SSIV_Reverse_Transcriptase_UG.pdf



2



Select the correct protocol on the liquid handler (Perkin Elmer Sciclone) and load the Sciclone as instructed.

This program will perform steps corresponding to steps 3 to 20 in the manual protocol. The code can be obtained from PerkinElmer.

The SuperScript RT IV protocol was automated and optimized by Perkin Elmer personnel following ThermoFisher's commercially available protocol linked below.

SuperScript RT IV



Fig. 01: Empty Sciclone deck.

Sciclone G3 NGSx iQ Workstation Automated Liquid Handling

Perkin Elmer CLS145321



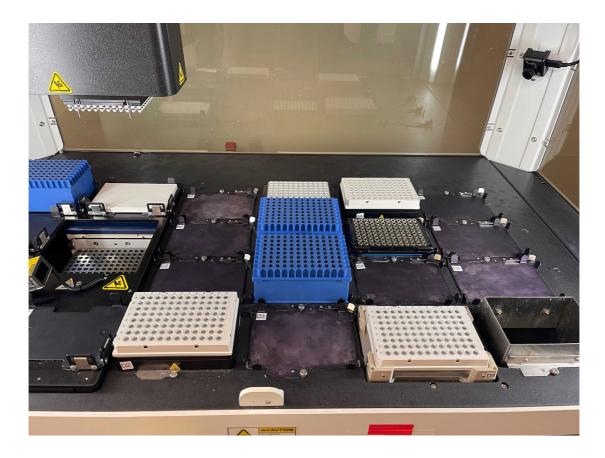


Fig. 02 Final (fully set-up) Sciclone desck set-up for sscDNA synthesis.

Sciclone movement: sscDNA synthesis

3



Moved Plate from D4 to B4 (magnet with no spacer)

- 4 Moved Plate from D2 to D4
- 5 Aspirate 11 μL of RNA from B4 and Dispense to plate on D4
- 6 Pipette mix and shake for © 00:01:00

1m

6.1 Plate on B4 moved to C0 concurrently

6.2 Plate on C0 discarded

- 7 Moved plate from D4 to onboard thermocycler
- 8 Move ODTC lid from A1 to thermocycler
- 9 Thermocycler door closes and runs © 00:07:40

7m 40s

- 10 Thermocycler door opens
- 11 ODTC lid moved from Thermocycler to A1
- 12 Moved plate from Thermocycler to D4
 - 12.1 Empty tip box discarded. Fresh tip box placed onto C3
- 13 Load new tips
- 14 Aspirate ¬7 μL of PCR master mix from A4 and dispense into plate on D4

- 15 Move plate from D4 to onboard thermocycler
- 16 Move ODTC lid from A1 to thermocycler
- 17 Thermocycler door closes and runs for © 00:53:04

53m 4s

- 18 Thermocycler door opens
- 19 Plate moved from thermocycler to D4
- 20 End product of **20 μL**

Swift SNAP v2 Reagent Set-up

21



Prepare reagents for required number of samples as per the manufacturer's instructions (including 10% overage) linked below.

Swift Normalase[™] Amplicon Panels (SNAP) SARS-CoV-2, Additional Genome Coverage, and SARS-CoV-2 S Gene Panel

Swift SNAP v2 Sciclone Set-up

22



Select the correct protocol on the liquid handler (Perkin Elmer Sciclone) and load the Sciclone as instructed.

The selected program will perform steps corresponding to steps 23 to 127. The code can be obtained from PerkinElmer.

The SWIFT SNAP v2 protocol was automated and optimized by Perkin Elmer personnel following SWIFT Bioscience's commercially available protocol linked below. SWIFT SNAP v2

Insert labeled photo of the deck without components (deck only)

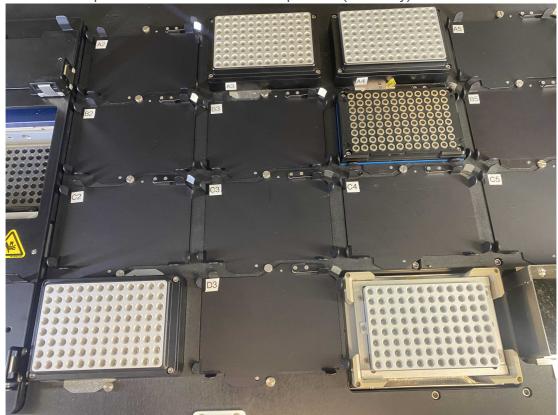


Fig. 03. Empty Sciclone deck.

Insert labeled photo of the deck with labeled components and volume information

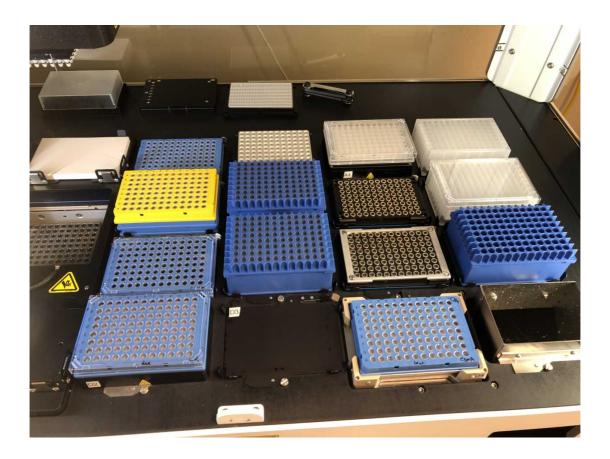


Fig. 04. Final (fully-loaded) deck set-up for SWIFT SNAP v2 protocol on the Sciclone.

Sciclone G3 NGSx iQ Workstation Automated Library Preparation

Perkin Elmer *******

G

Sciclone Movements: Swift SNAP v2 Multiplex PCR

1h 20m

2



Steps 23-127 include detailed Sciclone deck movements for troubleshooting.

24 Moved Lid from A4 to A2



25	Pipette multiplex master-mix (A4 to samples on D4)	
26	Dispose off tips	
27	Lid from A2 to A4	
28	Mix samples on D4	
29	Tips disposed	
30	Shake mix of D4 (Noted: 30uL)	
31	D4 to Thermocycler	
32	ODTC lid placed	
33	Thermo-cycler cover closed © 01:15:00	1h 15m
34	PCR door opened	

35	ODTC lid removed
----	------------------

- 36 The PCR plate was moved to D4
- 37 Moved plate from B2 to A3
- 38 Mix beads in B2, 10 X
- 39 Picked up 30 uL Beads from B2 to D4 **□30 µL Beads**
- 40 Mixed D4 wells, 10 X
- 41 Shake mix D4 for several minutes
- 42 Plate moved from D4 to B4
- 43 Allow beads to sediment © 00:05:00
- 44 Supernatant discarded
- 45 Tips discarded

5m

- 49 Dispensed to D2
- 50 Supernatant removed to Wash plate
- 51 Tips disposed
- 52 Fresh tips picked
- 53 Alcohol (150uL) dispensed to plate B4 \blacksquare 150 μ L
- 54 Removed supernatant from B4
- 55 Tips discarded

- 56 Fresh tips picked
- 57 Alcohol transferred to B4 **150 μL 80% EtOH**
- 58 Supernatant collected and disposed
- 59 Tips discarded
- 60 Alcohol lid replaced back
- 61 Sample plate moved to D4
- 62 New tips
- 63 TE buffer lid
- 64 TE buffer added **17.4 μL NEED TO CHECK**
- 65 Mixed
- 66 Tips discarded

- 67 Plate D4 shake mixed
- 68 Lid from A4 lifted
- 69 Transfer master mix A4 to A3 **28.9 μL Indexing PCR Reaction Mix**
- 70 Tips disposed

Sciclone Movements: SWIFT SNAP v2 Indexing PCR 1h 20m

71 🛠

Pipette Index from D2 to A3 (3.7 uL of indexes)

- 72 Pipette master mix plus indexes from A3 to D4
- 73 Shake mixed D4
- 74 Lid of A2 to A4
- 75 Fresh tips
- 76 D4 mixed by pipetting

77	Tips	disposa
----	------	---------

- 78 D4 shake mixed
- 79 D4 plate moved to thermocycler
- 80 21 min PCR run © **00:21:00**
- 81 Incubation 26 min © **00:26:00**
- 82 Thermocycler door open
- 83 OPTC lid removed
- 84 Sample plate moved to D4
- 85 Beads plate moved from B2 to D3
- 86 Pipetted PEG NaCl from plate on B2 to plate on D4 **32.5 μL PEG NaCl (ratio: 0.65)**
- 87 Mixing in D4 by pipetting (10 times)

21m

26m

88 Shake mixed D4, for 4 min 35 sec © **00:04:35**

4m 35s

89 Plate moved from D4 to C4 (on the magnet)

90 Incubation at room temperature © 00:05:00

5m

91 Supernatant from C4 discarded

92 Tips discarded

93 Alcohol lid removed

94 New tips picked

95 Picked up alcohol **150 μL 80% EtOH**

96 Transferred alcohol to C4

97 1 min wait

- 98 Supernatant discarded
- 99 Tips dropped off
- 100 New tips picked up
- 101 Picked up alcohol **150 μL 80% EtOH**
- 102 Dispensed to C4
- 103 1 min wait
- 104 Supernatant discarded
- 105 Tips discarded
- 106 Moved Tip container from C3 to D0
- 107 New tip container picked up from A0
- 108 New tips picked up

109	Alcohol	picked (up 🔲 1	50 μ	L 80%	EtOH
-----	---------	----------	--------	------	-------	-------------

- **111** 1 min wait
- 112 Supernatant from C4 pipetted up and discarded
- 113 Tips discarded
- 114 Alcohol cover replaced
- 115 2 min wait (incubate at room temperature until residual alcohol evaporates)
- 116 Moved Plate from C4 to D4
- 117 Picked up TE
- 118 Transfer TE to D4 **■20 µL TE**

119 Mixing 120 Tips discarded 121 Mix by pipetting and shaking in plate D4 122 Moved plate D4 to B4 123 Moved empty Plate A2 to D4 124 Tips picked up 125 Transfer from B4 to D4 **■20 µL Eluate** 126 Tips discarded 127 Plate lid from A4 to D4

Library Quality Control I

128 \wedge

Following the manufacturer's instructions for library quality control, quantify (Qubit or other fluorometric instruments) and determine the size (TapeStation or other electrophoretic

instruments) of the library.

4200 TapeStation System

Electrophoresis tool for DNA and RNA sample quality control.

TapeStation Instruments G2991AA



Invitrogen™ Qubit™ 3 Fluorometer Accurately measures DNA, RNA, and protein using the highly sensitive fluorescence-based Qubit quantitation assays

Invitrogen™ Q33216

Q33216

 \ominus

Normalase

129



Following manufacturer's instructions as linked below, normalize the libraries using SWIFT's proprietary enzymatic normalization.

Normalase Pages: 9-12


Pooled Library Quality Control II

130



Determine the library concentration for final library quality control before loading on the sequencer.

Quality Control Page: 12

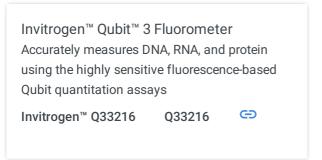

Sequencing

131



Following manufacturer's recommendation, load the pooled library for sequencing on the selected/preferred sequencer.

Sequencing Page: 13

NextSeq 500 System
Sequencer
Illumina *********

