



DEC 19, 2022

WORKS FOR ME

1

UPitt TriState SenNet TMC Human lung single cell suspension (test run)

COMMENTS 0

DOI

dx.doi.org/10.17504/protocols.io.eq2ly768plx9/v1

Marta Bueno¹, Lanping Guo¹, koenigshoffm¹,
Oliver Eickelberg¹

¹Division of Pulmonary and Critical Care Medicine.
Department of Medicine. School of Medicine

TriState SenNet

Cellular Senescence Network (SenNet) Method
Development Community



Marta Bueno

[Division of Pulmonary and Critical Care Medicine. Department...](#)

ABSTRACT

Human lung tissue dissociation protocol to produce single cell suspension that can be used for different purposes. Optimal incubation times to obtain cell suspension with great cell numbers and high cell viability from non-disease lungs.

DOI

dx.doi.org/10.17504/protocols.io.eq2ly768plx9/v1

PROTOCOL CITATION

Marta Bueno, Lanping Guo, koenigshoffm, Oliver Eickelberg 2022. UPitt TriState SenNet TMC Human lung single cell suspension (test run). **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.eq2ly768plx9/v1>

FUNDERS ACKNOWLEDGEMENT

TriState SenNET (Lung and Heart) Tissue Map and Atlas consortium – NIA
Grant ID: U54AG075931

LICENSE

_____ This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Dec 14, 2022

LAST MODIFIED


Dec 19, 2022

PROTOCOL INTEGER ID

73971

Before you start

1 Stock solutions

Can be prepared beforehand & stored  -20 °C :



1.1 Antibiotic-Antimycotic (100X) **Thermo Fisher Scientific Catalog #15240062** :

Aliquot in  5 mL in sterile conditions.

Store at  -20 °C



1.2 DNase I from bovine pancreas (Roche) **Sigma Aldrich Catalog #11284932001** 0 mg :

Stock of  100 mg in  10 mL of sterile ddH₂O


Aliquot in  300 µL in sterile conditions ( 3 mg per aliquot)

Store at  -20 °C

1.3 Elastase (3U/mg) 100mg **Worthington Biochemical Corporation Catalog #LS002292** :

Stock of  100 mg in  1 mL of sterile ddH₂O


The enzyme should be 0.22 microns filtered after reconstitution and prior to use

Aliquot in  500 µL in sterile conditions (150U per aliquot)

Store at  -20 °C

Collagenase IV (160U/mg) 1g **Worthington Biochemical Corporation Catalog #LS004188** :


Stock of  1 g in  10 mL of sterile ddH₂O

Aliquot in  500 µL in sterile conditions (8000U per aliquot)

Store at  -20 °C

1.4 Liberase TL **Roche Catalog #05 401 020 001** :

Stock of  5 mg in  2 mL of sterile ddH₂O

Aliquot in  500 µL in sterile conditions (1.25mg per aliquot)

Store at  -20 °C

2 Same-day preparation:

-  DMEM high glucose HEPES **Thermo Scientific Catalog #12-430-062**
-  RPMI-1640 **Thermo Fisher Scientific Catalog #22400089**
-  FBS **Invitrogen - Thermo Fisher**
-  1X PBS (Phosphate-buffered saline)
-  Red Blood Cell Lysis Buffer (Roche) **Sigma Aldrich Catalog #11814389001**

2.1 Base media
500 mL DMEM + 5 mL anti/anti (100x)

2.2 Digestion media (10 mL per digestion -> scale accordingly)

I. Modified base media with Final DNaseI 100 µL per 10 mL (final 0.1mg/ml)

II. Then add corresponding proteases:

II.a. Elastase / Collagenase (EC tubes)

Final elastase 250 µL per 10 mL (final 7.5U/ml) +

Final collagenase IV 250 µL per 10 mL (final 400U/ml)

II.b Liberase TL (TL tubes)

Final Liberase TL 1 mL per 10 mL (final 0.25mg/ml)

2h

Tissue dissociation

3 Cut the lung tissue with a scalpel into approximately 3 cm³ pieces (size of a thumb) upon receipt.

4 Place lung tissue into a 100mm petri dish. Place a ruler/control object next to it and take pictures for archives.

5 Transfer them to the lid of a fresh 100mm petri dish.


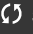





6 Each piece should be minced with scissors.









Locations: A-> Pleura; B-> bronchovascular bundle; C-> Parenchyma

(Label tubes accordingly)

A	B
EC-A	TL-A
EC-B	TL-B
EC-C	TL-C

Table 1. Samples generated.

- 7 Transfer to a 50 ml conical Falcon tube with  10 mL added digestion media. Use one tube per piece of tissue.
- 8 Repeat the procedure with the remaining pieces.
- 9 Incubate **30min at 37C**, shaking 500rpm.  500 rpm, 37°C, 00:30:00 30m
- 10 Add  5 mL of FBS to the cell suspension to stop the digestion (final volume ~15ml)
- 11 Filter through 100 µm cell strainer using the syringe plunger to press and smash the tissue, washing the tube and rinse the strainer with  5 mL of plain RPMI.
- 12 Filter through 70 µm cell strainer using the syringe plunger to press and smash the tissue, washing tube and rinse the strainer with  5 mL of plain RPMI.
- 13 Centrifuge. 300g, 15°C for 10 minutes.  300 x g, 15°C, 00:10:00 10m
- 14 Remove and discard the supernatant.
- 15 Resuspend the cell pellets with  10 mL of plain RPMI and transfer to a fresh 15 ml Falcon conical tubes.

- 15.1 Count cells here. (This count can be skipped -not informative because of RBC numbers)
- 15.2 To count: Take a 20 μ L aliquot to determine total cell number, and mix with 20 μ L of VitaStain. Run AO/PI program in Nexcelom cell counter.
- 16 Centrifuge. 300g, 15°C for 10 minutes. 10m
 300 x g, 15°C, 00:10:00
- 17 Remove supernatant and resuspend the cell pellets with  5 mL of Red Blood Cell Lysis Buffer. Mix tubes by inversion and incubate for  00:01:00 at  Room temperature . 1m
(Note: It can be done by one after another)
- 18 Fill up the tubes with  9 mL of plain RPMI.
- 19 Centrifuge. 300g, 15°C for 10 minutes. 10m
 300 x g, 15°C, 00:10:00
- 20 Remove and discard supernatant. Resuspend cell pellets with  10 mL of RPMI buffer kept on ice at all times.
- 20.1 Count cells here: Take a 20 μ L aliquot to determine total cell number, and mix with 20 μ L of VitaStain. Run AO/PI program in Nexcelom cell counter.
- 21 Centrifuge. 300g, 15°C for 10 minutes. 10m
 300 x g, 15°C, 00:10:00
- 22 Resuspend in the corresponding media for the next step

22.1 For scRNA-seq:

- Resuspend in 500µl of PBS+2%FBS in the 15ml tube and then take 200µl of that suspension to a new Eppendorf tube.
- Pellet the cells by a centrifugation pulse, remove the supernatant and then resuspend to prepare ~4Million cells in 200µl of PBS+2%FBS (minimum 2M in 100µl)

22.2 For storage:

- Resuspend in 500µl of PBS
- Pellet the cells by a centrifugation pulse, remove the supernatant and stored at -80C