



SEP 10, 2023

OPEN  ACCESS

Protocol Citation: leroi.bwh
2023. Basic Cell Culture
Maintenance: Splitting Cells.
protocols.io
<https://protocols.io/view/basic-cell-culture-maintenance-splitting-cells-czsmx6c6>

License: This is an open
access protocol distributed
under the terms of
the [Creative Commons
Attribution License](#), which
permits unrestricted use,
distribution, and reproduction
in any medium, provided the
original author and source
are credited

Protocol status: Working
We use this protocol and it's
working

Created: Sep 09, 2023

Last Modified: Sep 10,
2023

PROTOCOL integer ID:
87597


Basic Cell Culture Maintenance: Splitting Cells

 Forked from [Basic Cell Culture Maintenance: Splitting Cells](#)

leroi.bwh¹

¹Technion

Ellison Lab

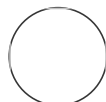


leroi.bwh


ABSTRACT


Basic protocol to split THP1

MATERIALS

MATERIALS


 HyClone Classical Liquid Media Dulbeccos Modified Eagles Medium (DMEM) Fisher Scientific Catalog #SH3024301

 Gibco™ (Phosphate Buffered Saline) Solution, pH 7.4 (PBS) Fisher Scientific Catalog # 10010-049

 Trypsin-EDTA (0.25%), phenol red Thermofisher Catalog #25200-056

Keywords: cell culture, hek293, cell, maintenance, split, cell splitting, splitting, maintaining cell lines, cell lines, maintain, cells

SAFETY WARNINGS

-  Human Embryonic Kidney (HEK293) cells are **biosafety level 2 (BSL-2)** and should be handled according to the CDC's [Biosafety in Microbiological and Biomedical Laboratories \(BMBL\)](#) guidelines. They are considered BSL-2 not because they are inherently hazardous or infectious, but because of their potential to be infected with pathogens and in turn infect their handlers. Due to the impossibility to regularly screen this cell line for every human pathogen, **HEK293 cells should always be handled as potentially infectious**. Other BSL-2 cell lines include those positive for *Legionella pneumophila*, HIV, and other disease-causing pathogens in humans.
- Dispose of ALL waste that comes into contact with cells such as pipettes, gloves, and materials as biohazardous waste.
- Bleach all direct cell waste thoroughly. In our lab, our vacuum line tube empties in to a sealed waste jug with bleach already added to the bottom of it, making up at least 10% of the total volume. This way, aspirated media and cells immediately come into contact with the bleach. Before disposing of glass pipettes, we aspirate a small amount of 10% bleach through to clean both the pipette and tubing, then dispose of the pipettes as biohazardous sharps.


BEFORE START INSTRUCTIONS

Make complete DMEM:

<i>Reagent</i>	<i>Volume</i>
DMEM	432.5 mL
FBS	50 mL
Pen/Strep	5 mL
HEPES (1M, pH 7.4)	12.5 mL

Preparation

- 1 Confirm that cells are at least 80% confluent by microscopy.

- 2 Warm complete RPMI in 37°C water bath.
- 3 UV light for 30 minutes then spray down the biosafety cabinet with 70% ethanol and use as a secondary decontaminant.
- 4 Prepare  20 µL of Trypan blue into a 200 µl Eppendorf.

Remove Media

- 5 Aspirate the media from the flask using a sterile autoclaved glass pipette. **Do not touch the cells with the pipette.**



Note

To avoid touching cells, is best to tilt the flask and gently remove media from a corner.

Transfer

- 6 Transfer ALL contents/cells to a  50 mL falcon tube.


Spin

- 7 Spin down 130 rcf for 7 minutes.  00:07:00  130 x g, 37°C 7m
- 8 While spinning, clean surfaces with EtOH and label new flasks, noting the +1 passage number and dilution.


Remove Media

- 9 Aspirate media from falcon tubes with cells; make sure to not disturb the pellet.


Resuspend Pellet

- 10 Add  5 mL media to pellet and pipette *violently* up and down.


Count the cells

- 11 Pipette  20 μL and add them to the Trypan blue


- 12 Count the cells using a Neubauer counting chamber

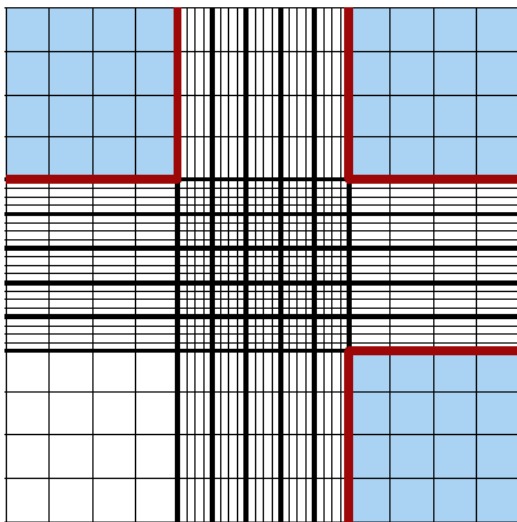

$$N = \frac{n_{\text{cells}}}{n_{\text{squares}}} \cdot D \cdot V_{\text{tot}} \cdot 10000$$

- 13 Calculate the required amount of medium to keep a concentration of  100 kC/ml

Prepare New Flask

- 14 Add the required amount of medium.
- 15 Distribute in flasks, no more than  50 mL per flask.

Incubate

- 16 Gently shuffle, ensure even dispersal, and return the fresh flask to incubator.