



Jul 24, 2020

# COVID-19 ARTIC v3 Illumina library construction and sequencing protocol - high throughput 384 format

Forked from [COVID-19 ARTIC v3 Illumina library construction and sequencing protocol](#)

DNA Pipelines R&D<sup>1</sup>, Benjamin Farr<sup>1</sup>, Diana Rajan<sup>1</sup>, Emma Betteridge<sup>1</sup>, Lesley Shirley<sup>1</sup>, Michael Quail<sup>1</sup>, Naomi Park<sup>1</sup>, Nicholas Redshaw<sup>1</sup>, Iraad F Bronner<sup>1</sup>, Louise Aigrain<sup>1</sup>, Scott Goodwin<sup>1</sup>, Scott Thurston<sup>1</sup>, Stefanie Lensing<sup>1</sup>, Carol Scott<sup>1</sup>, Nicholas Salmon<sup>1</sup>, Charlotte Beaver<sup>1</sup>, Rachel Nelson<sup>1</sup>, Alex Alderton<sup>1</sup>, Ian Johnston<sup>1</sup>

<sup>1</sup>Wellcome Sanger Institute

1

Works for me

[dx.doi.org/10.17504/protocols.io.biefkbnn](https://dx.doi.org/10.17504/protocols.io.biefkbnn)

Coronavirus Method Development Community

Diana Rajan

## ABSTRACT

This SOP describes the procedure for generating cDNA from SARS-CoV-2 viral nucleic acid extracts and subsequently producing 400nt amplicons tiling the viral genome using V3 nCov-2019 primers (ARTIC). This is followed by library construction, equivolume pooling of samples and quantitation, prior to sequencing on the Illumina NovaSeq.

It offers the benefit of higher density sample processing in 384 format, whilst matching the data quality achieved in 96 format described in the original protocol:

<https://www.protocols.io/view/covid-19-artic-v3-illumina-library-construction-an-bgxjxkn>

Both the above protocols were adapted from the nCov-2019 sequencing protocol:

<https://www.protocols.io/view/ncov-2019-sequencing-protocol-bbmuk6w>

## DOI

[dx.doi.org/10.17504/protocols.io.biefkbnn](https://dx.doi.org/10.17504/protocols.io.biefkbnn)

## PROTOCOL CITATION

DNA Pipelines R&D, Benjamin Farr, Diana Rajan, Emma Betteridge, Lesley Shirley, Michael Quail, Naomi Park, Nicholas Redshaw, Iraad F Bronner, Louise Aigrain, Scott Goodwin, Scott Thurston, Stefanie Lensing, Carol Scott, Nicholas Salmon, Charlotte Beaver, Rachel Nelson, Alex Alderton, Ian Johnston 2020. COVID-19 ARTIC v3 Illumina library construction and sequencing protocol - high throughput 384 format.

**protocols.io**

[dx.doi.org/10.17504/protocols.io.biefkbnn](https://dx.doi.org/10.17504/protocols.io.biefkbnn)

## FORK FROM

Forked from [COVID-19 ARTIC v3 Illumina library construction and sequencing protocol](#), Diana Rajan

## KEYWORDS

COVID-19, SARS-Cov-2, amplicon sequencing, ARTIC, Illumina library construction, coronavirus

## LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

## CREATED

Jul 09, 2020

LAST MODIFIED

Jul 24, 2020

PROTOCOL INTEGER ID

39079

#### GUIDELINES

It is vital cDNA setup is performed in a laboratory in which post-PCR COVID-19 amplicons are not present, to minimise any risk of sample contamination.

**Note:** Throughout the protocol we have indicated the liquid handling automation in use at Sanger for specific parts of the process. However, these steps could be performed on alternative liquid handlers or manually.

#### MATERIALS

NAME	CATALOG #	VENDOR
NEBNext Ultra II DNA Library Prep Kit for Illumina - 96 rxns	E7645L	New England Biolabs
2x Kapa HiFi Hotstart Readymix	KK2602	Kapa Biosystems
LunaScript RT SuperMix Kit	E3010L	New England Biolabs
Illumina Library Quantitation Complete kit (Universal)	KK4824	Kapa Biosystems
NEB Q5® Hot Start High-Fidelity 2X Master Mix	M0494L	New England Biolabs
AccuClear® Ultra High Sensitivity dsDNA Quantitation Kit with DNA Standards	#31028	Biotium

#### STEPS MATERIALS



NAME	CATALOG #	VENDOR
AccuClear® Ultra High Sensitivity dsDNA Quantitation Kit with DNA Standards	#31028	Biotium
NEBNext Ultra II DNA Library Prep Kit for Illumina - 96 rxns	E7645L	New England Biolabs
LunaScript RT SuperMix Kit	E3010L	New England Biolabs
NEB Q5® Hot Start High-Fidelity 2X Master Mix	M0494L	New England Biolabs
2x Kapa HiFi Hotstart Readymix	KK2602	Kapa Biosystems

#### MATERIALS TEXT

Primer pool sequences (v3) can be found here:

[https://github.com/joshquick/artic-ncov2019/blob/master/primer\\_schemes/nCoV-2019/V3/nCoV-2019.tsv](https://github.com/joshquick/artic-ncov2019/blob/master/primer_schemes/nCoV-2019/V3/nCoV-2019.tsv)

#### cDNA generation

- 1 Important!** This step must be performed in a RNase free, pre-PCR environment in which post PCR COVID-19 amplicons are not present, to minimise risk of sample contamination.  
  
Decontaminate bench surfaces, pipettes and gloves with RNase ZAP before starting work. Keep reagents and samples chilled throughout the process.
- 2** Defrost PCR plate containing  10 µl extracted RNA  On ice .
- 3**



### LunaScript RT SuperMix Kit

by New England Biolabs



Catalog #:

E3010L

Prepare RT mastermix in a dedicated UV treated pre-PCR area to minimise contamination risk.

RT Master Mix	Vol / RXN (µL)	Vol/384 RXN (µL) inc. excess
LunaScript Super Mix	4	1843
Nuclease-free water	6	2765
<b>Total</b>	<b>10</b>	<b>4608</b>

Mix thoroughly by vortexing.

- 4 Use the SPT Labtech Dragonfly Discovery to dispense  10 µl of RT mastermix into the PCR plate containing  10 µl extracted RNA.
- 5 Seal plate and place on a BioShake plate shaker for 30 seconds at 1500rpm to mix. Briefly centrifuge plate.
- 6 Place plate on a thermocycler and run the following program:

Temperature	Time
25°C	2 minutes
55°C	20 minutes
95°C	1 minute
4°C	∞
Lid temp: Tracking	

- 7 **PAUSE POINT** cDNA can be stored at 4°C (same day) or -20°C (up to a week).

## cDNA amplification

8



Primer pool sequences (v3) can be found here:

[https://github.com/joshquick/artic-ncov2019/blob/master/primer\\_schemes/nCoV-2019/V3/nCoV-2019.tsv](https://github.com/joshquick/artic-ncov2019/blob/master/primer_schemes/nCoV-2019/V3/nCoV-2019.tsv)



**NEB Q5® Hot Start High-Fidelity 2X Master Mix**  
by New England Biolabs  
Catalog #: M0494L

Prepare the following mastermixes:

PCR Primer Pool 1 Master Mix	Vol/PCR RXN (µl)	Vol/384 plate (µl) inc. excess
Q5 Hotstart 2X Master Mix	12.5	5760
Primer <b>Pool 1</b> (11µM total)	3.6	1659
Nuclease-free water	2.9	1336
<b>Total</b>	<b>19</b>	<b>8755</b>

Final concentration of each primer in the reaction is 0.015µM

PCR Primer Pool 2 Master Mix	Vol/PCR RXN (µl)	Vol/384 plate (µl) inc. excess
Q5 Hotstart 2X Master Mix	12.5	5760
Primer <b>Pool 2</b> (10.8µM total)	3.6	1659
Nuclease-free water	2.9	1336
<b>Total</b>	<b>19</b>	<b>8755</b>

Final concentration of each primer in the reaction is 0.015µM


Mix thoroughly by vortexing.

- 9 Use the SPT Labtech Dragonfly Discovery to dispense **19 µl** mastermix per well into 2x384 well plates.
  - 10 Use the Agilent Bravo to add **6 µl** of cDNA template to each primer pool reaction and mix.
  - 11 Heat seal and place the plates onto a thermocycler and run the following program.  
**Important!** Heat seal to minimise evaporation.
- Note:** Amplification should ideally be performed in a different lab to minimise the risk of contaminating other samples.



**Critical step:** We strongly recommend performing a gradient PCR to determine the optimal annealing temperature for your thermocycler. Subtle differences in thermocycler calibration can result in specific

amplicons dropping out. Reducing our annealing temperature from 65°C to 63°C for identical cDNA input recovered amplicon #64 as shown in the image below.

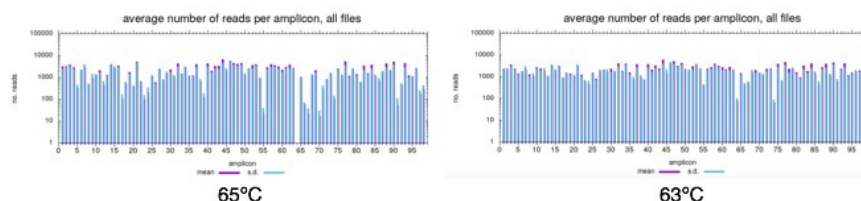

Step	Temperature	Time
1	98°C	30 seconds
2	95°C	15 seconds
3	63°C	5 minutes
4	Repeat steps 2 & 3 for a total of 35 cycles	
5	4°C	∞

**12 PAUSE POINT** Amplified cDNA can be stored at 4°C (overnight) or -20°C (up to a week).

#### Amplified cDNA SPRI

**13** Allow AMPure XP beads to equilibrate to room temperature (~30 minutes). Ensure solution is homogenous prior to use, mixing gently by inversion.

**14** Centrifuge amplified cDNA plates. **1000 x g 00:01:00**

**15** Use the SPT LabTech HV Mosquito to combine **5 µl** of primer 1 PCR and **5 µl** of primer 2 PCR reactions per sample into a new plate. Store the unused portion of primer 1 and primer 2 PCR plates at **-20 °C**. Proceed as follows with the recombined plate.

**16** Use the Hamilton STAR with a 384 well multichannel head to perform the following steps:

**16.1** Add **10 µl** nuclease-free water to each sample and mix well by pipetting.

- 16.2 Add **0.8X** volume of SPRI beads per sample ( **16 µl** SPRI : **20 µl** amplified cDNA), mix well by pipetting.
- 16.3 Incubate for **00:05:00** at **Room temperature**
- 16.4 Transfer the plate to the magnet, allow **00:02:00** for the beads to settle.
- 16.5 Carefully remove and discard the supernatant without disturbing the bead pellet.
- 16.6 Wash the beads with **45 µl** 75% freshly prepared ethanol for **00:00:30** , then remove ethanol and discard.  
(First wash)
- 16.7 Wash the beads with **45 µl** 75% freshly prepared ethanol for **00:00:30** , then remove ethanol and discard.  
(Second wash)
- 16.8 Allow beads to dry **00:05:00**
- 16.9 Remove plate from magnet, add **20 µl** nuclease-free water and resuspend by mixing well.
- 16.10 Incubate for **00:03:00** at **Room temperature**
- 16.11 Transfer the plate to the magnet, allow **00:05:00** for the beads to settle.
- 16.12 Carefully transfer supernatant into a new plate, taking care not to disturb the bead pellet.

**17 PAUSE POINT** Purified amplified cDNA can be stored at -20°C for several weeks prior to library preparation.

Amplified cDNA quantification

18



Purified amplified cDNA is quantified with a fluorescence based assay. We use the **AccuClear Ultra High Sensitivity dsDNA Quantitation kit with 7 DNA standards** (Biotium) according to manufacturer's instructions.

To streamline the workflow, we do not normalise sample input for library preparation. Instead we confirm samples are in the range of **50ng-1ug** per **20µl sample** and take the entire volume into library preparation.



AccuClear® Ultra High Sensitivity  
dsDNA Quantitation Kit with DNA  
Standards

by Biotium

Catalog #: #31028

- 19 Pipette **20 µl** of each DNA standard into wells A1 – G1 of a PCR plate. Add nuclease-free water to H1.
- 20 Dilute the AccuClear dye (100X) to working concentration by mixing **300 µl** dye with **30 mL** AccuClear buffer in a 50ml Falcon. Mix thoroughly by vortexing and transfer to a 384 well reservoir.
- 21 Use the SPT Labtech Mosquito LV to stamp 1µl of amplified cDNA and 1µl of known standards into a 384 assay plate. Immediately proceed to the next step.
- 22 Use the Agilent Bravo 384ST to add **50 µl** 1X AccuClear dye from the reservoir to the assay plate, mix thoroughly by pipetting.
- 23 Measure fluorescence values on a BMG FLUOstar Omega plate reader calibrated for use with AccuClear dye.
- 24 Confirm known standards are performing as expected.
- 25 Dilute any samples >125ng/µl with nuclease free water so they are in the range of 10 - 125ng/µl and repeat quantitation.
- 26 Ensure all samples (20µl total volume) are in the range of 2.5-50ng/µl prior to proceeding with library preparation.

27



We use the NEB NEBNext® Ultra™ II DNA Library Prep Kit for Illumina, which we have automated on the Agilent Bravo platform with some modifications. 200ng is our standard input for library preparation, an acceptable range is 50ng – 1ug per sample. We use a custom adapter set, however any TruSeq adapters are suitable.



**NEBNext Ultra II DNA Library Prep Kit  
for Illumina - 96 rxns**

by New England Biolabs

Catalog #: [E7645L](#)

28

Prepare end prep mastermix **On ice** :

NEBNext End Prep	Vol/PCR RXN (µl)	Vol/384 plate (µl) inc. excess
NEBNext Ultra II End Prep Enzyme Mix	1.2	552
NEBNext Ultra II End Prep Reaction Buffer	2.8	1288
<b>Total</b>	<b>4</b>	<b>1840</b>

Mix well by pipetting.

29

The Bravo will combine **4 µl** of end prep mastermix with **20 µl** amplified cDNA and mix by pipetting.

30

Seal and transfer the plate to a thermocycler and run the following program:

Temperature	Time
20°C	30 minutes
65°C	30 minutes
4°C	∞

31

Prepare adapter ligation mastermix **On ice** :



Adapter Ligation	Vol/PCR RXN (µl)	Vol/384 plate (µl) inc. excess
NEBNext Ultra II Ligation Master Mix	12	5520
NEBNext Ligation Enhancer	0.4	184



TruSeq adapter (10µM)	1	460
<b>Total</b>	<b>13.4</b>	<b>6164</b>

Mix well by pipetting.

32 The Bravo will add  **13.4 µl** adapter ligation mastermix to each sample and mix by pipetting.

33 The plate is incubated on deck at  **20 °C** for  **00:15:00** , however this step may also be performed on a thermocycler.

34



**Note:** We use alternative TruSeq compatible adapters, which do not require the USER enzyme incubation step. If using NEBNext adapters, follow the steps in the NEB protocol to add USER enzyme to the ligation reaction.

35 A **0.8X** SPRI is performed post-ligation.

Ensure AMPure XP beads have been equilibrated to room temperature (~30 minutes) and the solution is homogenous prior to use.



The Bravo will perform a **0.8X** SPRI clean-up and elute sample in 25 µl nuclease-free water as follows:



35.1 Add 0.8X volume of SPRI beads per sample, mix well by pipetting.

35.2 Incubate for  **00:05:00** at  **Room temperature** .

35.3 Transfer the plate to the magnet, allow  **00:02:00** for the beads to settle.

35.4 Carefully remove and discard the supernatant without disturbing the bead pellet.

35.5 Wash the beads with  **45 µl** 75% freshly prepared ethanol for  **00:00:30** , then remove ethanol and discard.  
(First wash)

35.6 Wash the beads with  **45 µl** 75% freshly prepared ethanol for  **00:00:30** , then remove ethanol and discard.

(Second wash)

35.7 Allow beads to dry for ⌚ 00:05:00

35.8 Remove plate from magnet, add 🧴 25 µl nuclease free water and resuspend by mixing well.

35.9 Incubate for ⌚ 00:03:00 at 🌡 Room temperature .

35.10 Transfer the plate to the magnet, allow ⌚ 00:05:00 for the beads to settle.

35.11 Carefully transfer supernatant into a new plate, taking care not to disturb the bead pellet.  
🧴 10 µl of this eluate is used as input for library PCR.

#### Library PCR

36



We use KAPA HiFi HotStart ReadyMix and unique dual indexed (UDI) tag plates for library PCR.

**Note:** this deviates from the standard NEB protocol which uses NEBNext Ultra II Q5 Master Mix and different cycling conditions.



2x Kapa HiFi Hotstart Readymix

by Kapa Biosystems

Catalog #: KK2602

🧴 6900 µl KAPA HiFi HotStart ReadyMix is required per 384 plate (including excess).

37 The Bravo will add 🧴 15 µl KAPA HiFi HotStart ReadyMix and 🧴 10 µl sample into a 🧴 5 µl plate of UDIs and mix thoroughly by pipetting. The final concentration of each UDI in the PCR reaction is 2µM.

38 Seal and transfer the plate to a thermocycler and run the following program:

Temperature	Time
95°C	5 minutes
98°C	30 seconds

65°C	30 seconds
72°C	2 minutes
Repeat 4 times	
72°C	5 minutes
4°C	∞

#### Construct equivolume pool


- 39 In a post-PCR lab, use the Hamilton STAR to combine **3 µl** of each sample per plate to form an equivolume pool of 384 samples.

#### Equivolume pool SPRI

- 40 Allow AMPure XP beads to equilibrate to room temperature (~30 minutes). Ensure solution is homogenous prior to use, mixing gently by inversion.

The Hamilton STAR will perform a **0.8X** SPRI clean-up and elute the final pool in **200 µl** elution buffer as follows:

- 40.1 Add **0.8X** volume of SPRI beads per pool tube, mix well by pipetting.
- 40.2 Incubate for **00:06:00** at **Room temperature**.
- 40.3 Transfer the tube to a magnet, allow **00:04:00** for the beads to form a pellet.
- 40.4 Carefully remove and discard the supernatant, taking care not to disturb the bead pellet.
- 40.5 Wash the beads with **500 µl** 75% ethanol for **00:00:15** then carefully remove ethanol and discard.  
(First wash)
- 40.6 Wash the beads with **500 µl** 75% ethanol for **00:00:15** then carefully remove ethanol and discard.  
(Second wash)
- 40.7 Wash the beads with **500 µl** 75% ethanol for **00:00:15** then carefully remove ethanol and discard.  
(Third wash)
- 40.8 Allow beads to dry for **00:05:00**.

40.9 Remove tube from magnet and resuspend beads in  **200 µl** elution buffer, mix well by pipetting.

40.10 Incubate for  **00:05:00** at  **Room temperature**

40.11 Transfer tube to magnet, allow  **00:00:45** for the beads to form a pellet.

40.12 Carefully transfer supernatant into a new tube, taking care not to disturb the bead pellet.

#### Equivolume pool quantification

41



Equivolume pools may be quantified either by qPCR or on an Agilent Bioanalyzer. Pools are then diluted to 1nM for sequencing.

#### qPCR

Quantify pools in triplicate using the KAPA Complete kit (Universal) for Illumina (KK4824) plus the KAPA Library Quantification Dilution Control (KK4906).

We use the SPT Labtech Mosquito LV to stamp library pools in triplicate into a 384 assay plate, and the Agilent Bravo to setup the qPCR reactions (1:1600 dilution).

qPCR is performed on the Roche LightCycler 480.

#### Agilent Bioanalyzer

Prepare 3 dilutions of the equivolume pool (1:10, 1:100, 1:1000). Run 1 µl of each dilution in triplicate using the High Sensitivity DNA assay kit.

Confirm size distribution is as expected, check there is no primer-dimer or adapter-dimer present.

#### Sequencing

42



We currently sequence samples on an Illumina NovaSeq SP flow cell, using the XP workflow.

Alternatively, samples may be sequenced on an Illumina MiSeq using either v2 (500 cycle) or v3 (600 cycle) reagent kits. We have plexed up to 96 samples per run, this could be increased further depending on coverage requirements. Loading concentration will need to be optimised for MiSeq.  
MiSeq run parameters: Read length 212 paired end + 16bp.

- 43 The following protocol is for loading a NovaSeq. We currently plex up to 384 samples per NovaSeq SP lane.
- 44 Steps must be performed within a given timeframe or data quality may be affected. Therefore, ensure the instrument is washed, waste containers emptied and ready for use prior to beginning step 46.
- 45 Defrost Illumina NovaSeq SP SBS and cluster reagent cartridges for 2-4 hours in a **Room temperature** water bath. Use a lint free tissue to blot any water present on the foil seal. Gently mix cartridges 10X by inversion. Gently tap the bottom of the cartridges on the bench to reduce air bubbles.
- 46 Defrost components DPX1, DPX2 and DPX3 from a NovaSeq XP-2 lane kit, then keep **On ice**
- 47 Bring flow cell to **Room temperature** (~10 minutes) prior to use.
- 48 **18 µl** of each **1 Nanomolar (nM)** pool is required per SP lane.  
Denature pools by adding **4 µl** 0.2N NaOH per 18µl. Vortex briefly to mix.
- 49 Incubate at **Room temperature** for **00:08:00**
- 50 Add **5 µl** 400mM Tris-HCl, pH8.0 to each tube to neutralise the reaction. Vortex briefly to mix, then keep **On ice**.





For the following steps, keep samples and mastermix **On ice** until ready for loading onto the flow cell.




- 51 **Important!** Use mastermix within **01:00:00** of preparation for optimal sequencing performance.

Prepare ExAmp mastermix on ice:

ExAmp Master Mix	Volume per SP flow cell (µl)
DPX1	126
DPX2	18

DPX3	66
<b>Total</b>	<b>210</b>

Vortex  **00:00:30** to mix, then centrifuge briefly up to  **280 x g**

- 52 Add  **63 µl** ExAmp mastermix to each denatured pool, mix well by pipetting.
- 53 Prepare the flowcell for sample loading by placing into the flow cell dock with the 2-lane manifold clamped in place.
- 54 Pipette  **80 µl** of library + ExAmp pool mix per manifold well. Wait for approximately 2 minutes to allow the solution to fill the lane.
- 55 **Important!** The sequencing run must be started within  **00:30:00** of libraries being loaded onto the flow cell.
  - 55.1 Unclamp the flow cell dock and discard the manifold. Load the flow cell onto the NovaSeq flow cell stage.
  - 55.2 Load the SBS and cluster reagent cartridges.
  - 55.3 Start sequencing run (250PE).