



Version 2 ▼


Sep 02, 2022

# Processing of Radial K-space DW-MRI Data V.2

Miguelrj<sup>1</sup>, Hee Kwon Song<sup>1</sup>, Stephen Pickup<sup>1</sup>, Mamtaaryagupta<sup>1</sup>, Rong Zhou<sup>1</sup>

<sup>1</sup>University of Pennsylvania, Department of Radiology

1 Works for me

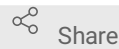

Share

[dx.doi.org/10.17504/protocols.io.j8nlkkwj115r/v2](https://dx.doi.org/10.17504/protocols.io.j8nlkkwj115r/v2)



Miguelrj

## ABSTRACT

The protocol includes reconstruction of diffusion weighted images from radial k-space data and using a 3-parameter fit to derive metrics of ADC and kurtosis index.

## DOI

[dx.doi.org/10.17504/protocols.io.j8nlkkwj115r/v2](https://dx.doi.org/10.17504/protocols.io.j8nlkkwj115r/v2)

## PROTOCOL CITATION

Miguelrj, Hee Kwon Song, Stephen Pickup, Mamtaaryagupta, Rong Zhou 2022. Processing of Radial K-space DW-MRI Data. **protocols.io**  
<https://protocols.io/view/processing-of-radial-k-space-dw-mri-data-cf5qtq5w>  
Version created by Miguelrj



## LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

## CREATED

Sep 02, 2022

## LAST MODIFIED

Sep 02, 2022

## PROTOCOL INTEGER ID

69520

## 1 Image reconstruction

Radially acquired diffusion-weighted images (DWIs) are reconstructed using the following steps:

1. Apply zero-order phase correction to each radial spoke using the average phase offset at the center slice of the lowest b-value image
2. Zerofill k-space by a factor of 2 to double field of view
3. Multiply signal of each point by its respective area on a Voronoi diagram of the points (including added zerofill points) in k-space
4. Re-grid each radially defined point to its nearest Cartesian coordinate using its Kaiser-Bessel index
5. Apply Fourier transform to now Cartesian-defined k-space

This process can be easily done using the **DWI Processing resources** [\[link\]](#).

To use this, place both **DWI\_shell.ipynb** and **DWI\_processing\_functions.py** in the same directory (your working directory). The **readme.txt** file outlines specific file formats, function inputs, and outputs.

Quickly, once all files are formatted according to **readme.txt**, the first cell defines all image and reconstruction parameters:

```
##### Parameters to change #####
input_file = 'fid' # k-space data
# Acquisition data
xres_ro = 128 # total readout points stored
views = 403 # number of radial views
# angl = (math.sqrt(5)-1)/2 * PI # golden angle in radians
angl = 2*PI/views # used by Steve Pickup
# Image data
xres = 96 # actual number of points collected
yres = 96
slices = 16
bvalues = 5
b_array = np.array([10,535,1070,1479,2141]) # b-values
```

Each cell should then be run consecutively.