

### **VERSION 1**

MAY 12, 2023

## OPEN BACCESS

#### DOI

dx.doi.org/10.17504/protocol s.io.e6nvwjwpwlmk/v1

**Protocol Citation:** sandra.cr ameri 2023. TS Procure 812 - primary fixation Karnovsky's - transwells (TM - 013).

#### protocols.io

https://dx.doi.org/10.17504/protocols.io.e6nvwjwpwlmk/v1

License: This is an open access protocol distributed under the terms of the Creative Commons
Attribution License, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working We use this protocol and it's working

Created: May 12, 2023

Last Modified: May 12,

2023

**PROTOCOL** integer ID:

81814

# TS Procure 812 - primary fixation Karnovsky's - transwells (TM - 013) V.1

sandra.crameri<sup>1</sup>

<sup>1</sup>CSIRO ACDP




sandra.crameri


#### **ABSTRACT**


The method is used to Karnovsky's EM fixed transwells to Procure 812 resin blocks. The presence of 4% paraformaldehyde is only used to satisfy microsecurity requirements.

#### **GUIDELINES**

Sample is not optimal for EM, best to dissect and use only the outer 2mm of tissue mass. Sample is usually derived from histology.

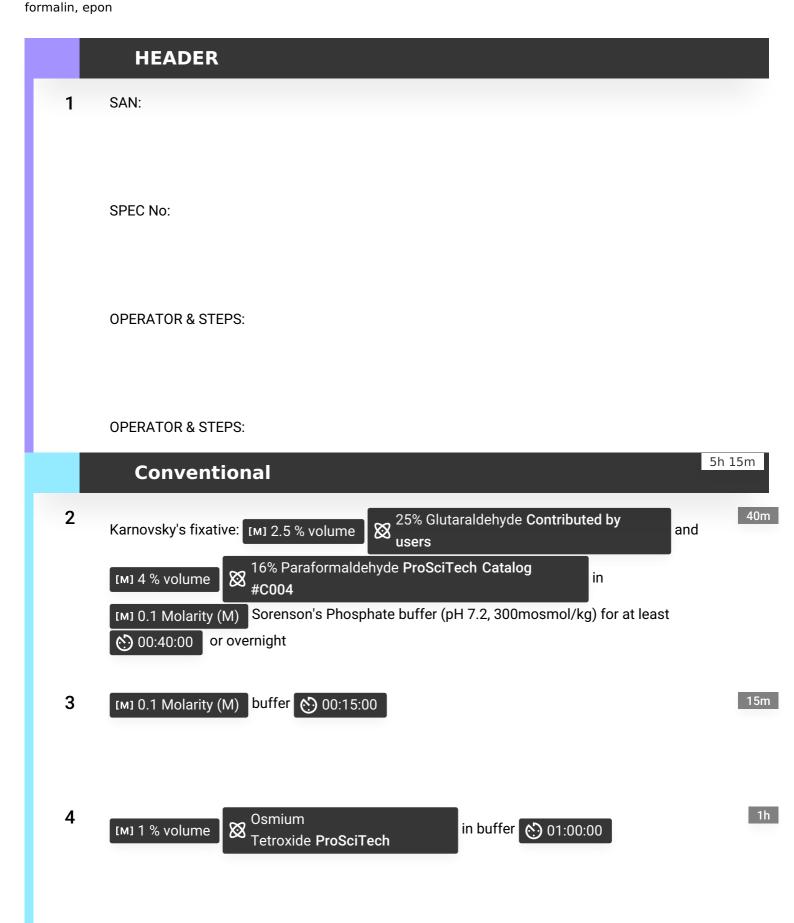

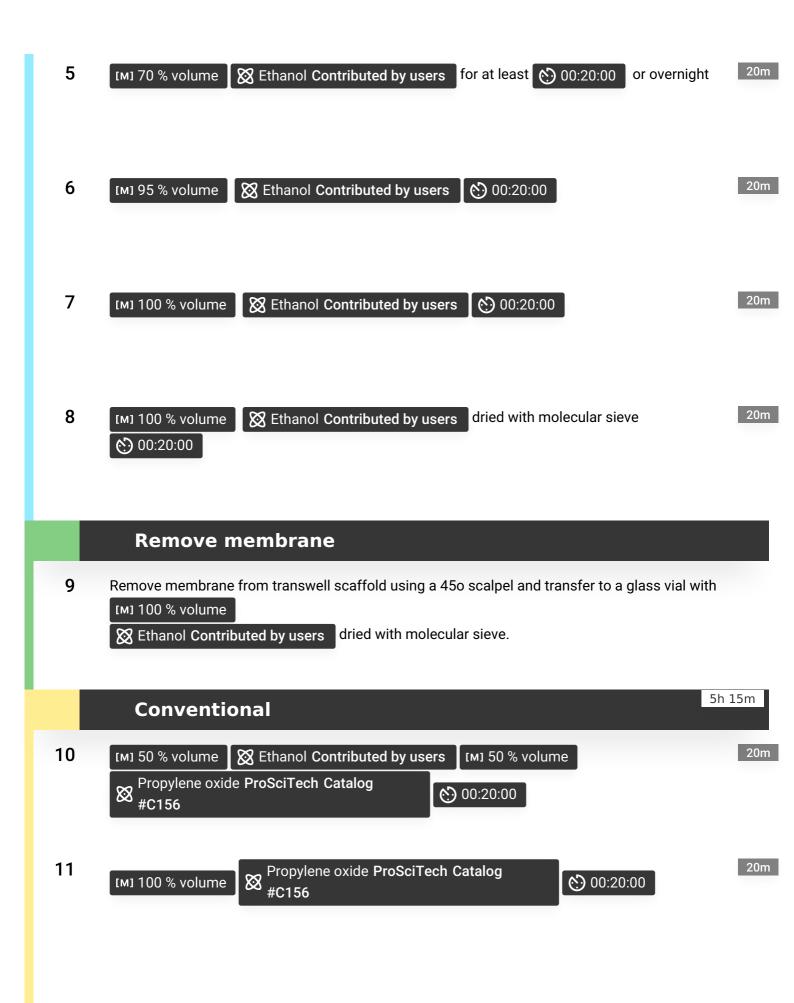

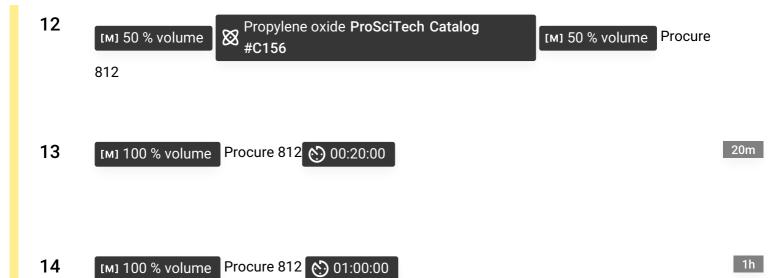

15 Embed & polymerise overnight in oven at  $[60\ ^{\circ}\text{C}]$  , section with diamond knife as usual.