



OCT 26, 2023

OPEN ACCESS



**DOI:**  
[dx.doi.org/10.17504/protocols.io.j8nlko9k5v5r/v1](https://dx.doi.org/10.17504/protocols.io.j8nlko9k5v5r/v1)

**Protocol Citation:** Jason D Limberis, Alina Nalyvayko, Michele Tomasicchio, Shameem Jaumdally, anil pooran, Keertan Dheda, john.metcalfe 2023. Sputum Sample Processing for Single Cell Isolation and Live Recovery for Single Cell RNA Sequencing. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.j8nlko9k5v5r/v1>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working  
 We use this protocol and it's working

## Sputum Sample Processing for Single Cell Isolation and Live Recovery for Single Cell RNA Sequencing

Jason D Limberis<sup>1</sup>, Alina Nalyvayko<sup>1</sup>, Michele Tomasicchio<sup>2</sup>, anil Keertan Dheda<sup>2</sup>, Shameem Jaumdally<sup>2</sup>, pooran<sup>2</sup>, john.metcalfe<sup>1</sup>

<sup>1</sup>University of California, San Francisco;

<sup>2</sup>University of Cape Town, South Africa



Jason D Limberis  
 University of California, San Francisco

### ABSTRACT

Processing of clinical sputum samples for single live cell isolation for downstream applications.

### MATERIALS

⊗ Zymo DNA/RNA Shield Fisher Scientific Catalog #50-125-1706

⊗ Phosphate buffer solution (PBS) Contributed by users

⊗ EDTA (0.5 M), pH 8.0 Life Technologies Catalog #AM9260G

⊗ Bovine Serum Albumin (BSA) Merck MilliporeSigma (Sigma-Aldrich) Catalog #A7906

⊗ Magnesium Chloride Fisher Scientific Catalog #AC223210010

⊗ DMSO Merck MilliporeSigma (Sigma-Aldrich) Catalog #472301

⊗ Fetal calf serum (FCS) Gemini Bio-Products Catalog # 900-108

⊗ DNase I, Bovine Pancreas, >500 Kunitz U/mg Gold Biotechnology Catalog #D-300

⊗ 100 g N-Acetyl-L-Cysteine biorbyt Catalog #orb320385

⊗ Collagenase, Type IV, powder Thermo Fisher Catalog #17104019

**Created:** Sep 11, 2023

## BEFORE START INSTRUCTIONS

**Last Modified:** Oct 26, 2023

Increasing the amount of sputum, by asking the patient to expectorate several time, will likely increase the final cell numbers.

**PROTOCOL integer ID:**  
87649

**Keywords:** sputum, single cell isolation, single cell RNA sequencing, scRNAseq, infectious disease

### Funders

### Acknowledgement:

UCB Catalyst

Grant ID: Chipotle

## Prepare Buffers

### 1 5mM EDTA in PBS

500ul 0.5M EDTA stock  
49.5ml PBS. Store at 4C

### 10% BSA stock


Add 1g BSA to 10ml PBS, mix gently till dissolved  
Filter sterilize & keep at 4C

### 50mM MgCl<sub>2</sub> solution

Add 238mg of MgCl<sub>2</sub> to 50ml of H<sub>2</sub>O

### 10% DMSO/FCS freezing mixture

 5 mL DMSO




 45 mL fetal calf serum (FCS)

 7.5 Mass Percent **NALC in PBS**

### DNase

Aliquoted into 100ul aliquots at 224ku/ml and stored at -20C  
Add 100ul to 40ml dissociation buffer to make a final conc of 0.56ku/ml.

### Collagenase IV

Add  500  $\mu$ L to  40 mL dissociation buffer to make a final conc of 0.05mg/ml  
Aliquoted  500  $\mu$ L at 4mg/ml and stored at -20C.

**Make fresh for each use**

### 1% BSA in PBS

🧴 5 mL 10% BSA

🧴 45 mL PBS

### Dissociation buffer

🧴 1.6 mL of 50 millimolar (mM)  $\text{MgCl}_2$

🧴 100  $\mu\text{L}$  DNase

🧴 500  $\mu\text{L}$  collagenase

🧴 37.8 mL PBS

## Sputum collection

2 Collect an induced sputum, preferably 🧴 5 mL or more in volume, and process as below as soon as possible.

2.1 Weigh the sputum collection cup before and after the sputum collection.



2.2 Estimate the volume of sputum collected using a graduated container of the same size.


2.3 Record the viscosity (salivary, mucosalivary, purulent, mucopurulent)


2.4 Record the time of sputum collection

## Bulk RNA sequencing


3 Record the time of sputum processing

4 Using a wide bore tip (or cut the bottom of a  1000  $\mu\text{L}$  ) draw up ~  200  $\mu\text{L}$  sputum. If the sample is too viscous for a pipette, use a syringe needle instead.


5 Add to a 1.5ml tube containing  500  $\mu\text{L}$  Zymo RNA Shield



6 Vortex tube for  00:00:10

10s


7 Immediately transfer to  -80 °C


## Make a single cell suspension (scRNAseq)

8 Add  40 mL dissociation buffer into the collection tube containing the sputum sample


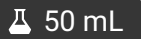

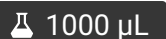

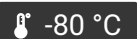


9 Secure cap, manually invert the tube several times. Secure to a hula mixer (low setting) and incubate for  00:10:00 at  Room temperature .

20m

Note: Big clumps should dissolve after  00:10:00 . Some small fibers may will remain.






10 Add DTT to a final concentration of 1-2mM.  
Invert the tube several times and incubate for  00:05:00

5m

- 11 Pass the dissociated sample through a  100 µm strainer on top of a new  50 mL tube using an automated pipettor (do not pour).  10m
- 12 Wash the strainer with  2 mL of PBS
- 13 Spin tube in a sealed bucket at  10m
- 14 Remove  1000 µL supernatant and place in a cryovial for Luminex Cytokine assays. Add  140 µL of protease inhibitor and store at  -80 °C
- 15 Aspirate the remainder of the supernatant and discard.
- 16 Gently resuspend cell pellet in  1 mL EDTA/PBS using a wide bore  1000 µL tip
- 17 Count cells using Turks protocol followed by paraformaldehyde fixation to decontaminate

## Long-term storage of cells (scRNAseq)


- 18 Pre-cool CellCool and cryovials at  4 °C for >  02:00:00 2h

- 19 Place cells  On ice
- 20 Gently mix the cells
- 21 Centrifuge at  300 rcf, 4°C, 00:05:00 5m
- 22 Discard supernatant and resuspend cell pellets in pre-chilled ( 4 °C) 10% DMSO/FCS solution (add dropwise)
- 23 Dispense cell suspension in  500 µL aliquots into the pre-cooled cryovials and place into CoolCell at  -80 °C


## Expected results

### 24 Pilot samples

A	B	C	D
Sample	S1	S2	S3
Viscosity	mucosa livery	mucosa livery	mucosal ivary
Volume used (ml)	10	4	6
Total cells before freezing (10 <sup>6</sup> )	2.2	0.9	2.9
Total cells recovered after thawing (2 weeks after frozen; 10 <sup>6</sup> )	0.88	0.48	1.33
Percentage cell recovery post freezing (%)	40	53	46

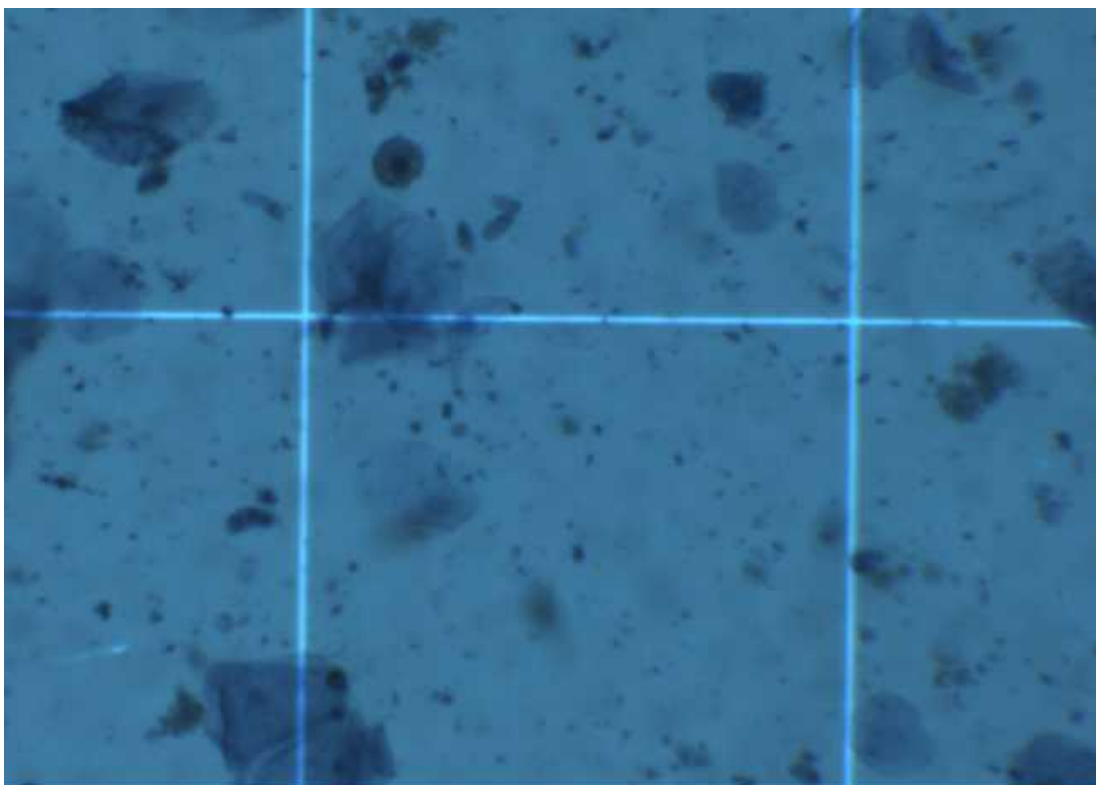

S2 400x

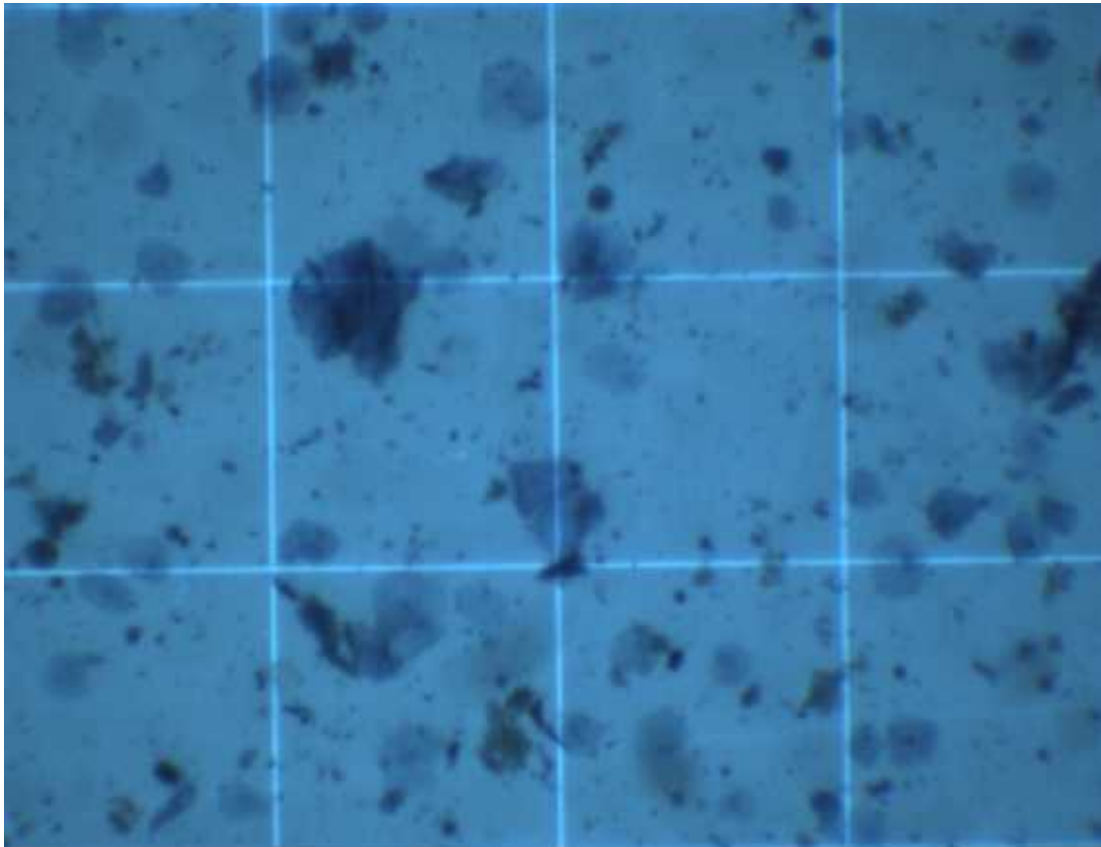

S3 100x

### Clinical samples

A	B	C	D	E	F	G	H	I	J	K	L	M	N
Sample ID	Time to sputum processing (h.m)	weight of sputum (g)	volume of sputum (ml)	viscosity of sputum	Time to thaw (days)	Total cell count pre freeze (x10 <sup>6</sup> )	Total cell count taken forward pre freeze (x10 <sup>6</sup> )	Total cell count post freeze (x10 <sup>6</sup> )	Total cell count difference (x10 <sup>6</sup> )	Live cell count	Dead cell count	Viability (%)	Viability of prefreeze total (%)
CC20	2.85	17.8	15	salivary	662	9.1	7	4.86	2.14	1.58	3.28	33	23
CC22	<6.00	2.1	2	mucoid	658	1.3	1.3	0.60	0.71	0.15	0.44	26	12
CC23	2.10	2.9	3	mucosalivary	657	0.43	0.43	0.37	0.06	0.08	0.29	21	18



	A	B	C	D	E	F	G	H	I	J	K	L	M	N
	CC2 4	1.10	8.7	8	muc osali vary	657	1.1	1.1	0.81	0.30	0.19	0.62	23	17
	CC2 5	4.90	0.7	0.9	muc osali vary	652	0.17	0.17	0.00	0.17	0.00	0.00	0	0
	CC2 6	2.00	1.5	1.4	muc osali vary	652	1.13	1.13	0.73	0.40	0.17	0.57	23	15
	CC2 7	2.20	3.2	3	muc osali vary	651	1.4	1.4	1.33	0.07	0.30	1.03	23	22
	CC2 8	2.00	12	10	muc oid	644	3.2	2	2.74	-0.74	0.29	2.45	11	15
	CC3 3	2.00	3	3	muc osali vary	627	0.8	0.8	0.67	0.13	0.25	0.42	37	31
	CC3 7	<6.0 0	4.1	4	muc osali vary	623	0.92	0.92	0.97	-0.05	0.34	0.63	35	37