



NOV 23, 2023

OPEN ACCESS



**Protocol Citation:** Allan JW Lui 2023. Lentivirus production & concentration. **protocols.io** <https://protocols.io/view/lentivirus-production-amp-concentration-byphpvj6>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working  
We use this protocol and it's working

**Created:** Sep 30, 2021

**Last Modified:** Nov 23, 2023


**PROTOCOL integer ID:** 53705

## Lentivirus production & concentration

Allan JW

Lui<sup>1</sup>

<sup>1</sup>Cancer Research UK, Cambridge Institute

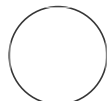

Allan JW Lui

### ABSTRACT


Lentivirus production protocol based on official protocol for Lipofectamine 3000 and TransIT-Lenti

### ATTACHMENTS

[Lipofectamine3000-LentiVirus-AppNote-Global-FHR.pdf](#)


### PROTOCOL MATERIALS


 LentiFuge Collecta Catalog #LFVC1 Step 17

 Opti-MEM<sup>®</sup>; I Reduced Serum Medium Thermo Fisher Catalog #31985047


Step 7

 GlutaMAX<sup>®</sup>; Supplement Thermo Fisher Catalog #35050038 Step 7

 Sodium Pyruvate (100 mM) Thermo Fisher Catalog #11360070 Step 7

 DMEM, low glucose, GlutaMAX<sup>®</sup>; Supplement, pyruvate Thermo Fisher Catalog #21885025


Before starting

 Fetal Bovine Serum, qualified, One Shot<sup>®</sup>; format, Brazil Thermo Fisher Catalog #A3160801


Before starting, Step 7

### BEFORE START INSTRUCTIONS

Prepare culture media for HEK293T:

 DMEM, low glucose, GlutaMAX<sup>®</sup>; Supplement, pyruvate Thermo Fisher Catalog #21885025

supplemented with 10%

 Fetal Bovine Serum, qualified, One Shot<sup>®</sup>; format, Brazil Thermo Fisher Catalog #A3160801

**Keywords:** Lentivirus  
production, lentivirus  
concentration, lentivirus

## HEK293T seeding density titration

1d

1 Detach, count & seed HEK293T cells into T75 flasks at:

1h

1. 4 x 10E6 cells
2. 6 x 10E6 cells
3. 8 x 10E6 cells
4. 10 x 10E6 cells

### Note

>28 x 10E6 cells required in total; recommend detaching cells from 1x T175 at 70-80% confluence.

2 Incubate plates at 37 °C incubator (5% CO<sub>2</sub>, humidified) for up to 24:00:00

1d

Observe confluence of cells under a microscope at 18 - 24 hours after plating.

The optimal plating density for transfection yield 80-90% confluence 18 - 24 hours after plating.

## Lentivirus production (Lipofectamine 3000)

3 Detach, count & seed HEK293T cells at a density that yields 80-90% confluence 18 - 24 hours after plating

1d

Incubate at 37 °C incubator (5% CO<sub>2</sub>, humidified) for up to 24:00:00

### Note

Lentivirus packaging medium (step 7 below) can be used when plating cells & during transfection the next day, but DMEM+10% FBS works well too

4 When HEK293T cells are at 80-90% confluence, transfect cells with lentiviral transfer and packaging plasmids using lipofectamine 3000

Culture plate	Scaling factor	Growth media (ml)	Lipofectamine mix		DNA mix		
			Opti-MEM (ml)	Lipofectamine 3000 (ul)	Opti-MEM (ml)	DNA (ug)	P3000 (ul)
6-well	1	2	0.25	7	0.25	3	6
100mm	6	10	1.5	42	1.5	18	36
T75	8	15	2.0	56	2.0	24	48
150mm	16	30	4.0	112	4.0	48	96
T175	18	35	4.5	126	4.5	54	108

Table 1.

#### Note

If transfecting multiple vessels, make a master mix of 1.1x the total volume required

Source document: [Lipofectamine3000-LentiVirus-AppNote-Global-FHR.pdf](#)


## 4.1 Calculate the amount of each plasmid required


Example using lentiCas9-Blast, psPAX2 and pMD2.G:

A	B	C	D	E	F	G	H
			Mass of plasmids per vessel (ug)				
Plasmid	Size (bp)	Equimolar ratio (ug)	6-well	10cm	T75	150mm	T175
lentiCas9-Blast	12859	2.209	1.313	7.88	10.50	21.00	23.63
psPAX2	10709	1.839	1.093	6.56	8.75	17.49	19.68
pMD2.G	5822	1.000	0.594	3.57	4.75	9.51	10.70
Total			3	18	24	48	54

### Note



Optimisation of plasmid ratios may be required for each transfer plasmid; an equimolar ratio is a good starting point. Zhang Lab uses a 4:3:2 transfer : psPAX2 : pMD2.G ratio that approximates an equimolar ratio for lentiCas9-blast, lentiGuide-Puro & lentiCRISPR v2 plasmids. With lentiCas9-blast, 4:2:1, 9:8:1, 7.5:1.5:4 ratios all result in similar titres (+/- 20% of 4:3:2), around 1-2e6 TU/ml.

- 4.2** Make up DNA mix by diluting transfer and packaging plasmids into Opti-MEM, then adding P3000 reagent (2ul / 1ug DNA)  
Mix well by vortexing  00:00:03 **3s**

- 4.3** Make up lipofectamine mix by diluting Lipofectamine 3000 (7ul / 3ug DNA) into Opti-MEM (volume specified in Table 1)  
Mix well by vortexing  00:00:03 **3s**

### Note

Lipofectamine 3000 reagent diluted in Opti-MEM medium should be used within 15 minutes of dilution. Longer times can result in a loss of transfection efficiency

- 5** Combine both lipofectamine and DNA mixes and incubate for  00:20:00  Room temperature **20m**  
allow DNA-lipid complexes to form

Remove 50% volume of media from HEK293T culture vessels intended for transfection

A	B	C	D	E
Volume of media to remove per vessel (ml)				
6-well	100mm dish	T75	150mm dish	T175
1	5	7.5	15	17.5

- 6** Gently add DNA-lipid complexes to cells **6h**

A	B	C	D	E
Volume of DNA-lipid complex to add per vessel (ml)				

A	B	C	D	E
6-well	10cm dish	T75	150mm dish	T175
0.5	3	4	8	9

Gently rock culture vessel back-and-forth and from side-to-side to evenly distribute

Incubate at 37 °C incubator (5% CO<sub>2</sub>, humidified) for 06:00:00

## 7 Prepare Lentivirus packaging medium (Opti-MEM + 1x GlutaMAX + 1mM sodium pyruvate + 5% FBS)

Supplement 500 mL

Opti-MEM<sup>®</sup>; I Reduced Serum Medium Thermo  
Fisher Catalog #31985047 with:

- 5 mL GlutaMAX<sup>®</sup>; Supplement Thermo Fisher Catalog #35050038

- 5 mL Sodium Pyruvate (100 mM) Thermo  
Fisher Catalog #11360070

- 25 mL Fetal Bovine Serum, qualified, One Shot<sup>®</sup>; format, Brazil Thermo  
Fisher Catalog #A3160801

## 8 5-6 hours after transfection, replace media with pre-warmed packaging medium

Incubate plates at 37 °C incubator (5% CO<sub>2</sub>, humidified)

### Note

Consider supplementation with 1 millimolar (mM) sodium butyrate or  
 2 millimolar (mM) caffeine to improve titer.

Supplements will be present in viral supernatant if not concentrated

## 9 24 hours post-transfection, collect supernatant from culture vessels and store at 4 °C

Replace media with pre-warmed packaging medium

Incubate plates at 37 °C incubator (5% CO<sub>2</sub>, humidified) for up to 28:00:00

## 10 48-52 hours post-transfection, collect supernatant from culture vessels and combine with supernatant harvested the day before.

Centrifuge at 300 x g, Room temperature, 00:05:00 to separate any detached cells, and collect the supernatant.

- 11 Filter the supernatant through a 45µm PES filter to remove any remaining cellular debris

**Note**

Smaller filters (22µm) improve purity but can lower the viral titer.

- 12 Prepare single-use aliquots in cryovials and snap freeze virus in a dry ice-ethanol bath, then store at -80°C.

Alternatively, proceed with virus concentration or infection of cells.

**Note**

It is recommended to avoid subjecting lentiviral preparations to multiple freeze-thaw cycles, since each cycle can result in a 10%–20% loss in functional titers.

## Lentivirus production (TransIT-Lenti)

- 13 Detach, count & seed HEK293T cells at a density that yields 80-90% confluence 18 - 24 hours after plating

1d

Incubate at  37 °C incubator (5% CO<sub>2</sub>, humidified) for up to  24:00:00

- 14 Change media for HEK293T cells

- 15 When HEK293T cells are at 80-90% confluence, transfect cells with lentiviral transfer and packaging plasmids using TransIT-Lenti

Culture vessel	Scaling factor	Growth media (ml)	Opti-MEM (ml)	DNA (ug)	TransIT-Lenti (ul)
6-well	1.0	2	0.2	2	6
10cm	6.0	10	1.2	12	36
T75	7.5	15	1.5	15	45
150mm	15.0	30	3.0	30	90
T175	17.5	35	3.5	35	105

#### Note

If transfecting multiple vessels, make a master mix of 1.1x the total volume required

### 15.1 Calculate the amount of each plasmid required


Example using lentiCas9-Blast, psPAX2 and pMD2.G:

A	B	C	D	E	F	G	H
			Mass of plasmids per vessel (ug)				
Plasmid	Size (bp)	Equimolar ratio (ug)	6-well	10cm	T75	150mm	T175
lentiCas9-Blast	12859	2.209	0.875	5.25	6.56	13.13	15.31
psPAX2	10709	1.839	0.729	4.37	5.47	10.93	12.75
pMD2.G	5822	1.000	0.396	2.38	2.97	5.94	6.93
Total			2	12	15	30	35

#### Note

Optimisation of plasmid ratios is required for each transfer plasmid but an equimolar ratio is a good starting point, Zhang Lab uses a 4:3:2 ratio that approximates an equimolar ratio for lentiCas9-blast, lentiGuide-Puro & lentiCRISPR v2 plasmids.

### 15.2 Warm TransIT-Lenti reagent to room temperature and vortex

 00:00:03

3s

15.3

Combine plasmids in quantities determined in step 14.1  
Add Opti-MEM (Volume specified in Table 2), mix by pipetting  
Add TransIT reagent (3ul / 1ug DNA), mix by pipetting

10m

Incubate at Room temperature 00:10:00 to allow transfection complexes to form

**Note**

Precipitation may be observed when excess DNA is used during complex formation. This may negatively impact transfection efficiency.

15.4

Add transfection complexes drop-wise to different areas of vessel  
Gently agitate culture vessel back-and-forth and from side-to-side to evenly distribute

2d

Incubate at 37 °C incubator (5% CO<sub>2</sub>, humidified) for 48:00:00

16

Harvest virus 48 hours post-transfection by following steps 10 - 12.

**Lentivirus concentration**

1h

17

Add 1 µL LentiFuge Collecta Catalog #LFVC1 per 1 mL lentiviral supernatant  
Incubate at 4 °C 01:00:00

1h

18

Centrifuge at 12000 rpm, 4°C, 01:00:00

1h


**Note**

Pellets generated from small supernatant volumes may be difficult to see. It is advised to mark the centrifuge tube or bottle with a marker at the site where you expect the virus pellet to be.

19

Resuspend the pelleted lentiviral preparation in 1/100 of original volume using sterile phosphate buffered saline (PBS).





**20** Prepare single-use aliquots and snap freeze virus in dry ice, then store at  $-80^{\circ}\text{C}$ .