


May 08, 2024 Version 1

🌐 Imaging single SYTOX Orange molecules on a PLL-coated cover glass V.1

DOI

dx.doi.org/10.17504/protocols.io.x54v9pbbmg3e/v1

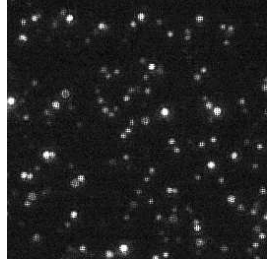

Ezra Bruggeman¹

¹University of Cambridge

ASAP Collaborative Rese...



Ezra Bruggeman

University of Cambridge

OPEN  ACCESS



DOI: dx.doi.org/10.17504/protocols.io.x54v9pbbmg3e/v1

Protocol Citation: Ezra Bruggeman 2024. Imaging single SYTOX Orange molecules on a PLL-coated cover glass. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.x54v9pbbmg3e/v1>

Manuscript citation:

POLCAM: Instant molecular orientation microscopy for the life sciences

Ezra Bruggeman, Oumeng Zhang, Lisa-Maria Needham, Markus Körbel, Sam Daly, Matthew Cheetham, Ruby Peters, Tingting Wu, Andrey S. Klymchenko, Simon J. Davis, Ewa K. Paluch, David Klenerman, Matthew D. Lew, Kevin O'Holleran, Steven F. Lee.

bioRxiv 2023.02.07.527479, doi: <https://doi.org/10.1101/2023.02.07.527479>

License: This is an open access protocol distributed under the terms of the **Creative Commons Attribution License**, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

We use this protocol and it's working.

Created: January 09, 2024

Last Modified: May 08, 2024

Protocol Integer ID: 93146



Keywords: ASAPCRN, Single-molecule, Single molecules, SYTOX Orange, Microscopy, Imaging, Fluorescence, Fluorescence microscopy, POLCAM

Funders Acknowledgement:

Aligning Science Across

Parkinson's

Abstract


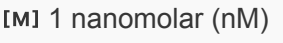
This is a protocol for the preparation of a microscopy sample of single SYTOX Orange molecules on a PLL-coated cover glass. This protocol was used to generate the data shown in **Figure 1a, 1b** and **1c** of the following publication:

- Bruggeman *et al.*, POLCAM: Instant molecular orientation microscopy for the life sciences. bioRxiv 2023.02.07.527479 (Feb **2023**), doi: <https://doi.org/10.1101/2023.02.07.527479>





Protocol

1h 30m

- 1 Argon plasma clean cover glass (VWR collection, 631-0124) for  00:30:00 in a plasma cleaner (Expanded Plasma Cleaner, PDC-002, Harrick Plasma). 30m
- 2 In the meantime:
 - Filter phosphate-buffered saline (PBS) using a 0.02 µm syringe filter (6809-1102, Whatman).
 - Dilute SYTOX Orange (S11368, Invitrogen) in filtered PBS to a concentration of  [M] 1 nanomolar (nM) .


Note

Always use a new aliquot of SYTOX Orange to prepare the 1 nM dilution, as dye doesn't store well at low concentrations.

- 3 Create a sample well on the cleaned cover glass by sticking a frame-seal slide chamber (9x9 mm, SLF0201, Bio-rad) on the cover glass.
- 4 Pipet  70 µL of 0.01% PLL (0.01% poly-L-lysine solution, P4707, Sigma-Aldrich) into the well and wait for  00:15:00 . The PLL will coat the surface of the cover glass. 15m

Note




Always use a freshly thawed aliquot of PLL. You can aliquot the PLL and store it in a -20 °C or -80 °C freezer.

- 5 Use a pipet to remove the excess PLL from the well and immediately replace it with  70 µL of filtered PBS.

Note

It is important to always have liquid on top of the PLL-coated glass and not let it dry out.



- 6 Use a pipet to remove the excess filtered PBS from the well and immediately replace with  70 μL filtered PBS. Gently pipet up and down in the corners of the well. Repeat this step 2 more times.
- 7 Use a pipet to remove the excess PBS from the well and immediately replace with  50 μL of  1 nanomolar (nM) SYTOX Orange (S11368, Invitrogen). The SYTOX Orange molecules will stick to the surface of the PLL-coated cover glass.
- 8 Image the sample straight away and make sure it doesn't dry out during imaging.