

Whole-cell proteomics and Analysis by Tandem Mass Tagging-based proteomics

AUG 04, 2023

OPEN ACCESS



DOI:
dx.doi.org/10.17504/protocols.io.kqdg3x4n1g25/v1

Protocol Citation: Felix Kraus, Sharan Swarup, Vinay V. Eapen, J. Wade Harper wade_harper@hms.harvard.edu 2023. Whole-cell proteomics and Analysis by Tandem Mass Tagging-based proteomics. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.kqdg3x4n1g25/v1>

MANUSCRIPT CITATION:
 - McAlister, G. C. et al. MultiNotch MS3 Enables Accurate, Sensitive, and Multiplexed Detection of Differential Expression across Cancer Cell Line Proteomes. Analytical chemistry 86, 7150-7158 (2014).
 - Paulo, J.A., et al. Quantitative mass spectrometry-based multiplexing compares the abundance of 5000 S. cerevisiae proteins across 10 carbon sources. J Proteomics 148, 85-93 (2016a).

Whole-cell proteomics and Analysis by Tandem Mass Tagging-based proteomics

Forked from [Whole-cell proteomics and Analysis by Tandem Mass Tagging-based proteomics](#)

In 1 collection


Felix Kraus^{1,2}, Sharan Swarup³, Vinay V. Eapen³,
 J. Wade Harper wade_harper@hms.harvard.edu^{1,2}

¹Department of Cell Biology, Blavatnik Institute, Harvard Medical School, 240 Longwood Ave, Boston MA 02115, USA;

²Aligning Science Across Parkinson's (ASAP) Collaborative Research Network, Chevy Chase, MD 20815, USA;

³Department of Cell Biology, Harvard Medical School, Boston MA 02115

Harper

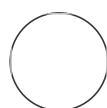

Melissa Hoyer

DISCLAIMER

DISCLAIMER – FOR INFORMATIONAL PURPOSES ONLY; USE AT YOUR OWN RISK

The protocol content here is for informational purposes only and does not constitute legal, medical, clinical, or safety advice, or otherwise; content added to [protocols.io](#) is not peer reviewed and may not have undergone a formal approval of any kind. Information presented in this protocol should not substitute for independent professional judgment, advice, diagnosis, or treatment. Any action you take or refrain from taking using or relying upon the information presented here is strictly at your own risk. You agree that neither the Company nor any of the authors, contributors, administrators, or anyone else associated with [protocols.io](#), can be held responsible for your use of the information contained in or linked to this protocol or any of our Sites/Apps and Services.

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

Created: Jun 23, 2023

Last Modified: Aug 04, 2023

PROTOCOL integer ID: 83944

Keywords: Tandem Mass Tagging-based proteomics, Whole-cell proteomics and Analysis, ASAPCRN

ABSTRACT

The analysis of relative protein abundance has emerged as an important tool in cell biology. Typically, it is possible to quantify >8000 proteins under standard conditions. Tandem Mass Tags (TMT) are isobaric reagents that contain a set of isotopically distinct reporter ions, which can be used to quantify individual peptides in distinct samples through multiplexing([McAlister et al., 2014](#)). Because the TMT analysis is performed in multiplexed format (up to 18 plex), it is possible to examine the effect of different perturbations (treatments, time courses, etc) on the total abundance of the proteome and include replicate samples as desired. This protocol is applicable to many different cell types, although the number of proteins quantified may differ, depending on the complexity of the proteomes in individual cell types. The small amount of protein needed (50-100 ug) makes application of this approach simple for many different types of cells

ATTACHMENTS

[Whole-cell proteomics and Analysis by Tandem Mass Tagging-based proteomics.pdf](#)

GUIDELINES

Mass spectrometry:

The analysis of TMT-labelled peptides by mass spectrometry will depend on the type of instrument/platform used. Typical instrument settings for analysis on a Thermo Fusion Lumos instrument are provided in the following section.

Inject 3 µl for each LC–MS/MS analysis using available mass spectrometer with a 120-minute online LC separation.

Search raw data against UniProt human protein database using any proteomic analysis software with the following parameters:

- Up to 3 missed cleavages allowed for trypsin/LysC digestion
- Carbamidomethyl (C), TMT (N-term peptide and K) set as a fixed modification
- Oxidation (M) and di-glycine (K) set as variable modifications

Extract signal to noise intensity values of each TMT reporter and identified proteins, and further calculate the ratio of each condition to the control sample's intensity.

Instrument settings:

Mass spectrometry data were collected using an Orbitrap Fusion Lumos mass spectrometer (Thermo Fisher Scientific, San Jose, CA) coupled to a Proxeon EASY-nLC1200 liquid chromatography (LC) pump (Thermo Fisher Scientific). Peptides were separated on a 100 µm inner diameter microcapillary column packed in house with ~35 cm of Accucore150 resin (2.6 µm, 150 Å, ThermoFisher Scientific, San Jose, CA) with a gradient consisting of 3%–26% (0–100 min), 26–32% (100–110min) (ACN, 0.1% FA) over a total 120 min run at ~400 nL/min. For analysis, we loaded 1/3 of each fraction onto the column. Each analysis used the Multi-Notch MS³-based TMT method([McAlister et al., 2014](#)).

The scan sequence began with an MS¹ spectrum (Orbitrap analysis; resolution 120,000 at 200 Th; mass range 400–1250 m/z; automatic gain control (AGC) target 1×10⁶; maximum injection time 100 ms). Precursors for MS² analysis were selected using a Top 4 sec method. MS² analysis consisted of collision-induced dissociation (quadrupole Orbitrap analysis; AGC 1×10⁵; isolation window 0.7 Th; normalized collision energy (NCE) 35; maximum injection time 300 ms resolution was 7,500 at 200 Th). Monoisotopic peak assignment was used, and previously interrogated precursors were excluded using a dynamic window (120 s ± 7 ppm). Following acquisition of each MS² spectrum, a synchronous-precursor-selection (SPS) MS³ scan was collected on the top 10 most intense ions in the MS² spectrum([McAlister et al., 2014](#)). MS³ precursors were fragmented by high energy collision-induced dissociation (HCD) and analyzed using the Orbitrap (NCE 65; AGC 2×10⁵; maximum injection time 500 ms, resolution was 50,000 at 200 Th).

MATERIALS

Chemicals, Peptides, and Recombinant Proteins

⊗ Protease Inhibitor Cocktail Roche Catalog #11873580001

⊗ PBS (10x) Santa Cruz Catalog #sc-24947

⊗ tris(2-carboxyethyl)phosphine (TCEP) Gold Biotechnology Catalog #TCEP2

⊗ Formic Acid (FA) Sigma – Aldrich Catalog #94318

⊗ Urea Sigma – Aldrich Catalog #U5378

⊗ Acetonitrile (ACN) Sigma – Aldrich Catalog #34851

⊗ Sodium Chloride Sigma – Aldrich Catalog #S9888

⊗ MOPS Sigma – Aldrich Catalog #M1254

⊗ Sequencing grade Trypsin Promega Catalog #V5111

⊗ Lys-C Wako Chemicals Catalog #129-02541

⊗ EPPS Sigma – Aldrich Catalog #E9502

⊗ 2-Chloroacetamide Sigma – Aldrich Catalog #C0267

⊗ Protein A Plus Ultralink resin Thermo Fisher Scientific Catalog #53142

⊗ Sodium metaborate Sigma – Aldrich Catalog #S0251

⊗ Aeris peptide XB-c18 column Contributed by users Catalog #00G-4507-E0

⊗ Dimethyl pimelimidate dihydrochloride (DMP) Sigma – Aldrich Catalog #D8388

Critical Commercial Assays

⊗ Tandem Mass Tags Thermo Fisher Scientific Catalog #90406

⊗ Bio-Rad Protein Assay Dye Reagent Concentrate BIO-RAD Catalog #5000006

Software

- SEQUEST Eng et al., 1994

- COMET Eng et al., 2013 <http://comet-ms.sourceforge.net/>

Other



Sep-Pak C18 1cc Vac Cartridge 50
mg Waters



Empore™ SPE Disks C18 Contributed by users Catalog
#2215

- Orbitrap Fusion Lumos Mass Spectrometer, ThermoFisher Scientific, Cat#IQLAAEGAAPFADBMBHQCat#IQLAAEGAAPFADBMBHQ
- high-pH reverse-phase HPLC - LC1260, Agilent
- Protein Lo-Bind Tubes

BUFFERS:

1. Urea lysis buffer:

A	B
Compound	[Compound]final
Urea	8 M
NaCl	75 mM
Tris pH 7.4	150 mM
Protease Inhibitors Phosphatase Inhibitors	1 x 1 x

2. EPPS buffer (50 mM EPPS, pH 8.5)
3. Phosphate buffered saline (pH 7.4)


SAFETY WARNINGS



Please refer to the Safety Data Sheets (SDS) for health and environmental hazards.

BEFORE START INSTRUCTIONS

WORKFLOW:

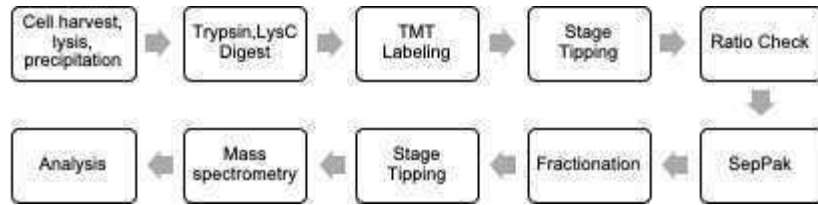

Harvest, precipitation and digestion

9h 10m

- 1 For whole proteome analysis, 50 µg of protein lysate is required for each replicate. Lyse cells in lysis buffer and pass them through a 21G needle 10 times. Alternatively, lyse cells by sonication as per manufactures instructions.
- 2 Centrifuge suspensions at 13000 rpm, 4°C , (high speed) for 00:10:00 and collect the supernatant. 10m
- 3 Transfer quantified protein lysate concentration and 50 µg of lysate for each sample to a clean 1.5 mL protein Lo-Bind Eppendorf tubes. Reduce lysates for 00:20:00 at Room temperature with 5 millimolar (mM) TCEP, and alkylate cysteine residues with 20 millimolar (mM) Chloroacetamide (Room temperature, 00:30:00). 50m
- 4 Extract protein content by methanol-chloroform precipitation and subsequent MeOH washes.
 - 4.1 Add 4x volumes of MeOH and vortex.
 - 4.2 Add 1x volume of chloroform and vortex.



4.3 Add 3x volume of water and vortex.

4.4 Spin down at  Room temperature for  00:05:00 at high speed.

5m

4.5 Aspirate and discard the upper aqueous phase. Do not disturb the protein disc at the interface of the aqueous phase (top) and organic phase (bottom).

4.6 Add 4x volumes of MeOH and vortex.

4.7 Spin down at  Room temperature for  00:05:00 at high speed.

5m



4.8 Aspirate and discard supernatant. Do not disturb the protein pellet at the bottom of the tube.

4.9 Repeat MeOH wash.



4.10 Air dry (or speed-vac) protein pellet down to remove all traces of MeOH.





5 Resuspend protein pellets  100 μL of  200 millimolar (mM) EPPS buffer ( 8.5).



6 Digest samples at  37 $^{\circ}\text{C}$ for  02:00:00 with endoproteinase Lys-C (Wako, Japan) at a 1/200 enzyme/protein ratio. 2h

7 Digest with Trypsin (1:100) for  06:00:00 at  37 $^{\circ}\text{C}$. 6h




TMT-Labeling of samples

1h 15m

8 Add  5 μL of the TMT reagent to each sample. Solubilize TMT reagents are in ACN as per manufacturer's instructions and  5 μL of TMT label is used for every  50 μg of protein lysate . Perform TMT labeling in a final concentration of 20-25% ACN. Add  20 μL of ACN to bring the reaction volume at 125 μL . The number of samples, and hence the number of individual TMT reagents, will depend upon the design of the experiment.



9 Incubate for  01:00:00 at  Room temperature to label the samples. 1h




10 Combine 1% of each labeled sample together in a tube, quench the reaction with  4 μL of 5% hydroxylamine for  00:15:00 at  Room temperature , and dry down using the speed-vac. This combined sample is used to perform the ratio check to test labeling efficiency. The remaining amount of each sample can be stored in the freezer. 15m


Stage tip for Ratio Check

15m

11 Resuspend the dried sample in  100 μ L of 5% FA/5% ACN . Check to ensure that the pH of the sample is ~  3 (or lower) using pH strip.

12 Make stage tip by placing 6-8 “cookies” of C-18 embedded membranes in  200 μ L Rainin tip .

13 Perform C-18 cleanup:

13.1 Equilibrate C-18 with  100 μ L of 100% methanol .


13.2 Wash C-18 with  50 μ L of 50% ACN/5% FA .



13.3 Wash C-18 with  100 μ L of 5% ACN/5% FA .



13.4 Load sample on to C-18 to bind peptides.

13.5 Wash bound peptides on C-18 with  50 μ L of 5% ACN/5% FA .



13.6

Elute peptides off C-18 with

50 μ L 75% ACN/ 5%
FA

into a mass spec vial.

14

Dry down eluted peptides in speed-vac.

15

Reconstitute peptides in

10 μ L 5% ACN/5%
FA

16

Perform ratio check by analyzing the total amount of reporter ions present, as measured by mass spectrometry, for each TMT reporter ion channel.

17

Quench the entire volume of each sample using

8 μ L of 5%
hydroxylamine

for

15m

00:15:00 at Room temperature. Combine samples in 2 ml Eppendorf tube to equal amounts, based on normalization values you obtain from the ratio check.

18

Dry down labeled, combined sample in speed-vac.

19

Re-constitute the sample in






750 μ L of 5% ACN/5%
FA


SepPak clean-up (C18 solid-phase extraction (SPE))

15m

20

Place SepPak column into vacuum slot on a vacuum manifold.



- 21 Fill SepPak with  1 mL 100% MeOH.
- 22 Start the vacuum, gently, try to ensure that the pressure gauge is below 10. This can be achieved by opening one of the valves gently. The fluid should pass through slowly, drop wise.
- 23 Fill SepPak with  1 mL 5% FA/50% ACN, Flow Liquid through.
- 24 Fill SepPak with  1 mL 5% FA/5% ACN, Flow Liquid through, repeat 2 more times.
- 25 Replace tube with 2 ml collection tube.
- 26 Add Peptides, Flow Liquid through (~ 750 µl).
- 27 Wash with  1 mL 5% FA/5% ACN, Flow Liquid through, 2 times.
-  28 Replace 15 ml collection tube with 2 ml collection tube.

29 Elute with  750 μ L 75% ACN/5% FA. Since there is residual liquid left in the SepPak, ensure that all the liquid flows through SepPak.


30 Dry down in SpeedVac  Room temperature  Overnight or  30 °C. 15m


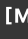




Fractionation

1h 30m


31 Resuspend sample in  100 μ L of 10 mM NH_4HCO_3  8.

32 Fractionate using pH reverse-phase HPLC:


32.1 Fractionate samples by high-pH reverse-phase HPLC (Agilent LC1260) into 96 fractions over  01:30:00 run. 1h 30m

32.2 Fractions are run through an Aeris peptide XB-c18 column (Phenomenex; 250 mm x 3.6 mm), with mobile phase A containing  5 Mass Percent ACN and  10 millimolar (mM) NH_4HCO_3  8 and phase B containing  90 Mass Percent ACN and  10 millimolar (mM) NH_4HCO_3  8 (all in LC-MS grade H_2O).

33 Combine fractionated samples (either 12 or 24 fractions) in a non-continuous manner into individual 1.5 ml Eppendorf tube (see outline below from [Paulo et al., 2016](https://doi.org/10.1002/prot.24016)).

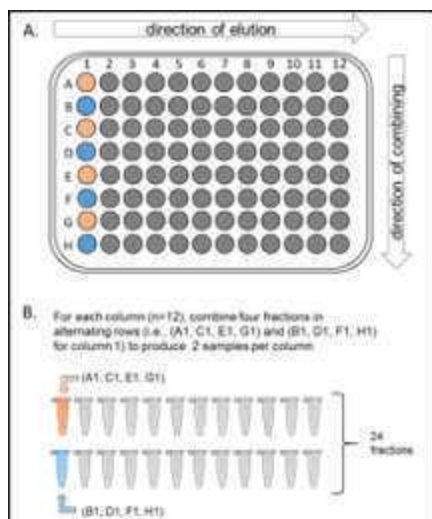

34 Dry down in SpeedVac.

35 Resuspend peptides in  100 μ L 5%FA /5% ACN

36 Check pH (\sim 3.5) with pH indicator strips.


Stage tip for proteomics sample

3m

37 Stage tip each fraction.

37.1

Make stage tips and equilibrate. Spin down at

 3000 rpm,
00:03:00


3m

37.2 Perform C-18 cleanup:

37.3




Wash with

 50 μ L 50% ACN /5% FA

37.4



Wash with

 100 μ L 5% ACN /5% FA


37.5 Load sample


37.6 Collect flow through and freeze.

37.7




Wash with

 50 μ L 50% ACN /5% FA

37.8 Elute with  50 μ L 75% ACN/1% FA in mass spec vial.

38 Dry down in SpeedVac.

39 Reconstitute pellet in  10 μ L 5% FA / 5% ACN .

40 Freeze sample at -  -20 °C until ready to run proteomics.