

Version 4

Sep 11, 2020

# covid 19 indirect detection thru rise of 100nm filter fluid resistance V.4

Ladislav Nevery

In Development

[dx.doi.org/10.17504/protocols.io.bmazk2f6](https://dx.doi.org/10.17504/protocols.io.bmazk2f6)

Coronavirus Method Development Community

XPRIZE Rapid Covid Testing

1 more workspace



Ladislav Nevery

## ABSTRACT

By passing saliva sample thru 150nm filter to remove particles larger than avg 125nm sized covid 19 virus particles and concentrating resulting fluid to very small 1 mm<sup>2</sup> area of 110nm filter. Pores in this small area should be clogged up fast raising filter air/fluid resistance. Due to high virus density 5.2 log<sub>10</sub> ml saliva of infected person should clogg filter much faster than healthy person. Test is setup in such way that we let gravity pass 2l of fluid thru clogged filter and record rate of drops from its bottom as audio on mobile phone placed bellow glass cup. rate of drops thru clogged filter should be measurably slower than clean filter.

## THIS PROTOCOL ACCOMPANIES THE FOLLOWING PUBLICATION

[dx.doi.org/10.17504/protocols.io.bkznkx5e](https://dx.doi.org/10.17504/protocols.io.bkznkx5e)

## ATTACHMENTS

[DropRateV2.zip](#)

## DOI

[dx.doi.org/10.17504/protocols.io.bmazk2f6](https://dx.doi.org/10.17504/protocols.io.bmazk2f6)

## PROTOCOL CITATION

Ladislav Nevery 2020. covid 19 indirect detection thru rise of 100nm filter fluid resistance. **protocols.io**  
<https://dx.doi.org/10.17504/protocols.io.bmazk2f6>  
 Version created by [Ladislav Nevery](#)

## MANUSCRIPT CITATION please remember to cite the following publication along with this protocol

[dx.doi.org/10.17504/protocols.io.bkznkx5e](https://dx.doi.org/10.17504/protocols.io.bkznkx5e)

## LICENSE

———— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

## CREATED

Sep 11, 2020

## LAST MODIFIED

Sep 11, 2020

## PROTOCOL INTEGER ID

42041

## MATERIALS

NAME

CATALOG #

VENDOR

NAME	CATALOG #	VENDOR
<a href="#">pragopor 9</a>		
<a href="#">pragopor 10</a>		
<a href="#">120nm nanoparticles</a>		
<a href="#">25mm x2 rubber gasket</a>		

#### STEPS MATERIALS


NAME	CATALOG #	VENDOR
<a href="#">3d print filament 8m</a>		
<a href="#">pragopor 10</a>		<a href="#">pragochema</a>
<a href="#">2l soda bottle</a>		
<a href="#">glass cup</a>		
<a href="#">mobile phone</a>		
<a href="#">120nm nanoparticles</a>		
<a href="#">pragopor 9</a>		
<a href="#">25mm x2 rubber gasket</a>		

#### DISCLAIMER:

DISCLAIMER – FOR INFORMATIONAL PURPOSES ONLY; USE AT YOUR OWN RISK

The protocol content here is for informational purposes only and does not constitute legal, medical, clinical, or safety advice, or otherwise; content added to [protocols.io](#) is not peer reviewed and may not have undergone a formal approval of any kind. Information presented in this protocol should not substitute for independent professional judgment, advice, diagnosis, or treatment. Any action you take or refrain from taking using or relying upon the information presented here is strictly at your own risk. You agree that neither the Company nor any of the authors, contributors, administrators, or anyone else associated with [protocols.io](#), can be held responsible for your use of the information contained in or linked to this protocol or any of our Sites/Apps and Services.

1

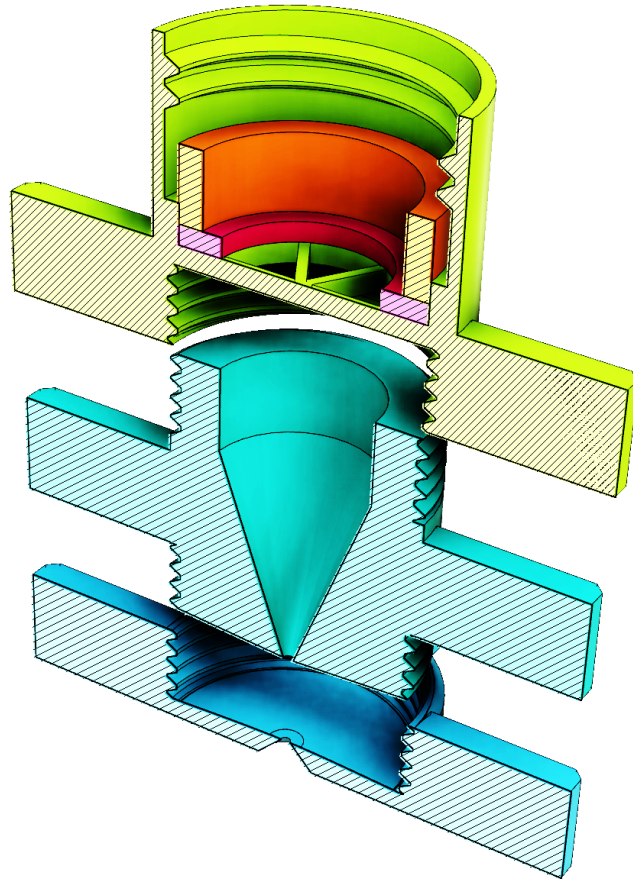

3d print all required parts



3d print filament 8m

2 place 25mm ideally 110nm pc or nylon membrane filter between blue parts. I used what I had



pragopor 10

by pragochema

[View](#)




0.24eur piece

3 place 25mm ideally 150nm membrane filter in green saliva holder bellow red rubber gasket. I used what I had




pragopor 9

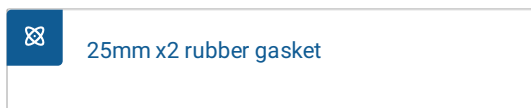

**0.24eur** piece

- 4 .attach empty 2L soda bottle and squeeze until all saliva passes thru both filters

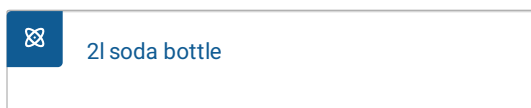


- 5 remove green part remove red parts and filter. clean green part

- 6 fill and reinstall bottle with 2l of clean fluid

- 7 place whole assembly with bottle on top over any glass cup capable of holding blue part



. place mobile phone bellow glass cup and start audio recording to count and record sounds of drops.

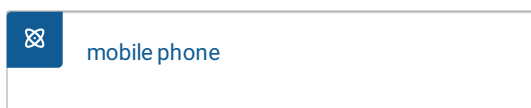

- 8 stop audio recording. fill sample holder in step 4 with



fluid containin 120nm test nanoparticles in concentration resembling covid19 as  
5.2log10ml and repeat whole process once more with new filters and clean parts → **go to step #2**

- 9 compare distance between drop sounds in all recordings.



drop frequency with 120nm nanoparticle clogged filter should be measurably different