



MAR 03, 2023

OPEN  ACCESS

**DOI:**  
[dx.doi.org/10.17504/protocols.io.b2c8qazw](https://dx.doi.org/10.17504/protocols.io.b2c8qazw)

**External link:**  
<https://doi.org/10.1016/j.neuron.2021.12.031>

**Protocol Citation:** Daehun Park, Pietro De Camilli 2023. Surface protein biotinylation. **protocols.io**  
<https://dx.doi.org/10.17504/protocols.io.b2c8qazw>

**MANUSCRIPT CITATION:**  
 Yang S, Park D, Manning L, Hill SE, Cao M, Xuan Z, Gonzalez I, Dong Y, Clark B, Shao L, Okeke I, Almoril-Porras A, Bai J, Camilli PD, Colón-Ramos DA, Presynaptic autophagy is coupled to the synaptic vesicle cycle via ATG-9. *Neuron* 110(5). doi: [10.1016/j.neuron.2021.12.031](https://doi.org/10.1016/j.neuron.2021.12.031)

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited


**Protocol status:** Working

## Surface protein biotinylation

Daehun Park<sup>1,2</sup>, Pietro De Camilli<sup>1,2</sup>

<sup>1</sup>Departments of Neuroscience and of Cell Biology, Howard Hughes Medical Institute, Program in Cellular Neuroscience, Neurodegeneration and Repair, Yale University School of Medicine, New Haven, Connecticut 06510, USA;

<sup>2</sup>Aligning Science Across Parkinson's (ASAP) Collaborative Research Network, Chevy Chase, MD, 20815

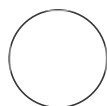

daehun.park

### ABSTRACT

This protocol describes surface protein labeling with biotin using EZ-Link Sulfo-NHSLC-Biotin. This chemical reacts with primary amines such as lysine but does not permeate cell membranes because of the charge. Thus, it only biotinylates surface proteins.

### ATTACHMENTS

[326 - 699.pdf](#)

**Created:** Nov 26, 2021

## MATERIALS

**Last Modified:** Mar 03, 2023

### Reagents

**PROTOCOL integer ID:**  
55424

**Keywords:** surface protein, biotinylation, EZ-Link Sulfo-NHS-LC-Biotin, ASAPCRN

#### 1. EZ-Link Sulfo-NHS-LC-Biotin

**[M]** 0.25 mg/mL EZ-Link Sulfo-NHS-LC-Biotin (Thermo) in ice-cold PBS. This solution should be made fresh just before use.

#### 2. Quenching solution

**[M]** 50 millimolar (mM) glycine in ice-cold PBS (keep at **4 °C**).

#### 3. 1% triton X-100 lysis buffer

A	B
Tris-HCl, pH 8	20 mM
Triton X-100	1%
Glycerol	10%
NaCl	137 mM
EDTA	2 mM

Keep at **4 °C** and add protease inhibitor cocktails just before use.

#### 4. 2x sample buffer

**[M]** 100 millimolar (mM) Tris, 4% SDS, 0.2% bromophenol blue, 20% glycerol in DW.

## Surface protein biotinylation

5h 5m

### 1 Wash cells with ice-cold PBS.

#### Note

**Note1.** There are many washing steps. Thus, cells may detach from the dishes. If this occurs, coat dishes with poly-D-lysine.

1.1 Wash cells with ice-cold PBS. (1/3)





1.2 Wash cells with ice-cold PBS. (2/3)



1.3 Wash cells with ice-cold PBS. (3/3)

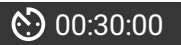


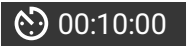
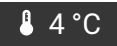
2 Incubate cells with EZ-Link Sulfo-NHS-Biotin for  00:30:00 at  4 °C to label surface proteins.

30m



#### Note

A rocking platform is recommended.

3 Discard biotin containing medium. Quench and remove unbound biotin using  50 millimolar (mM) glycine in ice-cold PBS for  00:10:00 at  4 °C .

10m

4 Wash 2-3 times with ice-cold PBS.

4.1 Wash with ice-cold PBS. (1/3)



4.2 Wash with ice-cold PBS. (2/3)



#### 4.3 Wash with ice-cold PBS. (3/3)



#### 5 Lysis cells with 1% triton X-100 lysis buffer and centrifuge the samples at

20m



14000 x g, 4°C, 00:20:00 .

#### 6 Collect supernatants and discard the pellets.

#### 7 Measure the protein concentrations using BCA Protein Assay Kit.

#### 8 Incubate the same amount of lysates ( 500 µg to 1000 µg ) with streptavidin or NeutraAvidin beads for 02:00:00 to Overnight at 4 °C to pull-down the biotinylated proteins.

4h



#### 9 Wash the beads with lysis buffer by cycles of suspension.



##### 9.1 Wash the beads with lysis buffer by cycles of suspension. (1/3)





##### 9.2 Wash the beads with lysis buffer by cycles of suspension. (2/3)



### 9.3 Wash the beads with lysis buffer by cycles of suspension. (3/3)



**10** Centrifugation and elute proteins from the packed beads by adding an equal volume of 2x sample buffer and boiling for  00:05:00  95 °C . 5m



**11** Run eluate samples on a SDS polyacrylamide gel and perform western blotting to visualize labeling (and thus evidence of surface expression) of the protein of interest. Total cell lysates can be used to determine the expression level.