



Apr 18, 2022

Purification of cytosolic fraction and quantification of mtDNA by qPCR

Will Hancock-Cerutti^{1,2,3}, Zheng Wu^{4,5}, Gerald S. Shadel⁵,
Pietro De Camilli^{1,3}

¹Departments of Neuroscience and of Cell Biology, Howard Hughes Medical Institute, Program in Cellular Neuroscience, Neurodegeneration and Repair, Yale University School of Medicine, New Haven, Connecticut 06510, USA;

²Interdisciplinary Neuroscience Program and MD-PhD Program, Yale University School of Medicine, New Haven, Connecticut 06510, USA;

³Aligning Science Across Parkinson's (ASAP) Collaborative Research Network, Chevy Chase, MD, 20815;

⁴Department of Genetics, Yale University School of Medicine, New Haven, Connecticut 06510, USA;

⁵Salk Institute for Biological Studies, La Jolla, CA, USA

1



dx.doi.org/10.17504/protocols.io.14egnz12yg5d/v1

William Hancock-Cerutti

This protocol describes the purification of a cytosolic fraction depleted of membrane from cultured cells, and the quantification of mitochondrial DNA (mtDNA) in this fraction by qPCR.

[dxwwbkjz7.pdf](#)

DOI

dx.doi.org/10.17504/protocols.io.14egnz12yg5d/v1

Will Hancock-Cerutti, Zheng Wu, Gerald S. Shadel, Pietro De Camilli 2022.
Purification of cytosolic fraction and quantification of mtDNA by qPCR.

protocols.io

<https://dx.doi.org/10.17504/protocols.io.14egnz12yg5d/v1>



Mitochondria DNA (mtDNA), Cytosol fractionation, qPCR, ASAPCRN

protocol ,

Jul 08, 2021

Apr 18, 2022

Jul 08, 2021

Urmilas

Jul 14, 2021

William Hancock-Cerutti

51399

Solutions to prepare:

DMEM solution:

A	B
FBS	10%
Penicillin	100 U/ml
Streptomycin	100 mg/ml
L-glutamine	2 mM


Cytoplasmic buffer:












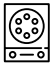












A	B
NaCl	150 mM
HEPES	50 mM
Digitonin, pH 7.4	1 mg/mL

Lysis buffer:

A	B
NaCl	150 mM
HEPES	50 mM
Digitonin, pH 7.4	1 mg/mL
SDS supplemented with Protease Inhibitor Cocktail (Roche)	1%

Cell culture and purification of cytosolic fraction

- 1 Plate HeLa cells in DMEM  15 cm plates (3.5×10^6 cells per plate).
- 2 The following day, prepare cytosolic buffer with fresh digitonin.
- 3 Prepare lysis buffer.

- 4  5m
- Trypsinize cells and centrifuge at  **1500 rpm** for  **00:05:00** at  **22 °C** .
- 5 Resuspend cells in PBS and count cells.
- 6  5m
- Collect the same number of cells from each genotype (5×10^6) and centrifuge at  **1500 rpm** for  **00:05:00** at  **22 °C** .
- 7 Resuspend cells in  **1 mL** PBS and transfer  **50 μ L** to a prechilled Eppendorf another (WCE). Keep  **On ice** .
- 8  5m
- Transfer the remaining  **950 μ L** to a prechilled Eppendorf and centrifuge at  **4500 rpm** for  **00:05:00** at  **22 °C** .
- 9 Remove supernatant and resuspend pellet in  **500 μ L** cytosolic buffer.
- 10 Rotate for  **00:10:00** at  **4 °C** . 10m
- 11  3m
- Centrifuge extract at  **980 x g** for  **00:03:00** at  **4 °C** . Transfer supernatant to new Eppendorf tube and save pellet for analysis (Pel).
-  10m

12 



Centrifuge supernatant at  **17000 x g** for  **00:10:00** at  **4 °C**.

13 Collect supernatant (Cyt).




14 Purify DNA from WCE and Cyt fractions using DNeasy Kit (Qiagen).

15 Measure DNA concentration and dilute samples 1:10.

qPCR

16 Combine  **10 µL** SYBR Green Master Mix (BioRad) with  **6.78 µL** Sterile Water (American Bio) per sample.

17 

Combine  **16.78 µL** diluted SYBR Green Master Mix with  **0.61 µL** each of  **10 Micromolar (µM)** forward and reverse primers per sample. Pipette this mixture into wells of 96-well qPCR plate. Perform at least two technical replicates for each sample.

18 

Pipette  **2 µL** of diluted DNA from step 15 in well with SYBR Green Master Mix.

19 Cover plate with Optical Adhesive Covers (Applied Biosystems).

20 

Spin down plate in table top centrifuge.

21

Run qPCR in CFX96 Real-Time System (BioRad) using the following protocol:

A	B	C
95 °C	3 min	
95 °C	10 sec	Repeat 39x
55 °C	10 sec	
72 °C	30 sec	
95 °C	10 sec	
65 °C	5 sec	
95 °C	5 sec	

Data analysis

- 22 Subtract the nuclear gene (hB2M) mean threshold cycle (Ct) values from WCE samples from mtDNA amplicon of interest mean Ct values from Cyt samples to calculate ΔC_t .
- 23 Subtract the ΔC_t of the control sample from each sample ΔC_t to calculate the $\Delta\Delta C_t$ value.
- 24 Calculate relative expression using the $2^{-\Delta\Delta C_t}$ method.
- 25 WT mtDNA abundance was given a value of 1.