



Jun 17, 2022

Soil Sample (Laboratory, English)

Hsin-Mao Wu¹¹Kaohsiung Medical College

1

dx.doi.org/10.17504/protocols.io.14egn7rezv5d/v1

Hsin-Mao Wu

This protocol is for 1 sample site

DOI

dx.doi.org/10.17504/protocols.io.14egn7rezv5d/v1

Hsin-Mao Wu 2022. Soil Sample (Laboratory, English). **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.14egn7rezv5d/v1>



Soil sample

protocol ,


Jun 17, 2022

Jun 17, 2022

64803

This protocol is for 1 sample site

- 1 Shovel
- 2 HDPE plastic bottles
- 1 Desiccant
- 1 Recode paper
- 1 Plastic zipper bag
- 1 Insulated bag

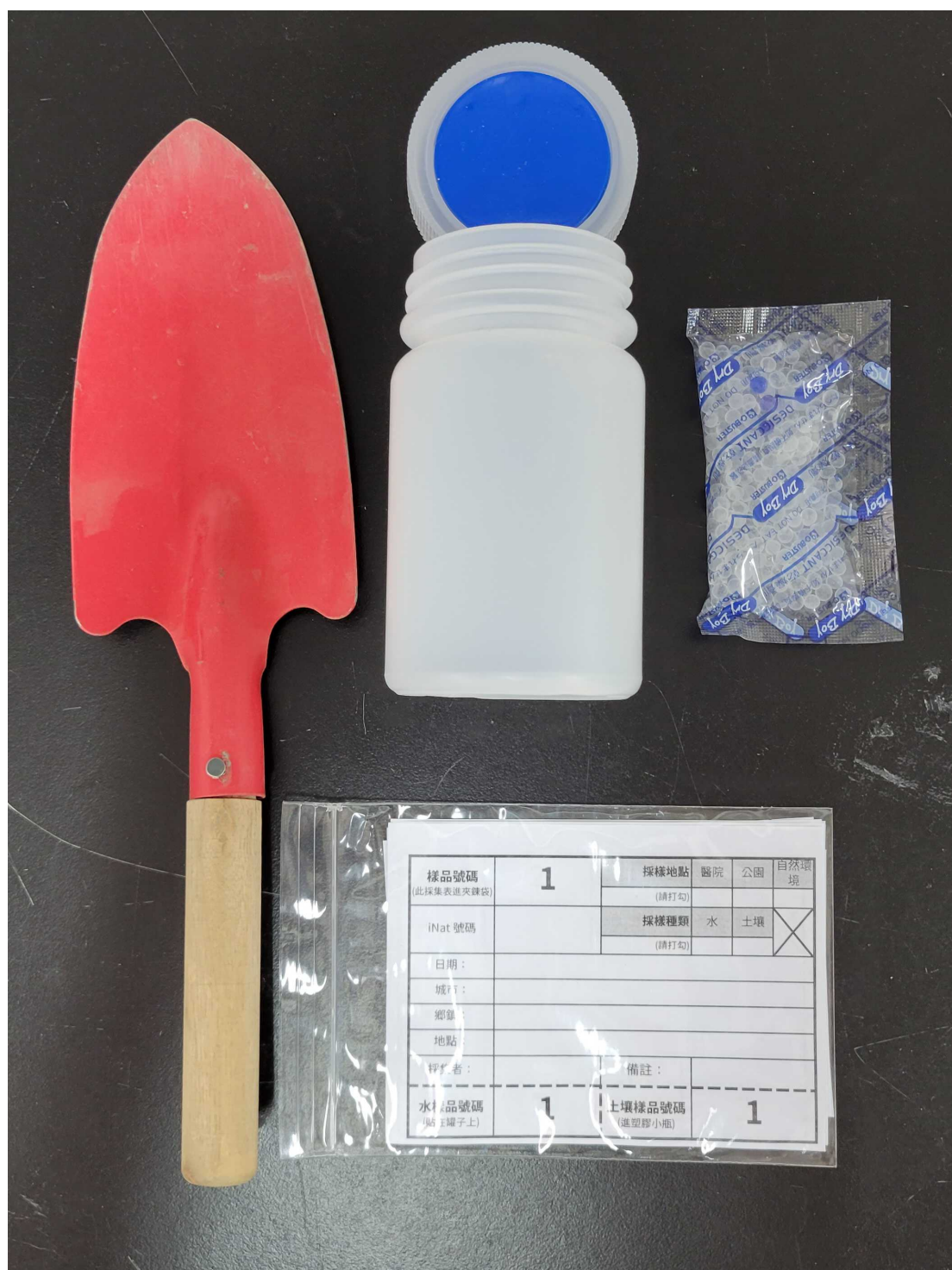

Materials (shovel, HDPE plastic bottles, desiccant, recode papers, plastic zipper bag)



Materials (insulated bag)

Soil sampling

- 1 Sample sites have 3 different types

A
Hospitals
Parks
Natural region

- 2 In one sampling site, choose 5 places as the average
- 3 Remove the topsoil, take the soil at a distance of 5-7 cm from the soil surface, and put it into a plastic zipper bag
- 4 Mix the soil evenly, and then add soil into 2 separate HDPE plastic bottles

1 bottle with desiccant called "RT" (meaning room temperature); 1 bottle called "Cold"

- 5 Put the "Cold" one into the insulated bag; put the "RT" into a backpack
- 6 Take 2 pictures:

A	B
1st	Take a picture at the sampling site
2nd	Take another picture at sampling site, and put the recode paper beside