

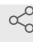


Aug 01, 2022

# Mito-Keima assay to assess mitophagy

Thanh Ngoc Nguyen<sup>1</sup><sup>1</sup>Laboratory of Michael Lazarou, Walter and Eliza Hall Institute of Medical Research, Parkville, Victoria, Australia.Thanh Ngoc Nguyen: [nguyen.tha@wehi.edu.au](mailto:nguyen.tha@wehi.edu.au)

1 Works for me

 Share[dx.doi.org/10.17504/protocols.io.q26g74e1qgwz/v1](https://dx.doi.org/10.17504/protocols.io.q26g74e1qgwz/v1) [nguyen.tha](#)

## ABSTRACT

This protocol details the procedure of mito-keima assay to assess mitophagy.

## ATTACHMENTS

[ihdmbr7hf.docx](#)

## DOI

[dx.doi.org/10.17504/protocols.io.q26g74e1qgwz/v1](https://dx.doi.org/10.17504/protocols.io.q26g74e1qgwz/v1)

## PROTOCOL CITATION

Thanh Ngoc Nguyen 2022. Mito-Keima assay to assess mitophagy. **protocols.io**  
<https://dx.doi.org/10.17504/protocols.io.q26g74e1qgwz/v1>



## KEYWORDS

Mito-Keima, Mitophagy, BFP-Parkin

## LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited


## CREATED

Jul 19, 2022

## LAST MODIFIED

Aug 01, 2022

## OWNERSHIP HISTORY

Jul 19, 2022  madhavi.d

Jul 29, 2022  nguyen.tha

## PROTOCOL INTEGER ID

67007

## MATERIALS TEXT

### Growth media:

#### DMEM:

A	B
FBS	10%
Glucose	4.5 g/l
GlutaMAX™	1x
MEM NEAA	1x
HEPES	25 mM

 45% D-(-)-

-  Glucose **Sigma Catalog #G8769**

 GlutaMAX™ Supplement **Thermo Fisher**

- **Scientific Catalog #35050061**

▪

 MEAA (MEM Non-Essential Amino Acids) **Gibco - Thermo**

**Fischer Catalog #11140050**

 Antimycin A from Streptomyces sp. **Sigma –**

- **Aldrich Catalog #A8674**

(made up in

100% Ethanol to 20 mg/ml)

- Oligomycin (Calbiochem, 495455; made up in DMSO to 10 mg/ml) and

-  qVD **MedChemExpress Catalog #HY-12305** (made up in DMSO to 10 mM)

- 1x PBS

- **FACS buffer:** 10% FBS, **[M]0.5 millimolar (mM)** EDTA in 1xPBS

## Procedure

2h 9m

- 1 Seed the HeLa cells the day before the treatment day in 24 well plates.

Each well contained 0.5 ml of growth media; 120,000 cells were seeded for penta KO

expressing BFP-Parkin, mito-Keima (mtKeima) and GFP-OPTN or -NDP52; the number of cells of other cell lines were adjusted so that the next day they are all in similar confluency with penta KO expressing BFP-Parkin, mtKeima and GFP-OPTN or -NDP52.

- 1.1 Control cells include unstained cells (without any fluorescence proteins), cells that only express GFP-OPTN or GFP-NDP52 (to set up compensation so that GFP signal doesn't bleach into mtKeima signal).

- 2 The next day, make sure the seeded cells are spreading out.

Not concentrated in the middle of the well because this can affect the results.



Aspirate off the old media and treat each well with **0.5 mL** of growth media containing **4 micromolar ( $\mu\text{M}$ )** Antimycin A, **10 micromolar ( $\mu\text{M}$ )** Oligomycin and **10 micromolar ( $\mu\text{M}$ )** QVD for indicated times.

Make sure all drugs are vortexed well, mix the media well after adding each drug.

- 3.1 Treat the longest time points first.



2 hours prior to harvesting, feed the untreated wells with **1 mL** of warm growth media.

- 5 After treatment, harvest the cells by trypsinisation.

- 5.1 Pre-chill eppies **On ice**.

5.2 Aspirate the media thoroughly from the wells.





5.3 

Wash the wells with  **0.5 mL** of 1x PBS.

Make sure swirl around after adding the PBS to wash the cells properly.

5.4   

7m

Aspirate 1x PBS and add  **150 µL** of trypsin and incubate at  **37 °C** for  **00:02:00** -  **00:05:00** .

Check under microscope to make sure all the cells were trypsinised properly.

5.5 

Add  **300 µL** of growth media to each well.

5.6 



Mix well with a P1000 and transfer the cells to cold eppies.

5.7 

2m

Centrifuge the eppies at  **1000 x g** for  **00:02:00** at  **4 °C** .

5.8 Aspirate off the liquid.

- 6 Resuspend the cell pellets in  **150 µL** of FACS buffer, transfer to pre-chilled FACS tubes and keep the samples  **On ice** in the dark for analysis.