

VERSION 2

FEB 23, 2023

OPEN ACCESS

DOI:

dx.doi.org/10.17504/protocols.io.kxyg9edzg8j/v2

Protocol Citation: Mirsada Causevic, Edin Begic 2023. Characterization of the VKORC1 and CYP2C9 genotypes. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.kxyg9edzg8j/v2> Version created by [mirsadacausevic](#)

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Feb 23, 2023

Last Modified: Feb 23, 2023


PROTOCOL integer ID:
77502

Characterization of the VKORC1 and CYP2C9 genotypes V.2

Mirsada Causevic¹, Edin Begic^{2,1}

¹Sarajevo Medical School, Sarajevo School of Science and Technology, Sarajevo, Bosnia and Herzegovina;

²General Hospital "Prim.Dr. Abdulah Nakas" Sarajevo, Bosnia and Herzegovina

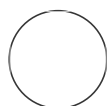

mirsadacausevic

ABSTRACT

Vitamin K antagonists are anticoagulants which represent widely prescribed drugs for prevention and treatment of thromboembolic disorders.

The molecular target of these drugs is vitamin K epoxide reductase enzyme and their metabolism is performed by the cytochrome P-450 2C9 enzyme, both of which are encoded by polymorphic genes - vitamin K epoxide reductase complex subunit 1 (*VKORC1*) and cytochrome P-450 2C9 (*CYP2C9*), respectively.

Identification of *VKORC1* and *CYP2C9* genotypes was recommended by the U.S. Food and Drug Administration in 2007, in order to guide the initial dosing of warfarin, the first oral vitamin K antagonist drug, to achieve optimum anticoagulation and prevent hemorrhagic events.

This protocol is designed for undergraduate students of pharmacology or medicine, or anyone receiving training in pharmacogenomics/pharmacogenetics and understanding a relationship between a gene and drug response for the first time.

The protocol represents a step-by-step guide to using polymerase chain reaction restriction fragment length polymorphism (PCR-RFLP) method in determining *VKORC1* -1639G>A, *CYP2C9**1/*1, *CYP2C9**1/*2, *CYP2C9**1/*3, and *CYP2C9**2/*3 genotypes.

Based on the combination of the *VKORC1* -1639G>A and *CYP2C9* genotypes, that is, the pattern of DNA fragments, which are visualized using agarose gel electrophoresis and DNA staining, a conclusion on an individual's sensitivity to vitamin K antagonists and appropriate drug dose is reached.

Keywords: Vitamin K Epoxide Reductase, Vitamin KO Reductase, VKORC1, Cytochrome P-450, CYP2C9

IMAGE ATTRIBUTION

Mirsada Causevic, Honours BSc, PhD, Sarajevo Medical School, Sarajevo School of Science and Technology, Sarajevo, Bosnia and Herzegovina

GUIDELINES

This protocol uses human whole blood as a starting material for extracting human, genomic DNA. Human whole blood could be a source of hepatitis B (HepB) virus infection (if taken from persons with HepB virus infection) to the individuals who work with human blood and, as such, represents health and safety or biological hazard.

MATERIALS

Agarose - Sigma A9539-100G;
Avall restriction endonuclease enzyme - New England BioLabs R0153S;
Bromophenol blue - Sigma B0126-25G;
Boric acid - Sigma B6768-500G;
DNA dye - Ethidium bromide - Sigma E1510-10ML;
DNA ladder (1) - BenchTop 100bp DNA Ladder - Promega G8291;
DNA ladder (2) - MiniSizer 50bp DNA Ladder - Norgen 11200;
dNTPs - dATP, dCTP, dGTP and dTTP - at a stock concentration of 100mM each - Promega U1330;
EDTA - Sigma EDS-500G;
Ethanol - Sigma 32205-1L;
Glycerol - Sigma G5516-500ml;
HCl - Sigma 30721-2.5L;
Isopropanol - 2-propanol - Sigma 59304-500ml;
KpnI restriction endonuclease enzyme - Promega R6341;
NaCl - Honeywell-Fluka 31434-1KG;
NsiI restriction endonuclease enzyme - New England BioLabs R0127S;
MspI restriction endonuclease enzyme - New England BioLabs R0106S;
Proteinase K - Sigma P2308;
SDS - Sigma L3771-100G;
GoTaq G2 DNA Polymerase - Promega M7845;
Trizma base - Sigma RDD008-1KG.

SAFETY WARNINGS

! Ethidium bromide, which is used during the preparation of agarose gels, is a DNA intercalating dye and it can, therefore, act as a mutagen. It should be handled with care and only when wearing gloves and safety glasses.

BEFORE START INSTRUCTIONS

In order to protect the individuals from the potential HepB virus infection while working with human blood, [the U.S.A. Centers for Disease Control and Prevention \(CDC\) issued a recommendation](#) for HepB vaccination to be carried out in individuals before they start working with human blood. Vaccine-induced antibody titer to the HepB virus of ≥ 10 mIU/mL should be present in individuals before commencing work with human blood.

Genomic DNA extraction

- 1 Patients' whole blood was collected in ethylenediaminetetraacetic acid (EDTA)-containing tubes and stored at -20°C until use. Genomic DNA extraction from the human whole blood, that is, leukocytes, was carried out according to the protocol described by [Subbarayan PR and colleagues](#) (doi: 10.2144/02336bm10), with modifications.

Upon thawing of the blood, for every patient's blood sample, five (5) sterile tubes were labeled and 300 μl of the whole blood was added. This was the amount of the starting material that yielded ample amount of human, genomic DNA at the end of the DNA extraction protocol.

- 2 Concentrations of human, genomic DNA was determined by using Qubit 4 Fluorometer (Q33226, Invitrogen and Thermo Fisher Scientific).

Subsequently, for polymerase chain reactions (PCRs) 100 ng of genomic DNA was used in every PCR reaction.

Identification of the *VKORC1* -1639G>A genotypes

- 3 Identification of the *VKORC1* -1639G>A genotypes was performed according to the polymerase chain reaction-restriction fragment length polymorphism (PCR-RFLP) method-based protocol published by [Sconce EA and colleagues](#) (doi: 10.1182/blood-2005-03-1108).

- 4 PCR reactions were assembled as follows: reactions of 25 μl contained 1) sterile, milliQ- H_2O ; 2) 1xPCR buffer (for GoTaq G2 DNA polymerase, Promega), from a stock of 5xPCR buffer; 3) 0.2 mM deoxyribonucleotide triphosphates (dNTPs) (U1330, Promega), from a stock of 2.5 mM; 4) 1 μM forward primer, from a stock of 25 μM ; 5) 1 μM reverse primer, from a stock of 25 μM ; 6) 100 ng of genomic DNA and 7) 2.5 U GoTaq DNA polymerase (M7845, Promega).

Sequences of the PCR primers were as follows: 1) Forward primer: 5'-GCCAGCAGGAGAGGGAAATA-3'; 2) Reverse primer: 5'-AGTTTGGACTACAGGTGCCT-3'.

PCR cycles used for the amplification of the *VKORC1* gene promoter, which were based on the touchdown PCR protocol published by [Korbie and Mattick in 2008](#), were as follows:

- 1) 95°C for 3 min (1 cycle);
- 2) 95°C for 2 s, 60°C for 2 s, 72°C for 30 s (2 cycles);
- 3) 95°C for 2 s, 59°C for 2 s, 72°C for 30 s (2 cycles);
- 4) 95°C for 2 s, 58°C for 2 s, 72°C for 30 s (2 cycles);
- 5) 95°C for 2s, 57°C for 2s, 72°C for 30 s (2 cycles);
- 6) 95°C for 2 s, 56°C for 2 s, 72°C for 30 s (2 cycles);
- 7) 95°C for 2 s, 55°C for 2 s, 72°C for 30 s (2 cycles);
- 8) 72°C for 7 min (1 cycle);
- 9) 4°C for 2 h (1 cycle).

Expected result

The expected size of the PCR product is 290 bp.

- 5** In order to characterise the *VKORC1* -1639G>A genotypes by using the PCR-RFLP method, the following reactions were set up with MspI restriction endonuclease enzyme:

- 1) 25 µl of the PCR product;
- 2) 3 µl of the 10xCutSmart Buffer (B7204S, New England Biolabs);
- 3) 1 µl of sterile, milliQ-H₂O;
- 4) 1 µl of MspI enzyme (R0106S, New England Biolabs). Total reaction volume: 30 µl.

Reactions were incubated in a 37°C water bath overnight.

DNA fragments were examined by using 2% agarose gel electrophoresis with 1xTris, boric acid, EDTA (TBE) buffer and agarose gels containing 0.5 mg/ml ethidium bromide (E1510-10ML, Sigma). The results were viewed with a UV lamp (UVstar, Biometra, Analytik Jena). Gel images were captured by using BioDocAnalyze (BDA) camera (Biometra, Analytik Jena).

Expected result

The following pattern of DNA fragments can be expected after the *MspI* restriction endonuclease enzyme reaction:

- 1) GG or *VKORC1* -1639GG genotype is characterized by two (2) DNA fragments of 168 and 122 base pairs (bp);
- 2) AG or *VKORC1* -1639AG genotype is characterized by three (3) DNA fragments of 290, 168 and 122 bp;
- 3) AA or *VKORC1* -1639AA genotype is characterized by one (1) DNA fragment of 290 bp.

Identification of the *CYP2C9* genotypes

- 6 Identification of the *CYP2C9* genotypes (**1/*1*, **1/*2*, **1/*3*, and **2/*3*) was performed according to the PCR-RFLP method-based protocol published by [Sullivan-Klose TH and colleagues](#) (doi: 10.1097/00008571-199608000-00007) and [Yasar U and colleagues](#) (doi: 10.1006/bbrc.1998.9992).
 - 7 PCR reactions were assembled as follows: reactions of 25 µl contained 1) sterile, milliQ-H₂O; 2) 1xPCR buffer (for GoTaq G2 DNA polymerase, Promega), from a stock of 5xPCR buffer; 3) 0.2 mM deoxyribonucleotide triphosphates (dNTPs) (U1330, Promega), from a stock of 2.5 mM; 4) 1 µM forward primer, from a stock of 25 µM; 5) 1 µM reverse primer, from a stock of 25 µM; 6) 100 ng of genomic DNA and 7) 2.5 U GoTaq DNA polymerase (M7845, Promega).
 - 7.1 Sequences of the PCR primers for *Avall* restriction endonuclease enzyme reactions were as follows: 1) *Avall*-Forward primer: 5'-TACAAATACAATGAAAATATCATG-3'; 2) *Avall*-Reverse primer: 5'-CTAACAACCAGACTCATAATG-3'.
- PCR cycles used with the above primers were as follows:
- 1) 94°C for 5 min (1 cycle);
 - 2) 94°C for 60 s; 55°C for 90 s; 72°C for 30 s (35 cycles);
 - 3) 72°C for 7 min;

4) 4°C for 2 h (1 cycle).

Expected result

The expected size of the PCR product is 691 bp.

7.2 In order to characterise the *CYP2C9* *1/*1 or *1/*2 genotypes by using the PCR-RFLP method, the following reactions were set up with Avall restriction endonuclease enzyme:

- 1) 10 µl of the PCR product;
- 2) 2 µl of the 10xCutSmart Buffer (B7204S, New England Biolabs);
- 3) 7 µl of sterile, milliQ-H₂O;
- 4) 1 µl of Avall enzyme (R0153S, New England BioLabs). Total reaction volume: 20 µl.

Reactions were incubated in a 37°C water bath overnight.

DNA fragments were examined by using 2% agarose gel electrophoresis with 1xTBE buffer and agarose gels containing 0.5 mg/ml ethidium bromide (E1510-10ML, Sigma). The results were viewed with a UV lamp (UVstar, Biometra, Analytik Jena). Gel images were captured by using BioDocAnalyze (BDA) camera (Biometra, Analytik Jena).

Expected result

The following pattern of DNA fragments can be expected after the *Avall* restriction endonuclease enzyme reaction:

1) **1/*1* or *CYP2C9 *1/*1* genotype is characterized by two (2) DNA fragments of 527 and 164 base pairs (bp).

In addition, the **1/*1* genotype is characterized by one (1) DNA fragment, running between 100 and 200 bp markers when compared to a DNA ladder, which is the result of the *Nsil* restriction endonuclease enzyme reaction (see sub-steps 7.3 and 7.4).

Furthermore, the **1/*1* genotype is characterized by one (1) DNA fragment, running between 100 and 200 bp markers when compared to a DNA ladder, which is the result of the *KpnI* restriction endonuclease enzyme reaction (see sub-steps 7.5 and 7.6).

2) **1/*2* or *CYP2C9 *1/*2* genotype is characterized by three (3) DNA fragments of 691, 527 and 164 bp.

In addition, the **1/*1* genotype is characterized by one (1) DNA fragment, running between 100 and 200 bp markers when compared to a DNA ladder, which is the result of the *Nsil* restriction endonuclease enzyme reaction (see sub-steps 7.3 and 7.4).

Furthermore, the **1/*1* genotype is characterized by one (1) DNA fragment, running between 100 and 200 bp markers when compared to a DNA ladder, which is the result of the *KpnI* restriction endonuclease enzyme reaction (see sub-steps 7.5 and 7.6).

7.3

Sequences of the PCR primers for *Nsil* restriction endonuclease enzyme reactions were as follows: 1) *Nsil*-Forward primer: 5'-AATAATAATATGCACGAGGTCCAGAGATGC-3'; 2) *Nsil*-Reverse primer: 5'-GATACTATGAATTTGGGACTTC-3'.

PCR cycles used with the above primers were as follows:

- 1) 94°C for 5 min (1 cycle);
- 2) 94°C for 60 s; 60°C for 90 s; 72°C for 30 s (35 cycles);
- 3) 72°C for 7 min;
- 4) 4°C for 2 h (1 cycle).

Expected result

The expected PCR product is running between 100 and 200 bp markers when compared to a DNA ladder.

- 7.4** In order to characterise the *CYP2C9* *1/*3 or *2/*3 genotypes by using the PCR-RFLP method, the following reactions were set up with Nsil restriction endonuclease enzyme:
- 1) 10 µl of the PCR product;
 - 2) 2 µl of the 10xNEBuffer 3.1 (B7203S, New England Biolabs);
 - 3) 7 µl of sterile, milliQ-H₂O;
 - 4) 1 µl of Nsil enzyme (R0127S, New England BioLabs). Total reaction volume: 20 µl.

Reactions were incubated in a 37°C water bath overnight.

DNA fragments were examined by using 2% agarose gel electrophoresis with 1xTBE buffer and agarose gels containing 0.5 mg/ml ethidium bromide (E1510-10ML, Sigma). The results were viewed with a UV lamp (UVstar, Biometra, Analytik Jena). Gel images were captured by using BioDocAnalyze (BDA) camera (Biometra, Analytik Jena).

Expected result

The following pattern of DNA fragments can be expected after the Nsil restriction endonuclease enzyme reaction:

1) *1/*3 or *CYP2C9* *1/*3 genotype is characterized by two (2) DNA fragments of 527 and 164 bp that are achieved by Avall restriction endonuclease enzyme.

In addition, the *1/*3 genotype is characterized by two (2) DNA fragments, running between 100 and 200 bp markers when compared to a DNA ladder, which is the result of the Nsil restriction endonuclease enzyme reaction.

Furthermore, the *1/*3 genotype is characterized by two (2) DNA fragments, running between 100 and 200 bp markers when compared to a DNA ladder, which is the result of the KpnI restriction endonuclease enzyme reaction (see sub-steps 7.5 and 7.6).

2) *2/*3 or *CYP2C9* *2/*3 genotype is characterized by three (3) DNA fragments of 691, 527 and 164 bp that are achieved by Avall restriction endonuclease enzyme.

In addition, the *2/*3 genotype is characterized by two (2) DNA fragments, running between 100 and 200 bp markers when compared to a DNA ladder, which is the result of the Nsil restriction endonuclease enzyme reaction.

Furthermore, the *2/*3 genotype is characterized by two (2) DNA fragments, running between 100 and 200 bp markers when compared to a DNA ladder, which is the result of the KpnI restriction endonuclease enzyme reaction (see sub-steps 7.5 and 7.6).

- 7.5** Sequences of the PCR primers for KpnI restriction endonuclease enzyme reactions were as follows: 1) KpnI-Forward primer: 5'-AATAATAATATGCACGAGGTCCAGAGGTAC-3'; 2) KpnI-Reverse primer: 5'-GATACTATGAATTTGGGACTTC-3'.

PCR cycles used with the above primers were as follows:

- 1) 94°C for 5 min (1 cycle);
- 2) 94°C for 60 s; 60°C for 90 s; 72°C for 30 s (35 cycles);
- 3) 72°C for 7 min;
- 4) 4°C for 2 h (1 cycle).

Expected result

The expected PCR product is running between 100 and 200 bp markers when compared to a DNA ladder.

7.6 In order to characterise the *CYP2C9* *1/*3 or *2/*3 genotypes by using the PCR-RFLP method, the following reactions were set up with KpnI restriction endonuclease enzyme:

- 1) 10 µl of the PCR product;
- 2) 2 µl of the 10xBuffer J (R009A, Promega);
- 3) 0.5 µl Bovine Serum Albumin Acetylated, 10mg/ml (R3960, Promega);
- 4) 6.5 µl of sterile, milliQ-H₂O;
- 5) 1 µl of KpnI enzyme (R6341, Promega). Total reaction volume: 20 µl.

Reactions were incubated in a 37°C water bath overnight.

DNA fragments were examined by using 2% agarose gel electrophoresis with 1xTBE buffer and agarose gels containing 0.5 mg/ml ethidium bromide (E1510-10ML, Sigma). The results were viewed with a UV lamp (UVstar, Biometra, Analytik Jena). Gel images were captured by using BioDocAnalyze (BDA) camera (Biometra, Analytik Jena).

Expected result

The following pattern of DNA fragments can be expected after the KpnI restriction endonuclease enzyme reaction:

1) **1/*3* or *CYP2C9 *1/*3* genotype is characterized by two (2) DNA fragments of 527 and 164 bp that are achieved by *Ava*I restriction endonuclease enzyme.

In addition, the **1/*3* genotype is characterized by two (2) DNA fragments, running between 100 and 200 bp markers when compared to a DNA ladder, which is the result of the *Nsi*I restriction endonuclease enzyme reaction (see sub-steps 7.3 and 7.4).

Furthermore, the **1/*3* genotype is characterized by two (2) DNA fragments, running between 100 and 200 bp markers when compared to a DNA ladder, which is the result of the KpnI restriction endonuclease enzyme reaction.

2) **2/*3* or *CYP2C9 *2/*3* genotype is characterized by three (3) DNA fragments of 691, 527 and 164 bp that are achieved by *Ava*I restriction endonuclease enzyme.

In addition, the **2/*3* genotype is characterized by two (2) DNA fragments, running between 100 and 200 bp markers when compared to a DNA ladder, which is the result of the *Nsi*I restriction endonuclease enzyme reaction (see sub-steps 7.3 and 7.4).

Furthermore, the **2/*3* genotype is characterized by two (2) DNA fragments, running between 100 and 200 bp markers when compared to a DNA ladder, which is the result of the KpnI restriction endonuclease enzyme reaction.