

VERSION 3

MAR 03, 2023



DOI:

dx.doi.org/10.17504/protocol s.io.yxmvmk6k6g3p/v3

Protocol Citation: Rory Kruithoff, Lei Zhou, Douglas Shepherd 2023. Multiplexed Iterative FISH Experimental Protocol SOP002.v3.12.

protocols.io

https://dx.doi.org/10.17504/p rotocols.io.yxmvmk6k6g3p/v3 Version created by Rory Kruithoff

Multiplexed Iterative FISH Experimental Protocol SOP002.v3.12 V.3

Rory Kruithoff¹, Lei Zhou¹, Douglas Shepherd¹

¹Arizona State University

Human Cell Atlas Method Development Community



Rory Kruithoff
Arizona State University

ABSTRACT

This document, SOP002 - Multiplexed Iterative FISH Experimental Protocol, describes the process for in-situ fluorescence labeling of RNA transcripts in cells and tissues using a layered probe design, which allows for identity barcoding (MERFISH or similar). This protocol also provides the option for signal amplification using Branched DNA [bDNA] amplification. Iterative rounds of labeling and imaging are possible through the use of a readout probe with a cleavable disulfide (S-S) reporter molecule, a method that allows for minimal disruption to sample integrity between rounds. This document also describes cell and tissue preparation for RNA FISH as well as a system of mRNA anchoring using a specialized poly-t (locked nucleic acid, LNA) probe with an acrydite linker to bind mRNAs to a polyacrylamide (gel) matrix. Clearing and digestion techniques are used to reduce cellular autofluorescence and increase the signal to noise ratio of the final data. This protocol is *strongly* derived from Moffitt 2016 (https://doi.org/10.1016/bs.mie.2016.03.020) with some modifications with credit for a majority of this protocol due to Moffitt et al 2016.

ATTACHMENTS

SOP002.v3.12 Multiplexed Iterative FISH_Paintshop_Saber.pd

GUIDELINES

v3.12 revision notes

i. Minor revisions to language and formatting to make compatible with current protocol changes per current fluidics programming (flow speeds, solution priming). Also adjusted some formatting to improve flow and readability. ii. Updated tissue re-hydration time.

MATERIALS

MANUSCRIPT CITATION:

- Hershberg, E. A., Close, J. L., Camplisson, C. K., Attar, S., Chern, R., Liu, Y., ... & Beliveau, B. J. (2020). PaintSHOP enables the interactive design of transcriptome-and genomescale oligonucleotide FISH experiments. bioRxiv. - Moffitt, J. R., Hao, J., Bambah-Mukku, D., Lu, T., Dulac, C., & Zhuang, X. (2016). High-performance multiplexed fluorescence in situ hybridization in culture and tissue with matrix imprinting and clearing. Proceedings of the National Academy of Sciences, 113(50), 14456-14461. https://doi.org/10.1073/pnas. 1617699113
- Moffitt, J. R., & Zhuang, X. (2016). RNA imaging with multiplexed error-robust fluorescence in situ hybridization (MERFISH). In Methods in enzymology (Vol. 572, pp. 1-49). Academic Press.

https://doi.org/10.1016/bs.mi e.2016.03.020.

- Stellaris RNA FISH protocol for frozen tissues: https://biosearchassets.blob.c ore.windows.net/assets/bti_st ellaris_protocol_frozen_tissue.

License: This is an open access protocol distributed under the terms of the Creative Commons Attribution License, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working

Created: Feb 27, 2023

Last Modified: Mar 03, 2023

PROTOCOL integer ID: 77717

Solution Preparation:

Wash Buffer A (40% Formamide Wash Buffer)

- 2x SSC (20X), RNase-free Thermo Fisher Catalog #AM9763
- 1% (vol/vol)

 Tween 20 Sigma Aldrich Catalog #P9416-100ML
- 40% (vol/vol) S Formamide Thermo Fisher Scientific Catalog #AM9342
- Nuclease-free water

Wash Buffer A Master Mix, 45 ml:

- 36.75 ml nuclease-free water
- 7.5 ml SSC (20X), RNase-free Thermo Fisher Catalog #AM9763
- 750 µl 🔀 Tween 20 Sigma Aldrich Catalog #P9416-100ML
- Add 40% Someonide Thermo Fisher Scientific Catalog #AM9342 to prepare on demand

Wash Buffer B

2x SSC (20X), RNase-free Thermo Fisher Catalog #AM9763 prepared in nuclease-free water

Wash Buffer C (10% Formamide Wash Buffer)

- Nuclease-free water
- 2x SSC (20X), RNase-free Thermo Fisher Catalog #AM9763
- 10% (vol/vol)

 Formamide Thermo Fisher Scientific Catalog #AM9342

Wash Buffer C Master Mix, 49.5 ml:

- 44 ml nuclease-free water
- 5.5 ml SSC (20X), RNase-free **Thermo Fisher Catalog #AM9763**
- Add 10% Someonide Thermo Fisher Scientific Catalog #AM9342 formamide to prepare on demand

Saber Encoding Hybridization Buffer

- Nuclease-free water
- 2x SSC (20X), RNase-free Thermo Fisher Catalog #AM9763
- 40% (vol/vol)

 Formamide Thermo Fisher Scientific Catalog #AM9342

 Formamide

Keywords: In-situ hybridization, FISH, fluorescence, RNA, iterative FISH, formamide, amplified probes, bDNA, branched DNA, fluidics, acrydite, linker probe, anchor probe, encoding probes, tissue, cells, thick tissue, clearing, digestion, polyacrylamide gel

- 0.1% (wt/vol)

 Yeast tRNA Thermo Fisher Catalog #15401011
- 1% (vol/vol)
 - RNase Inhibitor, Murine 15,000 units **New England Biolabs Catalog**#M0314L
- 1% (vol/vol)

 Tween 20 Sigma Aldrich Catalog #P9416-100ML
- 10% (wt/vol)

 Dextransulfate Sigma Aldrich Catalog #D8906-100g
- Add 5-200 μM encoding probes depending on the size of the pool
- Prepare on demand

Saber Encoding Hybridization Buffer Master Mix, 4.8 ml (for 8.0 ml prep with formamide added)

- Nuclease-free water
- 800 µl 🔀 SSC (20X), RNase-free **Thermo Fisher Catalog #AM9763**
- 320 µl

 Yeast tRNA Thermo Fisher Catalog #15401011 reconstituted to 25mg/ml
- 80 µl
 - RNase Inhibitor, Murine 15,000 units **New England Biolabs Catalog**#M0314L
- 80 µl 🔀 Tween 20 Sigma Aldrich Catalog #P9416-100ML
- 0.8 g 🛭 Dextransulfate Sigma Aldrich Catalog #D8906-100g
- Aliquot mix and store at -20°C
- To prepare on demand, add 40% (vol/vol)
 - Solution Formamide Thermo Fisher Scientific Catalog #AM9342 to master mix at time of use
- Add 5-200 μM encoding probes depending on the size of the pool

Encoding Buffer Rinse (SSC-tw)

- Nuclease-free water
- 2x SSC (20X), RNase-free Thermo Fisher Catalog #AM9763
- 0.1% (vol/vol) 🔀 Tween 20 Sigma Aldrich Catalog #P9416-100ML
- Store at Room Temperature

PA Solution

- Nuclease-free water
- 4% (vol/vol)
 - 8 40% Acrylamide/Bis Solution 19:1 **Contributed by users Catalog** #1610144

- 60 mM

 Tris (1 M), pH 8.0, RNase-free Thermo Fisher Catalog #AM9856
- 0.3 M
 - X NaCl (5 M) RNase-free Thermo Fisher Scientific Catalog #AM9759
- One of the following:
- 1. For four-color experiments: 1:500 dilution 0.1 µm-diameter light yellow beads
 - 🔀 Fluorescent Light Yellow Particles SpheroTech Catalog #FP-0245-2
- 2. For two-color experiments:1:200,000 dilution of 0.1µm-diameter carboxylatemodified orange fluorescent beads
 - RluoSpheres™ Carboxylate-Modified Microspheres **Thermo Fisher** Scientific Catalog #F-8800
- De-gas solution before use
- Prepare on demand

PA Gel

- PA Solution including polymerizing agents:
- 1. 0.03% (wt/vol) Ammonium Persulfate Sigma Catalog #A3678
- 2. 0.15% (vol/vol) TEMED
- 3. Prepare on demand. Polymerizing agents will act rapidly. Make gel in small quantities (1ml) and right before use
- Prepare on demand

Storage Buffer (SSC-SB)

- Wash Buffer B
- 0.1% (vol/vol)
 - RNase Inhibitor, Murine 15,000 units **New England Biolabs Catalog**#M0314L
- Store in aliquots at -20°C

Amplifier Hybridization Buffer

- Nuclease-free water
- 2x SSC (20X), RNase-free Thermo Fisher Catalog #AM9763
- 10% (vol/vol)

 Formamide Thermo Fisher Scientific Catalog #AM9342
- 0.1% (wt/vol)

 Yeast tRNA Thermo Fisher Catalog #15401011
- 1% (vol/vol)
 - RNase Inhibitor, Murine 15,000 units **New England Biolabs Catalog** #M0314L
- 10% (wt/vol)

 Dextransulfate Sigma Aldrich Catalog #D8906-100g
- Prepare on demand

Amplifier Hybridization Buffer Master Mix, 7.2 ml (for 8.0 ml prep with formamide added)

- Nuclease-free water
- 800 µl 🔀 SSC (20X), RNase-free **Thermo Fisher Catalog #AM9763**
- 320 µl

 Yeast tRNA Thermo Fisher Catalog #15401011 reconstituted to 25mg/ml
- 80 µl
 - RNase Inhibitor, Murine 15,000 units **New England Biolabs Catalog**#M0314L
- 0.8 g Dextransulfate Sigma Aldrich Catalog #D8906-100g
- Aliquot mix and store at -20°C
- To prepare on demand, add 10% (vol/vol)
 - Solution Figure 12 Formamide Thermo Fisher Scientific Catalog #AM9342 to master mix at time of use
- Add 5 nM amplifiers

Readout Hybridization Buffer

- 2x SSC (20X), RNase-free Thermo Fisher Catalog #AM9763
- 10% (vol/vol)

 Ethylencarbonat 98% Sigma Aldrich Catalog #E26258
- 0.1% (vol/vol)
 - RNase Inhibitor, Murine 15,000 units **New England Biolabs Catalog**#M0314L
- Nuclease-free water
- 3 nM readout probes
- Prepare on demand

Wash Buffer D (Readout Wash Buffer)

- 2x SSC (20X), RNase-free Thermo Fisher Catalog #AM9763
- 10% (vol/vol)

 Ethylencarbonat 98% Sigma Aldrich Catalog #E26258
- Store at Room Temperature

Imaging Buffer

- 2x SSC (20X), RNase-free Thermo Fisher Catalog #AM9763
- 50 mM X Tris (1 M), pH 8.0, RNase-free Thermo Fisher Catalog #AM9856
- 10% (wt/vol) ⊠ D-()-Glucose Sigma Aldrich Catalog #DX0145-1
- 2 mM 🔀 Trolox Sigma Aldrich Catalog # 238813

- 0.5 mg/ml

 Glucose oxidase Sigma Aldrich Catalog #G2133
- 40 µg/ml 🔀 Catalase Sigma Catalog #C30
- 0.1% (vol/vol)
 - RNase Inhibitor, Murine 15,000 units **New England Biolabs Catalog**#M0314L
- Nuclease-free water
- Prepare on demand
- Store under layer of Mineral Oil when using fluidics system

Cleavage Buffer

- 2x 🔀 SSC (20X), RNase-free Thermo Fisher Catalog #AM9763
- 50 mM
 - Tris(2-carboxyethyl)phosphine hydrochloride solution Sigma Aldrich Catalo #646547-10X1ML
- Prepare on demand

DAPI Staining Solution

■ 50 µg/ml

★ 4,6-Diamidino-2-Phenylindole, Dihydrochloride (DAPI) Thermo Fisher Scientific Catalog #D1306

in Wash Buffer B for thick (40 µm) tissue

- 1-10 µg/ml
 - 4,6-Diamidino-2-Phenylindole, Dihydrochloride (DAPI) **Thermo Fisher**Scientific Catalog #D1306

in Wash Buffer B for thin (10 µm) tissue

Prepare on demand

Permeabilization Buffer (PBS-t)

- Nuclease-free water
- 1x
 - PBS Phosphate-Buffered Saline (10X) pH 7.4 Invitrogen Thermo Fisher Catalog #AM9625
- 0.5% (v/v) 🔀 Triton[™] X-100 **Sigma Aldrich Catalog #T8787-100M**L
- Store at Room Temperature

Permeabilization Buffer Wash (PBS-tw)

Nuclease-free water

- 1x
 - PBS Phosphate-Buffered Saline (10X) pH 7.4 Invitrogen Thermo Fisher Catalog #AM9625
- 0.1% (v/v)

 Tween 20 Sigma Aldrich Catalog #P9416-100ML
- Store at Room Temperature

SDS Clearing Solution (SDS-CS)

- Nuclease Free Water
- 4% **SDS Sigma Catalog #75746**
- 200 mM ⊗ Boric Acid Sigma Catalog #B0394
- pH to 8.5
- Store at Room Temperature

SAFETY WARNINGS

This protocol uses formamide which is a teratogen and can cause developmental malformation. Always work in a fume hood with formamide to avoid inhalation and avoid physical contact.

BEFORE START INSTRUCTIONS

All reagents for this protocol should be prepared sterile and RNase-free.

All incubation periods should be done in the dark.

Find reagent/buffer preparation instructions at index "Materials".

Quick Overview:

Part 1 - Tissue or Cell-Based Experiment Preparation

- Step 1 Coverslip functionalization
- Step 2 Mount, SDS pretreat and permeabilize sample
- Step 3 Hybridize linker (optional)
- Step 4 Wash away residual linker
- Step 5 Gel embed sample (optional)
- Step 6 Clear & digest sample (optional)
- Step 7 Hybridize encoding probes
- Step 8 Wash away residual encoding probes
- Step 9 Hybridize amplifiers (optional)

Part 2a - Multiplexed Iterative FISH Imaging with Fluidics System

- Step 1 Prepare solutions for imaging
- Step 2 Assemble fluidics system
- Step 3 MULTIPLEXED ITERATIVE FISH imaging protocol

Part 2b - Alternate MULTIPLEXED ITERATIVE FISH Imaging without Fluidics System

Step 1a: MULTIPLEXED ITERATIVE FISH imaging protocol - using coverslip mounted sample.

Step 1b: MULTIPLEXED ITERATIVE FISH imaging protocol - multiple hybridizations using chamber slip.

- Step 2: Dapi stain the sample.
- Step 3: Proceed to imaging of the sample.

Document Summary

Document Summary: This document, SOP002 - Multiplexed Iterative FISH Experimental Protocol, describes the process for in-situ fluorescence labeling of RNA transcripts in cells and tissues using a layered probe design, which allows for identity barcoding (MERFISH or similar) and/or signal amplification (Branched DNA [bDNA] amplification), along with a cleavable disulfide (S-S) reporter molecule, attached to a readout oligo, to allow for iterative rounds of labeling and imaging of the same sample with minimal disruption to sample integrity between rounds. This document also describes cell and tissue handling for the labeling process, and the RNA labeling process which uses an mRNA binding using a specialized poly-t (locked nucleic acid, LNA) probe with an acrydite linker to bind mRNAs to a polyacrylamide matrix and clearing techniques used to

Part 1 - Tissue or Cell-Based Experiment Preparation

Part 1 of this protocol describes the steps to setup a multiplexed iterative FISH experiment for tissue or cell-based samples. These steps are focused on the biochemical requirements for tissue or cell preparation, probe hybridization and imaging. This protocol does not cover the requirements of the microscope for imaging. Additional detail can for the imaging setup can be found at https://doi.org/10.1016/bs.mie.2016.03.020.

Part 1 - Step 1: Coverslip Functionalization

Refer to current version of SOP003 for protocol on Coverslip Functionalization. PDL-coated coverslips are preferable as tissue can be post-fixed to the coating using 4% PFA.

Part 1 - Step 2: Mount, SDS Pretreat and Permeabilize Sam

- 4 If using [M] 4 % PFA-fixed tissue , follow step 4. For [M] 4 % (V/V) PFA-fixed cells , skip below to step 5.
- **Using 4% PFA-fixed tissue.** Note: For some tissue, it is simpler to mount directly to the coverslip after slicing. In this case, mount and post-fix first then follow the remaining steps in order.

1h 15m

- i. Slice tissue and place slices in 1xPBS for 00:05:00. Remove PBS and repeat this for a second wash.
- ii. Pretreat tissue by washing in 4% SDS Clearing Solution (SDS-CS), once for 👏 00:05:00
- iii. To permeabilize the tissue, immerse the slip mounted tissue in 70% (vol/vol) ethanol

 Overnight at 4°C (recommended) in a Pyrex 60mm petri dish (Fisher 08-747A) or similar. (For faster results, sample can be incubated in EtOH for 1 hour at RT).
- iv. Move tissue slices to functionalized (PDL-coated) coverslip, aspirate off the 70% ethanol and incubate in Room temperature PBS buffer for 00:30:00 to rehydrate the sample.
- v. To bring the sample in sufficient contact with the coverslip surface, aspirate the PBS buffer and

place in \$\ \begin{align*} 45 \cdot \cdot

vi. Post-fix tissue to the coverslip by incubating in 4% PFA at 00:10:00 • Room temperature for

- vii. Remove PFA from the sample and rinse with 1x PBS for 00:05:00 at RT, two times.
- 6 Using 4% PFA-fixed cells grown on coverslip (optionally, use 8-chamber well or similar)

1h 15m

- i. To permeabilize the cells, immerse the slip mounted sample in 70% (vol/vol) ethanol

 Overnight at 4 °C (recommended) in a Pyrex 60mm petri dish (Fisher 08-747A). (For faster results, sample can be incubated in EtOH for 01:00:00 at 8 Room temperature).
- ii. Alternatively, pipette $100\mu L$ permeabilization buffer (PBS-t) to each well and incubate at $0 \, ^{\circ}C$ for 00:10:00 with gentle rocking.
- iii. Rinse with 8 Room temperature permeabilization buffer rinse (PBS-tw).
- iv. Aspirate rinse from the sample and let dry.
- Using a hydrophobic pen, draw a barrier around your sample and let dry before hybridizations. You may want to add a very small volume of PBS during this process to tissue samples to prevent sample desiccation.

Part 1 - Step 3: Hybridize Linker (optional; use when gel em.

- 8 Wash & equilibrate sample by immersing slip-mounted sample in 37 °C pre-heated ~ 30m A 200 µL Wash Buffer A for 60 00:30:00 .
- Assemble humidified chamber (empty pipette box with lid or otherwise that can house the sample-mounted coverslip with a single, saturated and folded paper used to line the inner edge of the chamber to prevent evaporation of probe solution).

10 Remove slip from Wash Buffer A and carefully wipe away excess buffer surrounding sample. 11 Dispense Δ 125 μL Encoding Hybridization Buffer containing [M] 1 micromolar (μM) linker to your sample, replace the petri dish lid, parafilm the dish and place the dish with the sample in the humidified chamber. **Adjust concentration of the linker according to sample size. 12 Incubate at 37 °C in a humidified chamber for 5 18:00:00 1d 12h up to 24:00:00 36:00:00 1h 10m Part 1 - Step 4: Wash Away Residual Linker 13 Remove the hybridization buffer and carefully remove excess buffer surrounding sample. 14 Immerse slip in pre-heated 4 37 °C Wash Buffer A for 30 min , two times. 15 Wash two times in 37 °C pre-heated Encoding Wash Buffer (SSC-tw) for 5 min each. 16 Wash two times in | | | | | | | | Room temperature Part 1 - Step 5: Gel Embed Sample (optional - gel embed

1. Wash sample for 00:02:00 ** with de-gassed PA Solution. **adjust time based on sample

size. For 100µm tissue slices, increase this to 3 hours.

17



18



Wash sample for 00:02:00 with PA Gel Solution and then remove.

2m



- 19 1. Cast a thin PA film by adding Δ 50 μL gel solution - Δ 100 μL gel solution to the sample and invert a smaller (25 mm) gel-slick coated coverslip onto the gel solutions being careful to avoid air bubbles. Adjust the volume and make sure your gel film is thin. Aspirate any extra gel solution away.
- 20 Allow casting for 60 01:30:00 at Room temperature

1h 30m

21 After casting, carefully remove the smaller coverslip from your sample. If the coverslip is stuck, you can loosen the coverslip by immersing in SDS-CS at 1 37 °C


Part 1 - Step 6: Digest & Clear Sample (optional) (buffer

- 22 Note: For lung tissue start at step 22. Skip ahead to step 24 for brain tissue.
- 23 Incubate sample in A 3 mL PBS with [M] 10 % collagenase/elastase at [M] 20.000 U/mL (5) 03:00:00 at 8 37 °C



24 Wash the sample with a quick rinse of RT 1x PBS followed by two 5 min washes of 1x PBS at

Room temperature

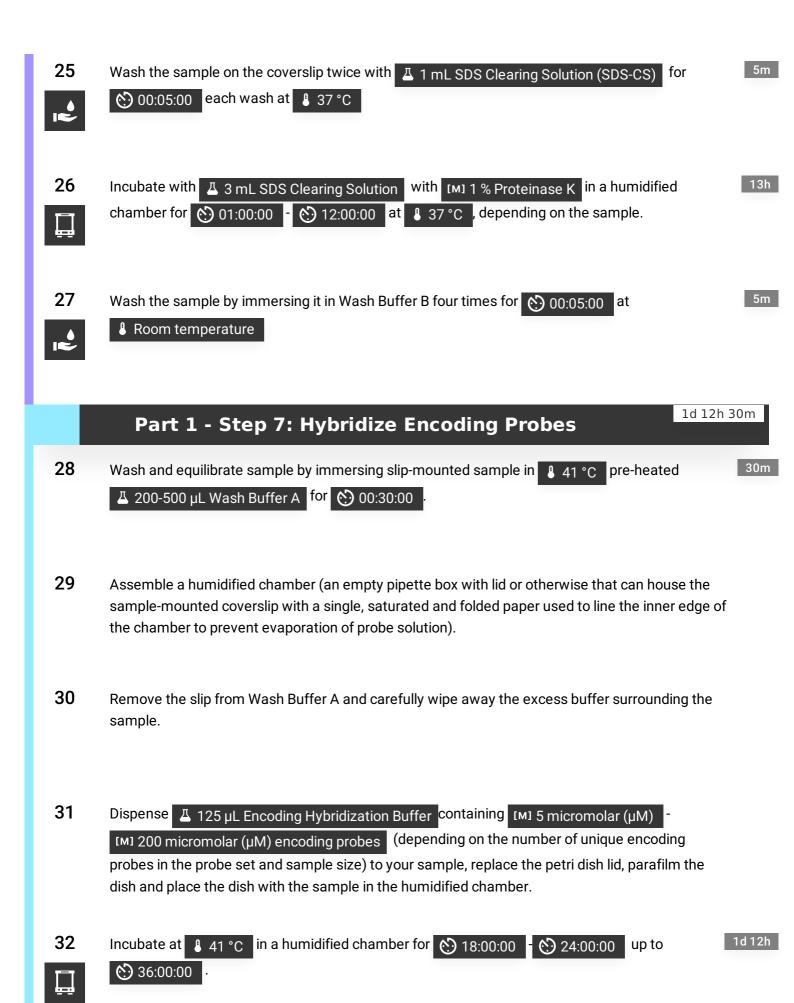

Part 1 - Step 8: Wash Away Residual Encoding Probes

1h 10m

- Remove the hybridization buffer and carefully remove the excess buffer surrounding the sample.
- Wash the sample in pre-heated 41 °C Wash Buffer A for 00:30:00 , two times.

30m

Wash two times in 41 °C pre-heated Encoding Wash Buffer (SSC-tw) for 00:05:00

5m

Wash two times in Room temperature 1x PBS

Part 1 - Step 9: Hybridize amplifiers (optional)

1h 25m

To label the gel embedded and cleared samples with primary and secondary amplifiers.

П

38

Incubate sample in Wash Buffer C at 8 37 °C for 00:30:00

30m

39 Aspirate to remove Wash Buffer C.

Hybridize primary amplifier. Incubate the sample in a A 125 µL droplet of

[M] 5 nanomolar (nM) primary amplifier * in amplifier hybridization buffer for humidity-controlled 37°C incubator. ** in amplifier hybridization buffer for humidity-controlled 37°C incubator.

- *Amplifier concentration may need to be increased based on the thickness of your sample.
- **Adjust incubation time based on sample size. For a 100µm tissue section, an overnight incubation is preferable for the primary amplifiers while the secondary amplifiers can be incubated for 5-6 hours, on the following day.
- Wash 3 times with Wash C for 5-10 min each at Room temperature
- Hybridize secondary amplifier. Incubate the sample in a Lambda 125 µL droplet

 [M] 5 nanomolar (nM) secondary amplifier in amplifier hybridization buffer for in humidity-controlled 37°C incubator.
 - *Amplifier concentration may need to be increased based on the thickness of your sample.
 - **Adjust incubation time based on sample size. For a 100µm tissue section, an overnight incubation is preferable for the primary amplifiers while the secondary amplifiers can be incubated for 5-6 hours, on the following day.
- Wash twice in Room temperature Wash C for 5 min each followed by a 00:15:00 00:30:00 wash in 37 °C Wash C.
- Perform MULTIPLEXED ITERATIVE FISH Imaging (Part 2) immediately or store sample for up to 24 hours in storage buffer at 4 °C.

Part 2a - MULTIPLEXED ITERATIVE FISH Imaging with Fluidi..

Note

**The following steps (Part 2a) are used for the Full MULTIPLEXED ITERATIVE FISH protocol. For experiments that don't use the fluidics system, move to Part 2B below.

45

Imaging for MULTIPLEXED ITERATIVE FISH involves multiple rounds of fluid exchange to hybridize, image, cleave and rinse samples. Automated fluid exchange and imaging approach is recommended. For setups lacking an automated fluidics exchange system, proceed to Part 2b.

Part 2a - Step 1: Prepare Solutions for Imaging

- 47 Prepare the following solutions with the corresponding volumes:
 - i. Readout Hybridization Buffer (RHB)
 - ii. Readout Wash Buffer (Wash D)
 - iii. Imaging Buffer (store under mineral oil) (IB)
 - iv. TCEP Cleavage Buffer (CB)
 - v. 2x SSC Wash Buffer (Wash B)
 - vi. DAPI Staining Solution

Part 2a - Step 2: Assemble Fluidics System

- 48 Make sure that all tubing is properly connected. MULTIPLEXED ITERATIVE FISH probes and preparation time are costly so leaks need to be avoided at all costs.
- 49 Ensure the system is fully assembled, plugged in and turned on.
- Double-check correctness of the details for the pump protocol for the MULTIPLEXED ITERATIVE FISH
 Fluidics for the current project.
- Load the sample to the flow cell and connect.
- **52** Carefully load all solutions to the proper reservoirs.

Part 2a - Step 3: MULTIPLEXED ITERATIVE FISH Imaging

1h 5m

Once the fluidics system is setup, solutions are prepped and loaded and the sample is in place in the chamber, an automated program should run the following cycle:

Readout hybridization buffer (with readout probes)

2h 22m

- a. Run 🗸 2 mL over 🚫 00:03:00 to prime buffer to the sample.
- b. Run additional A 2.5 mL over the sample for (5) 00:04:00
- c. Pause flow for 00:15:00 02:00:00 depending on sample size $10 \, \mu m = 15 \, min$, $30 \, \mu m = 60 \, min$, $100 \, \mu m = 120 \, min$).
- **55** Readout Wash Buffer (Wash D)

13m 30s

- a. Run <u>A</u> 2 mL over 00:03:30 to flush.
- b. Run 🚨 2 mL over 🚫 00:10:00 to wash.
- 56 Imaging Buffer

10m


- a. Run 2 mL over 00:10:00 to run remaining wash D over sample and prime imaging buffer to the sample.
- 1. 2x SSC Wash Buffer (Wash B)

3m

- a. Run 🗸 1 mL over 🔇 00:03:00 to move 1mL imaging buffer over the sample.
- 58 Imaging. Pause fluidics and proceed with imaging.
- TCEP Cleavage Buffer -

18m 30s

- b. Run 🚨 1 mL over sample for 🚫 00:05:00

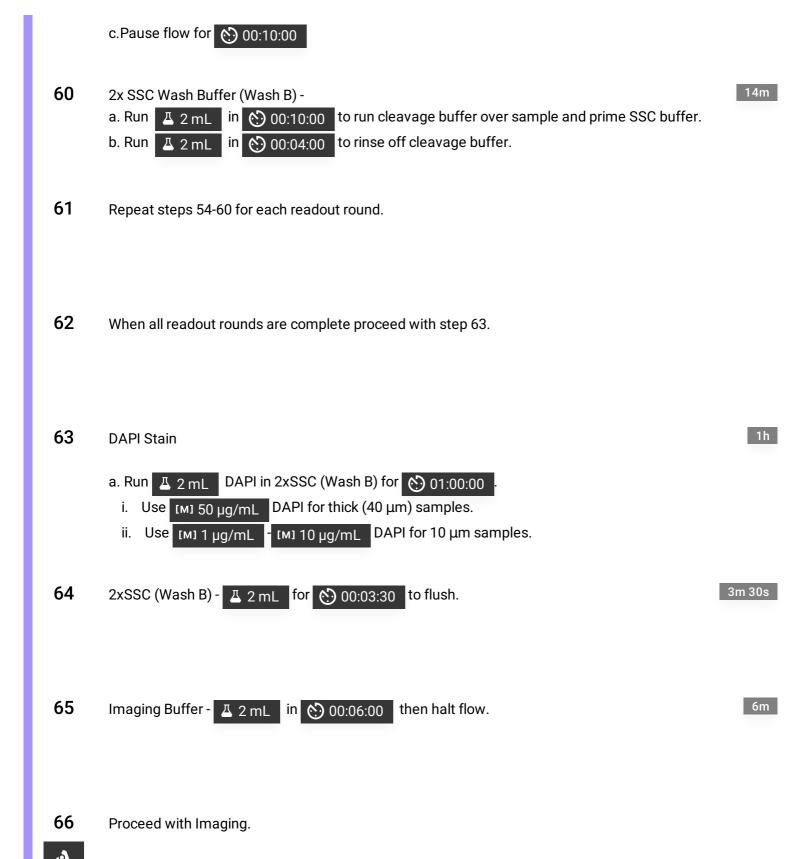

Part 2b - Alternate MULTIPLEXED ITERATIVE FISH Imaging ...

Note

The following steps are used for manual, iterative FISH without a fluidics system. For trial that uses the fluidics system, move to Part 2a (above).

For some MULTIPLEXED ITERATIVE FISH experiments, it may be simpler to proceed without the fluidics system for imaging. Once you have hybridized probes and amplifiers if desired, readout probes can be hybridized and imaged in a single round or in multiple rounds if necessary. If you are hybridizing more than one round of readouts, proceed to Steps 1b-3.

Part 2b - Step 1a: MULTIPLEXED ITERATIVE FISH Imaging P.

69 Readout Probe Hybridization.





- a. Pipette A 200 µL 3nM readout probes in Readout Hybridization Buffer to sample and incubate at Room temperature for 00:10:00.
- b. Aspirate Readout Hybridization Buffer from the sample.
- 71 Dapi Stain. Add Δ 200 μL Wash Buffer B with DAPI nuclear stain (at 1μg/mL) to sample and incubate for 00:30:00 at 3 37 °C.

30m



Remove the Dapi stain and wash with Wash Buffer B for 5 min, two times.



Add \perp 100 μ L - 200 μ L Imaging buffer to sample and mount to glass plate with clear nail polish.

74



Part 2b - Step 1b: MULTIPLEXED ITERATIVE FISH Imaging - .

75 Readout Probe Hybridization.

10m


- a. Pipette 200 µL 3nM readout probes in Readout Hybridization Buffer to sample and incubate at Room temperature for 00:10:00.
- b. Aspirate Readout Hybridization Buffer from the chambers.
- Wash away unbound probe by adding Δ 200 μL RT Readout Wash Buffer D to sample for 5 min, two times. Additional washes may improve the result.
- 77 Add \perp 100 μ L \perp 200 μ L Imaging buffer to sample.
- **78** Proceed with imaging of the round.



79 TCEP Cleavage Buffer − ♣ 100 µL for ♦ 00:15:00

15m

- 2x SSC Wash Buffer (Wash B) \pm 250 μ L each well, three times.
- **81** Repeat steps 81-85 for each probe set round.

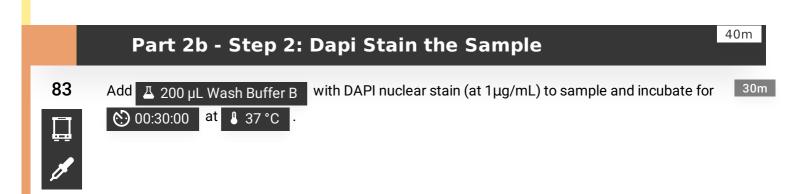

Wash sample in Wash Buffer B for 5 min two times.

Part 2b - Step 3:

85 Proceed to Imaging of the Sample.