




Oct 24, 2022

Automated 96-well PCR Purification

Ariel Rabines¹, Rob Lampe², Andrew E Allen²¹J. Craig Venter Institute; ²Scripps Institution of Oceanography

1 Works for me

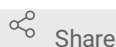

Share

dx.doi.org/10.17504/protocols.io.dm6gpr19jvzp/v1

Rob Lampe

Scripps Institution of Oceanography, J. Craig Venter Institu...

ABSTRACT

This protocol describes a fully-automated PCR cleanup in a 96-well PCR plate format with AMPure XP beads on an Eppendorf epMotion. Multichannel pipetting with the 300 uL and 50 uL pipettes is used. It is designed to use the thermomixer module; however, this can be replaced by modifying the program to mix by pipetting. The bead magnet used in this program is a Alpaqua Low Elution Magnet Plate. We also use Eppendorf twin.tec PCR plates (LoBind, skirted, 150uL).

IMPORTANT: The program is currently set to add **22 µL** of water to each well followed by **65 µL** beads. This results in an approximately 1:1.5 bead ratio with a starting volume of **22 µL** PCR product. Modify the first step accordingly with your starting volume and desired bead ratio.

The manufacturer's instructions are found here:

<https://www.beckmancoulter.com/wsrportal/techdocs?docname=B37419>

The epMotion program is attached here.

ATTACHMENTS

Application_1x bead
cleanup_210519_140208.
export6

DOI

dx.doi.org/10.17504/protocols.io.dm6gpr19jvzp/v1

PROTOCOL CITATION

Ariel Rabines, Rob Lampe, Andrew E Allen 2022. Automated 96-well PCR Purification. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.dm6gpr19jvzp/v1>



KEYWORDS

PCR purification, Eppendorf epMotion, bead cleanup, automatic

LICENSE

_____ This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

May 31, 2020

LAST MODIFIED

Oct 24, 2022

PROTOCOL INTEGER ID

37637

MATERIALS TEXT

☒ [Agencourt AMPure XP beads](#) **Contributed by users**

☒ [Eppendorf twin.tec® PCR 96-well plate,](#)
[skirted Eppendorf Catalog #951020401](#)

☒ [Molecular grade water nuclease-free](#) **Contributed by users**

☒ [100% Molecular grade ethanol](#) **Contributed by users**

LE Magnet Plate

Alpaqua A000350

- 1 Take aliquot of beads out of refrigerator and allow to warm to room temperature.
- 2 Start program and input the number of samples you have. It will then tell you the minimum volumes for each reagent. The positions in the reservoir rack are:
 1. Molecular grade H₂O
 3. Beads
 5. 70% EtOH
 7. Elution buffer (e.g. H₂O, TRIS-Acetate (10 mM pH 8.0), or TE Buffer (10 mM Tris-Acetate pH

8.0, 1 mM EDTA))

You can use the level sensor to detect volumes or measure them with a serological pipette. Do not start the program until all of your items are in the correct positions.


3 Mix beads well until they appear homogenous.

4 Prepare fresh 70% EtOH.

From the Agencourt XP beads manual: "Fresh 70% ethanol should be prepared for optimal results. There is also miscibility involved with ethanol and water. For example, measuring out 70 mL of ethanol and topping off to 100 mL with water will generate ~65% ethanol. Measuring 70 mL ethanol and 30 mL water separately, then combining them will generate ~95 mL of 70% ethanol."

5 Place reagents, your PCR plate, and a new clean plate on the epMotion.

Materials layout:

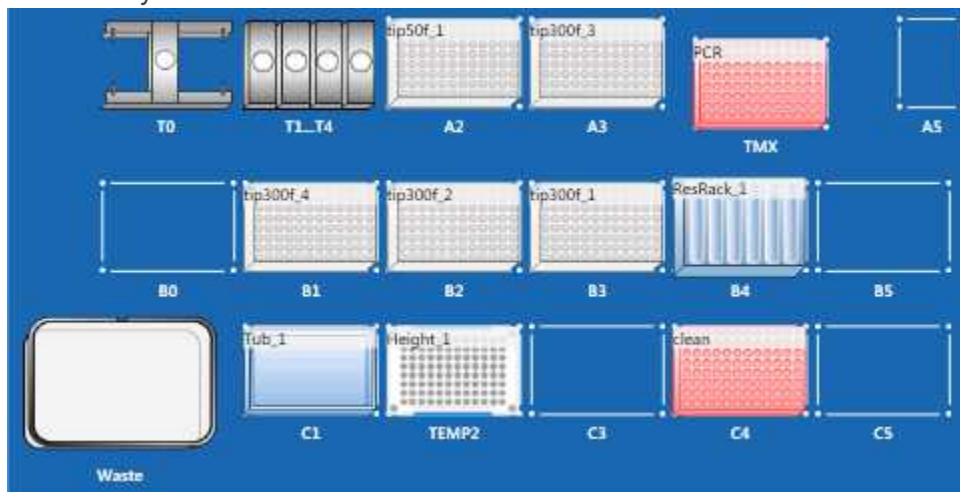

Run the program. It will take approximately 1 hour for a full plate.