


Dec 17, 2021

Modified Illumina DNA Prep (M) Tagmentation Library Preparation for cDNA amplicons from wastewater

 **Illumina DNA Prep (M) Tagmentation Library Preparation for use on an Illumina MiSeq Sequencer**

Padmini Ramachandran¹, Tamara Walsky¹, Amanda Windsor¹,
Maria Hoffmann, Chris Grim¹

¹FDA CFSAN



dx.doi.org/10.17504/protocols.io.b2yzqfx6

 Chris Grim

Please note that this protocol is public domain, which supersedes the CC-BY license default used by protocols.io.

This procedure outlines the protocol for sequencing of VarSkip Short SARS-CoV-2 cDNA amplicons using the Illumina DNA Prep library preparation kit for sequencing on an Illumina MiSeq sequencer.

DOI

dx.doi.org/10.17504/protocols.io.b2yzqfx6

Padmini Ramachandran, Tamara Walsky, Amanda Windsor, Maria Hoffmann, Chris Grim 2021. Modified Illumina DNA Prep (M) Tagmentation Library Preparation for cDNA amplicons from wastewater. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.b2yzqfx6>



Illumina DNA Prep (M) Tagmentation Library Preparation for use on an Illumina MiSeq Sequencer, Julie Haendiges

Illumina DNA Prep, WGS Library Preparation, GenomeTrakr, Whole Genome Sequencing

_____ protocol ,

Dec 16, 2021

Dec 17, 2021

56057

Illumina DNA Prep (M) Tagmentation Kit contain 3 components:

Box 1 of 3:

- SPB (Store at 2-8°C)
- TSB (Store at room temperature)
- TWB (Store at room temperature)

Box 2 of 3: (Store at -25 to -15°C)

- RSB
- TB1
- EPM

Box 3 of 3: (Store at 2-8°C)

- BLT

Abbreviations:

BLT: Bead-Linked Transposome

dsDNA: Double-Stranded DNA

EPM: Enhanced PCR Mix

HT1: Hybridization Buffer

PCR: Polymerase Chain Reaction

PR2: Incorporation Buffer

RSB: Resuspension Buffer

SPB: Sample Purification Beads

TB1: Tagmentation Buffer 1

TSB: Tagment Stop Buffer

TWB: Tagment Wash Buffer

MATERIALS

 [Ethanol \(100%, Molecular Biology Grade\)](#) **Fisher**

Scientific Catalog #BP2818500

 [Qubit® dsDNA HS assay kit, 100 reactions](#) **Life**

Technologies Catalog #Q32851

 [Nextera DNA CD Indexes \(96 samples\)](#) **Illumina,**

Inc. Catalog #20018708

 [Sodium Hydroxide 1N](#) **Sigma**

Aldrich Catalog #S2770-100ml

 [Molecular grade water nuclease-free](#) **Contributed by users**

 Illumina DNA Prep (M) Tagmentation (96 Samples) **Illumina,**

Inc. Catalog #20018705

 Illumina DNA Prep (M) Tagmentation (24 Samples) **Illumina,**

Inc. Catalog #20018704

 IDT for Illumina - DNA/RNA UD Indexes Set A (96 Samples) **Illumina,**

Inc. Catalog #20027213

 IDT for Illumina - DNA/RNA UD Indexes Set B (96 Samples) **Illumina,**

Inc. Catalog #20027214

 IDT for Illumina - DNA/RNA UD Indexes Set C (96 Samples) **Illumina,**

Inc. Catalog #20027215

 IDT for Illumina - DNA/RNA UD Indexes Set D (96 Samples) **Illumina,**

Inc. Catalog #20027216

Supplies:

- Qubit Assay Tubes (ThermoFisher cat# Q32856)
- Pipette Tips, sterile, filtered (assorted volumes)
- Conical Tubes, 10ml and/or 15ml (FisherSci cat# 14-959-53A or equivalent)
- Solution basins, sterile (FisherSci cat# 13-681-504 or equivalent)
- 96-well PCR Plates, semi-skirted, flat deck (FisherSci cat# AB-1400L or equivalent)
- Microcentrifuge tubes, 1.5 ml, sterile (ThermoFisher cat# AM12400 or equivalent)
- Plate Seals (FisherSci cat# AB-0558 or equivalent)

Equipment:

- Qubit 2.0 or 3.0 Fluorometer
- Thermocycler
- Microplate centrifuge
- Vortex
- Magnetic Stand-96 (ThermoFisher cat# AM10027) (If possible, have two; one for pre-PCR and one for post-PCR)
- Micropipettes (Single and Multichannel)
- Ice bucket
- Microcentrifuge



Chemical Safety Warning: Take proper precautions, and wear appropriate

PPE when handling potentially hazardous chemicals. Ensure that chemicals, spent containers, and unused contents are disposed of in accordance with governmental safety standards.

Illumina DNA Flex Library Preparation Kit: See Illumina SDSs for additional information. Take proper precautions and wear appropriate PPE when handling reagents.

TSB: GHS Category 1 for eye damage/irritant and is harmful to aquatic life.

TB1: GHS Category 4 for acute toxicity (dust/mist), Category 2A for eye irritant and Category 1B for reproductive toxicity. Contains N,N-Dimethylformamide.

EPM: GHS Category 4 for acute oral toxicity and Category 1 for specific organ toxicity. Contains tetramethylammonium chloride

Please note that this protocol is public domain, which supersedes the CC-BY license default used by protocols.io.

Run tape station or Bioanalyzer on the cDNA amplicon to determine the clean up of cDNA or no clean up as the next step. Quantify the cDNA using Qubit and determine the input cDNA concentration for the Illumina DNA prep.

Dilute and Tagment Input DNA

1

Based on the viral load in the sample, you will see a spectrum of fragment profiles. If there is a single peak around 400-600 bp, the clean-up step can be omitted and you can proceed directly to library preparation using this protocol.

Quantify the cDNA amplicon using Qubit HS kit. We recommend starting the input cDNA concentration at 100-200NG. Add 2-30ul DNA, so that the total input amount is 200NG. If DNA volume is less than 30ul, add nuclease free water to the DNA samples to bring the total volume to 30ul.

The minimum input cDNA recommended is 25NG. The number of PCR index cycle will vary depending on the starting input.

2 Bring BLT (stored in refrigerator) and TB1 (stored in freezer) to room temperature.

Ensure that BLT is stored upright at all times, so that the beads remain submerged in

the
buffer

- 3 Label a 96-well PCR plate with the Run ID.
- 4 Add molecular-grade water to the each sample well (from Column K of the workbook)
- 5 Add gDNA to the molecular-grade water (per volume in Column J) and mix well by gently pipetting 5 - 10 times.
- 6 Vortex BLT vigorously for 10 seconds, visually check the beads for complete resuspension and repeat vortexing if necessary.

Do not spin down the BLT tube, the beads must be resuspended

- 7 Vortex the TB1 and spin down the tube.
- 8 *****Scale up this step according to the number of reactions plus 3-4 for dead space volume/error*****

Prepare the Tagmentation Master Mix:

Combine  **10 µL** of TB1 with  **10 µL** of BLT

Reagent	Volume per Sample
TB1	10 µl
BLT	10 µl

Reagent volumes per sample for tagmentation master mix

Note: The 96 sample kit comes with 4 tubes of each reagent, each tube contains enough reagents for 24 samples

- 9 Vortex the tagmentation master thoroughly to make sure the BLT beads are evenly resuspended in the buffer.





10 Using fresh tips, transfer  **20 µL** of tagmentation master mix to each sample well.

Note: The master mix can be added to a reagent basin and distributed using a multichannel pipet

11 Pipette up and down 10 times mix the 50 µl reaction to resuspend the beads.

12 Apply an adhesive PCR plate seal to the plate.

13 Place the plate into the thermocycler and run the tagmentation program.


Program thermocycler to incubate at  **55 °C** for  **00:15:00** followed by a  **10 °C** hold with the lid heated at  **100 °C**

14 Check TSB for precipitate (if present, warm at 37 °C for up to 10 minutes and vortex) and ensure it is at room temperature prior to use.





15 Upon completion of the incubation, remove the plate from the thermocycler. Proceed to the Post Tagmentation Cleanup step.

Post Tagmentation Clean Up

16 Remove the plate seal.

17 Add  **10 µL** of TSB to each sample (A multi-channel pipette can be used) Gently pipette up and down 10 times to mix and fully resuspend the beads in the 50 µl reaction.

18 Apply an adhesive PCR plate seal to the plate.

19 Place the plate into the thermocycler and incubate at  **37 °C** for  **00:15:00** followed by a  **10 °C** hold with the lid heated at  **100 °C**

- 20 While samples are incubating, thaw EPM (stored in freezer) on ice and thaw indices at room temperature.
- 21 Remove the plate from the thermocycler, quick spin the plate and remove the seal.
- 22 Place the plate on a magnet for ⌚00:03:00 or until solution is clear.

Note: The DNA is tagged with adapters and bound to the beads.
- 23 Using a multichannel pipette, remove the supernatant and discard.
- 24 Remove the plate from the magnet and add 🧴100 µL of TWB directly to the pellet. Gently pipette to mix until beads are fully resuspended, try to avoid creation of foam from TWB.
- 25 Place the plate on the magnet for ⌚00:03:00 or until solution is clear.
- 26 Remove the supernatant and discard.
- 27 Remove the plate from the magnet and add 🧴100 µL of TWB directly to the pellet. Gently pipette to mix until beads are fully resuspended.
- 28 Place the plate on the magnet for ⌚00:03:00 or until solution is clear.
- 29 Remove the supernatant and discard.
- 30 Remove the plate from the magnet and add 🧴100 µL of TWB directly to the pellet. Gently

pipette to mix until beads are fully resuspended.

- 31 Place the plate with TWB on the magnet and allow to incubate until ready to proceed with adding the PCR master mix in the Amplify Tagmented DNA step. The plate should incubate for at least 3 minutes. It is important to keep the pellet in TWB to prevent overdrying of the beads.

Amplification and Index Addition of Tagmented DNA

- 32 Invert the EPM to mix, then briefly centrifuge.

- 33 Briefly centrifuge the Index plate.

- 34 *****Scale up this step according to the number of reactions plus 3-4 for dead space volume/error*****

Prepare the PCR master mix:

Combine  20 µL of EPM with  20 µL of Molecular grade water

Reagent	Volume per Sample
EPM	20 µl
Molecular grade water	20 µl


Reagent volumes per sample for PCR master mix

Note: The 96 sample kit comes with 4 tubes of each reagent, each tube contains enough reagents for 24 samples


- 35 Vortex and spin down the PCR master mix.
- 36 Remove the third TWB wash from the samples while on the magnet. Remove any excess liquid from the plate using a small volume pipette.

Note: Removal of TWB is crucial, as it can impede PCR. Any foam remaining on the wells will not negatively impact the library.

- 37 Remove the plate from the magnet and immediately proceed to adding the master mix.

- 38 Add  **40 µL** of PCR master mix to each sample well. Gently pipette to mix to ensure beads are resuspended.

Note: The master mix can be added to a reagent basin and distributed using a multichannel pipet

- 39 Add  **10 µL** of the index primer pair from the appropriate index wells in accordance with the sample sheet. The plate has a foil seal on it, P20 tips are sufficient to pierce the seal to pipette. The indexes are single-use only.

- 40 Use a pipette to gently mix a minimum of 10 times to ensure thorough mixing.

- 41 Apply an adhesive PCR plate seal to the plate.


42

NOTE: Depending on the initial concentration of cDNA amplicon that was used as the input, the number of cycles for index PCR WILL vary.

Aim for at least 100-200 NG of input cDNA amplicon in 30ul as recommended by Illumina.

If the input cDNA is around 100-200 NG, the index PCR cycle has to be set at 7/8 cycles.

If the cDNA amplicon input is from 25-99 NG, the index PCR cycle has to be set at 9/10 cycles.

- 43 Place the plate into the thermocycler and run the following pre-programmed settings with a heated lid at  **100 °C**

Thermocycler protocol (for use with DNA inputs above 100ng)

Step 1:  **68 °C** for  **00:03:00**

Step 2:  **98 °C** for  **00:03:00**

Step 3: 5 cycles of:

⌄ 98 °C for ⌚ 00:00:45

⌄ 62 °C for ⌚ 00:00:30

⌄ 68 °C for ⌚ 00:02:00

Step 4: ⌄ 68 °C for ⌚ 00:01:00

Step 5: Hold at ⌄ 10 °C

44 Centrifuge plate at ⌚ 280 x g for ⌚ 00:01:00

45 ⏸

This is a safe stopping point. The plate may be sealed and stored at ⌄ 2 °C to ⌄ 8 °C for up to 3 days.

Clean up Libraries

46 **NOTE:** *The steps listed below are critical for efficient size selection, product recovery and thus cluster generation and sequencing. Always check pipette tips for correct volumes and ensure that no beads have accidentally been aspirated. If beads have been aspirated or the bead pellet is disturbed, allow the pellet to reform (3-5 minutes on the magnet) and repeat the step.*

47 Before starting, prepare reagents:

47.1 Prepare fresh 80% ethanol sufficient for all samples.

47.2 Bring RSB to room temperature (from freezer) and vortex to mix.

47.3 Bring SPB to room temperature (at least 30 minutes) from refrigerator. Vortex and invert SPB several times to full resuspend the beads.


48 Prepare the SPB Master mix:


Reagent	Volume/Reaction (μl)	Total Volume (μl) for 24-32 Samples
SPB	45 μl	1350 μl
Nuclease-Free Water	40 μl	1200 μl

Note: The ratio of SPB:water has been validated by FDA-CFSAN for size-selection. CDC PulseNet uses a different ratio of SPB:water to select for large insert sizes.


49 If plate was retrieved from 4°C storage, centrifuge plate at 280 x g for 1 minute.

Remove the seal.

50 Place the sample plate on the magnet for  00:05:00

51 Transfer  45 μL of supernatant (now containing the DNA) to a set of new wells on the sample plate.

52 Remove sample plate from the magnet.

53 Vortex SPB master mix thoroughly and add  85 μL to each PCR product.

Note: The master mix can be added to a reagent basin and distributed using a multichannel pipet

54 Pipette to mix a minimum of 10 times or until thoroughly mixed.


55 Incubate at room temperature for  00:05:00

56 Place the 96-well plate on the magnet for  00:05:00 or until supernatant is clear.

During incubation, vortex the stock SPB to resuspend the beads.


57 With the plate still on the magnet, transfer  125 μL of supernatant (containing the DNA) to a

new set of wells.

58 Remove the plate from the magnet and add  15 µL of stock SPB to the supernatant.


59 Gently pipet at least 10 times to mix.

60 Incubate at room temperature for  00:05:00

61 Place on the magnet for  00:05:00 or until clear.

Remove and discard the supernatant (DNA is now bound to the beads) without disrupting the beads

62 Perform the steps below twice (for a total of two washes)

62.1 While the plate is on the magnet, add  180 µL of the prepared 80% ethanol

Note: Do not add directly to the bead and do not mix.

62.2 Incubate for  00:00:30





62.3 Remove and discard ethanol.


63 Use a pipette to remove any excess liquid from the plate.

Use a small volume pipette to get out any residual if necessary.

64 Allow beads to air dry for up to  00:05:00 (minimum of 3 minutes)


Note: Do not allow beads to over-dry. If bead pellet appears to be cracking, immediately resuspend beads regardless of drying time.

- 65 Remove the plate from the magnet and add  **32 µL** of RSB. Pipet thoroughly to mix
- 66 Incubate at room temperature for  **00:02:00**
- 67 Place the plate on the magnet for  **00:02:00**
- 68 Prepare a new 96-well plate as in below that is labelled with run date and initials. Transfer  **25 µL** of the supernatant to a new 96-well plate. This is the final library.

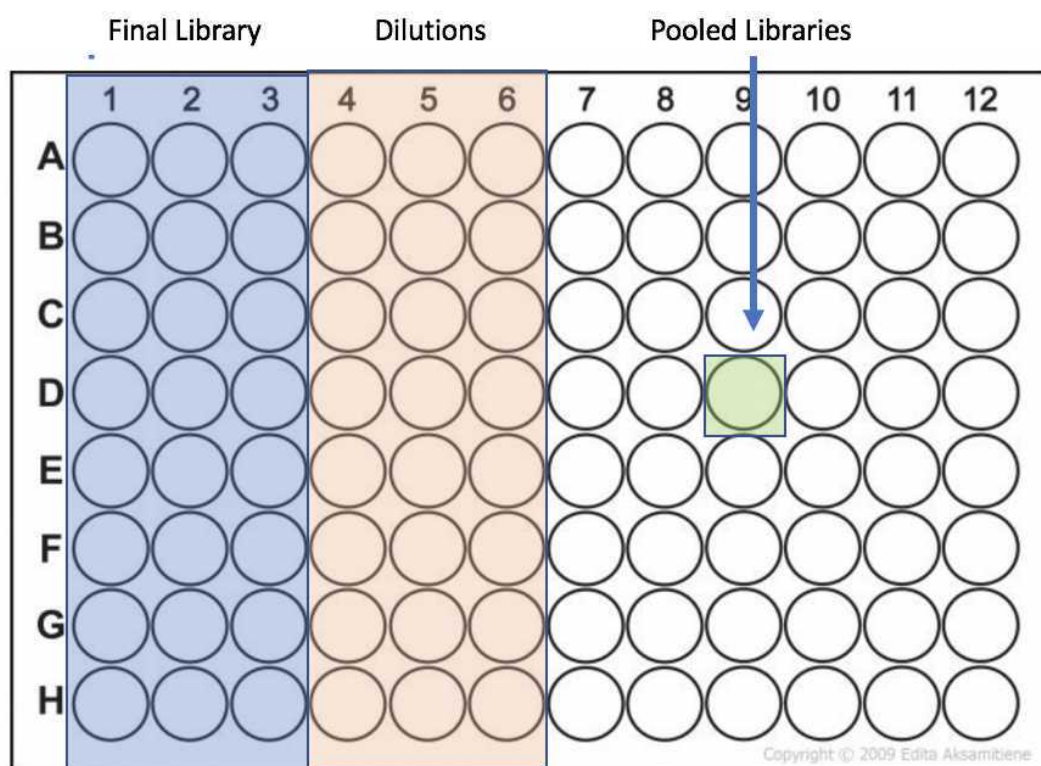



Example 96-well plate set up

69 




If ready to proceed go to the Pooling Libraries step. Otherwise, this is a safe stopping point. The plate may be sealed and stored at  **-20 °C** for up to 30 days.

Quantification, Normalization and Pooling of Libraries

- 70 Quantify each sample using the Qubit dsDNA High Sensitivity kit. See SOP titled "DNA Quantification using the Qubit Fluorometer" for more detailed information on performing DNA quantification.
- 71 Enter the Qubit values into the "Prep Workbook" tab on the DNA Prep Library worksheet. The DNA Prep worksheet can be found in section 8 of this SOP. The dilution values will be calculated automatically.
- 72 For v2 sequencing chemistry, it is recommended to dilute to 2nm; and for v3 sequencing chemistry to dilute to 4nm.
- 73 Dilute each sample according to the values on the worksheet in the designated dilution wells on the 96-well plate.
- 74 Pool  **5 µL** of each diluted library into the specified well and pipet to mix. **Note: if sequencing *Campylobacter*, pool  **2.5 µL** of diluted sample.**

Denaturing Pooled Library

- 75 Prepare a fresh aliquot of 0.2N NaOH. (This should be made fresh for each run)

Note: It is recommended to make aliquots of 1N NaOH and store in the freezer.
- 76 Transfer  **5 µL** of the diluted library to a new Eppendorf LoBind tube.
- 77 Add  **5 µL** of 0.2N NaOH and pipette to mix.
- 78 Incubate at room temperature for  **00:05:00** to denature the dsDNA.

79 Immediately add **990 µL** of HT1 and pipette to mix. The concentration for a 2 nm start is 10 pM and for 4 nm is 20 pM.

80 Dilute the denatured library to the final desired loading concentration.

	Denatured Pool Concentration = 10 pM		Denatured Pool Concentration = 20 pM	
Final Loading Concentration	Required Volume of HT1	Required Volume of Denatured Pool	Required Volume of HT1	Required Volume of Denatured Pool
8 pM	200 µl	800 µl	600 µl	400 µl
9 pM	100 µl	900 µl	550 µl	450 µl
10 pM	NA	NA	500 µl	500 µl
12 pM	NA	NA	400 µl	600 µl
14 pM	NA	NA	300 µl	700 µl
15 pM	NA	NA	250 µl	750 µl
18 pM	NA	NA	100 µl	900 µl
20 pM	NA	NA	NA	NA

Note: The final loading concentration for optimal Cluster Density may need to be adjusted based on data from previous runs.

81 Mix by repeated inversion of sample tube.

Optional: Denatured PhiX control can be spiked in at this point.

82 Heat the denatured DNA library to **96 °C** for **00:02:00** in a heat block to ensure complete denaturation of all dsDNA in the sample.

83 Immediately cool in an ice-water bath for at least **00:05:00** prior to loading.

The DNA library may sit on ice or at **2 °C - 8 °C** until ready for loading (<30 minutes)

84 Proceed to loading **600 µL** of the denatured library to the thawed cartridge.

- 85 Illumina MiSeq Operation and Maintenance: See “**Procedure for Operation and Maintenance of the Illumina MiSeq for Whole Genome Sequencing**” SOP for further instructions.

Appendix

86 Appendix 1

DNA Prep Library Workbook

 [GT_DNA Prep Library Worksheet.xlsx](#)

87 Appendix 2

Nextera DNA CD Indexes (96 indexes, plated)

	Run 1			Run 2			Run 3			Run 4		
	1	2	3	4	5	6	7	8	9	10	11	12
A	H505- H701	H506- H702	H517- H703	H505- H705	H506- H707	H517- H723	H505- H706	H506- H712	H517- H720	H505- H710	H506- H711	H517- H714
B	H517- H702	H505- H703	H506- H701	H517- H707	H505- H723	H506- H705	H517- H712	H505- H720	H506- H706	H517- H711	H505- H714	H506- H710
C	H506- H703	H517- H701	H505- H702	H506- H723	H517- H705	H505- H707	H506- H720	H517- H706	H505- H712	H506- H714	H517- H710	H505- H711
D	H503- H705	H503- H707	H503- H723	H503- H706	H503- H712	H503- H720	H503- H710	H503- H711	H503- H714	H503- H701	H503- H702	H503- H703
E	H516- H706	H516- H712	H516- H720	H516- H710	H516- H711	H516- H714	H516- H701	H516- H702	H516- H703	H516- H705	H516- H707	H516- H723
F	H522- H710	H510- H711	H513- H714	H522- H701	H510- H702	H513- H703	H522- H705	H510- H707	H513- H723	H522- H706	H510- H712	H513- H720
G	H513- H711	H522- H714	H510- H710	H513- H702	H522- H703	H510- H701	H513- H707	H522- H723	H510- H705	H513- H712	H522- H720	H510- H706
H	H510- H714	H513- H710	H522- H711	H510- H703	H513- H701	H522- H702	H510- H723	H513- H705	H522- H707	H510- H720	H513- H706	H522- H712