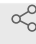




Sep 12, 2022

# Immunochemistry on paraffin sections

[miquel.vila](#)<sup>1</sup><sup>1</sup>Vall d'Hebron Research Institute1 *Works for me* Share[dx.doi.org/10.17504/protocols.io.kxygx9my4g8j/v1](https://dx.doi.org/10.17504/protocols.io.kxygx9my4g8j/v1)[joan.compte](#)

## ABSTRACT

Immunochemistry protocol on paraffin-embedded rat brain sections

## DOI

[dx.doi.org/10.17504/protocols.io.kxygx9my4g8j/v1](https://dx.doi.org/10.17504/protocols.io.kxygx9my4g8j/v1)

## PROTOCOL CITATION

miquel.vila 2022. Immunochemistry on paraffin sections. **protocols.io**  
<https://protocols.io/view/immunochemistry-on-paraffin-sections-cggyttxw>



## LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

## CREATED

Sep 12, 2022

## LAST MODIFIED

Sep 12, 2022

## PROTOCOL INTEGER ID

69880

## MATERIALS TEXT

### **Reagents :**

- TBS 10X : Tris base 121.1g + NaCl 90g in 1L H<sub>2</sub>O.pH 7.4.
- TBS 1X-Triton 0,5%
- Xilen
- Ethanol : 100%, 95%, 70%
- Unmasking buffer epitopes : Citrate solution 10mM pH6.0
- Blocking Buffer : TBS 1X + 5% NGS
- 1st Ab : Diluted in1X PBS +2%NGS
- 2nd Ab : Diluted in1X PBS +2%NGS
- Endogenous peroxidase blocking solution : TBS 1x + 3% H<sub>2</sub>O<sub>2</sub>(30%) + 10% methanol

### 1. Deparaffinization and hydratation :

- 1 Put the slides 30 min in the incubator at 60°C.
- 2 Wash 3x3 min in Xilen.
- 3 Wash 1x10 min in ethanol 100%.
- 4 Wash 1x10 min in ethanol 95%.
- 5 Wash 1x5 min in ethanol 70%.
- 6 Wash 2x5 min in TBS 1X.

### 2. Blocking endogenous peroxidase:

- 7 Put the slides 10 min in endogenous peroxidase blocking solution: TBS 1x + 3% H<sub>2</sub>O<sub>2</sub> + 10% methanol

#### Washing

- 8 Wash 3x5 min in buffer TBS 1X.

#### Antigen retrieval

- 9 Put sections in 200 mL of citrate buffer 10mM, pH 6

- 9.1 Put sections in a boiling water bath for 20 min

- 10 Let sections cool down for 20 min at RT

#### Washing

- 11 Wash 3x5 min in TBS 1X.

#### Blocking

- 12 Put wet paper in a black box for immunohistochemistry
- 13 Gently dry and circle sections with the hydrophobic pen ImmEdge Vector H 4000 without touching the tissue
- 14 Blocking in TBS 1X + 5% NGS (200uL/slides) 1h at RT

#### Primary antibody

- 15 Gently wipe off the water of the slides

16 Put 200ul/slide of TBS 1X + 2% NGS + primary Ab 48/72h (it depends of the Ab) at 4°C (cold room)

#### Washing

17 Wash 3x5min in TBS 1X Buffer.

#### 2ary antibody

18 Gently wipe off the water of the slides

19 Put 200ul/slide TBS 1X + 2% NGS + Secondary Ab 1h at RT.

19.1 At this step it is important to prepare ABC solution and let it, at least, 30 min on the shaker

#### Washing

20 3x5min in TBS 1X

#### Peroxidase staining

21 Incubate 1 hour at RT with ABC solution (4 drops A/B in 10 mL of TBS 1x)

#### Washing

22 3x5min in TBS 1X

#### Developing preparation

23 In aluminium foil: DAB Standard Kit: 1 drop of reagent B in 1 mL of reagent A (gives rise to brown staining)

- 23.1 In aluminium foil: Vector SG: 3 drops Chromogen + 3 drops Hydrogen Peroxide in 5mL PBS (gives rise to blue staining)

#### Developing

- 24 Cover tray plates with aluminium foil

- 24.1 Put 200 uL on each section and put a cardboard box on it to keep darkness for a time ranging from 3-15 minutes depending on the antibody used

- 24.2 Remove with an air-pump equipment and clean the material with bleach

#### Washing

- 25 3x5min in TBS 1X

#### Dehydration

- 26 1 min in ethanol 70%


- 27 1 min in ethanol 95%

- 28 1 min in ethanol 100%.

- 29 2x5 min in xylene

#### Mounting

- 30 Put a line of mounting medium(DPX) by slide. Put the coverslip (washed with ethanol previously) on the slide. Remove bubbles



31 Let dry the slides on the tray in the hood overnight