

VERSION 2

MAY 22, 2023

OPEN ACCESS

DOI:
dx.doi.org/10.17504/protocols.io.6qpvr4n3ogmk/v2

Protocol Citation: jconsi, Ikram Isa, Alexandra Naba 2023. Desalting of Peptides to Prepare for Mass Spectrometry Analysis. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.6qpvr4n3ogmk/v2>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
 We use this protocol and it's working!

Created: May 17, 2023

Last Modified: May 22, 2023

PROTOCOL integer ID:
 82035

Keywords: Peptides, Desalting, Proteomics, Mass Spectrometry

Desalting of Peptides to Prepare for Mass Spectrometry Analysis V.2

jconsi¹, Ikram Isa¹, Alexandra Naba¹

¹University of Illinois Chicago

Human BioMolecular Atlas Program (HuBMAP) Method Development Community

The Matrisom



Alexandra Naba
 University of Illinois Chicago

ABSTRACT

Prior to proteomic analysis, peptide samples are desalted and eluted with freshly prepared 50% acetonitrile, 0.1% trifluoroacetic acid, followed by concentration in a vacuum concentrator. Peptides are then resuspended in freshly prepared 5% acetonitrile, 0.1% formic acid.

Note

- The last step can be conducted at a mass spectrometry facility according to their own preferred methods.
- After desalting, the concentration of the peptide solution can be measured by spectrophotometry.

MATERIALS

Materials:

1. HPLC-grade water Thermo Fisher Scientific Catalog #51140
2. Trifluoro-acetic Acid Thermo Fisher Scientific Catalog #85183
3. Acetonitrile mass spectrometry grade Thermo Scientific Catalog #51101
4. Formic Acid Thermo Scientific Catalog #28905
- 5.

Equipment

Peptide Desalting Spin Columns

Thermo Scientific

89851

<https://www.thermofisher.com/order/catalog/product/89851>

NAME

BRAND

SKU

LINK

6.

Equipment

Low Retention Tubes and Tips

Brand

0000000

NAME

BRAND

SKU

Note

- The maximum volume for spin columns is 300 μ L. If columns ever become unpacked, repeat the step that caused this by reloading the flowthrough and spinning at the recommended speed.

Reagents

Table 1: Reagent Preparation

Table Reagent Preparation	
Reagents	Amount
Priming buffer	100% Acetonitrile
Washing buffer	0.1% TFA in HPLC H_2O
Elution buffer	50% solution with 0.1% TFA in HPLC-grade water
Reconstitution buffer	5% solution with 0.1% Formic acid in HPLC-grade water

Column Preparation

4m

1 Column Preparation

- 1.1 Take a Pierce peptide desalting spin column and remove the white tip (do not remove the screw cap of the tube).

Place in a 2mL tube and spin column at $5000 \times g$ for 00:01:00.

1m

- 1.2 Add 300 μ L of **acetonitrile**. Spin at $5000 \times g$ for 00:01:00 and discard flow-through.

1m

1.3 Repeat this step once

2m

Note




- Note that if columns ever become unpacked, repeat that step as the columns will not work properly if unpacked.

Sample Loading

2m

2 Sample Loading

2.1 Place the spin column in a new low-retention 2 mL tube labeled "flowthrough".


2.2 Load  300 µL of peptide sample into the tube and spin at  3000 x g for  00:01:00 .

1m


Note

- You can save the flow-through to ensure it does not contain any unbound peptides and that peptides are binding to the columns.

2.3 Based on the total sample volume, if:

2.3a) More than  300 µL of sample -> place into a new 2 mL tube and load the remaining volume

2.3b) Less than  300 µL -> reload the flow-through.



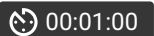
2.4 Spin the sample at  3000 x g for  00:01:00 . Store "FT" at  -80 °C for troubleshooting purpose.

1m

Wash

3m

3 Wash Sample

3.1 Place the spin column in a new low-retention 2mL-tube and load  300 µL of **0.1% TFA in HPLC-grade H₂O**.
Centrifuge at  for . Discard wash flow-through.

1m

3.2 Repeat step 3.1 2 more times.

2m

Note

- Note that if columns ever become unpacked, repeat that step as the columns will not work properly if unpacked.

Sample Elution

1m

4 Elute Samples

4.1 Place the spin column in a new 2mL low-retention tube labeled with the sample name.

4.2 Load  300 µL of **0.1% TFA, 50% acetonitrile in HPLC-grade H₂O**. Spin at  for .


1m

Note

- Note that if columns ever become unpacked, repeat that step as the columns will not work properly if unpacked.

4.3 Transfer the spin column to another 2mL-low retention tube and repeat the step.


4.4 Pool the two elution samples from 4.2 and 4.3. These are the desalted peptides.

4.5 Store at  -20 °C .


Lyophilization and Reconstitution

5 Lyophilization

5.1 To remove reagents incompatible with mass spectrometry place tubes in SpeedVac™ until completely dry.

5.2 Depending on the size of the peptide pellet, resuspend samples with  20 -75 µL of **0.1% formic acid, 5% acetonitrile in HPLC-grade H₂O**.

5.3 Vortex until completely resuspended.

5.4 Peptide concentration can be measured using a spectrophotometer or using peptide concentration measurement kits such as:  Pierce Quantitative Colorimetric Peptide Assay Thermo Fisher Scientific Catalog #23275)