



May 24, 2024

## DNA Extraction from Sterivex Filters - Phenol-Chloroform

 Forked from [DNA EXTRACTION Protocol Template](#)

This protocol is a draft, published without a DOI.



Colleen Kellogg<sup>1</sup>

<sup>1</sup>Hakai Institute

Better Biomolecular Ocea...



Andreas Novotny

University of British Columbia, Hakai Institute

OPEN  ACCESS



External link: <https://hakai.org>

**Protocol Citation:** Colleen Kellogg 2024. DNA Extraction from Sterivex Filters - Phenol-Chloroform. protocols.io <https://protocols.io/view/dna-extraction-from-sterivex-filters-phenol-chloro-dc5w2y7e>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** In development

**We are still developing and optimizing this protocol**

**Created:** November 07, 2023

**Last Modified:** May 24, 2024

**Protocol Integer ID:** 99222

### Disclaimer

Draft!

### Abstract

This protocol is used to extract genomic DNA from 0.22µl sterivex filters using Phenol-Chlorophorm. The protocol is part of the Hakai Institute's pipeline to analyze microbial and environmental DNA from seawater samples and is implemented as a standard procedure for ongoing sampling programs.

Guidelines

MIOP: Minimum Information about an Omics Protocol

MIOP Term	Value
analyses	Nucleic Acid Water Filtration
audience	scientists
broad-scale environmental context	marine biome ENVO_00000447
creator	Colleen Kellogg
environmental medium	sea water [ENVO:00002149]
geographic location	North Pacific Ocean[GAZ:00002410]
hasVersion	1
issued	2017
language	en
license	CC BY 4.0
local environmental context	oceanic epipelagic zone biome [ENVO:01000033]
materials required	Peristaltic Pump
maturity level	Mature
methodology category	Sample collection
personnel required	1
project	Hakai Institutes Marine Biodiversity
publisher	Hakai Institute, Genomics Lab
purpose	Sea water filtration [CHMO:0001640]
skills required	sterile technique   pipetting skills
target	DNA
time required	30

AUTHORS

PREPARED BY All authors known to have contributed to the preparation of this protocol, including those who filled in the template.	AFFILIATION	ORCID (visit <a href="https://orcid.org/">https://orcid.org/</a> to register)	DAT
Colleen Kellogg	Hakai Institute	<a href="https://orcid.org/0000-0003-4048-5316">https://orcid.org/0000-0003-4048-5316</a>	2017

RELATED PROTOCOLS

PROTOCOL NAME AND LINK	ISSUER / AUTHOR	RELEASE / ACCESS DATE
Seawater filtration	Hakai Institute	

This is a list of other protocols which should be known to users of this protocol. Please include the link to each related protocol.

ACRONYMS AND ABBREVIATIONS

ACRONYM / ABBREVIATION	DEFINITION

GLOSSARY

SPECIALISED TERM	DEFINITION

BACKGROUND

This document describes the required protocol to to filter seawater onto a 0.22 micrometer Sterivex filters using paristaltic pump setup.

*Method description and rationale*

Thiswater filtration is part of the standard best - practice method for analysing microbial and environmental DNA from seawater samples at the Hakai Institutes Genome Lab. The method is part of a pipeline that includes seawater filtration, DNA extraction, and amplicon sequencing.

*Spatial coverage and environments of relevance*

The protocol has been used to analyze water samples (commonly 2L) from all parts of the water columns, nearshore and offshore stations on the coast of British Columbia, Canad. The protocol is also designed to target a broad range of organisms and has been used for the downstream analysis of Bacteria, Protists, Zooplankton, and Fish - using 16S rRNA, 18S rRNA, COI and 12S respectively.

*Personnel Required*

1 Technician

*Safety*

Identify hazards associated with the procedure and specify protective equipment and safety training required to safely execute the procedure!

*Training requirements*

Sterile technique, pipetting skills. Work-safe laboratory practices.

*Time needed to execute the procedure*

The set-up, filtration, and clean-up steps each take about 20-30 minutes. If using two peristaltic pumps simultaneously, up to 8 samples may be filtered at the same time.

Materials

DESCRIPTION e.g. filter	PRODUCT NAME AND MODEL Provide the official name of the product	MANUFACTURER Provide the name of the manufacturer of the product.	QUANTITY Provide quantity
Durable equipment			
Content Cell	Content Cell	Content Cell	Content Cell
Content Cell	Content Cell	Content Cell	Content Cell
Consumable equipment			
Content Cell	Content Cell	Content Cell	Content Cell
Content Cell	Content Cell	Content Cell	Content Cell
Chemicals			
Content Cell	Content Cell	Content Cell	Content Cell
Content Cell	Content Cell	Content Cell	Content Cell

Safety warnings

⚠ Several of the chemicals are harmful to humans and the environment. Make a proper risk assessment before starting this procedure. Fumes from the extraction are toxic and waste must be treated with extreme care!

Before start

Read Minimum Information about an Omics Protocol (MIOP) and other recommendations under the "Guidelines" tab.



## Preparations

1



### Safety information

Several of the chemicals are harmful to humans and the environment. Make a proper risk assessment before starting this procedure. Fumes from the extraction are toxic and waste must be treated with extreme care!

Thaw filters on ice, if applicable.

Replace Parafilm on each sample if needed.

Prepare extraction negative - add 1800  $\mu\text{L}$  of SLB (filter-sterilized and UVed) into a new Sterivex.

Seal the bottom with parafilm and top with cap or parafilm, treat the same as the rest of samples.

## Lysis and Incubation

- 2 Thaw filters on ice, if applicable.  
Replace Parafilm on each sample if needed.  
Prepare extraction negative - add 1800  $\mu\text{L}$  of SLB or Longmire's buffer (filter-sterilized and UVed) into a new Sterivex.  
Seal bottom with parafilm and top with cap or parafilm, treat same as rest of samples.

- 3 Add 100  $\mu\text{L}$  lysozyme (125 mg fully dissolved in 1000  $\mu\text{L}$  1 x TE) and 20  $\mu\text{L}$  RNase A (10  $\mu\text{g}/\text{mL}$ : 1  $\mu\text{L}$  in 999  $\mu\text{L}$  1 x TE) to each filter.  
Reseal the top with Parafilm or a luer cap.  
Incubate Sterivex in a rotating incubator at 37  $^{\circ}\text{C}$  for 01:00:00 (in labelled 50 ml falcon tubes).

1h



- 4 Add 100  $\mu\text{L}$  Proteinase K and 100  $\mu\text{L}$  [w] 20 Mass / % volume SDS to each filter.  
Reseal using Parafilm or a luer cap.  
Incubate at 55  $^{\circ}\text{C}$  for 1-2 hours in a rotating incubator (in 50ml falcon tubes; be sure to keep the same falcon tube for the same sterivex throughout these incubation steps).



- 5 Remove lysate from sterivex filter using a 10cc syringe into a labeled 15 ml falcon tube, use a new syringe for each sample.  
*Do not pull up on the plunger while sterivex is attached. You will likely need to do multiple plunges to completely empty the Sterivex, so unscrew the syringe from the filter, pull back the plunger, reattach to Sterivex, and push and hold to plunge liquid out of Sterivex.*  
Repeat as needed until (nearly) all liquid and foam are out of Sterivex.  
Add 1 mL sucrose lysis buffer (SLB), rotate Sterivex with your hands several times to rinse the filter out and then plunge out liquid as described above.  
Pool this with the lysate in the 15 ml tube.

## EXTRACTION & CENTRIFUGATION

6

6m

### Safety information

Phenol:Chloroform is toxic when inhaled and has to be handled in a fume hood.
















In the fume hood, add an equal volume (about 3 mL ) of Phenol:Chloroform:IAA (25:24:1), pH 8.0 to the lysate tube.

Invert for 10 seconds by hand to mix.



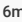

Spin at 2500 g for 00:06:00 or until the aqueous layer is clear. Wait at least 10 minutes before opening the centrifuge.

Transfer aqueous (top) layer into a new 15 ml falcon tube.



- 7 Add an equal volume (approx  3 mL ) of Chloroform:IAA (24:1) to the tube containing the aqueous layer.  
Invert for 10 seconds by hand.  
Spin at 2500 g for  00:06:00 or until aqueous layer is clear with no debris. Wait at least 10 minutes before opening the centrifuge.  
Transfer aqueous layer into a labeled Amicon Ultra centrifuge tube (UFC801096, EMD Millipore).  
- 8 Top up Amicon with  1-2 mL of 1 x TE buffer.   
Spin at 3500 g for 10-15 minutes.  3500 x g, Room temperature  
Check to make sure there is less than 1ml liquid left in Amicon at the end of this (if not, fill up with  2 mL 1xTE and spin again).
- 9 Add  2 mL TE buffer to Amicon and spin at 3500 g for  00:06:00 .   
Remove filtrate. 
- 10 Repeat Step 8 twice more (total of 3 washes with  2 mL TE).  
For the last wash, spin until 200 – 500 µl remain in Amicon (typically 8-10 minutes).  
Note the final volume and transfer to a labeled 1.5 µl Eppendorf tube.  
Rinse Amicon sides with 50 µL of 1xTE and pool with therest of the sample in the a labelled 1.5mL Eppendorf tube.

## DNA QUANTIFICATION AND STORING

- 11 If desired, aliquot  50 µL from the final sample volume into a 1.5ml eppendorf tube to use as working stock and place in  -20 °C freezer.  
Place the remaining DNA stock in the  -70 °C freezer for long term storage.
- 12 Quantify DNA stock using Qubit (following manufacturer's instructions).  
Use  2 µL of stock DNA when quantifying.

## DNA NORMALISATION

13

### Note

Optional, but does tend to yield good success during amplification

- 14 Make a normalized 2.5 ng/µl DNA stock for PCR :  
**Calculate volume of DNA for 2.5 ng/µL stock** (there is a spreadsheet for this)

$$\left(2.5 \frac{\text{ng}}{\mu\text{L}}\right) (25 \mu\text{L}) = \left([DNA] \frac{\text{ng}}{\mu\text{L}}\right) (x \mu\text{L})$$

Where **[DNA]** is the DNA concentration from Qubit and **x** is the number of uL of the DNA stock that you need to add for a 25 uL volume of a 2.5 ng/uL DNA stock.

**Calculate volume of water for 2.5 ng/µL stock**

$$25 \mu\text{L} - x \mu\text{L of DNA} = y \mu\text{L of water}$$

- 15 Place the 2.5 ng/µL DNA stock in the  -20 °C freezer.



## Protocol references

This protocol is modified from **DNA extraction from 0.22 microM Sterivex filters and cesium chloride density gradient centrifugation**. Jody J. Wright, Elena Zaikova, David A. Walsh and Steven J. Hallam

<https://doi.org/10.3791/1352>