



Mar 30, 2022

Chemically Competent Cell Preparation [Calcium Chloride Method]

Harry Akligoh¹, Jenny Molloy²¹Open Bioeconomy Lab; ²University of Cambridge

1

dx.doi.org/10.17504/protocols.io.261ge4k8dv47/v1[Reclone.org \(The Reagent Collaboration Network\)](https://reclone.org)[Open Bioeconomy Lab](#)Harry Akligoh
Open Bioeconomy Lab

This protocol is a low-cost approach for preparing competent cells using 0.1M CaCl₂ solution which is useful in cases where commercial competent cells are unaffordable, it is challenging to receive products on dry ice or you do not have an ultra low temperature freezer.

DOI

dx.doi.org/10.17504/protocols.io.261ge4k8dv47/v1

Harry Akligoh, Jenny Molloy 2022. Chemically Competent Cell Preparation [Calcium Chloride Method]. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.261ge4k8dv47/v1>



Competent Cells, Molecular Biology

_____ protocol ,

Jun 25, 2021

Mar 30, 2022

51054

Maintain aseptic technique at all point during this experiment

1. 0.1M ice cold CaCl_2 Solution
2. Sterile microcentrifuge tubes
3. Pipettes
4. LB broth
5. LB agar
6. Cold centrifuge (4 deg Celsius)
7. E. coli

Keep materials for this experiment cold as much as possible in situations where you're not using a cold centrifuge

Disinfect your bench and pipettes using bleach or 70% isopropanol solution

Overnight plating of E. coli

1d


- 1 Plate an inoculum of *Escherichia coli* (e.g. BL21(DE3), JM109) glycerol stock on LB agar ^{1d} without antibiotics.

Addgene. Streaking and Isolating Bacteria on an LB Agar Plate.

<https://www.addgene.org/protocols/streak-plate/>

There are no antibiotics in the culture, so it will be prone to contamination. Use good aseptic technique at all times!



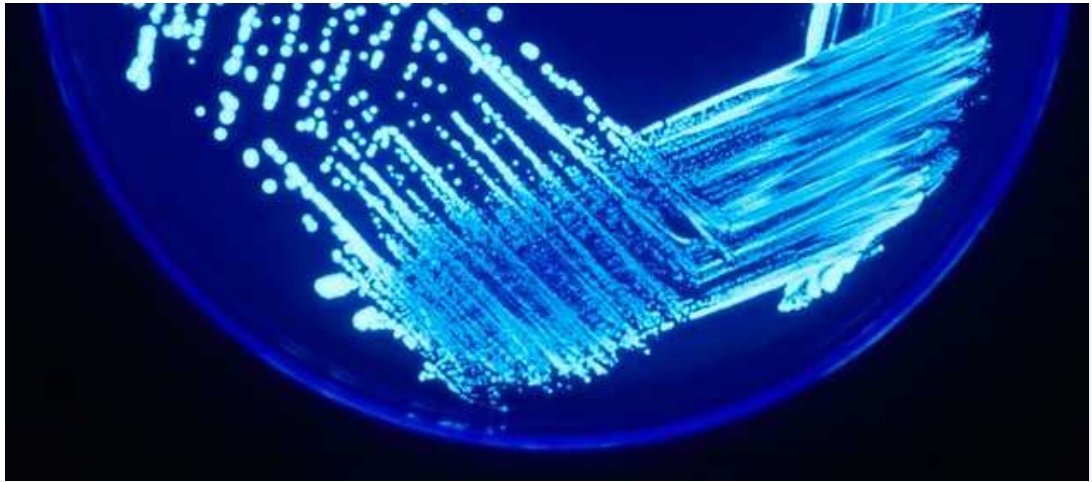

Source: CDC Public Health Image Library via [Wikimedia Commons](https://commons.wikimedia.org/wiki/File:Petri_dish_with_bacterial_culture.jpg), in the Public Domain

2 Incubate the plate 🌡️ **37 °C** ⌚ **Overnight**

2h

Invert the plate(s) downward when incubating overnight. This is to prevent condensed water in the plates flooding the LB agar and smearing your growing bacteria.

Overnight culture 14h

3 Pick and transfer a discrete single colony of *E. coli* into 10 ml LB broth using a sterile inoculation loop or toothpick.

14h

4 Incubate the culture overnight in a shaker incubator at 🌀 **140 rpm, 37°C, 12:00:00**

Harvesting cells and Washing of Cells 1h 8m

5 1. Sub-culture the overnight culture by adding 🧴 **500 µL** into 🧴 **50 mL** LB broth (i.e. a 1:100 dilution) and shake at 🌀 **140 rpm, 37°C, 02:00:00**

1h 8m

6 Monitor OD₆₀₀ every 30 min using a photometer or [McFarland Standards](https://www.fishbase.org/identification/identification-mcfarland-standards) until OD₆₀₀ = 0.3

Take particular care when aspirating the culture that it does not get contaminated as your media does not contain antibiotics. Aspirate over a Bunsen burner or in a microbiological safety cabinet or use other measures to ensure excellent aseptic technique.

- 7 Keep the culture **On ice** and only mix gently while harvesting the cells, avoid vigorous shaking.

- 7.1 Pipette **1 mL** culture into a **1.5 mL** microcentrifuge tube and centrifuge at **8000 rpm, 4°C, 00:05:00** . ^{5m}

If you have a centrifuge that can fit larger tubes, this will speed up the process and avoid multiple steps of pipetting and centrifugation.

It is important to keep the cells close to **0 °C** during this step and all subsequent steps. If you do not have a refrigerated centrifuge, you can try:

1. putting the rotor in the fridge or freezer prior to use
2. putting the whole centrifuge in a cold room or walk-in fridge
3. keeping centrifugation time to the minimum needed to form a pellet that is compact enough not to break apart while you discard the supernatant.

- 8 After each centrifugation, discard the supernatant making sure not to dislodge the cell pellets ^{1m} and repeat previous step as needed to harvest more cells. If using **1.5 mL** tubes, we recommend repeating 3-5 times.


- 9 Resuspend the harvested cells in **1 mL** ice cold **0.1 Molarity (M)** calcium chloride (CaCl₂) solution.

- 10 Incubate **On ice** for **00:30:00**




30m

11 Centrifuge to pellet the cells at  **8000 rpm, 4°C, 00:05:00** 5m

12 Finally, gently resuspend the cells in  **500 µL** ice cold  **0.1 Molarity (M)** calcium chloride (CaCl₂) solution. 1m

You can also supplement the CaCl₂ solution with 15 % v/v glycerol if you wish to store the cells at  **-80 °C** . They should then remain competent for at least 1 year. If you do not add glycerol, you must use the cells immediately.

Competent cells are not generally stable at  **-20 °C**

13 Pipette  **50 µL** into sterile, ice-cold  **1.5 mL** microcentrifuge tubes to use immediately 1m
for transformation or for storage at  **-80 °C** .

It is good practice to quantify the transformation efficiency of your competent cells each time you make a batch by transforming them with a known amount (say 1 µL of 50 ng/µL) of plasmid that contains a positive selection marker and counting the number of transformant colonies. This will help track any discrepancies and ensure you get consistent results.