



Halo assay to assess mitophagy

COMMENTS 0

DOI

dx.doi.org/10.17504/protocols.io.dm6gpjzm8gzp/v1

nguyen.tha¹

¹Laboratory of Michael Lazarou, Walter and Eliza Ha Il Institute of Medical Research, Parkville, Victoria, A ustralia

WORKS FOR ME

2

ASAP workspace



nguyen.tha

ABSTRACT

This protocol describes how to assess mitophagy using Halo assay developed by Mizushima lab (DOI: 10.7554/eLife.78923).

ATTACHMENTS

585-1231.docx

DOI

dx.doi.org/10.17504/protocols.io.dm6gpjzm8gzp/v1

PROTOCOL CITATION

nguyen.tha 2022. Halo assay to assess mitophagy. **protocols.io** https://dx.doi.org/10.17504/protocols.io.dm6gpjzm8gzp/v1

KEYWORDS

Halo assay, mitophagy

LICENSE

This is an open access protocol distributed under the terms of the <u>Creative Commons</u>

<u>Attribution License</u>, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Dec 05, 2022

LAST MODIFIED

Dec 08, 2022

OWNERSHIP HISTORY

Dec 05, 2022 renuka.s

Dec 07, 2022 nguyen.tha

👸 protocols.io

1

Citation: nguyen.tha Halo assay to assess mitophagy https://dx.doi.org/10.17504/protocols.io.dm6gpjzm8gzp/v1

PROTOCOL INTEGER ID

73537

MATERIALS TEXT

Buffers and reagents:

Growth media:

А	В
DMEM with 10% FBS	
Glucose (Sigma, G8769)	4.5 g/l
1x GlutaMAXTM (ThermoFisher, 350500	
1x MEM NEAA (ThermoFisher, 11140-05	
HEPES (1688449)	25 mM

- Antimycin-A (Sigma, A8674; made up in 100% Ethanol to 20 mg/ml),
- Oligomycin (Calbiochem, 495455; made up in DMSO to 10 mg/ml)
- qVD (MedChemExpress, HY-12305; made up in DMSO to 10 mM)
- TMR-conjugated Halo ligand (Promega, G8251)
- Lysis buffer. (diluted from 4x LDS (NP007; ThermoFisher)

A
1x LDS
0.1 M DTT

can be aliquoted and stored at -20 or -80°C.

- 4-12% Bis-Tris NuPAGE gels (ThermoFisher)
- NuPAGETM Antioxidant (NP0005, ThermoFisher; use 0.5 ml/ 200 ml of gel running buffer)
- 20x NuPAGETM MOPS SDS running buffer (NP001, ThermoFisher)
- 20x NuPAGE transfer buffer (NP00061, ThermoFisher)

PVDF destain.



- 1x PBS
- 1x PBS/0.1% Tween20 (PBS/Tween)
- Blocking buffer. 5% skim milk in PBS/Tween (make fresh)
- VCP (Cell Signaling, 2649), HALO (Promega, G9211)

- **⋈** 45% D-()-Glucose **Sigma Catalog #G8769**
- **⊠** GlutaMAX™ Supplement **Thermo Fisher Catalog #35050061**
- X MEM Non-Essential Amino Acids Solution (100X) Thermo Fisher Scientific Catalog #11140050
- 🔀 HEPES Buffer 1M Solution Cell Culture Grade MP Biomedicals Fisher Scientific Catalog #ICN1688449
- X Antimycin A from Streptomyces sp. Sigma Aldrich Catalog #A8674
- **☒** Oligomycin **Millipore Sigma Catalog #495455**
- **☒** Q-VD-OPh **MedChemExpress Catalog #HY-12305**
- ★ HaloTag(R) TMR Ligand, 30ul Promega Catalog #G8251
- X NuPAGE™ LDS Sample Buffer Thermo Fisher Scientific Catalog #NP0007
- NuPAGE Antioxidant Thermo Fisher Scientific Catalog #NP0005
- X NuPAGE™ MOPS SDS Running Buffer (20X) Thermo Fisher Scientific Catalog #NP0001
- X Tris-Glycine Transfer Buffer Invitrogen Thermo Fisher Catalog #NP00061
- X VCP (7F3) Rabbit mAb Cell Signaling Technology Catalog #2649
- X Anti-HaloTag(R) Monoclonal Antibody Promega Catalog #G9211

ATTACHMENTS

<u>585-</u> 1231.docx

4h 57m 10s

Procedures

- Generating cells expressing mitochondrially targeted Halo-GFP using pSu9-Halo-mGFP from Mizushima lab (Addgene #184905; DOI: 10.7554/eLife.78923).
- 2 Seed HeLa cells the day before the treatment day in 6 well plates.
 - Each well contained <u>A</u> 2 mL of growth media;
 - Seed 350,000 cells for penta KO expressing BFP-Parkin and GFP-OPTN or -NDP52 and 380,000 cells for other knockout lines such as ATG13 KO/penta KO expressing GFP-NDP52;
 - Adjust the number of cells of other cell lines, so that the next day they are all in similar confluency with penta KO expressing BFP-Parkin and GFP-OPTN or -NDP52.

Note			

3 The next day, make sure the seeded cells are spreading out (not concentrated in the middle of the well because this can affect the results).

protocols.io

3

Citation: nguyen.tha Halo assay to assess mitophagy https://dx.doi.org/10.17504/protocols.io.dm6gpjzm8gzp/v1

	4	Aspirate off the old media and treat each well with 2 1 mL growth media containing TMR-conjugated Halo ligand.	
		Note	
ı	5	Incubate in a normal TC incubator for 00:20:00	20m
l	6	Aspirate off the media and wash thoroughly.	
	6.1	Aspirate off the media and wash thoroughly with 1x PBS. (1/2)	
	6.2	Aspirate off the media and wash thoroughly with 1x PBS. (2/2)	
	6.3	Harvest the non-mitophagy-induced samples immediately by scraping (see step 8).	
	7 1	For mitophagy-induced samples, treat each well with Δ 2 mL of growth media containing [M] 4 micromolar (μM) Antimycin A, [M] 10 micromolar (μM) Oligomycin and [M] 10 micromolar (μM) QVD for desired period.	
		Note	
	protoco	ols.io 4	


Citation: nguyen.tha Halo assay to assess mitophagy https://dx.doi.org/10.17504/protocols.io.dm6gpjzm8gzp/v1

8	After the treatment, harvest the cells	
	Note	
8.1	Pre-chill eppies and 1x PBS On ice	
	Note	
8.2	Aspirate the media thoroughly from the wells, wash the wells with 1 mL _ of cold 1x PBS*, aspirate off	
8	the PBS and add 🔼 1 mL of cold 1x PBS.	
	Note	
8.3	After that, use a plastic cell scraper to scrape all the cells off the wells (I use one scraper for each well. You can wash and reuse them again). Transfer the cells-containing PBS to eppies.	
8.4	Centrifuge the eppies at 3000 x g, 4°C, 00:02:00 . Aspirate off PBS.	2m
₩		

protocols.io

5

 $\textbf{Citation:} \ \ nguyen. tha \ \ Halo \ assay \ to \ assess \ mitophagy \ \underline{https://dx.doi.org/10.17504/protocols.io.dm6gpjzm8gzp/v1}$

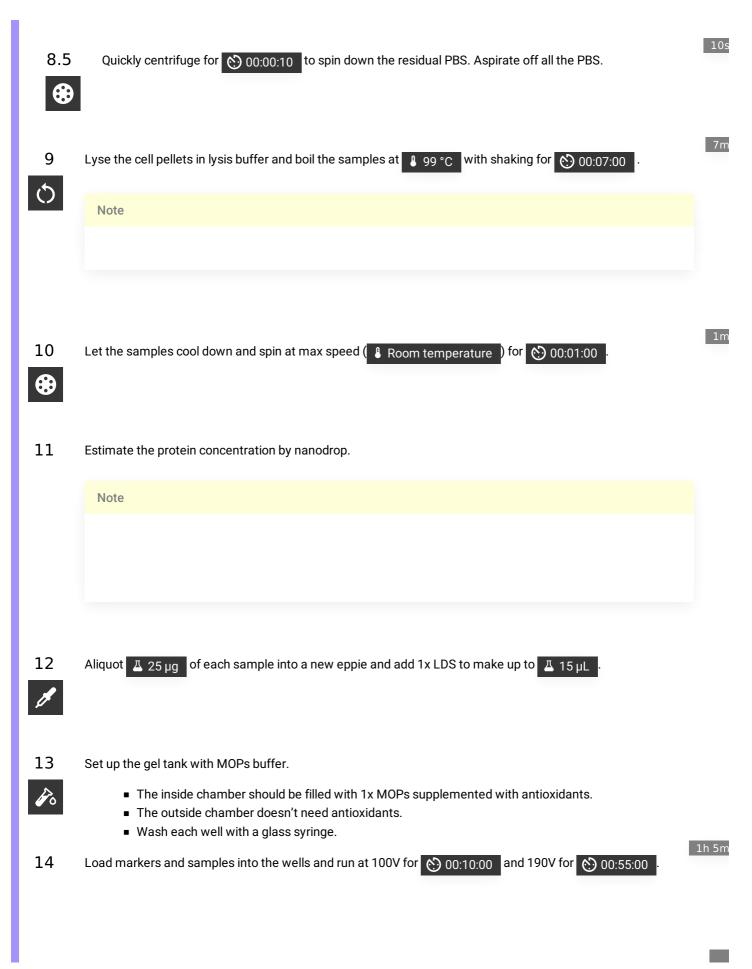



Citation: nguyen.tha Halo assay to assess mitophagy https://dx.doi.org/10.17504/protocols.io.dm6gpjzm8gzp/v1

