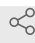




Sep 23, 2022

# 🌐 Assay for Dual Rab GTPase binding to the LRRK2 Armadillo Domain

Claire Y Chiang<sup>1</sup>, Suzanne R Pfeffer<sup>1</sup><sup>1</sup>Department of Biochemistry, Stanford University, Stanford, CA 94305-53071 *Works for me* Share[dx.doi.org/10.17504/protocols.io.81wgbyzvvpk/v1](https://dx.doi.org/10.17504/protocols.io.81wgbyzvvpk/v1) Suzanne R Pfeffer

## ABSTRACT

The LRRK2 Armadillo domain contains multiple Rab GTPase binding sites. To show that the sites can be occupied simultaneously, we use this assay. The idea is to immobilize Rab8A, bind Armadillo domain, and test if phosphoRab10 can bind to Rab8A-immobilized Armadillo domain.

## DOI

[dx.doi.org/10.17504/protocols.io.81wgbyzvvpk/v1](https://dx.doi.org/10.17504/protocols.io.81wgbyzvvpk/v1)

## PROTOCOL CITATION

Claire Y Chiang, Suzanne R Pfeffer 2022. Assay for Dual Rab GTPase binding to the LRRK2 Armadillo Domain. **protocols.io**  
<https://protocols.io/view/assay-for-dual-rab-gtpase-binding-to-the-lrrk2-arm-cg2ztyf6>



## FUNDERS ACKNOWLEDGEMENT

Aligning Science Across Parkinson's  
Grant ID: 000463

## MANUSCRIPT CITATION

 please remember to cite the following publication along with this protocol<https://doi.org/10.7554/eLife.79771>

## LICENSE

\_\_\_\_\_ This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Sep 23, 2022

LAST MODIFIED

Sep 23, 2022

PROTOCOL INTEGER ID

70457

MATERIALS TEXT

His-MST3 protein (pET15b 6HIS MST3 TV1; MRC-PPU DU62980)

His-Rab8A Q67L full length

His-GFP-Rab10 Q68L 1-181

Reaction buffer: 50 mM HEPES **pH 8**, 150 mM NaCl, 5 mM MgCl<sub>2</sub>, 0.2 mM TCEP, 100 μM GTP, 2 mM ATP, 5% (v/v) glycerol

## Dual Rab GTPase binding to LRRK2 Armadillo Domain

3h 36m 30s

- 1 Phosphorylate His-Rab10 Q68L 1-181 with His-MST3 kinase at a molar ratio of 1:3 (kinase:substrate) at **30 °C** **02:00:00** in reaction buffer. See below for details.

2h

Axel Knebel, Kerryen Berndsen, Pawel Lis, Paul Davies, Dario R Alessi.  
Expression and purification of Rab8A (1-181) stoichiometrically phosphorylated at pThr72 (the LRRK2 site).  
<http://dx.doi.org/10.17504/protocols.io.butinwke>

- 2 Pellet **50 μL** glutathione agarose slurry at **3000 rpm, 4°C, 00:05:00**.

5m

- 3 Add GST-Rab8A Q67L to glutathione beads to achieve a concentration of **6 micromolar (μM)** in a total volume of **50 μL** reaction buffer. Incubate at **Room temperature** for **00:30:00** on a rotator.

30m

- 4 Pellet beads by spinning at **3200 x g, Room temperature, 00:00:30** and discard

30s

supernatant.

- 5 Add His-LRRK2 Armadillo domain 1-552 ([M]**10 micromolar (μM)** final in **50 μL**) or buffer<sup>30m</sup> alone and incubate at **Room temperature** for **00:30:00** on a rotator.
- 6 Pellet beads by spinning at **3200 x g, 00:00:30** and discard supernatant. <sup>30s</sup>
- 7 Add phosphorylated His-Rab10 Q68L 1-181 ([M]**4 micromolar (μM)** in **50 μL**) and<sup>30m</sup> incubate at **Room temperature** for **00:30:00** on a rotator.
- 8 Wash beads 2X with **500 μL** reaction buffer.
- 9 Elute protein from beads using **50 μL** elution buffer (reaction buffer + **50 millimolar (mM)** reduced glutathione).
- 10 Pellet beads by spinning at **3200 x g, 00:00:30** and collect supernatant. <sup>30s</sup>
- 11 Analyze eluate by SDS-PAGE and immunoblot for phosphoRab10; image blots with Li-COR, and quantify bands using ImageJ (see below for details).

Francesca Tonelli, Dario Alessi. Quantitative Immunoblotting  
Analysis of LRRK2 Signalling Pathway.  
<http://dx.doi.org/10.17504/protocols.io.bsgnrv6>