



Sep 06, 2022

# Lentivirus production for primary neuron transduction

Miguel Da Silva Padilha<sup>1</sup>, Irina Dudanova<sup>1</sup>, Itika Saha<sup>2</sup>, F. Ulrich Hartl<sup>2</sup>, Mark S. Hipp<sup>3,4</sup>


<sup>1</sup>Molecular Neurodegeneration Group, Max Planck Institute of Neurobiology, 82152 Martinsried, Germany;

<sup>2</sup>Department of Cellular Biochemistry, Max Planck Institute of Biochemistry, Am Klopferspitz 18, 82152 Martinsried, Germany;

<sup>3</sup>Department of Biomedical Sciences of Cells and Systems, University Medical Center Groningen, University of Groningen, Antonius Deusinglaan, 1, 9713 AV Groningen, The Netherlands;

<sup>4</sup>School of Medicine and Health Sciences, Carl von Ossietzky University Oldenburg, 26129 Oldenburg, Germany

1 Works for me

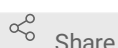

Share

[dx.doi.org/10.17504/protocols.io.kqdg39eo1g25/v1](https://dx.doi.org/10.17504/protocols.io.kqdg39eo1g25/v1)

Felix Kraus

## ABSTRACT

This protocol describes the production of lentiviruses to transduce mouse primary neurons and has to be performed in a biosafety level 2 laboratory

## DOI

[dx.doi.org/10.17504/protocols.io.kqdg39eo1g25/v1](https://dx.doi.org/10.17504/protocols.io.kqdg39eo1g25/v1)

## EXTERNAL LINK

<https://www.biorxiv.org/content/10.1101/2022.02.18.481043v1.full>

## PROTOCOL CITATION

Miguel Da Silva Padilha, Irina Dudanova, Itika Saha, F. Ulrich Hartl, Mark S. Hipp 2022. Lentivirus production for primary neuron transduction. **protocols.io** <https://protocols.io/view/lentivirus-production-for-primary-neuron-transduct-cf7ytrpw>

## MANUSCRIPT CITATION please remember to cite the following publication along with this protocol

Itika Saha, Patricia Yuste-Checa, Miguel Da Silva Padilha, Qiang Guo, Roman Körner, Hauke Holthusen, Victoria A. Trinkaus, Irina Dudanova, Rubén Fernández-Busnadiego, Wolfgang Baumeister, David W. Sanders, Saurabh Gautam, Marc I. Diamond, F. Ulrich Hartl, Mark S. Hipp *bioRxiv* 2022.02.18.481043; doi: <https://doi.org/10.1101/2022.02.18.481043>

#### LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

#### CREATED

Sep 05, 2022

#### LAST MODIFIED

Sep 06, 2022

#### PROTOCOL INTEGER ID

69592

#### SAFETY WARNINGS

This protocol describes the production of lentiviruses to transduce mouse primary neurons and has to be performed in a biosafety level 2 laboratory

- 1 Expand HEK293T cells (LentiX 293T cell line, Takara) for lentiviral packaging to 70-85% confluency in DMEM Glutamax (+4.5 g/L D-glucose, - pyruvate) supplemented with 10% Fetal Bovine Serum (FBS)(Sigma), 1% G418 (Gibco), 1% Non-Essential Amino Acids (Thermo Fisher) and 1% HEPES (Biomol).

NOTE: Only low passage cells should be used.

- 2 For lentiviral production, plate cells ( $\sim 4.8 \times 10^6$ ) in a three-layered 525 cm<sup>2</sup> flask (Falcon).
- 3 On the following day, transfect cells with 59.52  $\mu$ g expression plasmid, 35.2  $\mu$ g packaging plasmid psPAX2 (RRID:Addgene\_12260) and 20.48  $\mu$ g envelope plasmid pVsVg (gift from Dieter Edbauer) using 345.6  $\mu$ L TransIT-Lenti transfection reagent (Mirus) in 9.6 mL DMEM without FBS. <sup>1d</sup>
- 4 Incubate transfection mix for 20 min at room temperature and exchange the cell medium in <sup>20m</sup> the meantime.
- 5 Add 10 mL of transfection mix to the flask, followed by incubation overnight.

- 6 Exchange the medium on the following day. 1d
- 7 After 48-52 h, collect culture medium containing the viral particles and centrifuge for 10 min at 1,200 x g. 2d
- 8 Filter the supernatant through 0.45 µm pore size filters using 50 mL syringes and add 20 mL Lenti-X concentrator (Takara) to filtered supernatant.
- 9 Incubate overnight at 4 °C and centrifuge samples at 1,500 x g for 45 min at 4 °C. 1d
- 10 Remove the supernatant and resuspend the lentivirus pellet in 150 µL TBS-5 buffer (50 mM Tris-HCl pH 7.8, 130 mM NaCl, 10 mM KCl, 5 mM MgCl<sub>2</sub>).
- 11 After aliquoting, store lentivirus at -80 °C.
- 12 Thaw virus preparation immediately before adding to freshly prepared neuronal culture medium.
- 13 Remove a fifth of the medium from cultured neurons and add the equivalent volume of virus-containing medium.

NOTE: Volume of concentrated virus to be added to the neurons and the length of transduction should be determined empirically.