




Sep 24, 2021

Protocol for Mouse IP Anesthesia (Ketamine-Xylazine)

Haowei Li¹, Markus Holzl¹, Mohsen Khosravi-Maharlooeei¹, Nichole Danzl¹, Austin Chen¹, Megan Sykes¹

¹Columbia Center for Translational Immunology, Columbia University, New York

1 Works for me

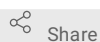

Share

dx.doi.org/10.17504/protocols.io.bvxtn7nn

Human Islet Research Network



Lili Liang

ABSTRACT

This protocol details the constitution of our anesthesia cocktail that we use for sedating our mice when performing surgeries. The surgeries we typically perform are thymectomies and implantation of material under the kidney capsule (through open abdomen or minimally invasive approaches).

The changes of this protocol was made by Austin Chen on Markus Holzl original protocol.

Corresponding Authors

Mohsen Khosravi-Maharlooeei

Email: mkm2182@cumc.columbia.edu

Austin Chen

Email: ac4274@cumc.columbia.edu

Tel: 425-283-6900

DOI

dx.doi.org/10.17504/protocols.io.bvxtn7nn

PROTOCOL CITATION

Haowei Li, Markus Holzl, Mohsen Khosravi-Maharlooeei, Nichole Danzl, Austin Chen, Megan Sykes 2021. Protocol for Mouse IP Anesthesia (Ketamine-Xylazine). **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.bvxtn7nn>

LICENSE

This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Jun 18, 2021

LAST MODIFIED

Sep 24, 2021

PROTOCOL INTEGER ID

50899

MATERIALS TEXT

Materials:

Anesthesia is a mixture of ketamine and xylazine in 1x PBS.

- **Ketamine:** Ketathesia (100 mg/mL; cat: 056344; Henry Schein)
- **Xylazine:** Xylazine sterile solution (20 mg/mL; Columbia ICM)
- **1x PBS:** Phosphate Buffered Solution (1X) ([cat: BW17516F24; Lonza](#))
- **0.22 um filter:** Syringe Filter Unit, 0.22 µm, polyethersulfone ([Millex-GP; cat: SLGPR33RS](#))

Preparation:

Ketamine stock: 100 mg/mL

Xylazine stock: 20 mg/mL

Desired ketamine concentration: 0.08 mg/gram of mouse

Desired xylazine concentration: 0.012 mg/gram of mouse

- We will make anesthesia solution to dose **10 µl** anesthesia/gram of mouse for IP injections

For 10 mL of anesthesia:

1. **0.8 mL** Ketamine (100 mg/mL) + **0.6 mL** Xylazine (20 mg/mL) + **8.6 mL** 1x PBS
2. Filter solution with 0.22 um Millex-GP filter

Calculations:

- Ketamine: 0.8 mL * 100 mg/mL in 10 mL → 8 mg/mL Ketamine → 0.08 mg/10 uL
- Xylazine 0.6 mL * 20 mg/mL in 10 mL → 1.2 mg/mL Xylazine → 0.012 mg/10 uL

For Intraperitoneal (IP) Approach:

- 1 Use 10 uL/gram of anesthesia

Comments:

1. It takes time until the mice fall asleep

2. Very often, mice do not fall asleep but move around like zombies. This requires another small dose of IP ketamine-xylazine injection or a small dose of IV injection

- a. When doing IP injection again, we have injected IP with a volume of about 1/3 of original volume

- b. For IV approach after IP, it takes very little anesthesia for them to fall asleep. We usually use around 10-20 uL. Inject this slowly into the mouse. The mouse should fall asleep rapidly. It is very easy to overdose with IV injection so be careful.