




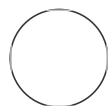
APR 20, 2023

## 2-step PCR mixture and conditions (Barcoded-head primers for seqs pooling)

 Forked from [2-step PCR mixture and conditions \(Barcoded-head primers for seqs pooling\)](#)

Yin-Tse Huang<sup>1</sup>, Tsu-Chun Hung<sup>1</sup>

<sup>1</sup>Kaohsiung Medical University



Tsu-Chun Hung

### ABSTRACT

#### PCR mixture and condition (2X SUPERGREEN PCR MASTER MIX)

### OPEN ACCESS

#### DOI:

[dx.doi.org/10.17504/protocols.io.yxmvm263og3p/v1](https://dx.doi.org/10.17504/protocols.io.yxmvm263og3p/v1)

**Protocol Citation:** Yin-Tse Huang, Tsu-Chun Hung 2023. 2-step PCR mixture and conditions (Barcoded-head primers for seqs pooling).

**protocols.io**

<https://dx.doi.org/10.17504/protocols.io.yxmvm263og3p/v1>

#### MANUSCRIPT CITATION:

Herbold CW, Pelikan C, Kuzyk O, Hausmann B, Angel R, Berry D, Loy A. 2015. A flexible and economical barcoding approach for highly multiplexed amplicon sequencing of diverse target genes. *Front. Microbiol.* [Internet] 6:731. Available from: <http://dx.doi.org/10.3389/fmicb.2015.00731>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** Working  
We use this protocol and it's working

**Created:** Apr 19, 2023

**Last Modified:** Apr 20, 2023

**PROTOCOL integer ID:**  
80784

- 1 Wear glove, clean up the working bench w. 1% bleach

## For 1' PCR head-primers


- 2 Prepare 1' PCR master mixutre for **head-primers** (prepare 1.2X of solutions for pipetting error if needed)



PCR mixture for head-primers for each reaction

A	B	C	D
Component	Volume	Volume (1.2X)	Final conc.
Forward Primer (10 µM)	1 µl	1.2 µl	0.4 µM
Reverse Primer (10 µM)	1 µl	1.2 µl	0.4 µM
ZEJU PCR Master Mix	12.5 µl	15 µl	-
ddH2O	9.25 µl	11.1 µl	-
Total volume	23.75 µl	28.5 µl	-

### Note

Negative control ALWAYS NEEDED! For example, if you have 5 PCR reactions to run, prepare master mixture for 6 reactions (5 DNA template + 1 negative control).

- 3 Mix the 1' PCR master mixture gently by pippeting. Quick spin the tube.
- 4 Transfer  23.75 µL 1' PCR master mixutre in 8-strip PCR tubes.

- 5 Add  1.25  $\mu\text{L}$  DNA template in 8-strip PCR tubes, resulting in a  25  $\mu\text{L}$  reaction mixture for 1' PCR.



#### Note

**Negative control** contains only  23.75  $\mu\text{L}$  master mixture but not DNA template

- 6 Mix the reaction mixture gently by tapping the tubes. Quick spin the tubes.

- 7 Carry out PCR using the following condition:

1' PCR condition for **head-primers**

A	B	C	D
Step	Temp	Sec	Cycle
Initial denaturation	95 °C	180	
Denaturation	95 °C	30	20-35 cycles
Annealing	60-66 °C varied (b)	30	
Extension	72 °C	180	
Final extension	72 °C	420	
Preservation	Preservation	4 °C	$\infty$

b. Annealing varied, **60-66°C** is working; Refer to 1' PCR primers for annealing temperature


c. 1kb ~ 1min extension; enough time allow full extension of sequence

### 7.1 1' hear-primers used in Huang lab

A	B	C	D
Name	Sequence	Tm°C	CG%


A	B	C	D
NS1B1ngs_H1	GCTATGCGCGAGCTGCcctngttgatyctgccag t	71.7	60
LR5_H1	GCTATGCGCGAGCTGCTcctgagggaaacttcg	70.2	60.6
EF1-526F_H1	GCTATGCGCGAGCTGCgtcgytygytygghca ygt	71	59.3
EF1-2218R_H1	GCTATGCGCGAGCTGCatgacaccracrgcracr gtytg	72.2	60.3

## 8 Carry out **electrophoresis** for inspection of DNA products

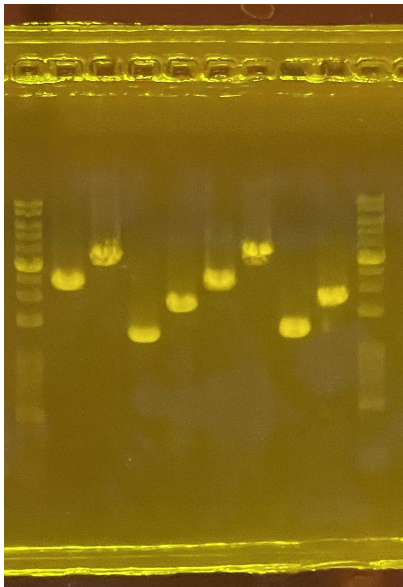

Gel before markdown

## 9 Markdown wells and upload the pictures to the Lab Google drive

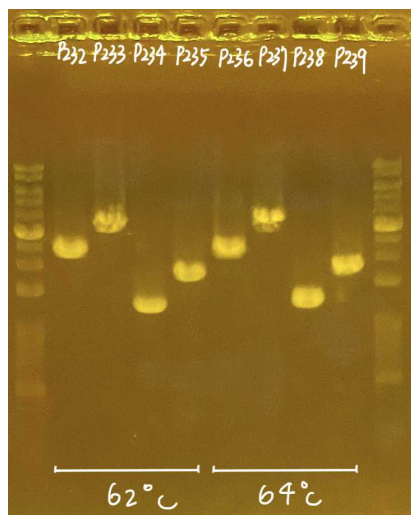

Marked gel picture go to the Lab Google drive

## For 2' PCR barcoded-head primers

- 10 Prepare 2' PCR master mixutre for **barcoded-primers** (prepare 1.2X of solutions for **pipetting error if needed**)


PCR mixture for barcoded-primers for each reaction (**NO PRIMERS at this point!!**)

A	B	C	D
Component	Volume	Volume (1.2X)	Final conc.
ZEJU PCR Master Mix	7.5 $\mu$ L	9 $\mu$ L	-
ddH2O	5.55 $\mu$ L	6.66 $\mu$ L	-
Total volume	13.05 $\mu$ L	15.66 $\mu$ L	-

### Note

Negative control ALWAYS NEEDED! For example, if you have 5 PCR reactions to run, prepare master mixture for 6 reactions (5 DNA template + 1 negative control).

- 11 Mix the 2' PCR master mixture gently by pippeting. Quick spin the tube.

12 Transfer  13.05 µL of the 2' PCR master mixture to 8-strip PCR tubes.

13 Add  1.2 µL **pre-mixed barcoded-head primers** (Forward + Reverse) to each PCR tubes.

14 Add  0.75 µL of **1' PCR product as template**, resulting in  15 µL reaction mixture for 2' PCR.



**Negative control** contains only  14.25 µL master mixture and premixed barcoded-head primers but not DNA template

15 Mix gently by tapping the tubes. Quick spin the tubes.

16 Carry out 2' PCR using the following condition:


2' PCR condition for **barcoded-head primers**

A	B	C	D
Step	Temp	Sec	Cycle
Initial denaturation	98 °C	30	
Denaturation	98 °C	15	10-15 cycles
Annealing	64-68 °C varied (a)	15	
Extension	72 °C	20 (b)	
Final extension	72 °C	210	
Preservation	Preservation	4 °C	∞

a. Annealing varied, **65 °C** is working based on test on 220531; Refer 2' PCR primers for annealing temperature


b. 1kb ~ 1min extension; enough time allow full extension of sequence

- 17 Carry out **electrophoresis** for inspection of DNA products

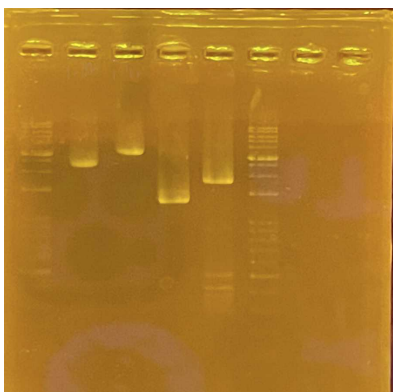

Gel before markdown

- 18 Markdown wells and upload the pictures to the Lab Google drive

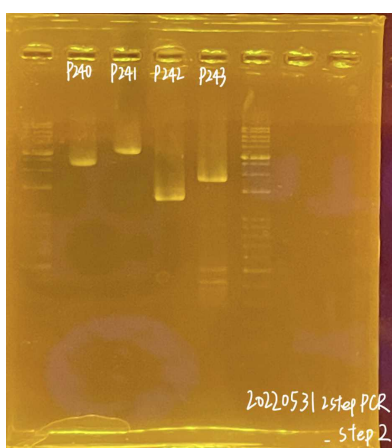

Marked gel picture go to the Lab Google drive