

Nov 21, 2020

## © CELL COUNT- 03 - Automated cell count with Trypan Blue Solution by Cellometer Auto T4 Cell Counter

Marco Cosentino<sup>1</sup>, Elisa Storelli<sup>1</sup>, Alessandra Luini<sup>1</sup>, Massimiliano LM Legnaro<sup>1</sup>, Emanuela Rasini<sup>1</sup>, Marco Ferrari<sup>1</sup>, Franca Marino<sup>1</sup>

<sup>1</sup>Center for Research in Medical Pharmacology, University of Insubria (Varese)

1 Works for me

dx.doi.org/10.17504/protocols.io.bpxtmpnn



Farmacologia Medica

ABSTRACT

Published work using this protocol:

- Kustrimovic N., Comi C., Magistrelli L., Rasini E., Legnaro M., Bombelli R., Aleksic I., Blandini F., Minafra B., Riboldazzi G., Struchio A., Mauri M., Bono G., Marino F., Cosentino M. Parkinson's disease patients have a complex phenotypic and functional Th1 bias: cross-sectional studies of CD4+ Th1/Th2/T17 and Treg in drug-naïve and drug-treated patients (2018). Journal of neuroinflammation, 15(1), 205. <a href="https://doi.org/10.1186/s12974-018-1248-8">https://doi.org/10.1186/s12974-018-1248-8</a>
- Kustrimovic, N., Rasini, E., Legnaro, M., Bombelli, R., Aleksic, I., Blandini, F., Comi, C., Mauri, M., Minafra, B., Riboldazzi, G., Sanchez-Guajardo, V., Marino, F., & Cosentino, M. (2016). Dopaminergic Receptors on CD4+ T Naive and Memory Lymphocytes Correlate with Motor Impairment in Patients with Parkinson's Disease. Scientific reports, 6, 33738. <a href="https://doi.org/10.1038/srep33738">https://doi.org/10.1038/srep33738</a>
- Cosentino M., Ferrari M., Kustrimovic N., Rasini E., Marino F. (2015). Influence of dopamine receptor gene polymorphisms on circulating T lymphocytes: A pilot study in healthy subjects. Human immunology, 76, 10, 747-752. https://doi.org/10.1016/j.humimm.2015.09.032
- Boydum A.Isolation of mononuclear cells and granulocytes from human blood. Scand.J.Clin.Lab. Invest. 21 (Suppl.97): 77-89, 1968
- Alex Mabou Tagne, Franca Marino, Massimiliano Legnaro, Alessandra Luini, Barbara Pacchetti and Marco Cosentino. A Novel Standardized Cannabis sativa L. Extract and Its Constituent Cannabidiol Inhibit Human Polymorphonuclear Leukocyte Functions. Int J Mol Sci2019 Apr; 20(8): 1833. Published online 2019 Apr 13. doi: 10.3390/ijms20081833.
- Angela Scanzano, Laura Schembri, Emanuela Rasini, Alessandra Luini, Jessica Dallatorre, Massimiliano Legnaro, Raffaella Bombelli, Terenzio Congiu, Marco Cosentino, Franca Marino. Adrenergic Modulation of Migration, CD11b and CD18 Expression, ROS and interleukin-8 Production by Human Polymorphonuclear Leukocytes. Inflamm Res. 2015 Feb;64(2):127-35. doi:10.1007/s00011-014-0791-8. Epub 2015 Jan 6.

DOI

dx.doi.org/10.17504/protocols.io.bpxtmpnn

🐼 protocols.io

11/21/2020

Marco Cosentino, Elisa Storelli, Alessandra Luini, Massimiliano LM Legnaro, Emanuela Rasini, Marco Ferrari, Franca Marino 2020. CELL COUNT- 03 - Automated cell count with Trypan Blue Solution by Cellometer Auto T4 Cell Counter. **protocols.io** 

https://dx.doi.org/10.17504/protocols.io.bpxtmpnn

LICENSE

This is an open access protocol distributed under the terms of the Creative Commons Attribution License, which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**CREATED** 

Nov 21, 2020

LAST MODIFIED

Nov 21, 2020

PROTOCOL INTEGER ID

44755

MATERIALS TEXT

Instruments needed:

- Plastic labware

Materials:

Trypan blue catalog number: T8154-100mL Sigma, Italy

SAFETY WARNINGS

Cancerogenic solution, BE CAREFUL WHEN HANDLING IT!

**ABSTRACT** 

Published work using this protocol:

- Kustrimovic N., Comi C., Magistrelli L., Rasini E., Legnaro M., Bombelli R., Aleksic I., Blandini F., Minafra B., Riboldazzi G., Struchio A., Mauri M., Bono G., Marino F., Cosentino M. Parkinson's disease patients have a complex phenotypic and functional Th1 bias: cross-sectional studies of CD4+ Th1/Th2/T17 and Treg in drug-naïve and drug-treated patients (2018). Journal of neuroinflammation, 15(1), 205. <a href="https://doi.org/10.1186/s12974-018-1248-8">https://doi.org/10.1186/s12974-018-1248-8</a>
- Kustrimovic, N., Rasini, E., Legnaro, M., Bombelli, R., Aleksic, I., Blandini, F., Comi, C., Mauri, M., Minafra, B., Riboldazzi, G., Sanchez-Guajardo, V., Marino, F., & Cosentino, M. (2016). Dopaminergic Receptors on CD4+ T Naive and Memory Lymphocytes Correlate with Motor Impairment in Patients with Parkinson's Disease. Scientific reports, 6, 33738. <a href="https://doi.org/10.1038/srep33738">https://doi.org/10.1038/srep33738</a>
- Cosentino M., Ferrari M., Kustrimovic N., Rasini E., Marino F. (2015). Influence of dopamine receptor gene polymorphisms on circulating T lymphocytes: A pilot study in healthy subjects. Human immunology, 76, 10, 747-752. https://doi.org/10.1016/j.humimm.2015.09.032
- Boydum A.Isolation of mononuclear cells and granulocytes from human blood. Scand.J.Clin.Lab. Invest. 21 (Suppl.97): 77-89, 1968
- Alex Mabou Tagne, Franca Marino, Massimiliano Legnaro, Alessandra Luini, Barbara Pacchetti and Marco

protocols.io
2
11/21/2020


Citation: Marco Cosentino, Elisa Storelli, Alessandra Luini, Massimiliano LM Legnaro, Emanuela Rasini, Marco Ferrari, Franca Marino (11/21/2020). CELL COUNT- 03 - Automated cell count with Trypan Blue Solution by Cellometer Auto T4 Cell Counter. https://dx.doi.org/10.17504/protocols.io.bpxtmpnn

Cosentino. A Novel Standardized Cannabis sativa L. Extract and Its Constituent Cannabidiol Inhibit Human Polymorphonuclear Leukocyte Functions. Int J Mol Sci2019 Apr; 20(8): 1833. Published online 2019 Apr 13. doi: 10.3390/ijms20081833.

- Angela Scanzano, Laura Schembri, Emanuela Rasini, Alessandra Luini, Jessica Dallatorre, Massimiliano Legnaro, Raffaella Bombelli, Terenzio Congiu, Marco Cosentino, Franca Marino. Adrenergic Modulation of Migration, CD11b and CD18 Expression, ROS and interleukin-8 Production by Human Polymorphonuclear Leukocytes. Inflamm Res. 2015 Feb;64(2):127-35. doi: 10.1007/s00011-014-0791-8. Epub 2015 Jan 6.
- 1 For automatic cell count with Cellometer machine use Trypan Blue.

The machine will calculate the n°of cells/ml and the % of viability.

Take  $\blacksquare 10 \ \mu l$  of cell suspention and add an equal amount of **Trypan Blue**. Use all the volume to place it in the given counting chamber. Place the chamber inside Cellometer and count.

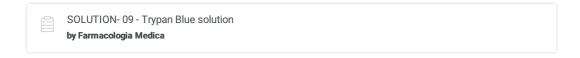

Cellometer Auto T4 Automated cell counter

Nexcelom Bioscience EuroClone