



AUG 01, 2023

OPEN ACCESS



DOI:
dx.doi.org/10.17504/protocols.io.6qpvr36yovmk/v1

Protocol Citation: Tae-Un Han, Ellen Sidransky 2023. Immunohistochemistry (using floating section). **protocols.io** <https://dx.doi.org/10.17504/protocols.io.6qpvr36yovmk/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

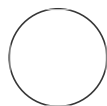
Protocol status: Working
 We use this protocol and it's working

Created: Jul 28, 2023

🌐 Immunohistochemistry (using floating section)

Tae-Un Han¹, Ellen Sidransky¹

¹National Institute of Health



Tae-Un Han

ABSTRACT

- This is an immunohistochemistry protocol used for evaluation of pathology in mouse brain.

- We successfully detected Gfap, Iba1, CD68, TH and alpha-Syn and phospho S129 alpha synuclein protein expression using this protocol.

- In this protocol, floating sections are mounted on slides first before starting blocking step. In this way, large number of sections can be stained more efficiently.

MATERIALS

30 um thickness of floating tissue sections in PBS with 0.1% sodium azide.

- Washing: Glass Coplin jar

- PBST: PBS 1X + 0.05% Tween20

- Normal Goat Serum Blocking Solution:
 2% goat serum (blocking)
 1%BSA (stabilizer)
 0.1% cold fish skin gelatin (blocking)
 0.1% Triton X-100 (penetration enhancer)
 0.05% Tween 20 (detergent and surface tension reducer)
 0.05% sodium azide (preservative)
 0.01M PBS, pH 7.2
 Mix well and store at 4 °C for up to 6 months.

Last Modified: Aug 01, 2023

PROTOCOL integer ID: 85658

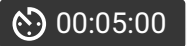

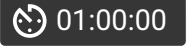



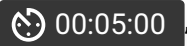

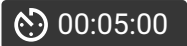
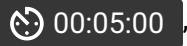
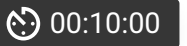
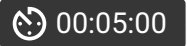
Keywords: ASAPCRN

- Primary Antibody Dilution Buffer
1%BSA (blocking & stabilizer)
0.1% cold fish skin gelatin (blocking)
0.5% Triton X-100 (penetration enhancer)
0.05% sodium azide (preservative)
0.01M PBS, pH 7.2-7.4
Mix well and store at 4 °C. up to 6 months.

- Secondary Antibody Dilution Buffer:
0.01M PBS, pH 7.2
0.05% Tween 20

- List of antibodies used for neuropathology evaluation with IHC

A	B	C	D
Antibodies	Cat#	Company	dilution
Anti-Synuclein α Antibody	AB5038	Sigma- Aldrich	1:1000
Recombinant Anti-Alpha-synuclein (phospho S129) antibody [EP1536Y]	ab51253	Abcam	1:1000
GFAP Antibody , Rabbit Polyclonal	16825-1- AP	Proteintech	1:4000
Recombinant Anti-Iba1 antibody [EPR16588]	ab178846	Abcam	1:1000
Anti-Tyrosine Hydroxylase	AB152	Sigma- Aldrich	1:1000
Anti-CD68 antibody	ab125212	Abcam	1:1000
Goat Anti-Rabbit IgG H&L (Alexa Fluor® 488) : 2nd antibody	ab150077	Abcam	1:1000

- 1 Mount floating sections onto slides and dry for  00:30:00 (<45 min), While section drying, PAP pen drawing line and dry  00:05:00 . 35m
- 2 PBS rehydrate  00:10:00 , Wash with PBS for  00:05:00 three times. 15m
- 3 Blocking sol  500 μ L , incubate  01:00:00 at  Room temperature . 1h
- 4 Replace buffer to primary antibody sol, incubate  Overnight at  4 °C . 10m
- 5 Wash in PBST three times for  00:05:00 ,  00:10:00 and  00:05:00 . 20m
- 6 Secondary antibody sol, incubate for  01:30:00 at  Room temperature . 1h 30m
- 7 Wash in PBS three times for  00:05:00 ,  00:10:00 and  00:05:00 . 20m
- 8 Mounting with Vectorshield mounting media

