

Thawing adherent cancer cell lines V.1

Emily Souster¹, Verity Goodwin¹, Charlotte Beaver¹, Adam Hunter¹, Fiona Behan¹, Rizwan Ansari¹, Mathew Garnett¹

¹Wellcome Sanger Institute

1 Works for me This protocol is published without a DOI.


Cellular Generation and Phenotyping

Emily Souster

ABSTRACT

This SOP is for thawing an adherent cancer cell line from a frozen cryovial.

Process diagram:

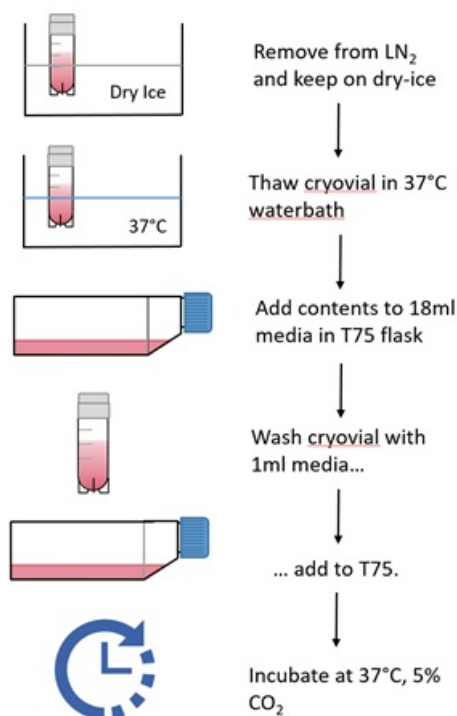

PROTOCOL CITATION

Emily Souster, Verity Goodwin, Charlotte Beaver, Adam Hunter, Fiona Behan, Rizwan Ansari, Mathew Garnett 2020. Thawing adherent cancer cell lines. **protocols.io**
<https://protocols.io/view/thawing-adherent-cancer-cell-lines-bgs8jwhw>

COLLECTIONS

Whole genome CRISPR screening in cancer cell lines

LICENSE

This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

May 24, 2020

LAST MODIFIED

Jul 06, 2020

PROTOCOL INTEGER ID

37440

PARENT PROTOCOLS

Part of collection

[Whole genome CRISPR screening in cancer cell lines](#)

MATERIALS

NAME	CATALOG #	VENDOR
Cell culture treated T75 flasks	430641	Corning

MATERIALS TEXT

Select an appropriate culture media for your cell line. Common culture medias used for cancer cell lines are serum supplemented Advanced DMEM F-12 or RPMI, in the presence of pen-strep.

Equipment

Pipette boy

25ml stripette

P1000 pipette and tips

Microbiology Safety Cabinet (MSC)

Light Microscope

 **37 °C** waterbath

 **37 °C** , 5% CO₂ incubator

SAFETY WARNINGS

Retrieving vials from liquid nitrogen tanks is a significant ergonomic risk due to potential weight of filled racks. To reduce the risk of manual handling injury, lifting racks out of dewers should be kept to a minimum.

BEFORE STARTING

- Pre-warm complete culture media to room-temperature.
- Add 18ml complete culture media to a T75 flask.


1 Remove the cryovial from liquid nitrogen storage and place on dry-ice for transfer to cell culture lab.




2 Take the cryovial from the dry-ice and hold in a  **37 °C** waterbath until thawed.



Be careful here to **avoid submerging** the lid or rim of the cryovial in the waterbath.

Dry the cryovial thoroughly, spray with 70% ethanol and transfer to cell culture hood.

3 Use a P1000 pipette to transfer the contents of the cryovial to the T75 flask containing  **18 mL** media.

- 4 Use  1 mL complete culture media to wash the cryovial and add this to the flask.
- 5 Place the flask in a  37 °C , 5% CO₂ incubator. Agitate the flask gently back and forth and side to side to ensure even distribution of cells across the flask.
- 6 After 24 hours, inspect adherence and confluency. Remove the media using a sterile aspirating pipette and replace with  12 mL complete culture media.