

JAN 30, 2023

## OPEN ACCESS

**Protocol Citation:** Carlos Goller, Carly Sjogren 2023. Nanopore Rapid PCR Barcoding for Genomic Samples. [protocols.io](https://protocols.io/view/nanopore-rapid-pcr-barcoding-for-genomic-samples-chxzt7p6) <https://protocols.io/view/nanopore-rapid-pcr-barcoding-for-genomic-samples-chxzt7p6>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** In development  
We are still developing and optimizing this protocol

**Created:** Oct 15, 2022

**Last Modified:** Jan 30, 2023

**PROTOCOL integer ID:**  
71385

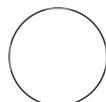
# Nanopore Rapid PCR Barcoding for Genomic Samples

Carlos Goller<sup>1</sup>, Carly Sjogren<sup>1</sup>

<sup>1</sup>North Carolina State University

Delftia and SCoOP

Tech. support phone: +91 95134-135 email: [ccgoller@ncsu.edu](mailto:ccgoller@ncsu.edu)



anmcleo

## ABSTRACT

This protocol is for Rapid sequencing of DNA using PCR Barcoding (SQK-RPB004) and has been adapted from Oxford Nanopore Technologies.

## MATERIALS

### Equipment:

- Microfuge
- Timer
- Thermal Cycler
- Pipette and tips P2, P10, P20, P100, P200, P1000.

### Consumables:

- 1.5 ml Eppendorf DNA LoBind tubes
- 0.2 ml thin-walled PCR tubes
- Nuclease-free water
- Agencourt AMPure XP Beads
- LongAmp Taq 2X Master Mix
- Fresh 70% ethanol in nuclease-free water
- 10 mM Tris-HCl pH 8.0 with 50 mM NaCl

### Materials:

- 1-5 ng high molecular weight genomic DNA
- Rapid PCR Barcoding Kit
- Flow Cell Priming Kit

## BEFORE START INSTRUCTIONS

Overall time is about 30 to 40 minutes.

### COLLECT THESE:

#### Equipment:

- Microfuge
- Timer
- Thermal Cycler
- Pipette and tips P2, P10, P20, P100, P200, P1000.

#### Consumables:

- 1.5 ml Eppendorf DNA LoBind tubes
- 0.2 ml thin-walled PCR tubes
- Nuclease-free water
- Agencourt AMPure XP Beads
- LongAmp Taq 2X Master Mix
- Fresh 70% ethanol in nuclease-free water
- 10 mM Tris-HCl pH 8.0 with 50 mM NaCl

#### Materials:

- 1-5 ng high molecular weight genomic DNA
- Rapid PCR Barcoding Kit
- Flow Cell Priming Kit

## Checklist

### 1 Gather Materials, Consumables and Equipment.

*These are listed in the Materials tab also*

#### **Equipment:**

- Microfuge
- Timer
- Thermal Cycler
- Pipette and tips: P2, P10, P20, P100, P200, P1000

#### **Consumables:**

- 1.5 ml Eppendorf DNA LoBind Tubes
- 0.2 ml thin-walled PCR tubes
- Nuclease-free water
- Agencourt AMPure XP Beads
- LongAmp Taq 2x Master Mix
- Freshly prepped 70% ethanol in Nuclease-free water
- 10 mM Tris-HCl pH 8.0 with 50 mM NaCl

### Materials:

- 1-5 ng high molecular weight Genomic DNA
- Rapid PCR Barcoding Kit
- Flow Cell Priming Kit


## Library Preparation

22m 30s



### 2 Thaw and prepare the reagents as follows:

- Barcodes (RLB 01-12A) at Room Temp
- Fragmentation Mix (FRM) on ice
- Rapid Adaptor (RAP) on ice

#### 2.1 Prepare the DNA in Nuclease-free water:


- Transfer 1-5 ng genomic DNA into a DNA LoBind tube
- Adjust volume to  3 µL with nuclease-free water
- Mix thoroughly by flicking (Avoid unwanted shearing)
- Spin down briefly in Microfuge

#### 2.2 In a thin 0.2 ml thin-walled PCR tube, Mix the following:

-  3 µL 3 1-5 ng template DNA
-  1 µL 1 Fragmentation Mix (FRM)

Mix Gently by flicking tube and spin down.





#### 2.3 Incubate tube in a thermal cycler at 30C for 1 minute. 00:01:00 30C

Then for 1 minute at 80C  00:01:00 80C

2m

Briefly put tube on ice to cool down




#### 2.4 Set up a PCR reaction as follows in 0.2 ml thin-walled PCR tube:



-  20 µL 2 Nuclease-free water
-  4 µL Tagmented DNA
-  1 µL RLB (01,12A, at 10µM)
-  25 µL LongAmp Taq 2X Master Mix

Mix gently by flicking tube and spin down

### 3 Amplifying using following cycling conditions:

15m 30s

- Initial Denaturation  00:03:00 95C (1 Cycle)
- Denaturation  00:00:15 95C (14 Cycles)
- Annealing  00:00:15 56C (14 Cycles)


- Extension  00:06:00 65C (14 Cycles)
- Final extension  00:06:00 65C (1 Cycle)
- HOLD AT 4C

**3.1** Transfer sample to clean 1.5ml Eppendorf DNA LoBind tube and Resuspend AMPure XP Beads by vortexing

-Add  30  $\mu\text{L}$  of resuspended AMPure XP Beads to reaction, mix by Pipetting

**3.2** Incubate on Hula Mixer  00:05:00 Room Temp

5m

Prepare  500  $\mu\text{L}$  of fresh 70% ethanol in Nuclease-Free water during incubation

**3.3** Take tube from Hula Mixer and put on a magnet. Once solution looks clear pipette off supernatant

While on magnet wash beads with  200  $\mu\text{L}$  of 70% ethanol.

#### Note

Take care to not disturb pellet when washing beads with the 70% ethanol


Remove ethanol using a pipette and discard.

**3.4** Repeat above washing and removal step.


**3.5** Spin down and place tube on magnetic rack. Pipette off any residual ethanol in the tube and keep the top open to allow the pellet to dry for about 30 seconds.


#### Note

DO NOT DRY PELLETT TO POINT OF CRACKING. DO NOT EXCEED 30 SECONDS.


**3.6** Take tube out of magnetic rack and resuspend pellet in  10  $\mu\text{L}$  of 10 mM Tris-HCl pH 8.0 with 50 mM NaCl.


2m


Incubate  00:02:00 at Room Temp

- 3.7** Put tube back on magnetic rack and wait until elute is clear and colorless. Once clear, remove and retain  10  $\mu\text{L}$  of eluate into a clean 1.5 ml Eppendorf DNA LoBind tube. Dispose the pelleted beads and previous tube.

**Note**

You may want to spend down the tube after getting  10  $\mu\text{L}$  of eluate and put back on magnetic rack to get the possible remaining eluate out to use in the next step of quantifying on Qubit Fluorometer. You can also use remaining eluate if present after retaining the initial 10  $\mu\text{L}$ .

- 4** Quantify  1  $\mu\text{L}$  of eluted sample using a Qubit Fluorometer.

- 5** Pool all barcoded libraries in desired ratios to a total of 50-100 fmoles in  10  $\mu\text{L}$  of 10 mM Tris-HCl pH 8.0 with 50 mM NaCl.

**Note**

PLEASE BE SURE TO POOL ALL BARCODED LIBRARIES TOGETHER BEFORE PROCEEDING TO NEXT STEP.

- 5.1** Add  1  $\mu\text{L}$  of RAP to barcoded DNA. Mix by gently flicking tube and spin down.

- 5.2** Incubate reaction for  00:05:00 at Room Temp

5m

**Note**

After incubating store this library on ice until ready to load into MinION flow cell.

6

Thaw the following at Room Temperature:

- Sequencing Buffer (SQB)
- Loading Beads (LB)
- Flush Tether (FLT)
- One tube of Flush Tether (FB)


Once thawed, mix reagents by vortexing then spin down at room temperature.

## Note

Sequencing Tether (SQT) tube will NOT be used in protocol.

6.1


Prepare the flow cell priming mix:

- add  30 µL of thawed and mixed Flush Tether (FLT) directly to tube of thawed and mixed Flush Buffer (FB).
- Mix by vortexing at room temperature.

6.2




Open the MinION device lid and slide the flow cell under the clip. Slide the priming port cover clockwise to open the priming port.

## Note

Be careful when drawing back buffer from flow cell. DO NOT remove more than  20-30 µL, and make sure that the array of pores are covered by buffer at all times. DO NOT introduce air bubbles into array; this can irreversibly damage pores.

6.3



After opening priming port, check for small air bubbles under the cover. Draw back small volume to remove any bubbles (a few µl):

- Set a P1000 to  200 µL
- Insert tip into priming port
- Turn wheel on pipette until dial shows  220-230 µL (This allows you to draw back  20-30 µL), or until you can see a SMALL amount of buffer entering pipette tip.

## Note

Visually check that there is continuous buffer from the priming port across the sensor array

5m

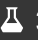
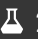
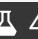

- Load  800 µL of the priming mix into the flow cell via the priming port. Avoid introduction of air bubbles. Wait  00:05:00 . During this 5 minutes, proceed to the next sub-step

#### 6.4 Thoroughly mix the contents of the Loading Beads (LB) tubes by vortexing.

##### Note


These loading beads (SB) will need to be mixed in immediately before use.

In a new tube, prepare the library for loading as follows:


-  34 µL of Sequencing Buffer (SQB)
-  25.5 µL Loading Beads (LB), mixed IMMEDIATELY before use.
-  4.5 µL Nuclease-free water
-  11 µL DNA library

#### 6.5 Complete the flow cell priming:

-Gently lift the SpotON sample port cover to make the SpotOn sample port accessible.

-Load  200 µL of priming mix into flow cell via the priming port (NOT the SpotON sample port), avoid introduction of air bubbles

#### 6.6 Mix prepared library gently by pipetting up and down prior to loading

Add  75 µL of sample to flow cell via SpotON sample port in dropwise fashion. Ensure each drop flows into port before adding the next drop

Gently replace the SpotON sample port cover, making sure the bung enters the SpotON port, close the priming port and replace the MinION device lid.