

DEC 22, 2023

Human iPSCs culture and cardiomyocyte subtype differentiation in fully chemically defined conditions

Natalia

Joanna Delimata-Raczek¹, Koralewska¹,


Marek

Figlerowicz¹,

Ireneusz Stolarek¹

¹Institute of Bioorganic Chemistry Polish Academy of Sciences

MOSAIC

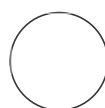

Joanna Delimata-Raczek

ABSTRACT

Cardiomyocyte in vitro differentiation represents a pivotal avenue in stem cell research, offering unprecedented insights into cardiac development and disease understanding. Over the years, protocols for inducing cardiomyogenesis have evolved substantially, particularly with the advent of induced pluripotent stem (iPS) cell technology. The transition from intricate and uncontrolled systems to more refined and robust protocols has been paramount, driven by a cumulative understanding of the molecular and cellular cues governing the process. This protocol addresses the intricacies of cultivating iPS cells and orchestrating their differentiation into distinct cardiomyocyte subtypes. The protocol encompasses not only the fundamental aspects of human induced pluripotent stem cell (hiPSC) culture but also delves into the nuanced stages of cardiomyocyte differentiation, sub-type specification, and the critical maturation phase. The refined methodologies outlined in this protocol contribute to advancing capabilities in disease modeling, drug discovery, developmental biology, and regenerative medicine.

PROTOCOL MATERIALS



Essential 8[®]; Flex Medium Kit Thermo
Fisher Catalog #A2858501

In [2 steps](#)



Penicillin/Streptomycin Thermo Fisher Scientific Catalog #Invitrogen 15140-122

Step 2.1



RPMI 1640 Medium, GlutaMAX[®]; Supplement, HEPES Thermo
Fisher Catalog #72400021

In [3 steps](#)



Pluripotent Stem Cell 4-Marker Immunocytochemistry Kit Thermo
Fisher Catalog #A24881

OPEN ACCESS



Protocol Citation: Joanna Delimata-Raczek, Natalia Koralewska, Marek Figlerowicz, Ireneusz Stolarek 2023. Human iPSCs culture and cardiomyocyte subtype differentiation in fully chemically defined conditions . **protocols.io** <https://protocols.io/view/human-ipscs-culture-and-cardiomyocyte-subtype-diff-c5my7u6>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Dec 06, 2023

Last Modified: Dec 22, 2023

PROTOCOL integer ID:
91885

Keywords: cardiomyocyte,
differentiation, in vitro, cell
culture, iPSC


Step 4.1

 Penicillin-Streptomycin (5,000 U/mL) Thermo
Fisher Catalog #15070063


In 2 steps

 StemMACS™ CHIR99021 Miltenyi
Biotec Catalog #130-103-926


Step 11.2

 Geltrex LDEV Free hESC Quality 5 ml Thermo Fisher
Scientific Catalog #A1413302

Step 1.1

 ROCK inhibitor (Ri) Y-27632 STEMCELL Technologies
Inc. Catalog #72308


In 3 steps

 Lookout® One-Step Mycoplasma PCR Detection Kit Merck MilliporeSigma
(Sigma-Aldrich) Catalog #MP0050-100TST


Step 5.1

 Retinoic acid Merck MilliporeSigma (Sigma-
Aldrich) Catalog #R2625-50MG


Step 11.2

 Fetal Bovine Serum, qualified, One Shot® format, United States Thermo
Fisher Catalog #A3160501


Step 18.1

 PSC Cardiomyocyte Differentiation Kit Thermo
Fisher Catalog #A2921201


Step 9

 dimethylsulfoxide (DMSO) Merck MilliporeSigma (Sigma-
Aldrich)


In 2 steps

 B-27™ Supplement, minus insulin Gibco - Thermo
Fisher Catalog #A1895601


Step 11.1

 B-27™ Supplement (50X), serum free Gibco - Thermo
Fisher Catalog #17504044

Step 11.1

 IWP-2 1 mg STEMCELL Technologies
Inc. Catalog #72122


Step 11.2

 Human Cardiomyocyte Immunocytochemistry Kit Thermo
Fisher Catalog #A25973

Step 16.2


 DMEM/F-12 Thermo Fisher Catalog #11320033

Step 1.2

 Trypan Blue Solution 0.4% Sterile-filtered Merck MilliporeSigma (Sigma-
Aldrich) Catalog #T8154

Step 6.6

SAFETY WARNINGS

 Appropriate safety procedures are recommended to be used when handling all cell lines, working with iPSCs and induced cardiomyocytes handle them under at least BSL-1 containment.

BEFORE START INSTRUCTIONS


Keep in mind that while working with cells everything should be kept sterile, every medium or any reagent should be prepared under the sterile hood, autoclaved, or filtered through 0,2um filter.

iPSC culture and maintenance

1 Coating culture plates with ECM proteins for enhanced attachment of iPSCs.

1.1

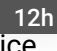

Thaw  Geltrex LDEV Free hESC Quality 5 ml Thermo Fisher Scientific Catalog #A1413302  on ice

 4 °C

 Overnight

It is prohibited to thaw Geltrex at room temperature or in a water bath.

Remember to keep it on ice during the whole aliquoting procedure and use only ice-cold pipette tips and tubes.

Note

Making small aliquot portions of Geltrex is recommended to limit the freeze-thaw cycles.




*Storage of aliquots in  -20 °C freezer for up to 1 year.

1.2

Resuspend single Geltrex aliquot in cold  DMEM/F-12 Thermo Fisher Catalog #11320033 medium to obtain desired dilution. Store diluted coating medium in the fridge.

Note

1:60-1:100 dilution usually works fine but it should be first tested on a particular cell line.

- 1.3 Coat 6-well culture plate by adding  1 mL of coating medium to each well and cover the entire  30m surface of each well. Incubate the plate in the incubator for at least  00:30:00 before proceeding further.


Note




Coated plates can be used immediately after removal of the coating medium or stored in the fridge for up to 2 weeks with added PBS and wrapped in parafilm to prevent from drying.


2 iPSC thawing




2.1


Prepare culture medium using the following kit:

 Essential 8[®]; Flex Medium Kit Thermo : add  48.5 mL of
Fisher Catalog #A2858501

Essential 8 flex medium to the  50 mL centrifuge tube and  1 mL of 10x supplement. It is not necessary, but it is advised to also add  500 µL of





 Penicillin/Streptomycin Thermo Fisher Scientific Catalog #Invitrogen 15140-122 to prevent cell culture from contamination.








- 2.2 Prepare the previously coated 6-well plate by adding  2 mL of warm culture medium to each well and  2 µL from  10 millimolar (mM) stock solution of

 ROCK inhibitor (Ri) Y-27632 STEMCELL Technologies to minimize the cell
Inc. Catalog #72308 death.

Note

Preparing 10mM stock solution of rock inhibitor Y-27632

1. Dissolve  5 mg of rock inhibitor in  1.56 mL DMSO.
2. Mix well until it dissolves completely.
3. Prepare small aliquot portions, for example  10 µL in PCR tubes and keep them in the  -20 °C freezer .


- 2.3 Take out the desired cryovial containing frozen cells from the liquid nitrogen and quickly transfer it to  37 °C water bath. Be careful not to submerge the cap of the cryovial. You can swirl the tube a little to help it thaw quicker but pay attention to save a little ice cube inside the cryovial before taking it out. Disinfect the cryovial with 70% ethanol and put it under the hood.
- 2.4 Transfer  1 mL of thawed cell suspension to the  15 mL centrifuge tube and add  9 mL of warm culture medium very gently drop by drop, swirling the tube to minimize osmotic shock to the cells.
- 2.5 Centrifuge the cell suspension  1400 rpm, Room temperature, 00:05:00 . Carefully remove and  5m discard the supernatant. Gently add  1 mL of medium to the remaining cell pellet. Gently resuspend cells in medium and seed into 2 wells of a previously prepared 6-well plate slowly drop by drop.
- 2.6 Gently rock plate side to side, and back and forth to spread the cells evenly across the well. Transfer the plate to the incubator. The next day in the morning change the medium to fresh culture medium without rock inhibitor.

3 iPSC culture

- 3.1 Incubate the cells at 37 °C/ 5% CO₂ incubator and observe their growth and confluence under the microscope daily. Perfect undifferentiated iPSC colonies should look compact and present distinct borders. Cells present large nuclei and smaller cytoplasm region. It is essential to check iPSCs for the presence of spontaneous differentiation. If any signs of differentiation occur, colonies should be discarded immediately.


Expected result

Morphological differences between correct and incorrect iPSC colonies

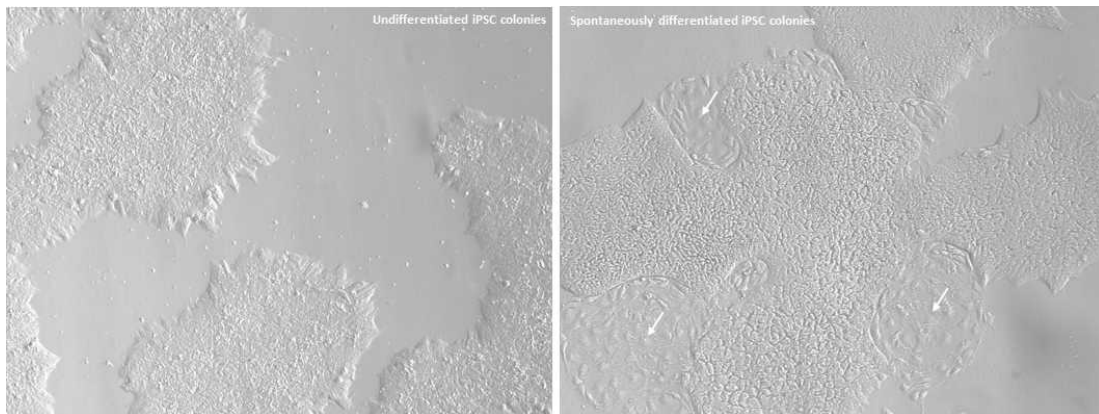

white arrows present differentiated regions

3.2 Medium should be changed every day. When using

 Essential 8[™]; Flex Medium Kit Thermo
Fisher Catalog #A2858501

it is possible to skip

weekend medium change by adding 2x more medium on Friday afternoon, but it is only recommended when cells are at low confluency.

4 Pluripotency characterization


4.1 To guarantee the reliability of the research iPSC lines should be verified for pluripotency markers either by qPCR or immunocytochemistry. The most common pluripotency markers are OCT4, SOX2, SSEA4, TRA-1-60, and alkaline phosphatase.

Expected result

iPSC colonies present positive staining for pluripotency markers OCT4 and SSEA4 using



Pluripotent Stem Cell 4-Marker Immunocytochemistry Kit Thermo
Fisher Catalog #A24881

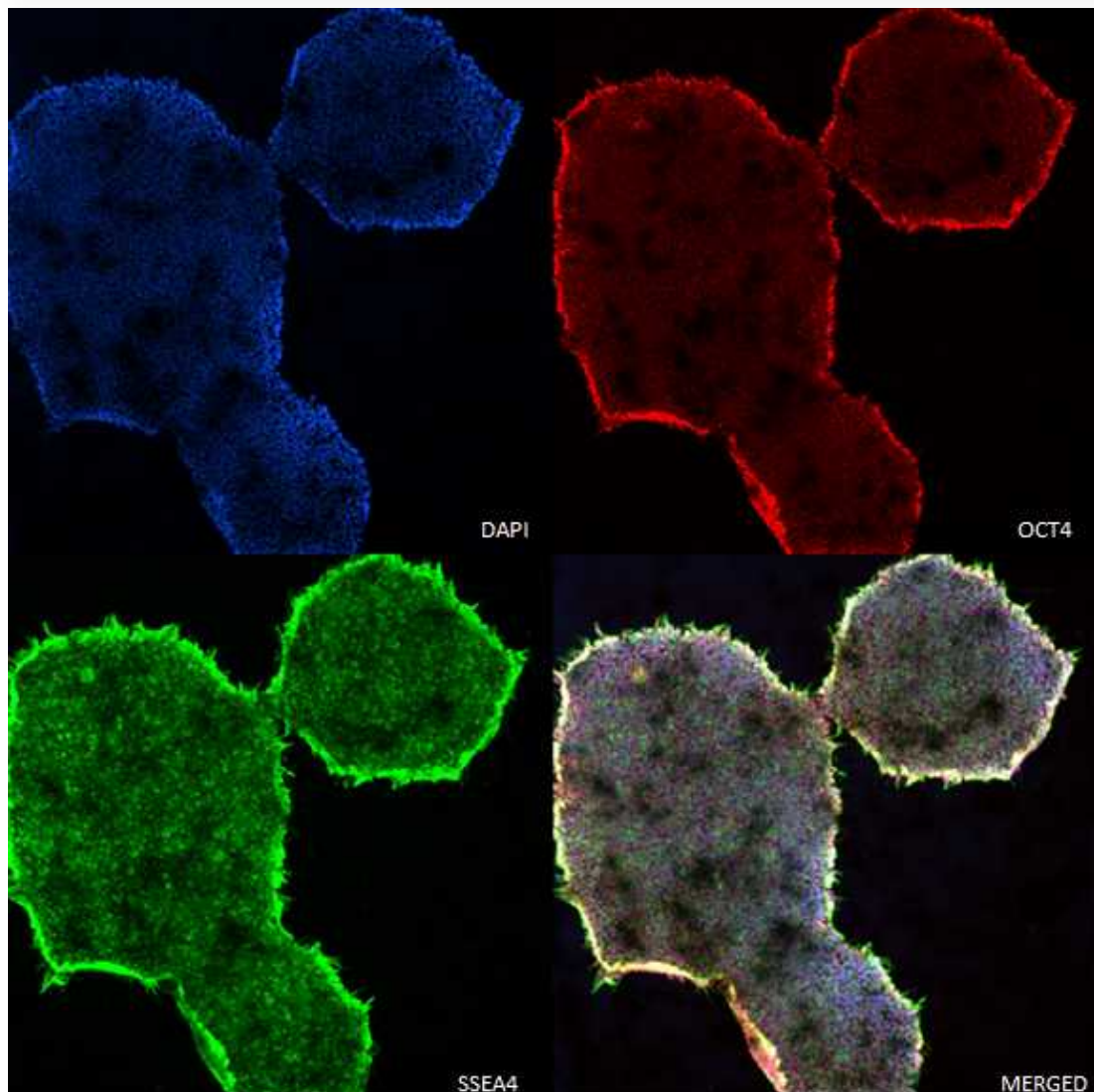

5 Mycoplasma testing


5.1 It is very important to maintain contamination-free cultures. A commonly known good laboratory

practice is to regularly check cell cultures for mycoplasma. Mycoplasma contamination is very hard to notice without special tests. This contamination usually does not cause cell death, but it negatively affects the cell metabolism and growth making the research not reliable. The best way to assess mycoplasma contamination in culture is to use available qPCR or immunocytochemistry ready-to-use kits, for example:





Lookout® One-Step Mycoplasma PCR Detection Kit Merck MilliporeSigma (Sigma-Aldrich) Catalog #MP0050-100TST


6 iPSC passage

- 6.1 The optimal time for harvest or passage for iPSCs is typically when cells reach approximately 70-80% confluency.
- 6.2 To start with prepare mild cell detachment solution, for example PBS with **[M] 0.5 millimolar (mM)** EDTA. Warm it to **🔥 Room temperature** before adding it to the cells.
- 6.3 Prepare a new Geltrex-coated plate as in steps **⇒ go to step #1.3 and #2.2**
- 6.4 Wash cells with warm PBS first, then add **🧴 1 mL** of cell detachment solution to each well in a **6 3m** well plate format. Transfer the plate to the incubator **🔥 37 °C**. After **🕒 00:03:00** take out the plate and observe the cells under the microscope. iPSCs colonies should look loose and be ready to detach, if not please transfer the plate back to the incubator for one minute or two more, but do not exceed the incubation time and do not let the iPSC colonies detach completely and float in the solution.
- 6.5 Remove the cell detachment solution and rinse the cells with **🧴 1 mL** of warm culture medium.  Pipette the cells a couple of times, but be careful not to over-pipette the cells. The cells should stay in clumps and easily detach from the well. Transfer the appropriate amount to new wells.
- 6.6 You can split the cells using ratio for example 1:6 to 1:30 or take a small amount for counting first to have an exact number of cells in a new well.


Note

Cell counting

1. Take  10 μL of cell suspension from the well and transfer it to a small PCR tube.
2. Pipette cell suspension a couple of times to make a single-cell suspension.
3. Add  10 μL of






 Trypan Blue Solution 0.4% Sterile-filtered Merck MilliporeSigma (Sigma-Aldrich) Catalog #T8154



and mix well by pipetting

4. Take  10 μL for counting using a cell counter or Neubauer chamber.


- 6.7 Agitate the plate gently from side to side and forwards and backwards to ensure that the cell aggregates are evenly distributed across the surface of the well containing culture medium with the addition of rock inhibitor.
- 6.8 The next day in the morning, exchange the medium to fresh culture medium without rock inhibitor. Continue to change the medium daily until iPSCs reach confluency again.

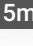7 iPSC cryopreservation

- 7.1 Prepare cryopreservation medium: culture medium without antibiotics with 5%-10%  dimethylsulfoxide (DMSO) Merck MilliporeSigma (Sigma-Aldrich). The cryopreservation medium should be stored in the fridge and keep cold during the procedure.
- 7.2 Firstly perform cell detachment as described in steps  go to step #6.5 and #6.6 but instead of adding the cell suspension to the new well transfer it to  15 mL centrifuge tube.
- 7.3 Centrifuge the tube  1400 rpm, Room temperature, 00:05:00 

7.4 Remove the supernatant with care not to move the cell debris. Add an appropriate amount of cold cryopreservation medium to the tube, and resuspend the cells. Transfer cryotubes to a freezing container and put them in  -80 °C freezer  Overnight . The next day, you can transfer cryotubes to the liquid nitrogen for longer storage.

Note

For example, it is recommended to use  1 mL of cell suspension for 1 cryotube and freeze around 1x10⁶ live cells per vial.

Equipment

Mr. Frosty Freezing Container	NAME
freezing container	TYPE
ThermoFisher Scientific	BRAND
5100-0001	SKU
https://www.thermofisher.com/order/catalog/product/5100-0001	LINK

Differentiation of hiPSCs into cardiomyocytes

8 Standard iPSC culture before differentiation.

8.1 When cells become confluent passage them to prepared geltrex-coated 12-well plate.

8.2 View cells under microscope every day to estimate their confluence. When cells reach desired confluence between 60%-80% they are ready for differentiation.

9 Differentiation of hiPSCs into cardiomyocytes using commercially available kit



PSC Cardiomyocyte Differentiation Kit Thermo
Fisher Catalog #A2921201

9.1 On day 0 of the differentiation remove iPSC culture medium from well and replace it with Cardiomyocyte Differentiation Medium A from the kit.

9.2 On day 2 replace the medium with Cardiomyocyte Differentiation Medium B from the kit.

Note

On this step of protocol cells may start to look different, they become more opaque and many dead cells are floating in the medium.

9.3 On day 4 replace the medium with Cardiomyocyte Maintenance Medium. Change the medium every 2-3 days with Cardiomyocyte Maintenance Medium. Starting from day 8 cells should begin to contract.

Subtype differentiation of hiPSCs into cardiomyocytes

25m

10 Standard iPSC culture before differentiation.

10.1 When cells become confluent passage them to prepared geltrex-coated 12-well plate.

11 Differentiation to cardiomyocyte subtypes

Note


When iPSCs reach 70-90% confluence they are ready to start differentiation


11.1 In the beginning, prepare two types of medium:

1. Cardiomyocytes induction medium
2. Cardiomyocytes culture medium


Note

1. Cardiomyocytes induction medium:


 RPMI 1640 Medium, GlutaMAX[®] Supplement, HEPES Thermo
Fisher Catalog #72400021


with 1x  B-27[™] Supplement, minus insulin Gibco - Thermo
Fisher Catalog #A1895601 from

50x stock and 1% of

 Penicillin-Streptomycin (5,000 U/mL) Thermo
Fisher Catalog #15070063

2. Cardiomyocytes culture medium:

 RPMI 1640 Medium, GlutaMAX[®] Supplement, HEPES Thermo
Fisher Catalog #72400021













with 1x  B-27[™] Supplement (50X), serum free Gibco - Thermo
Fisher Catalog #17504044

from 50x stock and 1% of

 Penicillin-Streptomycin (5,000 U/mL) Thermo
Fisher Catalog #15070063




11.2 Prepare small molecules working solutions of CHIR99021, IWP2 and retinoic acid (RA)

Note

1. Dissolve  StemMACS™ CHIR99021 Miltenyi Biotec Catalog #130-103-926
[M] 4.66 mg/mL in DMSO  37 °C  00:05:00 to obtain [M] 10 millimolar (mM) working solution, aliquot small portions, and store in the freezer  -20 °C
2. Dissolve  IWP-2 1 mg STEMCELL Technologies Inc. Catalog #72122 [M] 2.34 mg/mL
in DMSO  37 °C  00:20:00 to obtain [M] 5 millimolar (mM) working solution, aliquot small portions, and store in the freezer  -20 °C
3. Dissolve  Retinoic acid Merck MilliporeSigma (Sigma-Aldrich) Catalog #R2625-50MG
[M] 29.94 mg/mL in DMSO to obtain [M] 100 millimolar (mM) master stock, aliquot  10 μ L portions, and store in the freezer  -80 °C. To obtain [M] 1 millimolar (mM) working solution thaw and dilute RA aliquot in  990 μ L EtOH and store in the fridge for 2 weeks, protecting from light.

12 Atrial subtype differentiation

1w 1d

- 12.1 On day 0 change the iPSC culture medium to cardiomyocytes induction medium with an addition of [M] 5 micromolar (μ M) CHIR99021 ( 0.5 μ L /ml from [M] 10 millimolar (mM) working stock solution).
- 12.2 After 48h on day 2 change the medium to fresh cardiomyocytes induction medium with an addition of [M] 5 micromolar (μ M) IWP2 ( 1 μ L /ml from [M] 5 millimolar (mM) working stock solution).
- 12.3 On day 3, without changing the medium add  1 μ L /ml RA [M] 1 micromolar (μ M) final from [M] 1 millimolar (mM) working stock solution .
- 12.4 On day 4 change the medium to fresh cardiomyocytes induction medium with an addition of [M] 1 micromolar (μ M) RA.

12.5 On day 6 change the medium to fresh cardiomyocytes induction medium.

12.6 From day 8 change the medium to cardiomyocytes culture medium and continue media change every 2-3 days.

13 Ventricular subtype differentiation

1w 1d

13.1 On day 0 change the iPSC culture medium to cardiomyocytes induction medium with an addition of **5 micromolar (μM) CHIR99021** (**0.5 μL /ml** from **10 millimolar (mM) working stock solution**).

13.2 After 48h on day 2 change the medium to fresh cardiomyocytes induction medium with an addition of **5 micromolar (μM) IWP2** (**1 μL /ml** from **5 millimolar (mM) working stock solution**).

13.3 On day 4 change the medium to fresh cardiomyocytes induction medium.

13.4 On day 6 change the medium to fresh cardiomyocytes induction medium.

13.5 From day 8 change the medium to cardiomyocytes culture medium and continue media change every 2-3 days.

14 Nodal subtype differentiation

1w

- 14.1** On day 0 change the iPSC culture medium to cardiomyocytes induction medium with an addition of [M] 12 micromolar (μM) CHIR99021 (A 1.2 μL /ml from [M] 10 millimolar (mM) working stock solution).
- 14.2** After 24h on day 1 change the medium to fresh cardiomyocytes induction medium.
- 14.3** On day 3 change the medium to fresh cardiomyocytes induction medium with an addition of of [M] 5 micromolar (μM) final IWP2 (A 1 μL /ml from [M] 5 millimolar (mM) working stock solution .)
- 14.4** On day 5 change the medium to fresh cardiomyocytes induction medium with an addition of [M] 3 micromolar (μM) CHIR99021 (A 0.3 μL /ml from [M] 10 millimolar (mM) working stock solution).
- 14.5** From day 7 change the medium to cardiomyocytes culture medium and continue media change every 2-3 days.

Cardiomyocytes culture, characterization, passaging and cryo... 25m

15 Longer culture

- 15.1** Cardiomyocytes and their specific subtypes usually start contracting after 8-10 days from the beginning of the differentiation. After 12-15 days they become defined enough to be safely passaged or cryopreserved. Cells can be kept in a culture for at least 1 month or even longer if such experiments are desired.

Note

With time cardiomyocytes become more mature and gain more cardiomyocyte-like characteristics.

16 Characterization of cardiomyocytes


16.1 Successfully differentiated cells display completely different morphology, with characteristic net-like structure. The most important feature is their ability to contract. It is also recommended to perform PCR analysis or immunocytochemistry to confirm the expression of key cardiac markers such as TNNT2 and NKX2.5.

16.2 Immunocytochemistry confirmation of key cardiac markers expression in differentiated cardiomyocyte-like cells using commercially available kit



Human Cardiomyocyte Immunocytochemistry Kit Thermo
Fisher Catalog #A25973

Expected result

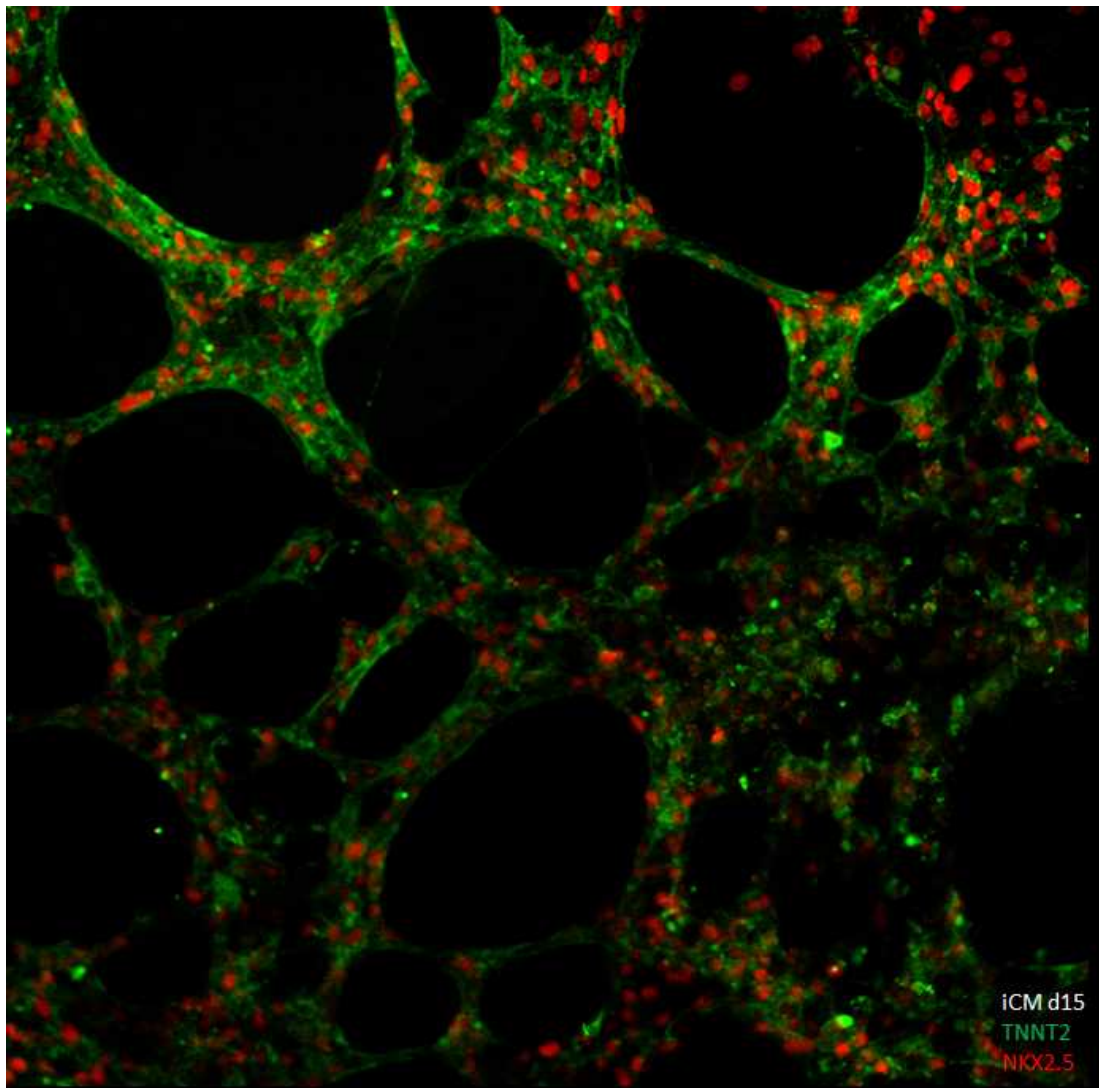

17 Passaging cardiomyocytes

17.1 Prepare a new Geltrex-coated plate [⇒ go to step #1.3](#)

- 17.2 For better cell survival post-passaging it is recommended to add **10 micromolar (μM)** **ROCK inhibitor (Ri) Y-27632 STEMCELL Technologies Inc. Catalog #72308** to the culture medium approximately 1 hour before passage.
- 17.3 Aspirate the medium and wash cells with warm PBS. Add warm **StemPro™ Accutase™ Cell Dissociation Reagent Thermo Fisher Scientific Catalog #A1110501** to the well (**1 mL per well in 6-well plate**) and incubate **37 °C** **00:05:00** - **00:10:00** **15m**
- 17.4 After 5 minutes observe the cells under the microscope to see if they are ready to detach, if not continue to incubate up to 10 minutes maximum. When cells look ready to be detached pipette the solution a couple of times to detach the cells from the well, be careful not to overpipette the cells, because it may disrupt them.
- 17.5 Transfer them to **15 mL** centrifuge tube, add an appropriate amount of culture medium (1:1), **5m** centrifuge the tube **1400 rpm, Room temperature, 00:05:00**
- 17.6 After centrifugation aspirate the supernatant and resuspend the cells in cardiomyocytes culture medium. You can take a small portion for cell counting and place desired amount of cells on a new geltrex-coated plate with addition of **10 micromolar (μM)** rock inhibitor for better cell recovery.
- 17.7 Agitate the plate gently from side to side and forwards and backwards to ensure that the cell aggregates are evenly distributed across the surface of the well. Put the plate in the incubator. The next day in the morning change the medium to fresh culture medium and continue media change every 2-3 days.

18 Cryopreservation of cardiomyocytes

Note

if cardiomyocytes and their subtypes are not required for any ongoing experiment, they can be cryopreserved from day 12-15 onwards and recovered for culture when needed. It is recommended to freeze around 2×10^6 cells per cryovial.

- 18.1 After centrifugation from step **go to step #17.5** you can aspirate the supernatant and resuspend cells in cold cryopreservation medium.

Note

Cardiomyocytes cryopreservation medium:

[M] 50 % volume



Fetal Bovine Serum, qualified, One Shot[®] format, United States Thermo
Fisher Catalog #A3160501

[M] 40 % volume



RPMI 1640 Medium, GlutaMAX[®] Supplement, HEPES Thermo
Fisher Catalog #72400021

[M] 10 % volume



dimethylsulfoxide (DMSO) Merck MilliporeSigma (Sigma-
Aldrich)

18.2 Add an appropriate amount of cold cryopreservation medium to the tube, and resuspend the cells. Transfer cryotubes to a freezing container and put them in -80 °C freezer Overnight . The next day, you can transfer cryotubes to the liquid nitrogen for longer storage.

18.3 Post-thaw recovery

18.4 Prepare the previously coated 6-well plate by adding 2 mL of warm culture medium to each well and 2 μ L from [M] 10 millimolar (mM) stock solution ROCK inhibitor (Ri) Y-27632 STEMCELL Technologies Inc. Catalog #72308 to minimize the cell death.

18.5 Take out the desired cryovial from the liquid nitrogen and quickly transfer it to 37 °C water bath. Be careful not to submerge the cap of the cryovial and keep a little ice cube in the cryovial before taking it out. Disinfect the cryovial with 70% ethanol.

18.6 Add 9 mL of warm cardiomyocyte culture medium to 15ml centrifuge tube. Centrifuge 1400 rpm, Room temperature, 00:05:00 **5m**

18.7 After centrifugation aspirate the supernatant and resuspend the cells in culture medium. You can take small portion for cell counting.

Note

It is recommended to seed around $1-3 \times 10^6$ cells per well in 6-well plate.

- 18.8** Transfer the cells to a new plate, agitate the plate gently from side to side and forwards and backwards to ensure that the cell aggregates are evenly distributed across the surface of the well. Put the plate in the incubator. Next day change the medium without rock inhibitor and continue media change every 2-3 days.