

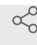


Jul 11, 2022

6-plex IF Protocol on the BOND RX (Leica Biosystems)

Maura A Sticco-Ivins¹, Abdallah Flaifel^{1,2}, Miriam Ficial^{1,2}¹Department of Pathology, Brigham and Women's Hospital, Boston, MA, 02115, USA;²Harvard Medical School, Boston, MA, 02215, USA

1 Works for me

 Sharedx.doi.org/10.17504/protocols.io.bjbzkip6 Miriam Ficial

ABSTRACT

A 6-plex immunofluorescence protocol was optimized on Leica Bond RX Autostainer.

DOI

dx.doi.org/10.17504/protocols.io.bjbzkip6

PROTOCOL CITATION

Maura A Sticco-Ivins, Abdallah Flaifel, Miriam Ficial 2022. 6-plex IF Protocol on the BOND RX (Leica Biosystems). **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.bjbzkip6>



LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Aug 04, 2020

LAST MODIFIED

Jul 11, 2022

PROTOCOL INTEGER ID

40025

Preparation

1

Multiplex IF assay on the Bond RX Autostainer (Leica Biosystems) using the Opal multiplex IHC system (PerkinElmer/Akoya Biosciences Cat# NEL871001KT) and

the BOND Polymer Refine Detection Kit (Leica Biosystems Cat# DS9800).

- Prepare the antibody and the fluorophore dilutions in 6mL Leica BOND tubes (see Ficial M et al for details about reagents, dilutions and order of antibodies).

- Fill an empty 7ml Leica BOND container with DAPI.

- Fill an empty 30ml Leica BOND container with Dako Protein Block.

2 Bake the unstained FFPE slides in the Isotemp Oven for 30 mins at 60°C.

Leica BOND RX Protocol

3 Deparaffinization and Rehydration

- BOND Dewax Solution (Leica Biosystems Cat# AR9222)

Time. (min): 30:00, Temperature: 72°C, Dispense type: open.

- Rehydration.

4 Bond ER Solution 1 (Leica Biosystems Cat# AR9961; citrate, pH 6.0)

Time. (min): 10:00, Temperature: 98°C, Dispense type: 150 µL.

5 Deionized Water

Time. (min): 2:00; Temperature: Ambient; Dispense type: 150 µL.

6 Bond Wash Solution (Leica Biosystems Cat# AR9590)

Time. (min): 2:00; Temperature: Ambient; Dispense type: 150 µL.

7 Peroxide Block (Leica Biosystems Cat# DS9800)

Time. (min): 10:00; Temperature: Ambient; Dispense type: 150 µL.

8 Bond Wash Solution (Leica Biosystems Cat# AR9590)

Time. (min): 4:00; Temperature: Ambient; Dispense type: 150 µL.

9 Dako Protein Block (Agilent Cat# X090930-2)

Time. (min): 15:00; Temperature: Ambient; Dispense type: 150 µL.

- 10 Primary Antibody (diluted in either Leica Biosystems BOND Primary Antibody Diluent Cat# AR9352 or Biocare Medical DaVinci Green Diluent Cat# PD900 M)
Time. (min): 40:00; Temperature: Ambient; Dispense type: 150 µL.

- 11 Bond Wash Solution (Leica Biosystems Cat# AR9590)
Time. (min): 4:00; Temperature: Ambient; Dispense type: 150 µL.

- 12 Post-Primary (Leica Biosystems Cat# DS9800) *
Time. (min): 15:00; Temperature: Ambient; Dispense type: 150 µL.

*Post-Primary should be used after mouse anti-human primary antibodies only; when using rabbit anti-human primary antibodies, please go directly to step 14.

For TIM-3 this step was replaced by:

Rabbit anti Goat Link

Time. (min): 15:00; Temperature: Ambient; Dispense type: 150 µL

- 13 Bond Wash Solution (Leica Biosystems Cat# AR9590)
Time. (min): 4:00; Temperature: Ambient; Dispense type: 150 µL.

- 14 Polymer (Leica Biosystems Cat# DS9800)
Time. (min): 10:00; Temperature: Ambient; Dispense type: 150 µL.

- 15 Bond Wash Solution (Leica Biosystems Cat# AR9590)
Time. (min): 6:00; Temperature: Ambient; Dispense type: 150 µL.

- 16 Opal Fluorophore (PerkinElmer/Akoya Biosciences Cat# NEL871001KT)
Time. (min): 5:00; Temperature: Ambient; Dispense type: 150 µL.

- 17 Bond Wash Solution (Leica Biosystems Cat# AR9590)
Time. (min): 6:00; Temperature: Ambient; Dispense type: 150 µL.

18 Bond ER Solution 1(Leica Biosystems Cat# AR9961; citrate, pH 6.0)
Time. (min): 10:00, Temperature: 98°C, Dispense type: 150 µL.

19 **Start a new cycle of staining from step 5 OR go to step 20**

20 Last cycle of staining:

20.1 Deionized Water
Time. (min): 2:00; Temperature: Ambient; Dispense type: 150 µL.

20.2 DAPI (PerkinElmer/Akoya Biosciences Cat# NEL871001KT)
Time. (min): 10:00; Temperature: Ambient; Dispense type: 150 µL.

20.3 Manually coverslip using Invitrogen ProLong Diamond Antifade Mountant (Thermo Fisher Scientific Cat# P36962).