



Version 2 ▼

Jul 07, 2022

Extraction of bacterial DNA using MagMAX™ CORE Nucleic Acid Purification Kit on KingFisher™ Flex Instrument V.2

Leyi Wang¹, Carol Maddox¹, Melanie Prarat², Yan Zhang², Lifang Yan³, Akhilesh Ramachandran⁴, Sai Sankara Narayanan⁴, Girish Patil⁴


¹University of Illinois Veterinary Diagnostic Laboratory;

²Ohio Department of Agriculture Animal Disease Diagnostic Laboratory;

³Mississippi Veterinary Research and Diagnostic Lab;

⁴Oklahoma Animal Disease Diagnostic Laboratory

1 Works for me

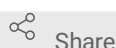

Share

dx.doi.org/10.17504/protocols.io.81wgb781ovpk/v2

Vet LIRN



Sarah Nemser

DISCLAIMER

Reference to any commercial materials, equipment, or process does not in any way constitute approval, endorsement, or recommendation by the Food and Drug Administration.

ABSTRACT

This procedure is used to extract genome DNA of bacteria isolates using the MagMAX™ CORE Nucleic Acid Purification Kit on the KingFisher Flex Robot. The process consists of Sample processing, Prepare plates for Robot, and Operate on KingFisher Flex machine.

DOI

dx.doi.org/10.17504/protocols.io.81wgb781ovpk/v2

PROTOCOL CITATION

Leyi Wang, Carol Maddox, Melanie Prarat, Yan Zhang, Lifang Yan, Akhilesh Ramachandran, Sai Sankara Narayanan, Girish Patil 2022. Extraction of bacterial DNA using MagMAX™ CORE Nucleic Acid Purification Kit on KingFisher™ Flex Instrument. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.81wgb781ovpk/v2>
Version created by Sarah Nemser

FUNDERS ACKNOWLEDGEMENT



Vet-LIRN

Grant ID: 1U18FD06983-01

KEYWORDS

Robot, KingFisher Flex, MagMAX CORE, DNA Extraction

LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Jul 07, 2022

LAST MODIFIED

Jul 07, 2022

PROTOCOL INTEGER ID

66208

GUIDELINES

In the whole procedure, need wear PPE

MATERIALS TEXT

Bacterial culture isolates

MagMAX™ CORE Nucleic Acid Purification Kit (A32702 or A32700)

SAFETY WARNINGS

Need work in the biosafety level 2 hood to scoop the bacteria in the buffer tube

DISCLAIMER:

Reference to any commercial materials, equipment, or process does not in any way constitute approval, endorsement, or recommendation by the Food and Drug Administration.

1 Prepare plates for Robot

1.1 Prepare Wash Plate 2 by adding 500 µL of MagMAX™ CORE Wash Solution 2

1.2 Prepare Wash Plate 1 by adding 500 µL of MagMAX™ CORE Wash Solution 1

1.3 Prepare Elution Plate by adding 100 µL of MagMAX™ CORE Elution Buffer

1.4 Set Tip Comb in a 0.5 ml 96-Well Plate

A	B	C	D
Plate setup of Processing Plates: KingFisher™ Flex instrument			
Plate ID	Plate Type	Reagent	Volume Per Well
Wash Plate 1	Deep Well	MagMAX™ CORE Wash Solution 1	500 µL
Wash Plate 2	Deep Well	MagMAX™ CORE Wash Solution 2	500 µL
Elution Plate	Standard	MagMAX™ CORE Elution Buffer	100 µL
Tip Comb Plate	Standard	Place a tip comb in the plate	Place a tip comb in the plate

1.5 Prepare Bead/PK mix by adding 20 µL of MagMAX™ CORE Magnetic Beads and 10 µL of MagMAX™ CORE Proteinase K for the required number of samples plus 10% overage.

1.6 Prepare Lysis/Binding solution by adding 350 µL of MagMAX™ CORE Lysis Solution and 350 µL of MagMAX™ CORE Binding Solution for the required number of samples plus 10% overage. Mix by inverting the tube or bottle at least 10 times.

1.7 Invert the tube of Bead/PK Mix several times to resuspend the beads, then add 30 µL of the Bead/PK Mix to the required wells in the 2.4 mL 96- deep well sample plate.

- 1.8 Transfer 200 µL sample to a well with 30 µL Bead/PK mix in Sample plate. Shake vigorously for 2 minutes on a plate shaker at room temperature, or pipette mix them and incubate for 2 mins at room temperature.
- 1.9 Add 700 µL of Lysis/Binding Solution to each sample-containing well.

2 Sample processing

From plate, scoop a full loop of 1-µl loop to 200 µL of PBS or molecular grade water solution. Mix the sample by vortexing.

3 Operate on KingFisher Flex machine

- 3.1 Select the right script on the instrument or load the script if it was not already in the instrument.

Download kit script KingFisher Flex heated script: MagMAX CORE_Flex.bdz at <https://www.thermofisher.com/order/catalog/product/A32700#/A32700>, and then use BindIt software to transfer this script from a computer to the machine with the connection line.

BindIt software can be requested and downloaded at the website <https://www.thermofisher.com/us/en/home/global/forms/life-science/download-bindit-software-kingfisher-instruments.html>

- 3.2 Start the run, then load the prepared plates in the appropriate positions when prompted by the instrument. The instrument will direct the user as to which plate to load in what position.
- 3.3 Load the Tip Comb Plate
- 3.4 Load the Elution Plate
- 3.5 Load the Wash Plate 1

3.6 Load Wash Plate 2

3.7 Lastly, load the sample plate

3.8 Click Start Button, and the machine will process and complete in around 20 mins

4 Laboratory specific variations from protocol

Ohio ADDL: use molecular grade water to resuspend the colonies.

Illinois VDL, Mississippi VRDL, and Oklahoma ADDL: use PBS to resuspend the colonies.

Oklahoma ADDL: uses Dry Garnet 0.7mm beads (component of QiAMP fecal kit) to lyse the colony solution on Qiagen TissueLyzer II (Cat. No. / ID: 85300) for 90seconds at 30Hz.