



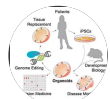
SEP 11, 2023

## Passaging of organoids (organoids split into cell clusters)

Gabriela Vallejo  
Flores<sup>1</sup>,

Annika  
Fendler<sup>1</sup>

<sup>1</sup>Charite



Gabriela Vallejo Flores

### ABSTRACT

The goal of this experiment is to split the organoids into cell clusters and induce organoids formation and organoids expansion.

### MATERIALS

Corning® Matrigel® Growth Factor Reduced (GFR) Basement Membrane Matrix Corning Catalog #356231

TrypLE™ Express Enzyme Thermo Fisher Scientific Catalog #12604013

Poly(2-hydroxyethyl methacrylate) Merck MilliporeSigma (Sigma-Aldrich) Catalog #P3932

Advanced DMEM/F-12 Thermo Fisher Catalog #12634028

FBS Superior stabil BIOSELL Catalog #FBS. S 0615

- 48 multiwell plate
- 24 multiwell plate
- 6 multiwell plate
- 15 ml tubes

OPEN ACCESS



**Protocol Citation:** Gabriela Vallejo Flores, Annika Fendler 2023. Passaging of organoids (organoids split into cell clusters). [protocols.io](https://protocols.io/view/passaging-of-organoids-organoids-split-into-cell-clusters) <https://protocols.io/view/passaging-of-organoids-organoids-split-into-cell-clusters>

**License:** This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

**Protocol status:** In development  
We are still developing and optimizing this protocol

**Created:** Jul 24, 2023

**Last Modified:** Sep 11, 2023

**PROTOCOL integer ID:**  
85419

**Keywords:** Organoids,  
Splitting, Clusters

## BEFORE START INSTRUCTIONS

- Prepare the organoids medium one day before organoids splitting
- Pre-warm seeding plate
- work with 1% BSA/PBS coated tubes and pipettes
- use pre-cooled pipette tips for matrigel
- Thaw the matrigel 1hr before starting the experiment and prepare the multiwell plate for each protocol;



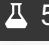

Protocol	Matrix	Preparation time
Minoli et al., 2023	Poly-HEMA	1 day before the experiment
Lee et al., 2018	Bed-Matrigel	1hr before the experiment
Mullenders et al., 2019	Dome-Matrigel	At the moment

## Organoids Splitting

1



Date	Sample ID	Splitting wells/passage/protocol	New wells/passage/protocol

Table 1. Sample description

- Remove the supernatant from the well and add  1 mL cool DMEM+1%FCS to release the organoids from the matrigel,  400 x g, 4°C, 00:05:00, and aspirate the supernatant and Matrigel above the pellet. Note: For the minoli protocol "organoids grow in cell suspension", collect the supernatant and do the centrifugation. 5m
- Resuspend the organoids in  500 µL of TrypLE Express (double the amount of Tryple Express when the matrigel is still visible after removing the supernatant, and incubate for  00:15:00 at RT with frequent perturbation by pipetting, check organoid dispersion under the bright-field microscope every 5 min. 15m

Quantity of TripLE exprees per well, Table 1

A	B	C
24 well plate	500µl	per well
48 well plate	250µl	per well

- 4 Add  1 mL cool DMEM+1%FCS to dilute the TripLE express and inactivate the enzyme,  400 x g, 4°C, 00:05:00, and aspirate the supernatant.

5m

- 5 Homogenize the pellet in appropriate volume of Matrigel and pipet into a pre warmed plate: prewarmed

A	B	C
24 well plate	50µl	per well
48 well plate	25µl	per well

## Cell seeding Mullenders

- 6 Add matrigel in pre warmed plate, as described in Table

Plate type	Matrigel	No. cells	Medium	E
48	25µl	$5 \times 10^4$	250µl	
6	5 drop of 50µl	$2.5 \times 10^5$ per drop	2ml	
24	50µl	$2.5 \times 10^5$	500µl	

- 7 Flip the 48,24 or 6 wells plate and incubate 15 min at 37°C 5%CO<sub>2</sub>.

- 8 Add 250µl, 500µl or 2ml of medium per well.

## Cell seeding Lee

- 9 Add the Matrigel per  $20 \times 10^3$  cells, per well and incubate the plate during 15 min at  $37^\circ\text{C}$ , 5%CO<sub>2</sub> ;  
Add the matrigel, as described in the Table 4.

	Plate Type	Matrigel	No. Cells per drop	Medium
		25µl	$5 \times 10^4$	250µl
	24	50µl	$2.5 \times 10^5$	500µl
	6	5 drop of 50µl	$2.5 \times 10^5$ per drop	2ml

- 10 Add 250µl, 500µl or 2ml of medium per well.