




Jun 24, 2021

RNA Extraction Full Protocol

Onkar Nath¹, Stephen Fletcher¹, Alice Hayward¹, Agnelo Furtado¹, Robert J Henry¹, Neena Mitter¹

¹Queensland Alliance for Agriculture and Food Innovation, University of Queensland, Brisbane 4072, Australia

1 Works for me

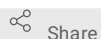

Share

This protocol is published without a DOI.

qaafi QAAFI, UQ

Onkar Nath

DISCLAIMER

DISCLAIMER – FOR INFORMATIONAL PURPOSES ONLY; USE AT YOUR OWN RISK

The protocol content here is for informational purposes only and does not constitute legal, medical, clinical, or safety advice, or otherwise; content added to protocols.io is not peer reviewed and may not have undergone a formal approval of any kind. Information presented in this protocol should not substitute for independent professional judgment, advice, diagnosis, or treatment. Any action you take or refrain from taking using or relying upon the information presented here is strictly at your own risk. You agree that neither the Company nor any of the authors, contributors, administrators, or anyone else associated with protocols.io, can be held responsible for your use of the information contained in or linked to this protocol or any of our Sites/Apps and Services.

ABSTRACT

RNA Extraction protocol suitable for all Avocado tissues. This protocol also suits extraction from many other plants like Capsicum, Finger lime and Cotton.

PROTOCOL CITATION

Onkar Nath, Stephen Fletcher, Alice Hayward, Agnelo Furtado, Robert J Henry, Neena Mitter 2021. RNA Extraction Full Protocol. **protocols.io**
<https://protocols.io/view/rna-extraction-full-protocol-bv2dn8a6>

KEYWORDS

avocado, rna extraction, high quality RNA, CTAB

LICENSE

This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Jun 23, 2021

LAST MODIFIED

Jun 24, 2021

PROTOCOL INTEGER ID

50981

GUIDELINES

for modifications required with different tissue types follow associate manuscript.

MATERIALS TEXT

Buffers and chemicals:

OB1:

CTAB Buffer (pH 7.6, filter sterilized with 33mm PES membrane filter unit): 2% CTAB, 1.4M NaCl, 10mM EDTA, 0.1M MES, 1% PEG

For DNA extraction 0.5% Na₂SO₃ (e.g. 0.05 g in 10 ml), 2% PVP (0.2 g in 10 ml), 2% β-mercaptoethanol (e.g. 20 µl in 10 ml) was added before use.

For RNA Extraction, 1/100th volume 10% SDS, 2% β-mercaptoethanol (e.g. 20 µl in 10 ml) was added before use.

For 100ml:

CTAB 2.0 g (hexadecyltrimethyl-ammonium bromide)

NaCl 8.18 g (sodium hydroxide)

EDTA 0.3723 g (ethylenediaminetetra acid disodium salt)

MES 1.96 g (2-(N-morpholino) ethanesulfonic acid hydrate)

PEG-8000 1.0 g (polyethylene glycol 8000)

OB2:

Buffer 2: 25mM sodium citrate, 0.5% lauroyl sarcosine, 1M sodium acetate, 2% PVP and 4.2M guanidine thiocyanate (optional)

For 100ml:

sodium citrate 0.535 g

lauroyl sarcosine 0.5 g

sodium acetate 13.605 g

PVP 2 g

guanidine thiocyanate 49.625 g

Optional Clean-up Buffer (OB3/SSTE) (Chang, Puryear, and Cairney 1993): 1M NaCl, 0.5% SDS, 10mM Tris pH 8.0 and 1mM EDTA pH 8.0.

chloroform:isoamyl alcohol (C:IAA) 24:1

Isopropanol 100%

Sodium dodecyl sulfate (SDS) 10% (To be used within two months as its efficiency reduces with time)

SAFETY WARNINGS

Use of various instruments and hazardous chemicals are required. Consult and follow GHS guidelines.

DISCLAIMER:

DISCLAIMER – FOR INFORMATIONAL PURPOSES ONLY; USE AT YOUR OWN RISK

The protocol content here is for informational purposes only and does not constitute legal, medical, clinical, or safety advice, or otherwise; content added to protocols.io is not peer reviewed and may not have undergone a formal approval of any kind. Information presented in this protocol should not substitute for independent professional judgment, advice, diagnosis, or treatment. Any action you take or refrain from taking using or relying upon the information presented here is strictly at your own risk. You agree that neither the Company nor any of the authors, contributors, administrators, or anyone else associated with protocols.io, can be held responsible for your use of the information contained in or linked to this protocol or any of our Sites/Apps and Services.

BEFORE STARTING

Prepare buffers and organise chemicals

Tissue Lysis

4h 45m

- 1 Finely grind frozen tissue. Grinding in the buffer is not recommended.

2 Preparation on the Day of Use

- Add **10 µl 10% SDS per 2ml OB1 buffer** and **1mM 2% (v/v) β-mercaptoethanol** (Fume-hood should be used for working with 2ME)

- Heat the buffer to \uparrow **65 °C** .

- 3 Add approximately \square **10 mg** (0.01 g) of ground frozen leaf tissue to \square **650 μ l** hot OB1 buffer in 2 ml tube and incubate at \uparrow **65 °C** for \odot **00:05:00** . 5m

Over-incubation leads to RNA degradation and delicate tissues for e.g. roots cannot tolerate this. Tissue-specific recommendations are reported in RNA extraction section

- 4 Centrifuge at \odot **20817 x g** (14,000 rpm) for \odot **00:10:00** at \uparrow **Room temperature** , pour supernatant to another tube. This improves the cleaning process. 10m

- 5 Add an equal volume of **OB2 buffer** .

Removing Impurity 4h 30m

- 6 Add an equal volume of **chloroform isoamyl alcohol** (C:IAA); to be performed in fume-hood. Vortex tubes well.

- 7 Centrifuge at \odot **20817 x g** (14,000 rpm) for \odot **00:20:00** at \uparrow **Room temperature** . 20m

- 8 Collect the supernatant and repeat chloroform wash.

Precipitation 4h 30m

- 9 Add an equal volume of 100% **Isopropanol** (Cold).

Preferably precipitate with a half volume of \square **10 Molarity (M)** LiCl instead of isopropanol, if the ratio of high molecular weight and low molecular weight RNA does not matter for the experiment.

- 10 10. Incubate at \uparrow **-20 °C** for \odot **01:00:00** or longer as needed. Incubating longer increases yield but adds more salt in precipitant. 1h

- 11 11. Centrifuge at \odot **12857 x g, 4°C, 00:25:00** (10,000 rpm). 25m

- 12 12. Wash pellet with \square **1 mL** of \square **75 % (v/v) ethanol** .

13 13. Centrifuge at **7674 x g, 4°C, 00:10:00** (8,500 rpm)

10m

Elution 4h 30m

14 14. Air Dry the pellet. Do not over-dry but get rid of maximum ethanol possible.

15 15. Dissolve in 20 µl of ultrapure DNase RNase free water.

15.1 Pool tubes if required at this stage for higher RNA content.

DNA Degradation 4h 30m

16 17. If DNA must be removed from the RNA sample, add 1/20th volume of

[DNase I Reaction Buffer - 6.0 ml New England](#)

Biolabs Catalog #B0303S

and 1/20th volume of

[DNase I \(RNase-free\) New England](#)

Biolabs Catalog #E6316

17 18. Incubate at **37 °C** for **00:30:00**.

30m

DNase Enzyme removal 4h 30m

18 19. Add water to make up **500 µl** of total volume.

RNAclean XP Bead protocol from Beckman Coulter can also be followed to remove DNase enzyme. TRISure protocol has been followed because of cost factor.

19 20. Add an equal volume of [TRISure Bioline Catalog #BIO-38033](#)

Use OB3 if RNA degrades with the use of TRISure.

20 21. Incubate for **00:05:00**.

5m

21 22. Add **100 µl** chloroform.

- 22 23. Centrifuge at  **20817 x g, 4°C, 00:20:00** (14,000 rpm), collect supernatant. 20m
- 23 24. Add an equal volume of **isopropanol**
- 24 25. Incubate at  **-20 °C** for  **01:00:00** or longer. 1h
- 25 26. Centrifuge at  **20817 x g, 4°C, 00:40:00** (14,000 rpm). 40m
- 26 27. Ethanol wash pellet.
- 27 28. Dissolve in 20 µl of ultrapure DNase RNase free water or 1x TE buffer as per requirement.

Allow the pellet to dissolve and store at -20 or -80 as required. Do not leave at room temperature or 4°C for long.

Optional clean up and assessment

- 27.1 Optionally if the pellet is still not clean, follow RNAClean XP Bead protocol from Beckman Coulter. Briefly, add beads to solution, separate beads attached to RNA using magnetic rack, ethanol wash and elute in DNase RNase free water (use 100µl buffer instead of the recommended, 90µl per 50µl sample) (Follow manuscript's RNA Extraction Section for specific recommendations and discussion).
- 28 30. Quality-assure using Nanodrop and Quantify using Quantifluor RNA, or Bioanalyzer.