




JUL 25, 2023

Live-cell imaging for synaptic vesicle precursors in human iNeuron axons

Dan
Dou¹, Erika Holzbaur¹

¹University of Pennsylvania

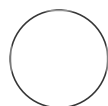

Dan Dou

OPEN ACCESS

DOI:
dx.doi.org/10.17504/protocols.io.5jyl8p9mdg2w/v1

Protocol Citation: Dan Dou, Erika Holzbaur 2023. Live-cell imaging for synaptic vesicle precursors in human iNeuron axons. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.5jyl8p9mdg2w/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Jul 20, 2023

Last Modified: Jul 25, 2023

PROTOCOL integer ID:
85304

Keywords: ASAPCRN

ABSTRACT

Here, we describe procedure and equipment used for live-imaging of synaptic vesicle precursors. This was performed using DIV21 human iPSC-derived excitatory glutamatergic neurons. Equipment and software used varied based on scheduled upgrades to microscopy equipment during the course of this study.

MATERIALS

1. GlutaMAX™ SupplementGibco - Thermo FisherCatalog #35050061
2. B-27™ Supplement (50X), serum freeGibco - Thermo FisherCatalog #17504044
3. Hibernate A low fluorescence media (BrainBits, Cat# HALF)
4. Recombinant Human NT-3peprotechCatalog #450-03
5. Recombinant Human/Murine/Rat BDNFpeprotechCatalog #450-02

SAFETY WARNINGS



Safety information

Investigators should be trained and familiar with the confocal microscope to avoid eye damage from lasers.

Note

Please refer to "Protocol: Culture and transfection of iPSC-derived neurons for live-imaging of axonal cargoes" for plating and transfection instructions.

Image human iNeurons on DIV21, 48-72 hours after transfection with PGK-mScarlet-synaptophysin.

2 Replace culture media with low fluorescence imaging media.


2.1 For iNeurons, use Hibernate A medium supplemented with:

A	B
BDNF	10 ng/mL
NT-3	10 ng/mL
B-27	2%

3 Image using spinning disk confocal microscope under 60x magnification (oil immersion objective). See "Materials and Methods" for specific microscopes and cameras used.

4 Identify axons of transfected neurons based on morphological parameters. (Boecker et al., 2020; Kaech and Banker, 2006). For example, axons can most reliably be identified by their length and should span over at least 500 μm .

5 Identify the neuronal soma and measure $\sim 100\text{-}150\text{ }\mu\text{m}$ from the soma. Create an ROI that includes a segment of the axon that is in the same Z plane. Image several frames at this field of view prior to photobleaching. Perform one photobleaching cycle with the 405 nm laser for 3 ms/pixel. Dramatically decreased mScarlet-SYP signal should be observed in the ROI following photobleaching, and entry of mScarlet-SYP+ vesicles into the bleached region should be readily observed.

6 Acquire time lapse recordings at a frame rate of 5 frames per second for  00:05:00 .

5m

Note

- Rapid framerate is preferable due to the high transport speed of SVPs
- Knowledge of the pixel/micron ratio for the specific objective and camera being used is necessary for accurately measuring distances.