



MAR 14, 2023

OPEN  ACCESS

DOI:
dx.doi.org/10.17504/protocols.io.36wgqjpkxvk5/v1

Protocol Citation: Presha Rajbhandari, Taruna Neelakantan, Brent R. Stockwell 2023. Immunohistochemistry of liver tissue sections. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.36wgqjpkxvk5/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Mar 14, 2023

Last Modified: Mar 14, 2023

PROTOCOL integer ID:
78790

Keywords: IHC,
immunohistochemistry

Immunohistochemistry of liver tissue sections


 Forked from [Immunohistochemistry of liver tissue sections](#)

Presha Rajbhandari¹, Taruna Neelakantan¹, Brent R. Stockwell¹

¹Columbia University

Human BioMolecular Atlas Program (HuBMAP) Method Development Community

Tech. support email: Jeff.spraggins@vanderbilt.edu

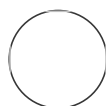


Presha Rajbhandari


ABSTRACT

This protocol outlines the steps used to perform standard immunohistochemistry for antibody validation in frozen human liver tissue samples, performed at Molecular Pathology Core facility at Columbia University.


Cryosection the frozen liver tissue at 5µm thickness and place it on charged slide


1


2 Air dry the sections for 3 minutes  03:00:00


3 Fix the tissue sections in cold acetone for 15 minutes  00:15:00 . Alternatively, fix with 1% paraformaldehyde at 4C for 15min followed by ice-cold methanol at -20C for 5min.



15m


4 Air dry the sections at room temperature for 2 minutes  00:02:00






5 Incubate the slides in 0.3% hydrogen peroxide in PBS for 5 minutes  00:05:00 to block peroxidase activity

6 Wash slides with PBS 3 times, for 5 minutes each  00:05:00 X3

7 1. Block the tissue sections in 10% normal goat serum or 5% Horse serum with 0.1% BSA for 20 minutes  00:20:00

8 Incubate the tissue sections with primary antibody diluted in DAKO antibody diluent at room temperature for 1.5-2 hours  01:30:00 -  02:00:00

9 Wash the slides with PBS 3 times, for 5 minutes each  00:05:00 each - 3 times

- 10** Incubate the tissue sections with biotinylated secondary antibody diluted in PBS at room temperature for 30-45 minutes  00:30:00 -  00:45:00
- 11** Wash the slides with PBS 3 times, for 5 minutes each  00:05:00 each - 3 times
- 12** Incubate the tissue sections with ABC (Avidin-Biotin complex) peroxidase solution at room temperature for 30 minutes  00:30:00
- 13** Wash the slides with PBS 3 times, 5 minutes each  00:05:00 each - 3 times
- 14** Incubate the tissue sections with DAB (3,3'-diaminobenzidine) peroxidase substrate solution until desired color intensity is reached and immerse slides in distilled water
- 15** Counterstain with hematoxylin and rinse with distilled water
- 16** Dehydrate the sections using 95% ethanol followed by 100% ethanol
- 17** Clear with xylene and mount coverslip using mounting medium

