



MAR 12, 2024

Blue Native PAGE of PINK1-TOM complex in mammalian cells

Hina Ojha¹, Miratul M. K. Muqit¹

¹Medical Research Council Protein Phosphorylation and Ubiquitylation Unit, School of Life Sciences, University of Dundee, Dow Street, Dundee DD1 5EH, UK

ASAP Collaborative Research Network



Francesca Tonelli
MRC-PPU at The University of Dundee

OPEN ACCESS



DOI:

dx.doi.org/10.17504/protocols.io.5qpvo3mezv4o/v1

Protocol Citation: Hina Ojha, Miratul M. K. Muqit 2024. Blue Native PAGE of PINK1-TOM complex in mammalian cells. **protocols.io** <https://dx.doi.org/10.17504/protocols.io.5qpvo3mezv4o/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](#), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

Protocol status: Working
We use this protocol and it's working

Created: Feb 19, 2024

ABSTRACT

Autosomal recessive mutations in PTEN-induced kinase 1 (PINK1) are associated with early-onset Parkinson's disease (PD) [1]. Upon inducing mitochondrial depolarization using uncouplers, the full-length human PINK1 stabilizes at mitochondria with the TOM (translocase of the outer mitochondrial membrane) complex, undergoing autophosphorylation and activation [2]. PINK1 plays a crucial role in phosphorylating ubiquitin and the E3 ubiquitin ligase, Parkin, initiating mitophagy for selectively removing damaged mitochondria [3]. To investigate PINK1 mutation impacts on PINK1-TOM complex formation, Blue Native Polyacrylamide Gel Electrophoresis (BN-PAGE) was done on stable cell lines in PINK1 knockout (KO) HeLa cells expressing doxycycline-induced PINK1/mutants. BN-PAGE substitutes SDS with Coomassie G250, preserving proteins in their native conformation, and utilizes non-ionic detergents for solubilization, ensuring complex integrity during analysis [4, 5]. Studying the PINK1-TOM complex has proved vital for our understanding of PINK1 stabilization on mitochondria. Here, we share our method for resolving protein complexes by blue native PAGE for PINK1-TOM complexes.

ATTACHMENTS

[1020-2633.pdf](#)

Last Modified: Mar 12, 2024

MATERIALS

PROTOCOL integer ID: 95588

Cells

Keywords: ASAPCRN

Funders Acknowledgement:

Aligning Science Across

Parkinson's





Grant ID: ASAP-000463

- Doxycycline induced WT-PINK1-3FLAG in PINK1 KO Hela Flip-in cells (Plasmid order no. DU43407)
- Doxycycline induced KI-PINK1-3FLAG in PINK1 KO Hela Flip-in cells ((Plasmid order no. DU46669)
- Doxycycline induced empty-3FLAG in PINK1 KO Hela Flip-in cells ((Plasmid order no. DU45919)
- Doxycycline induced L532A L539A L540A-PINK1-3FLAG in PINK1 KO Hela Flipin cells ((Plasmid order no. DU77629)
- Doxycycline induced R83E R88E R98E-PINK1-3FLAG in PINK1 KO Hela Flip-in cells ((Plasmid order no. DU77573)

Consumables

- Phosphate buffered saline (Invitrogen)
- 25G 1" (25mm) syringe needle (Orange)
- BD Plastipak Syringes 1ml x 120
- Benzonase® Nuclease, Merck, Purity >90%
- NativePAGETM 3 to 12% gels (Invitrogen/ThermoFisher)
- NativeMARKTM Unstained protein standard (Invitrogen)
- NativePAGETM 4x Sample Buffer (Invitrogen)
- NativePAGETM 20x Running buffer (Invitrogen)
- NativePAGETM 20x Cathode Buffer additive (Invitrogen)
- NativePAGETM 5% G-250 Sample additive
- 5% Digitonin (Invitrogen)
- 20x NuPAGE transfer buffer (Invitrogen)
- Pierce BCA protein assay kit (ThermoScientific)
- Clarity Western ECL subustrate (biorad)

Antibodies

-  PINK1 (D8G3) Rabbit mAb **Cell Signaling Technology Catalog #6946**
-  Recombinant Anti-TOMM40 antibody [EPR6932(2)] **Abcam Catalog #ab185543**
-  TOMM70 Antibody **Aviva Systems Biology Catalog #OAAN01138**
-  Recombinant Anti-TOMM20 antibody [EPR15581-54] **Abcam Catalog #ab186735**

Buffer and reagents

Mitochondrial fractionation buffer: Frozen stock (final Concentration)

A	B
HEPES pH 7.5	20 mM
EDTA	3 mM
Sodium β -glycerophosphate	5 mM
Sodium fluoride	50 mM
Sodium pyrophosphate	5 mM
Sucrose	250 mM
Added fresh before use:	
1x protease inhibitor cocktail tablet (Roche)	
Sodium orthovanadate	1mM

BN-PAGE buffer

A	B
NativePAGE 4x sample buffer	Make 1x buffer
1% digitonin	As required
milliQ water	As required
Benzonase (1-2 U/ul along with 2mM MgCl ₂)	If smears are observed

Dark Buffer

A	B
1x BN-PAGE buffer (1x)	200 ml
20x Cathode buffer additive	10 ml








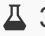
Light Buffer


A	B
1x BN-PAGE buffer (0.1x)	200 ml
20x Cathode buffer additive	1 ml

Equipment

- 150mm petri dishes for culturing cells
- VWR Micro Star 21R microcentrifuge for speed up to 17,000 g
- Eppendorf 5430R microcentrifuge for speed up to 20,000 g
- XCell SureLock Mini-Cell gel apparatus (Invitrogen)

Formation of mitochondrial enriched pellet

- 1 For making mitochondrial enriched fraction first harvest the cells after giving desired treatments (doxycycline for PINK1 induction and Antimycin/Oligomycin for mitochondrial damage).
- 2 For collection keep plates with cells  On ice covered with aluminium foil to provide even cool surface.
- 3 Wash the cells with PBS and collect the cells with cell scraper.
- 4 Collect the cells by centrifugation at  800 x g for  00:05:00 at  4 °C . 5m
- 5 Pellet down the cells at  800 x g for  00:05:00 at  4 °C . For 150 mm plate cell pellet, add  300 µL of mitochondria fractionation buffer. 5m

6 Disrupt cell membranes using a 25-gauge needle by passing through it for 25 times  On ice .

7 Clarify lysates by centrifugation at  800 x g for  00:10:00 at  4 °C .

10m



8 Discard the cytoplasmic membrane/nucleus/debris pellet.

9 Isolate supernatant and centrifuge at  17000 x g for  00:20:00 at  4 °C to collect mitochondrial enriched fraction.


20m



10 Try to remove as much buffer as possible from the pellets as it may interfere with the BN-PAGE.

11 Snap-freeze the mitochondrial enriched pellet for Blue native PAGE.



Mitochondrial enriched pellet Solubilization

12 Start with mitochondrial enriched pellets, snap-frozen and stored at  -80 °C . A single freeze-thaw cycle has negligible effects on TOMs/PINK1. Whole cells can also be used but the signal is clearer in mitochondrial enriched fraction.

Safety information



Do not freeze-thaw more than once. The lower PINK1 band near 480 kDa is more susceptible to freeze-thaw cycles.




13 Prepare Native-PAGE sample buffer per sample as described in buffer section.

14 Thaw mitochondrial pellets  On ice and very gently resuspend them in BNPAGE buffer  On ice with 1% digitonin.






Safety information

If you suspect traces of buffer in mitochondrial pellet, first make 1X buffer without digitonin and gently resuspend the pellets. Centrifuge samples at  17000 x g for  00:20:00 and proceed normally.

15 Store  On ice and incubate for  00:30:00, gently pipetting up and down 10-15 times every  00:10:00. Digitonin is chosen after optimizing with different reagents.

40m



16 Post-solubilization, centrifuge samples at  20000 x g for  00:30:00 at  4 °C. Transfer the supernatant to a cold Eppendorf, noting a small red pellet of insoluble material. If desired, solubilize this material in 1X SDS-PAGE buffer for further analysis.

30m





17 Quantify the lysate using BCA Protein Assay Kits.

Note


Do not use Bradford protein assay reagent.

- 18** Make samples of concentration around 1-2 ug/ul, typically loading 15-20 µg of solubilized mitochondrial pellets for one lane of a BN-PAGE gel. For diluting the lysates, make up the samples using 1x BN-PAGE lysis buffer with 1% digitonin. It is important to maintain the detergent percentage constant for proper Coomassie G250 and digitonin ratio.

- 19** Prepare gels in 1X Native PAGE running buffer (20x stock, make up  1 L in Millipore water). Also, prepare light and dark inner chamber buffers. Keep all buffers at  4 °C .

- 20** Add Coomassie G250 additive to samples (at 1/4 the detergent concentration, i.e., 0.25% if using 1% digitonin). Pipette thoroughly to distribute Coomassie in the sample; the sample will turn a deep blue.

Safety information

Do not heat the samples. Keep the samples  On ice the whole time.


- 21** Set up the gel apparatus in a cold room. Pour a small amount of Dark inner buffer into the inner chamber to check for leaks.

- 22** Remove any buffer from wells using an aspirator.




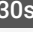
Note

Optional: You can give a wash to wells by filling them with clear running buffer and aspirate them.

- 23** Dry load samples onto the gel and use Native-Page ladder.


- 24 After loading, fill the inner chamber with  200 mL Dark inner buffer and 1x running buffer in the outer chamber.
- 25 Run the gel at 7mA constant current and allow the sample to enter the gel until it runs 1/3rd of the total gel. Subsequently, empty both inner and outer buffers, replacing the inner buffer with LIGHT buffer and filling the outer chamber with 1x running buffer again.
- 26 Continue running the gel at 7mA for 2-3 hours until the full gel is stained.

Transfer


- 27 Prepare 1x NuPAGE transfer buffer (20x stock) in a quantity sufficient for one transfer tank. Set up the transfer tank as usual.
- 28 For BN-PAGE, opt for PVDF membranes, as nitrocellulose binds Coomassie irreversibly. Activate the PVDF  30s membrane in 100% methanol for  00:00:30 before placing it over the gel in the transfer tank.
- 29 Upon completion, place the transfer tank in an ice bucket filled with ice. Run  Overnight at ~60mA,  30s limiting voltage to ~15 volts.

Staining and Probing



- 30 Disassemble the transfer tank. If successful, most or all of the Coomassie staining on the gel will transfer to the PVDF membrane.

31 Destain Coomassie with 100% methanol 2-3 times for  00:05:00 each. A slight Coomassie staining may persist (from the added samples; background staining will be removed). **5m**


32 Restain the gel with Ponceau to visualize protein ladders. Mark ladders and scan the blot.

33 Wash Ponceau with TBST (0.1% tween) until fully destained (3 times for  00:10:00 each). It might take longer to destain completely. If needed, incubate in 5% milk in TBST to destain. **10m**





34 Block in fresh 5% milk for  01:00:00 (do not use BSA for blocking) and probe  Overnight as for a standard Western blot. **2h**



35 Wash membranes three times in TBST for  00:10:00 each to remove residual primary antibody. The primary antibody may give different signals depending on how the antibody binding site is exposed in a complex. For visualization of PINK1 in PINK1-TOM complex, CST-PINK1 antibody works best. **10m**




36 Utilize HRP-conjugated secondary antibodies. Prepare HRP-conjugated AB (1:5000 in milk), incubate at  Room temperature covered for  01:00:00 . **1h**



37 Remove the secondary antibody and wash three times in TBST for  00:10:00 each and once in TBS. **10m**



38

Use ECL to visualize HRP-conjugated antibodies. Prepare the ECL reagent mix (1:1) just before scanning. Pour off TBS on one membrane and add  1 mL ECL mix; immediately scan, adjusting scanning times based on the primary antibody and protein abundance.