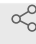




Sep 06, 2022

Culture and transfection of HEK293T cells

Hanqin Li¹, Dirk Hockemeyer¹¹University of California, Berkeley1 *Works for me* Sharedx.doi.org/10.17504/protocols.io.eq2lynkzpvx9/v1

Devin E Snyder

ABSTRACT

This protocol describes a standard procedure culturing and transfecting HEK293T cells

Protocol overview:

A. Culturing HEK293T cells

B. Transfection of HEK293T cells with Lipofectamine 2000

DOI

dx.doi.org/10.17504/protocols.io.eq2lynkzpvx9/v1

PROTOCOL CITATION

Hanqin Li, Dirk Hockemeyer 2022. Culture and transfection of HEK293T cells.
protocols.io
<https://protocols.io/view/culture-and-transfection-of-hek293t-cells-cd3bs8in>



FUNDERS ACKNOWLEDGEMENT

Aligning Science Across Parkinson's
Grant ID: ASAP-000486

LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

CREATED

Jul 23, 2022

LAST MODIFIED

Sep 06, 2022

A	B	C
Item	Vendor	Catalog #
DMEM, high glucose	Thermo Fisher	11965118
DPBS w/o Calcium and magnesium	Corning	MT21031CV
Fetal Bovine Serum (FBS)	Corning	35-011-CV
L-Glutamine	Sigma	G8540
MEM Non-Essential Amino Acids (100X)	Thermo Fisher	11140050
0.25% Trypsin with EDTA	Thermo Fisher	25200114
Opti-MEM	Thermo Fisher	31985062
Sodium Pyruvate 100 mM	Thermo Fisher	11360070
Lipofectamine 2000	Thermo Fisher	11668019

A. Culturing HEK293T cells

1m



- 1 HEK293T cells are cultured in HEK293T medium in 10 cm dishes.

1.1 HEK293T Medium

A	B
DMEM, high glucose	385 ml
Fetal Bovine Serum (FBS)	50 ml
L-Glutamine (100X)	5 ml
MEM Non-Essential Amino Acids (100X)	5 ml
Sodium Pyruvate 100 mM	5 ml

Final volume: 500 mL

- 2 Passage cells when the culture reaches 80% confluency.

- 3 Remove medium
- 4 Wash once with 5 ml DPBS
- 5 Add 1 ml 0.25% Trypsin with EDTA, tilt and shake the dish so that the Trypsin covers the entire dish.
- 6 Incubate at  **37 °C** for  **00:01:00** 1m
- 7 Add 5 ml fresh HEK293T medium to inactivate Trypsin
- 8 Pipet 10 times to dissociate the cells and mix
- 9 Transfer 1 ml of cell suspension to a new dish pre-added with 9 ml HEK293T medium. Shake to mix well. This is a 1:6 splitting.
- 10 HEK293T cells usually needs to be passaged every 2 days.

B. Transfection of HEK293T cells with lipofectamine 2000

45m

- 11 One day before transfection, dissociate HEK293T cells with Trypsin as described above
- 12 Seed 250,000 cells/1 well of 12-well plate

- 13 Change to fresh medium 1 h before transfection
- 14 Label two micro centrifuge tubes as I and II
- 15 In micro centrifuge tube I, add 125 µl Opti-MEM and 6 µl Lipofectamine 2000. Mix by gently^{5m} pipetting 3 times. Incubate at 🌡 **Room temperature** for ⌚ **00:05:00** .
- 16 While waiting, in micro centrifuge tube II, add 125 µl Opti-MEM and 250 ng total plasmid. Mix by pipetting.
- 17 Mix tube I and II by gently pipetting 3 times. Incubate at 🌡 **Room temperature** for 20m ⌚ **00:20:00** .
- 18 Mix one time by gently pipetting. Transfer all transfection reagents into one well of a 12-well plate, drop-wise.
- 19 After adding the reagent to all wells, shake the plate to mix.
- 20 Culture in 🌡 **37 °C** incubator ⌚ **Overnight** . 20m
- 21 Change to fresh medium and culture for another 2 days. Collect samples or passage once if longer culturing is needed.