



VERSION 1
OCT 29, 2023

OPEN ACCESS



DOI:
dx.doi.org/10.17504/protocols.io.n92ldmr7xl5b/v1

Protocol Citation: Erik.Hausner, micboe 2023. sgRNA library re-amplification in liquid culture. **protocols.io**
<https://dx.doi.org/10.17504/protocols.io.n92ldmr7xl5b/v1>

License: This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited


Protocol status: Working
We use this protocol and it's working

Created: Sep 15, 2023

🌐 sgRNA library re-amplification in liquid culture V.1

Erik.Hausner¹, micboe¹

¹Medical Faculty of the Martin Luther University Halle-Wittenberg

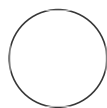

Erik.Hausner

ABSTRACT

In this protocol, we describe a stepwise procedure for the re-amplification of sgRNA libraries in liquid culture. In our hands, this protocol works reliably to amplify pre-cloned sgRNA libraries (e.g. order from Addgene) in a way that preserves the distribution of library elements.

MATERIALS

- 100 µL
- ElectroMAX[®]; Stbl4[®]; Competent Cells Thermo Fisher Catalog #11635018
- Up to 400 ng of Sample
- Electroporator and Electroporation cuvettes Biozym Catalog #748010

Last Modified: Oct 29, 2023

PROTOCOL MATERIALS

PROTOCOL integer ID:
87835

Keywords: sgRNA,
reamplification, Plasmid pool,
Library, Electroporation,
amplification

✕ 1 kb Plus DNA-Ladder Thermo Fisher
Scientific Catalog #10787018 In [2 steps](#)

✕ LB agar plates with the proper antibiotic (e.g. Kanamycin) Contributed by
users

In [3 steps](#)

✕ 1x TBE buffer Contributed by
users Step 13

✕ ElectroMAX[®] Stbl4[®] Competent Cells Thermo
Fisher Catalog #11635018

In Materials and [3 steps](#)

✕ NucleoBond Xtra Midi kit for transfection-grade plasmid DNA Macherey-
Nagel Catalog #REF 740410.50

Step 9

✕ GelRed[™] Nucleic Acid Gel Stain, 10,000X in Water Gold
Biotechnology Catalog #G-725

Step 14

✕ TriTrack DNA Loading Dye (6X) Thermo
Fisher Catalog #R1161 Step 16

✕ Agarose Low Melt Carl
Roth Catalog #6351.4 Step 13

✕ Electroporation
cuvettes Biozym Catalog #748010 In Materials and [2 steps](#)

✕ Liquid LB medium Contributed by
users In [4 steps](#)

✕ SOC Outgrowth Medium - 100 ml New England
Biolabs Catalog #B9020S

In [3 steps](#)

✕ 1.5 mL Eppendorf tubes Contributed by
users Step 5.1

✕ 1% Agarose gel Contributed by users Catalog #/ In [2 steps](#)


BEFORE START INSTRUCTIONS


Keep your original stock safe and aliquoted. For large plasmids with complementary sequences such as LTR sites, keep in mind that repeated reamplification from an already reamplified stock may lead to an accumulation of recombined plasmids and a poorer distribution of library elements.


Library transformation

25m


1 Prepare Sample 5m


 ElectroMAX[®]; Stbl4[®]; Competent Cells Thermo
Fisher Catalog #11635018

 Electroporation
cuvettes Biozym Catalog #748010 and

 SOC Outgrowth Medium - 100 ml New England
Biolabs Catalog #B9020S for
electroporation.

1.1 Thaw 5m

 ElectroMAX[®]; Stbl4[®]; Competent Cells Thermo
Fisher Catalog #11635018


 00:05:00 on ice.

1.2 Pre-cool Electroporation cuvettes Biozym Catalog #748010 by placing it on ice.

1.3 Pre-warm SOC Outgrowth Medium - 100 ml New England Biolabs Catalog #B9020S at 37 °C .

2 Add 100 ng Sample into 25 µL 10m




 ElectroMAX[®]; Stbl4[®]; Competent Cells Thermo
Fisher Catalog #11635018 ,

carefully mix by pipetting up and down.

3 Add 25 µL of the plasmid/cell mix into a cuvette, electroporate at 1.2 kV, 25 uF and 200 ohm or alternative setting (see note below). Directly after electroporation, add 1 mL of pre- 10m



warmed  SOC Outgrowth Medium - 100 ml New England
Biolabs Catalog #B9020S .

Note



The electroporator setting may vary from model to model and should be checked along with the test plasmids provided in the kit of the STBL4 cells.

Safety information

Make sure that any water or ice residue is removed from the cuvettes before inserting them into the electroporator to avoid arcing.

Library recovery

1h 15m

- 4 After electroporation, add the  1 mL resuspended cells in a 14 ml culture tube and incubate  600 rpm, 37°C, 01:00:00 . 1h 15m



Note

In general, an incubation temperature of 37°C is optimal for cell recovery. Since *E. coli* tend to recombine plasmids with complementary sequences (e.g. LTRs), recovery temperature can be reduced to 30°C. This may however, result in a lower total number of recovered cells.

Determination of transformation efficiency

16h 30m














- 5 Use a small fraction of your cells to determine the electroporation efficiency of the reamplification.

Note

In this step, much depends on the size of the particular plasmid and the number of elements in the library. Therefore, the dilution factor must be chosen based on properties of the library and the scale of the electroporation. Smaller plasmids yield significantly more colonies than large ones, and an upscaled plasmid input at the electroporation step may result in higher dilutions being required to achieve a countable number on the respective agar plates after plating.







5.1 For 1:10,000 dilution:

5m

Prepare  1.5 mL Eppendorf tubes Contributed by  users. Take  10 μL of recovery culture and dilute in  990 μL of  Liquid LB medium Contributed by  users (1:100 dilution). Take  100 μL of 1:100 dilution and dilute in  900 μL of  Liquid LB medium Contributed by  users (1:1,000) and plate  100 μL on  LB agar plates with the proper antibiotic (e.g. Kanamycin) Contributed by  users (1:10,000 dilution).



5.2 For 1:1,000,000 dilution:

5m

Take  10 μL of the 1:1,000 dilution and dilute in  990 μL of  Liquid LB medium Contributed by  users and plate 100 μL on a pre-warmed  LB agar plates with the proper antibiotic (e.g. Kanamycin) Contributed by  users (1:1,000,000 dilution).

Note

When preparing the dilution series, always mix stock solutions well by flicking the tube before diluting, to resuspend sedimented cells. Distribute the plated cells evenly over the plate by e.g. using glass beads.

6 Place the plates in an incubator at  37 °C  Overnight .

16h



Library extraction and quality control

16h

7 Use rest of recovery to inoculate up to 500 mL of 16h



Liquid LB medium Contributed by users

with an added selection marker specific antibiotic



like ampicillin in an Erlenmeyer flask for Overnight culture. 600 rpm, 30°C

Determination of transformation efficiency

15m

8 On the next day, check for overall coverage via colony counting on

15m



LB agar plates with the proper antibiotic (e.g. Kanamycin) Contributed by users

. The overall



colony count should be 1000x the element number of your library.

Note

Below we provide a simplified example for how to determine transformation coverage.

Example calculation of coverage: On the 1:10,000 dilution plate we count 100 colonies. This gives us $100 \times 10,000 = 1,000,000$ total colonies. This total colony number is divided by the number of elements (e.g. sgRNAs) in the respective library. For a library the size of 1,000 sgRNAs, the coverage would $1,000,000 / 1,000 = 1,000x$. For larger libraries, e.g. the size of the genome-wide Brunello library (80,000 sgRNAs), we would count the 1:1,000,000 dilution plate. In this case, 80 counted colonies would mean 80,000,000 total colonies which divided by the library size (80,000 sgRNAs) would again return a transformation coverage of 1,000x.

Library preparation and QC

1h

9 Follow the protocol instructions of the

25m

NucleoBond Xtra Midi kit for transfection-grade plasmid DNA Macherey-Nagel Catalog #REF 740410.50

for transfection-grade plasmid DNA for Midi Prep. Follow the protocol instructions of the for transfection-grade plasmid DNA for Midi Prep.

10 Determine your final Sample concentration via NanoDrop or Qubit measurement.



Equipment

| | |
|----------------------------------|----------------|
| | NAME |
| Qubit 2.0 Fluorometer instrument | BRAND |
| Q33226 | SKU |
| with Qubit RNA HS Assays | SPECIFICATIONS |


Equipment

| | |
|---|-------|
| NanoDrop™ One/OneC Microvolume UV-Vis Spectrophotometer | NAME |
| UV-Vis Spectrophotometer | TYPE |
| Thermo Scientific | BRAND |
| ND-ONE-W | SKU |

11

Send a sample of your reamplified  Sample for Sanger sequencing.

Below we show an example chromatogram of an expected sequencing result. We recommend using sequencing primers 50-100 nt upstream of the sgRNA region. You should see clean traces up- and downstream of the SPACER region, and a noisy 20 nt signal in the SPACER region, due to the sgRNA diversity in your library.

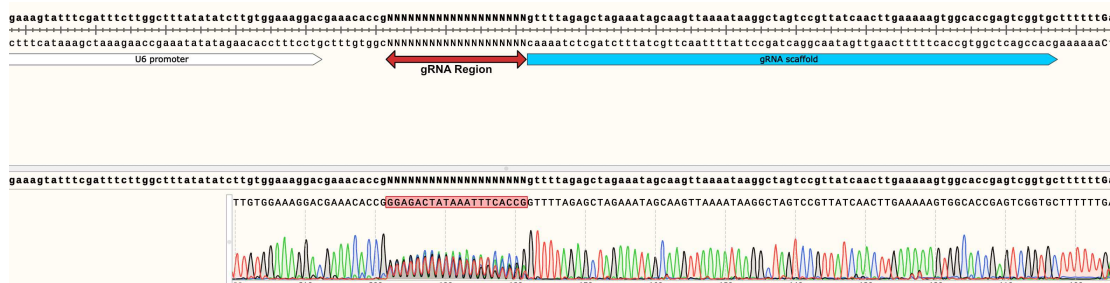



Example of a Sanger sequencing result.

Note

! To validate the distribution of elements in your library, we strongly recommend performing next generation sequencing of your plasmid pool before proceeding with downstream experiments. To do so, follow the NGS protocol provided with your library, using the plasmid pool as template, instead of the genomic DNA (as you would in CRISPR screens).

QC: Plasmid recombination check


3h 15m

12 Since sgRNA library plasmids can recombine during E.coli re-amplification, it is recommended to check for recombination via linearization of  200 ng of your reamplified  Sample via a restriction digest within the backbone of your library vector. 1h




13 Prepare a  1% Agarose gel Contributed by users Catalog #/ by melting  1 g of 20m



 Agarose Low Melt Carl
Roth Catalog #6351.4


in  100 mL

 1x TBE buffer Contributed by
users

14 Let the required amount for casting cool down till it is approximately  50 °C and add 30m



 1 µL of


 GelRed™ Nucleic Acid Gel Stain, 10,000X in Water Gold
Biotechnology Catalog #G-725 per

ml of melted  1% Agarose gel Contributed by users Catalog #/

15 Pour the warm, still liquid gel into an electroporation chamber and wait until it has cooled down. 15m

16 Mix your linearized reamplified  Sample with 5m



 TriTrack DNA Loading Dye (6X) Thermo
Fisher Catalog #R1161

and prepare



1 kb Plus DNA-Ladder Thermo Fisher
Scientific Catalog #10787018

for gel loading.

17


Add your linearized reamplified  Sample alongside with the prepared

1h



1 kb Plus DNA-Ladder Thermo Fisher
Scientific Catalog #10787018

onto the gel and run it

for  01:00:00 at 120 V.

18

Check the plasmid size on your gel using UV excitation.

5m



Expected result

In the best case, only one band will be present, corresponding to the size of the respective vector. In case of recombination, one or more additional bands will be present. For optimal downstream results, the band of the intact vector should be dominant.