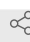


Version 1 ▼

Aug 21, 2022

# 🌐 Preparation of $\alpha$ -synuclein fibrils amplified from clinical material V.1

arpine.sokratian<sup>1</sup><sup>1</sup>Duke University*In Development* Share

This protocol is published without a DOI.

 arpine.sokratian

## ABSTRACT

This protocol describes the preparation of alpha-synuclein fibrils from clinical tissue with a quality control.

Optional conjugation step with phRodo STP ester dye is indicated

## PROTOCOL CITATION

arpine.sokratian 2022. Preparation of  $\alpha$ -synuclein fibrils amplified from clinical material.

**protocols.io**

<https://protocols.io/view/preparation-of-synuclein-fibrils-amplified-from-cl-b93gr8jw>



## LICENSE

————— This is an open access protocol distributed under the terms of the [Creative Commons Attribution License](https://creativecommons.org/licenses/by/4.0/), which permits unrestricted use, distribution, and reproduction in any medium, provided the original author and source are credited

## CREATED

May 26, 2022

## LAST MODIFIED

Aug 21, 2022

## PROTOCOL INTEGER ID

63304

## MATERIALS TEXT

[Protein LoBind Tubes, 1.5](#)

Eppendorf tube [mL Eppendorf Catalog #0030108116](#) Step 3

[PBS](#) Thermo Fisher

PBS [Scientific Catalog #28374](#)

[pHrodo™ iFL Red STP Ester \(amine-reactive\)](#) Thermo




Fisher Catalog #P36011

[DMSO](#) Sigma

Aldrich Catalog #472301

[Sodium bicarbonate](#) Sigma

Aldrich Catalog #S6014

- 1 Thaw down an aliquot of alpha-synuclein monomer stock (track down a batch number with identified EU number, concentration, A260/280 ratio) on ice
- 2 Spin down an aliquot of alpha-synuclein monomer stock solution (  **20000 x g** ,  **00:10:00**  **4 °C** ) <sup>10m</sup> and measure the concentration using nanodrop



Add  **3 µL** of 10x diluted aliquot in PBS onto nanodrop pedestal;


Parameters: other proteins; coefficient extinction: 5.98; MW: 14.4 kDa

Perform two measurements and confirm <10% standard error between two measurements

If necessary, prepare 20X and 30X dilutions to confirm findings.

NanoDrop™ One/OneC Microvolume UV-Vis Spectrophotometer  
UV-Vis Spectrophotometer  
Thermo Scientific ND-ONE-W

- 3 Calculate the reaction mix:  **5 mg/mL** concentration of alpha-synuclein in  **200 µL** in PBS (use safe-lock eppendorf tubes [mL Eppendorf Catalog #0030108116](#) or other [Thermo Scientific™ Low Protein Binding Collection Tubes \(1.5 mL\) Contributed by users Catalog #PI90411](#) )

- 3.1 Prepare a reaction mix containing 10% of clinical tissue and  **5 mg/mL** of alpha-synuclein in PBS

4



5d

Incubate the reaction at **1000 rpm, 37°C** for **120:00:00** using program settings: 1 min ON; 1 min OFF

Eppendorf Thermomixer C Model 5382  
Thermomixer C  
Eppendorf 5382000023

ThermoTop®  
Smart block  
Eppendorf 5308000003

5

Spin down the insoluble fraction at **15000 rpm, 10°C, 00:10:00** ; take **100 µL** of supernatant and <sup>10m</sup>  
add **1 mL** of fresh PBS to the pellet



Example of reaction mix after 5 days of incubation

6

Gently resuspend and spin down the fibrils at **15000 rpm, 10°C, 00:10:00** ; take **100 µL** of <sup>10m</sup>  
supernatant and add **1 mL** of fresh PBS to the pellet (repeat 3 times)

7

Step 7 includes a Step case.

### Conjugation of full-length fibrils with pHrodo dye

step case

#### Conjugation of full-length fibrils with pHrodo dye

1. Dissolve 100 ug of dye in 50 ul of sterile DMSO (2 mg/mL)
2. Dilute fibrils to 1mg/mL concentration in PBS containing 0.1M bicarbonate (total volume: 0.95 mL)
3. Add 50 ug of dissolved dye to each reaction tube
4. Incubate overnight in eppendorf tube (foil-wrapped) with continuous shaking
5. In the morning, spin down the dye-fibril solution at 10,000 g for 10 min at 10C
6. Transfer the supernatant and dissolve the pellet with 1 mL of PBS (repeat 5 times)

8 Take out 900 ul of PBS leaving 100 ul of fibril pellet in the tubes. Transfer the fibrils into 0.6 mL PCR tubes (thick wall).

9 Sonicate the fibrils using water bath sonicator at 10C, 30% amplitude and for 1 hour (no OFF ON cycles). Check the level of water in the water-bath is in a line with the tube content.

10 Measure size of sonicated particles using DLS



Dynamic Light Scattering measurements  
by **arpine.sokratian**

PREVIEW

RUN

